






Phantom QA results for / fbirn_out/ ses32ch_funcbold_acqcmrrstd_tasktsnr slice 34




bold_acqcmrrstd_tasktsnr_23.xml [percent fluct (trend removed), drift, driftfit] :


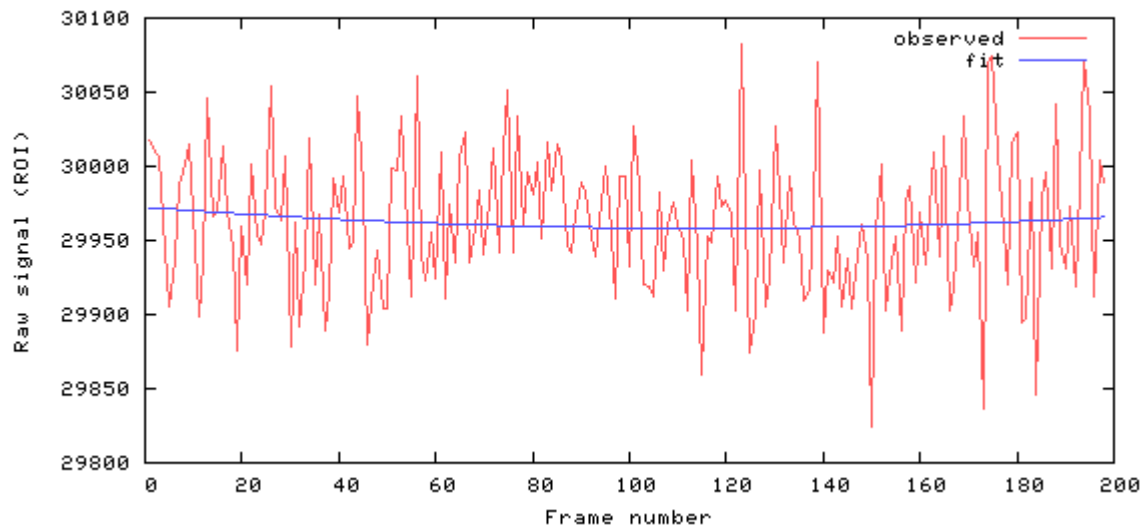


[mean, SNR, SFNR] = [29962.0 137.2 137.3]

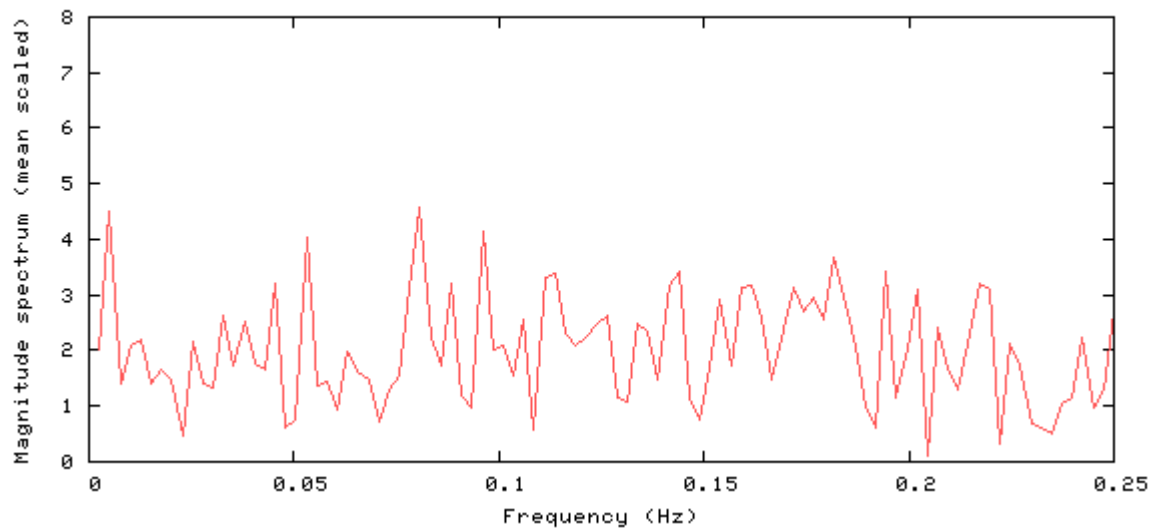

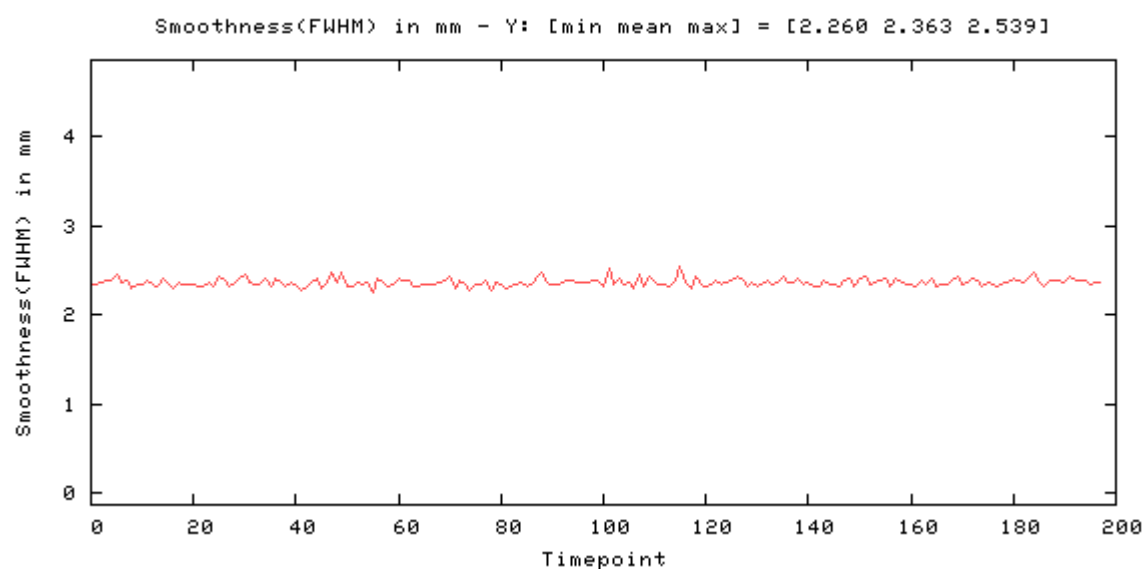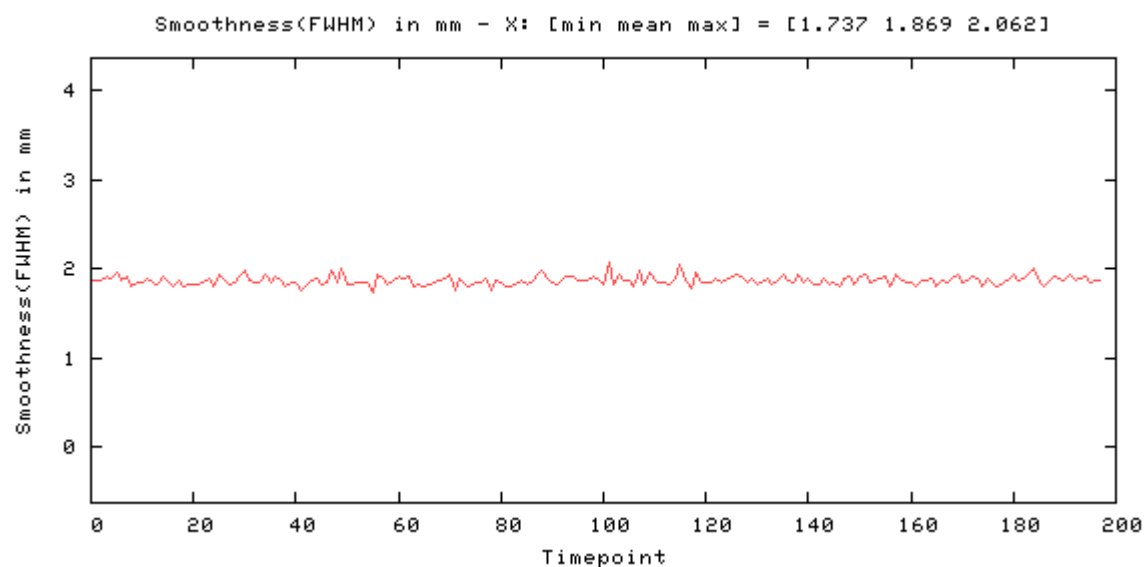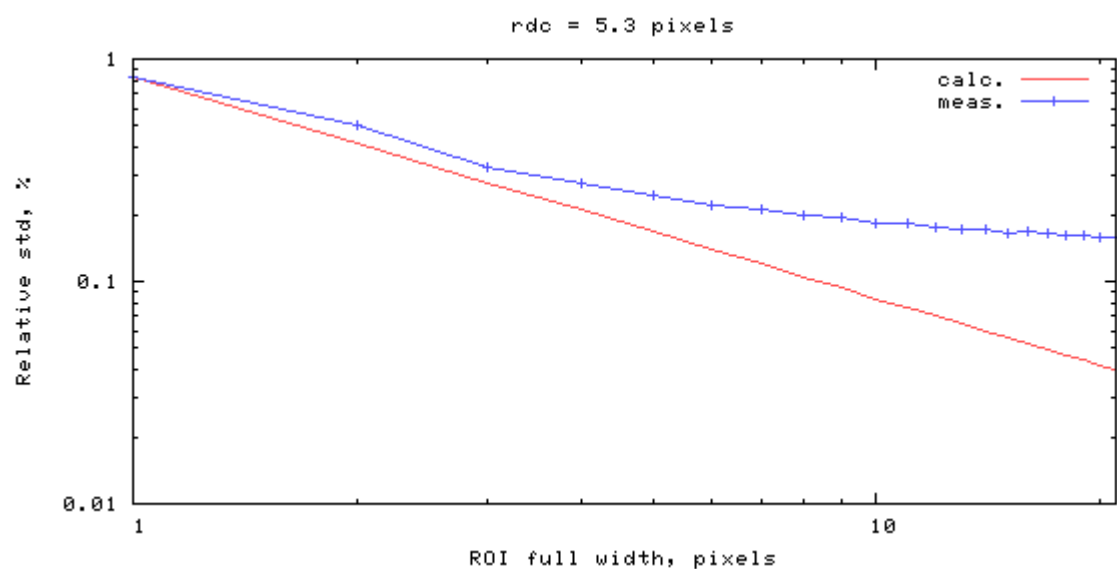

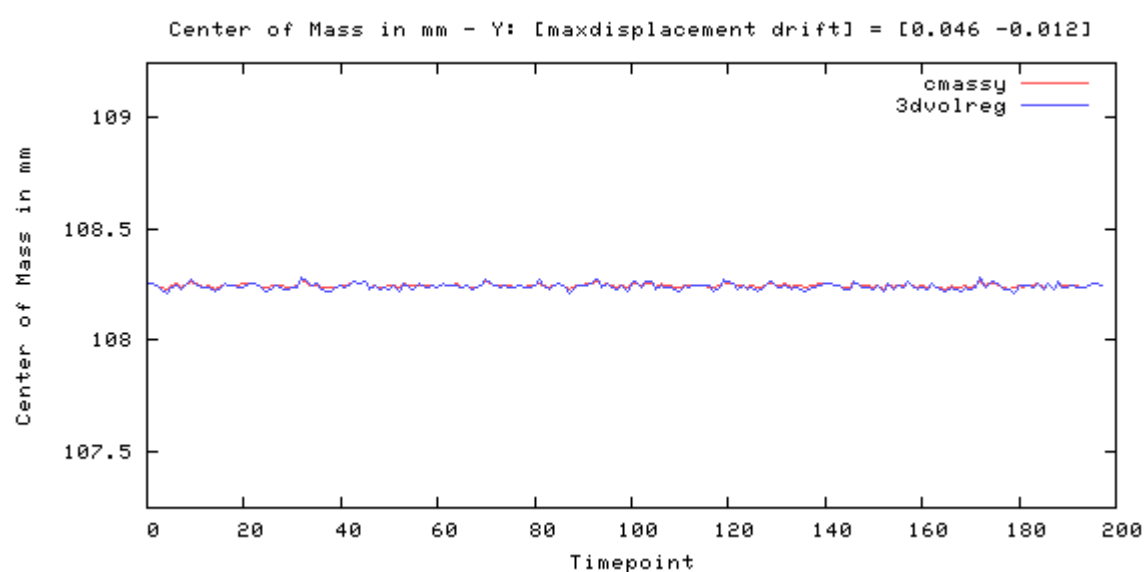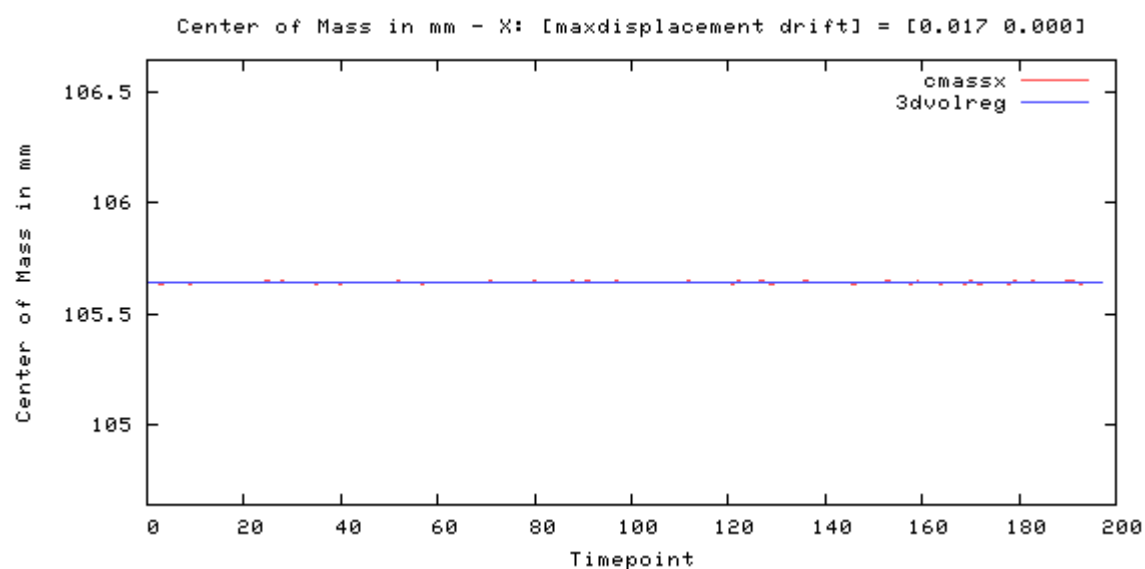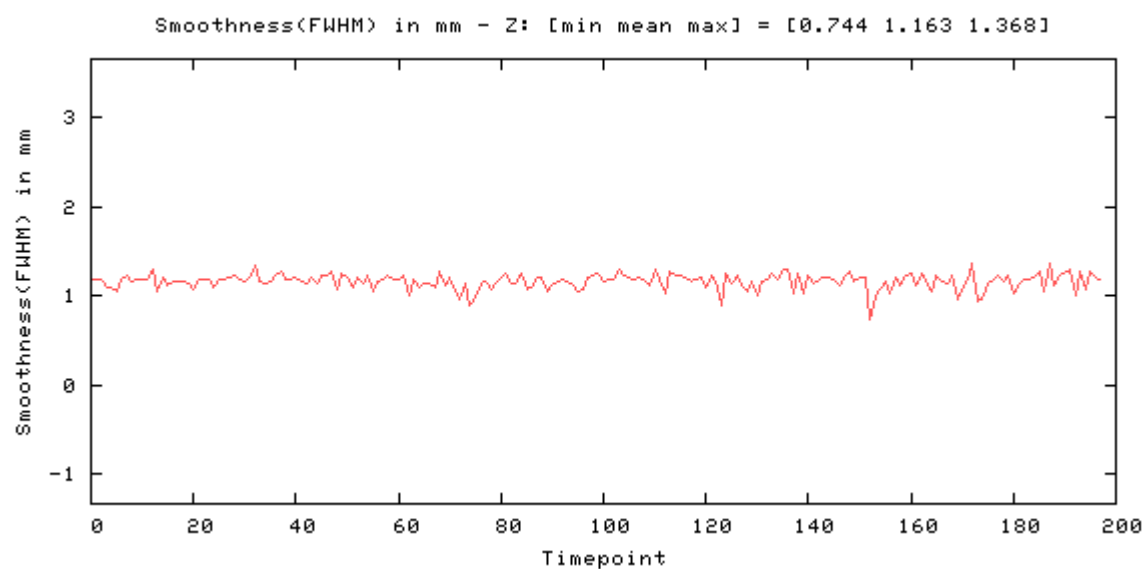

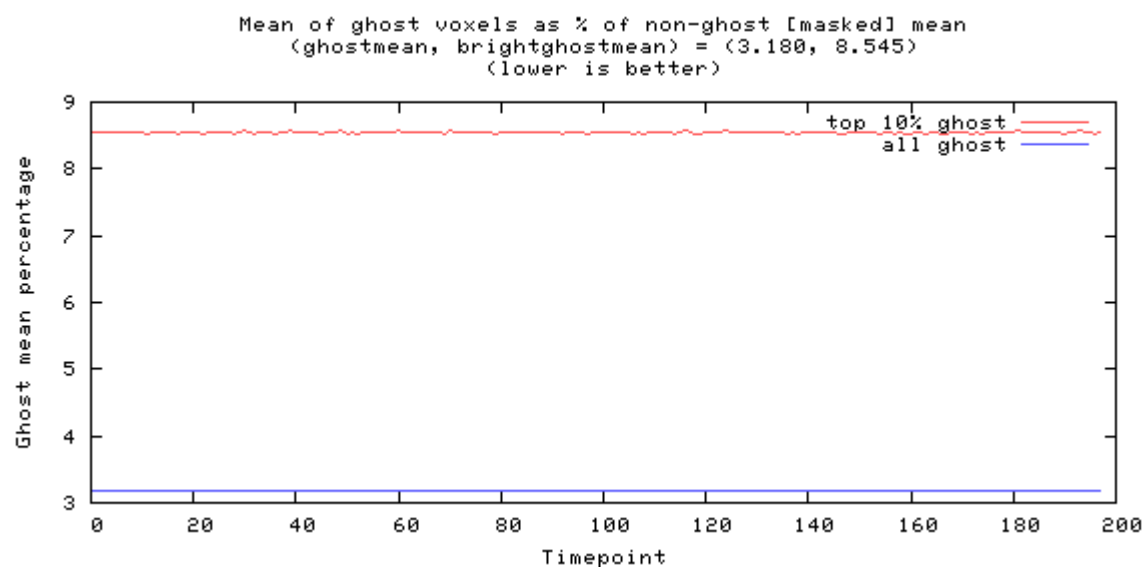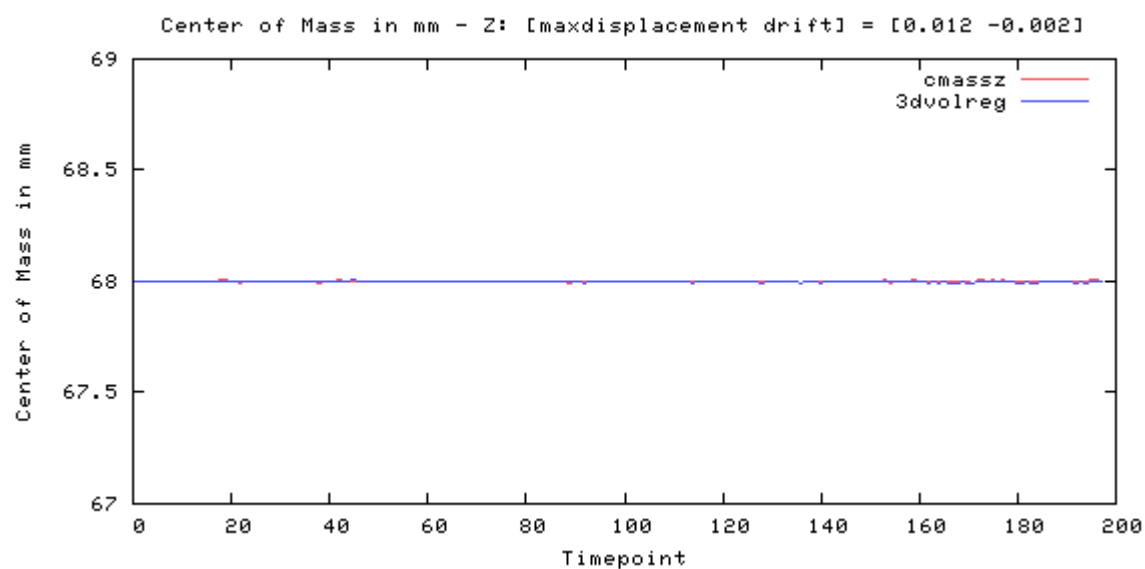


Odd-even difference image

-7966 11982




image min: -7966, image max: 11982

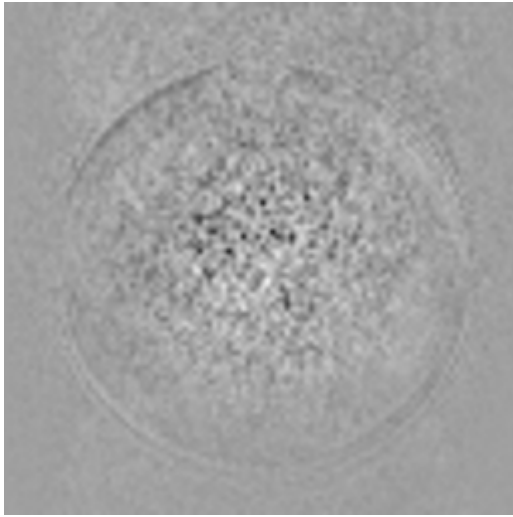

Mean image




image min: 0, image max: 41078.8

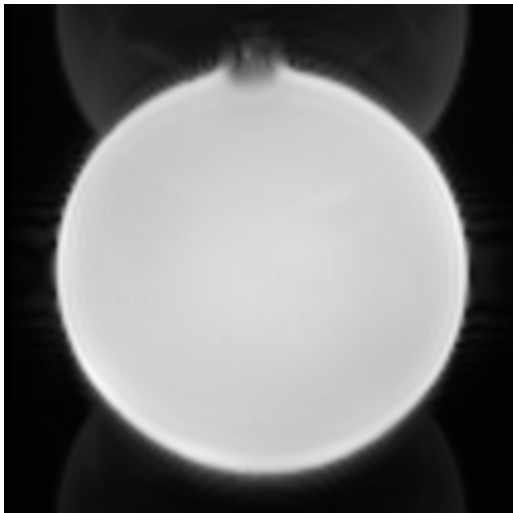

Standard Deviation image




image min: 0, image max: 297.55

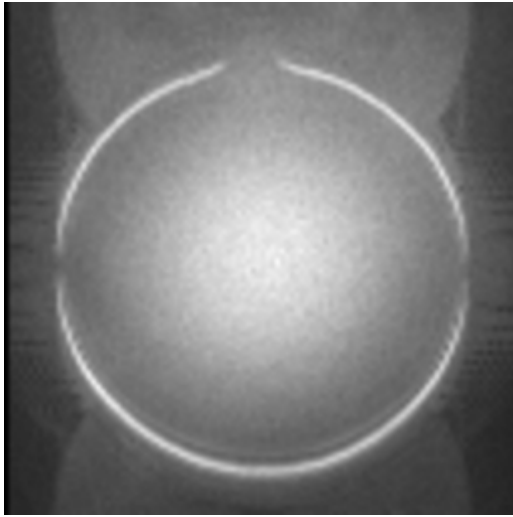

SFNR image



image min: 0, image max: 774.151

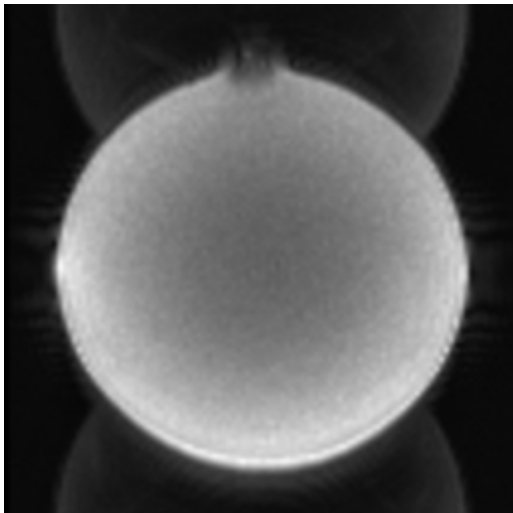

Acquisition parameters

dimensions	106x106x69x200
spacing	2mmx2mmx2mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	1300d8f2-96d2-4f
studyid	Stability_20260504
seriesnumber	23
runnumber	1
magneticfield	3

description	ses-32ch_func-bold_acq-cmrrstd_task-tsnr
scandate	2026-05-04
scantime	10:37:37
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
receivecoilname	Head_32
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-32ch_func-bold_acq-cmrrstd_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260504103735.100000
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	2
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	RESEARCH
acquisitionduration	417.89499999999998
resonantnucleus	1H

applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
imagecomments	Not for diagnostic use, Unaliased MB3/ PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	69
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00