






Phantom QA results for /fbirn_out/ses-20ch_func-bold_acq-cmrrstd_task-tsnrxacpc.xml, slice 34




old_acq-cmrrstd_task-tsnrxacpc.xml [percent fluct (trend removed), drift, driftfi-


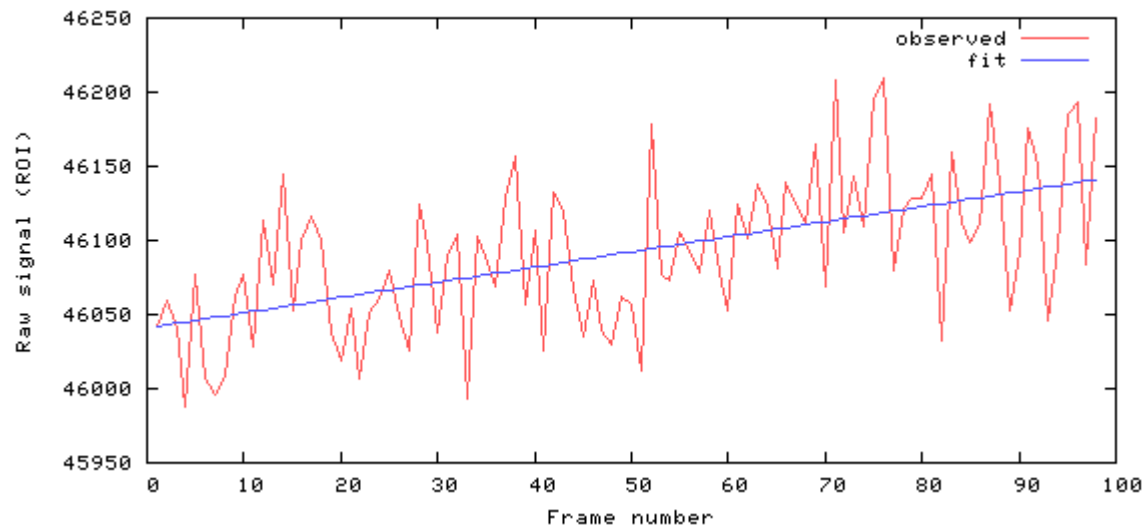


[mean, SNR, SFNR] = [46091.7 131.4 133.1]

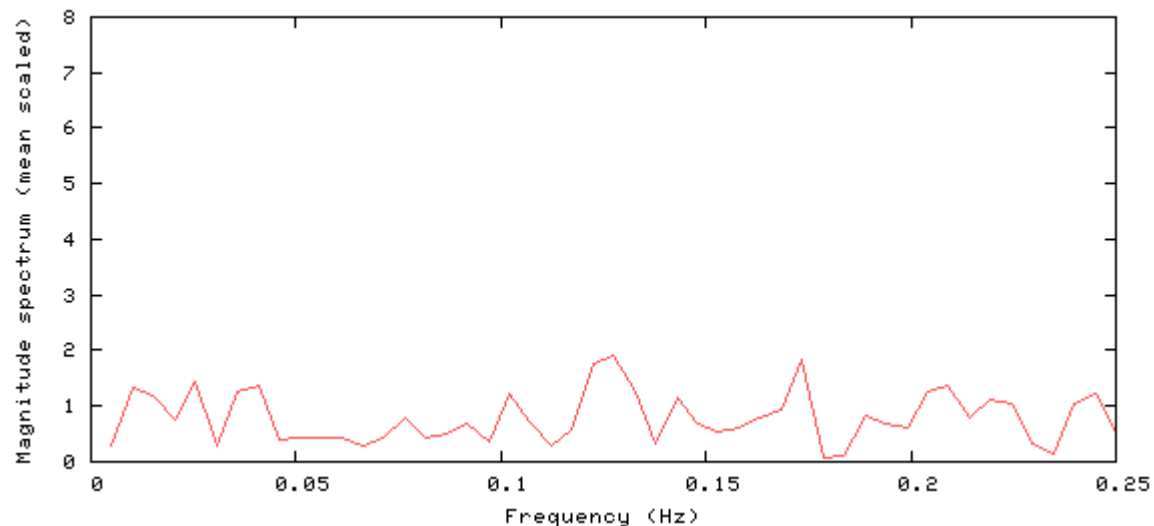

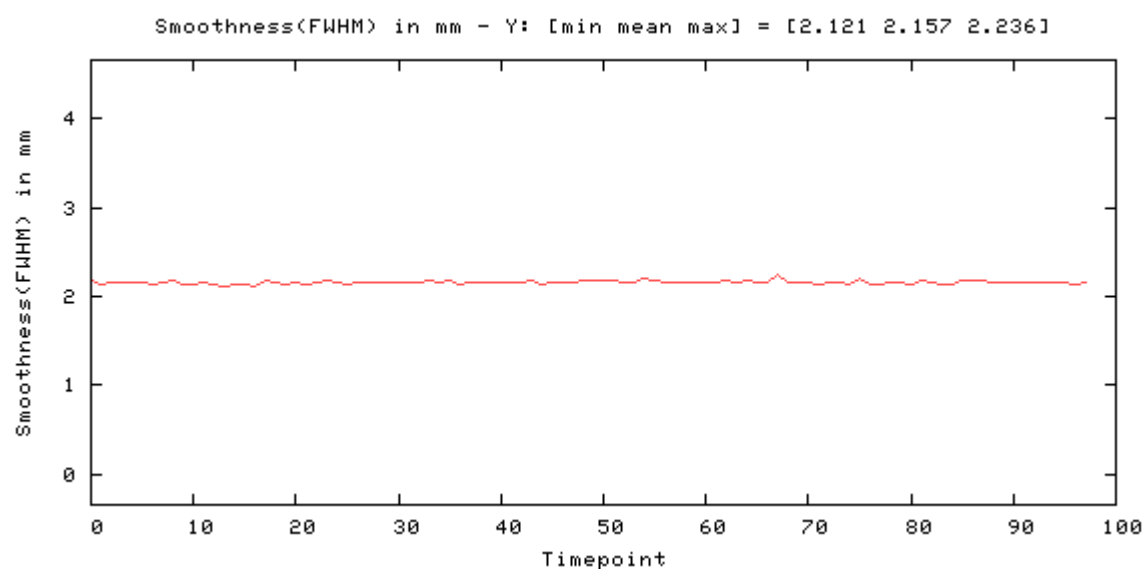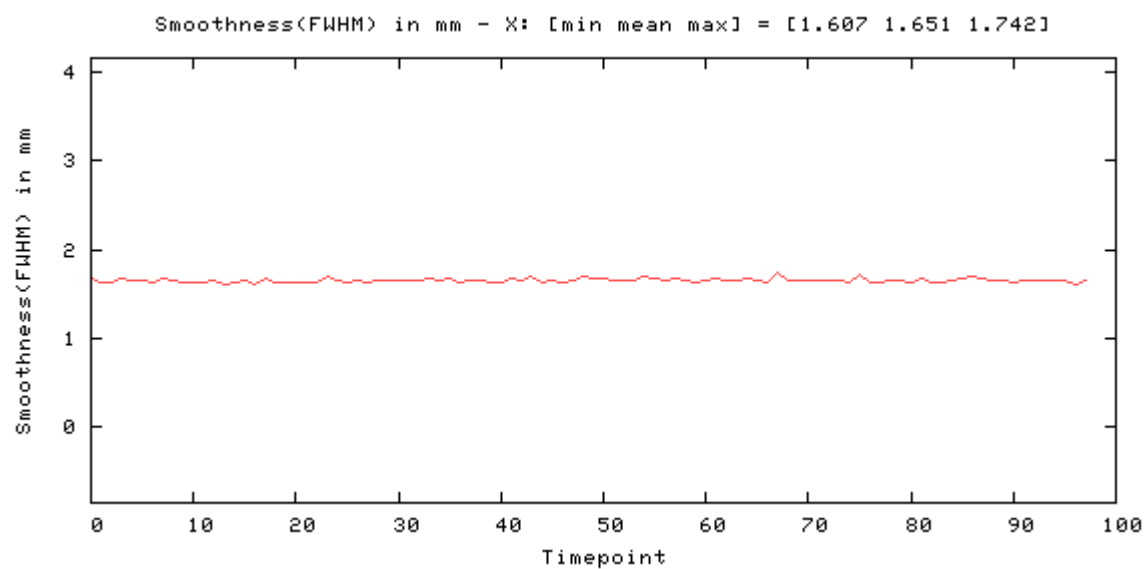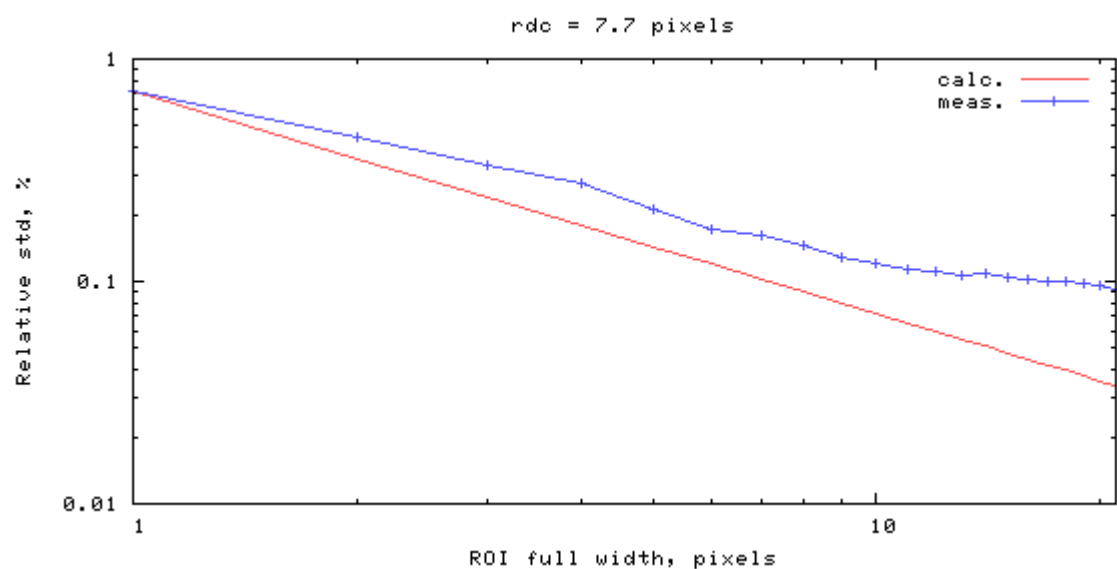

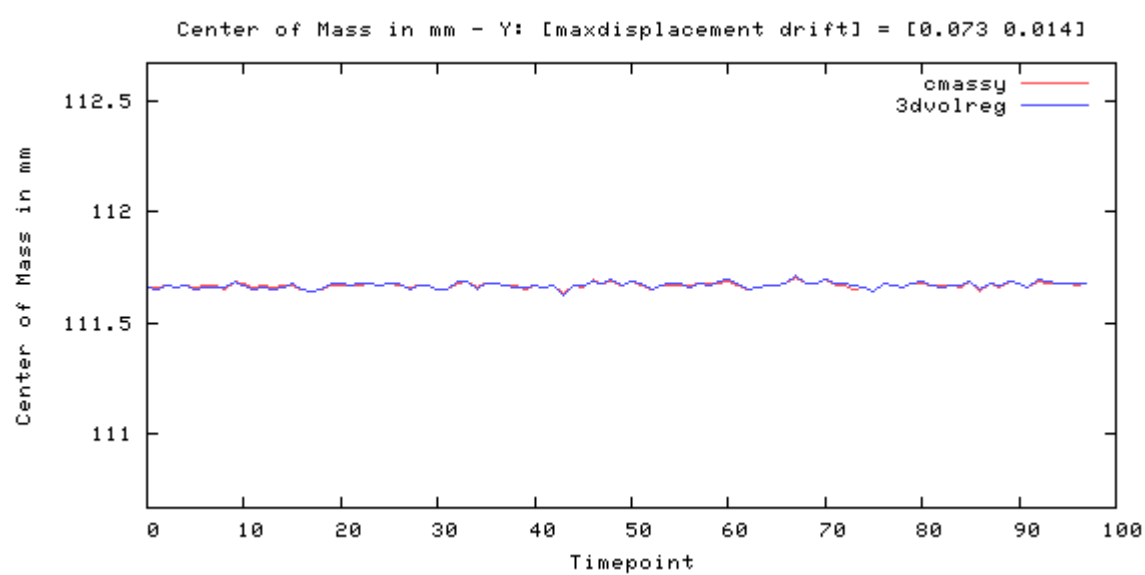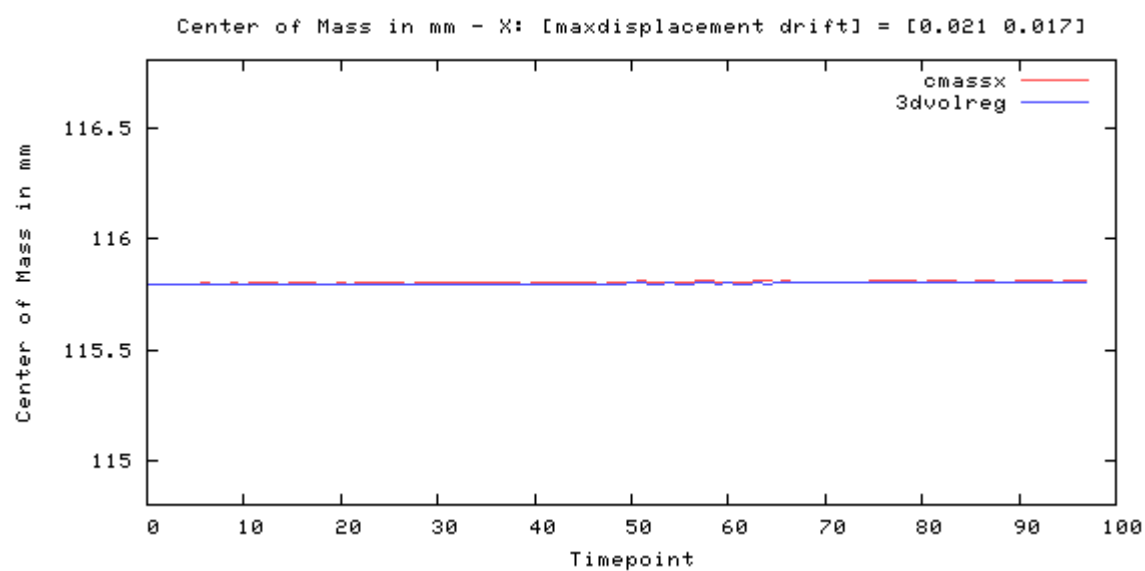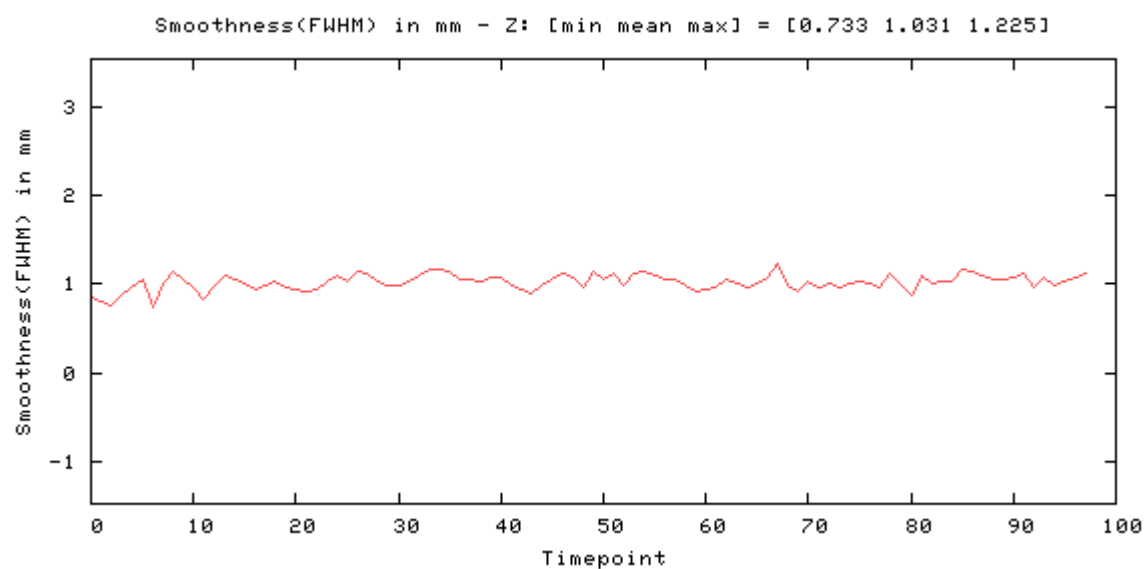

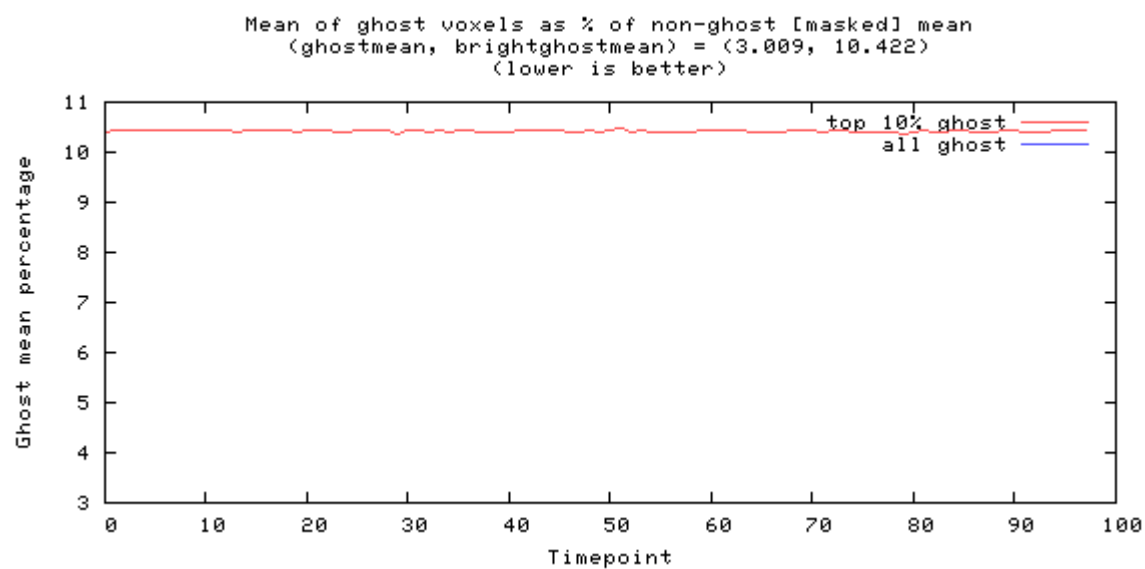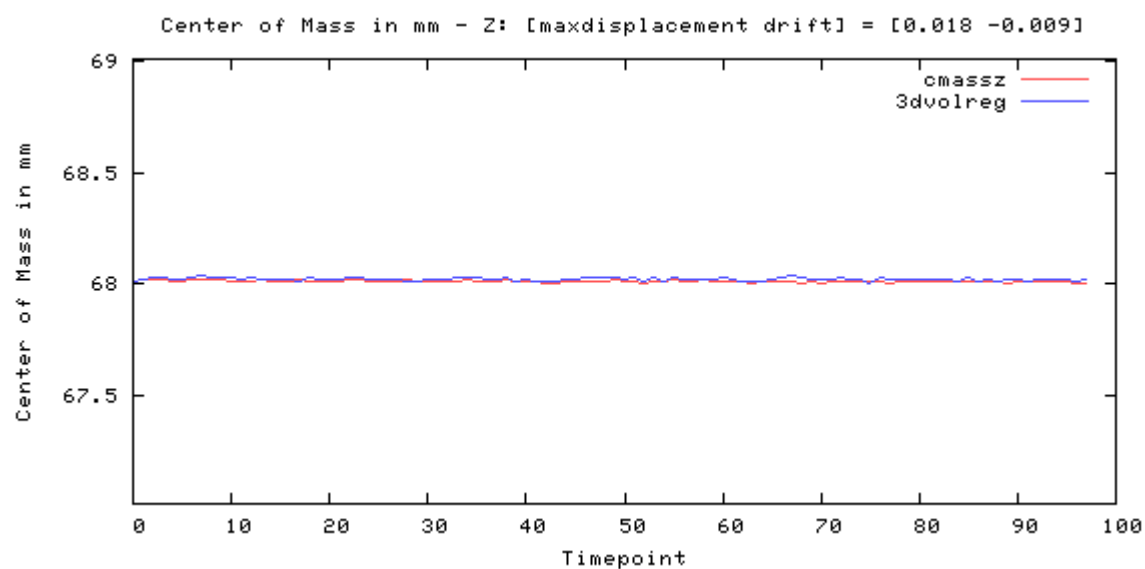


Odd-even difference image

-11109 13551




image min: -11109, image max: 13551

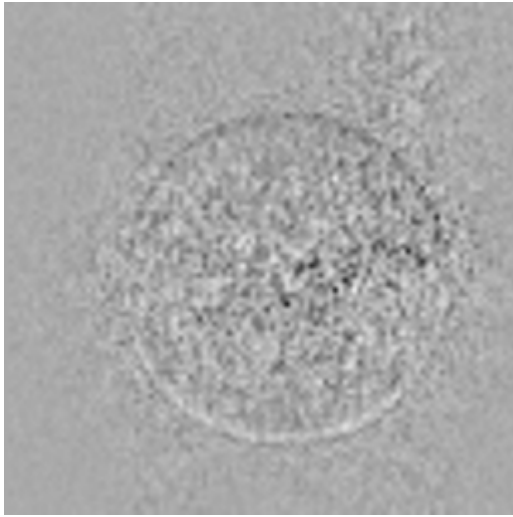

Mean image




image min: 0, image max: 52959.7

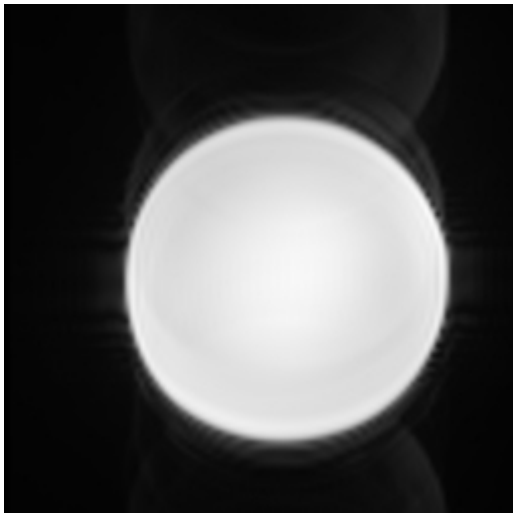

Standard Deviation image




image min: 0, image max: 506.016

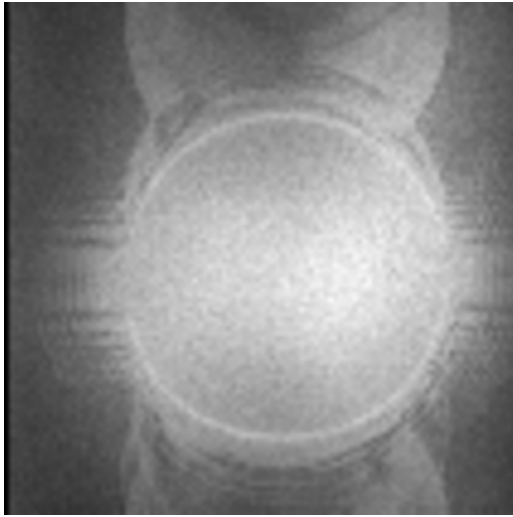

SFNR image



image min: 0, image max: 279.166

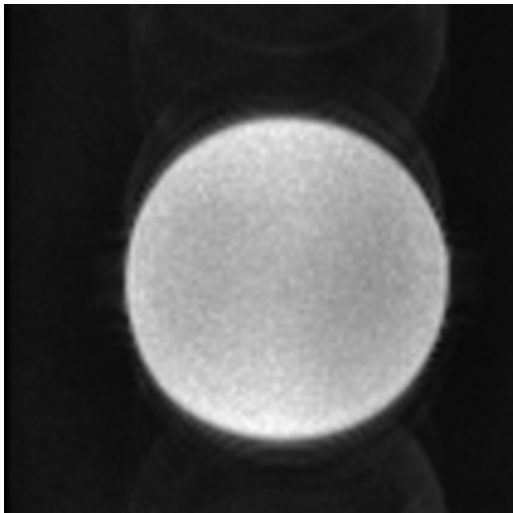

Acquisition parameters

dimensions	106x106x69x100
spacing	2mmx2mmx2mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	de5ccfd2-7e27-4f
studyid	Stability_20260930
seriesnumber	28
runnumber	1
magneticfield	3

description	ses-20ch_func-bold_acq-cmrrstd_task-tsnrxacpc
scandate	2026-09-30
scantime	14:08:17
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
receivecoilname	HeadNeck_20_TCS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-20ch_func-bold_acq-cmrrstd_task-tsnrxacpc
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260930140815.952500
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	2
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	RESEARCH
acquisitionduration	217.89500000000001

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	100
imagecomments	Not for diagnostic use, Unaliased MB3/ PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	69
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00