






Phantom QA results for /fbirn_out/ses-64ch_func-bold_acq-cmrrstd_task-tsnrxacpc.xml, slice 34




old acq-cmrrstd task-tsnrxacpc.xml [percent fluct (trend removed), drift, driftfi-


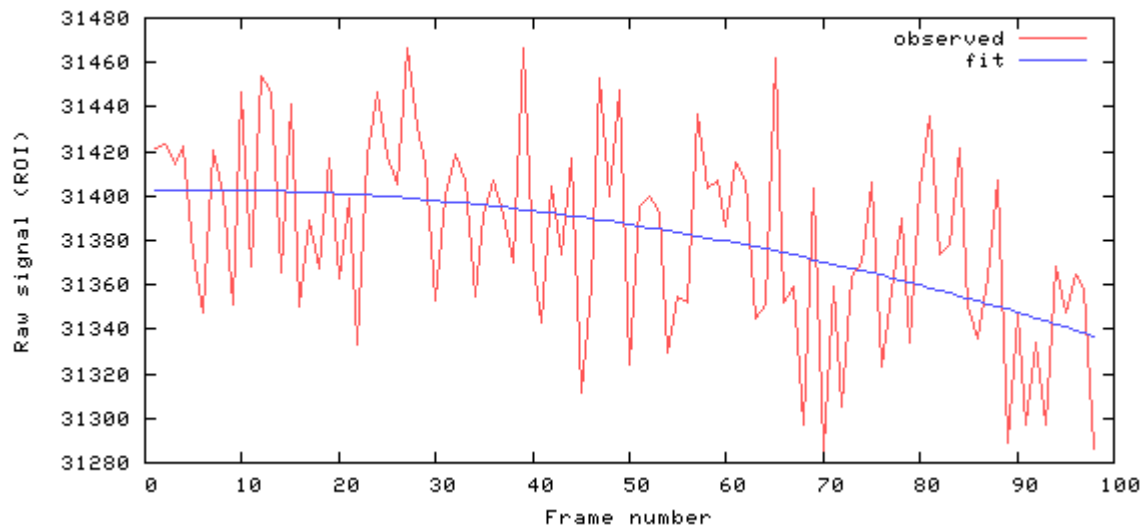


[mean, SNR, SFNR] = [31381.5 123.3 124.2]

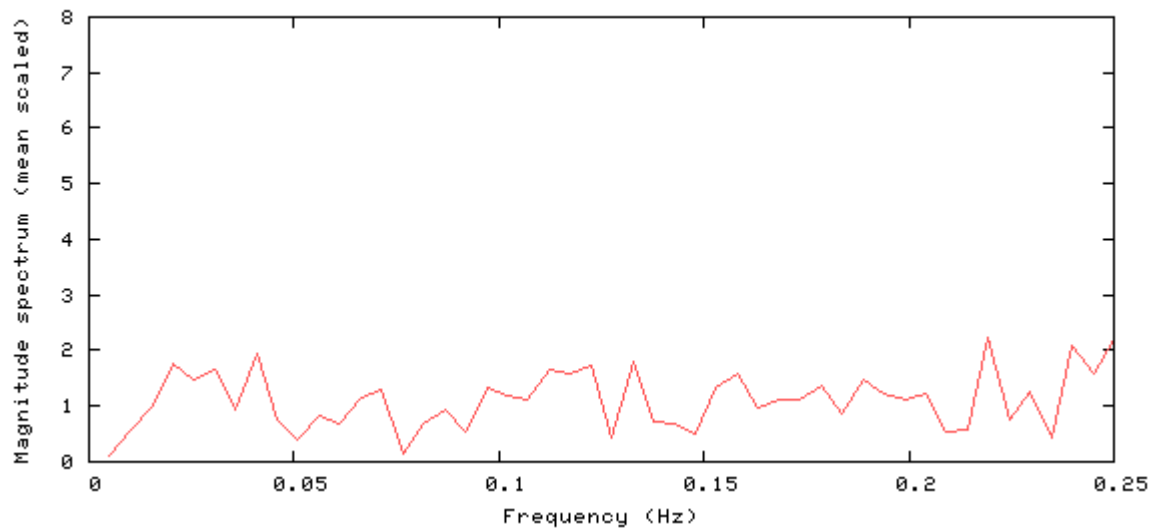

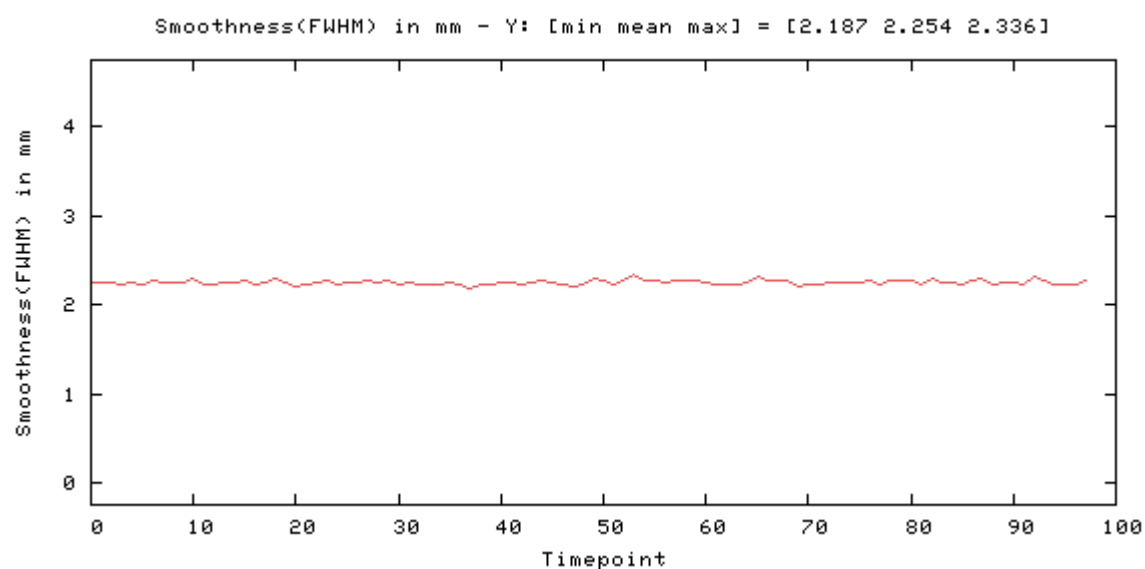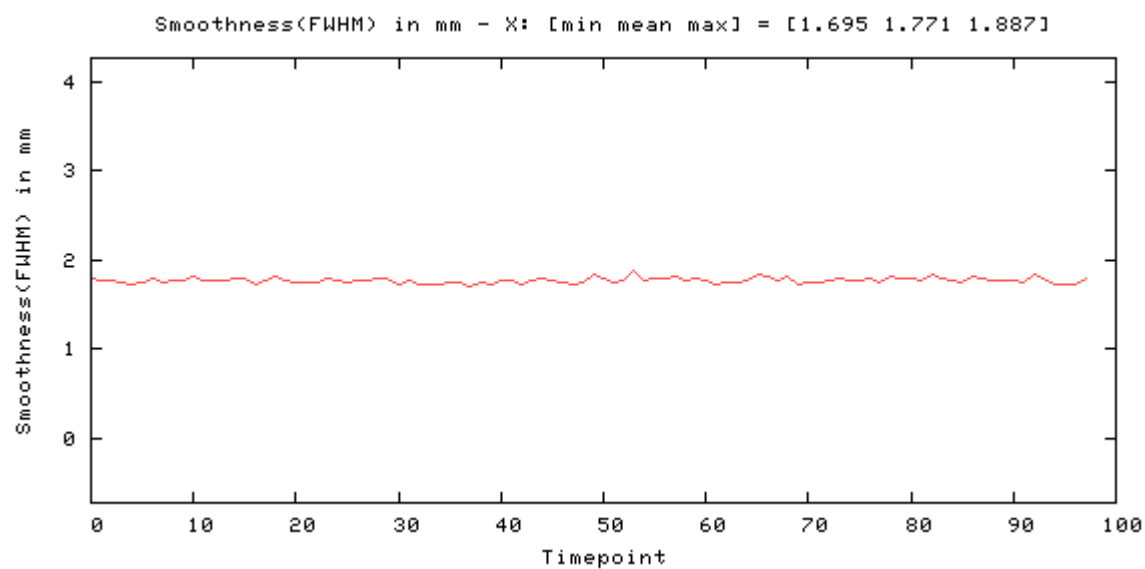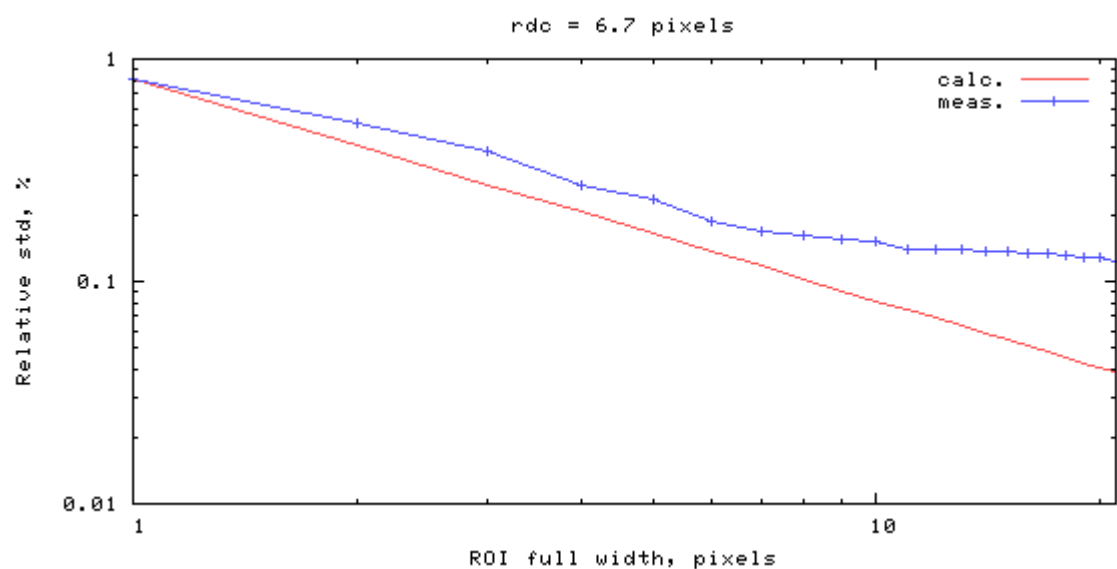

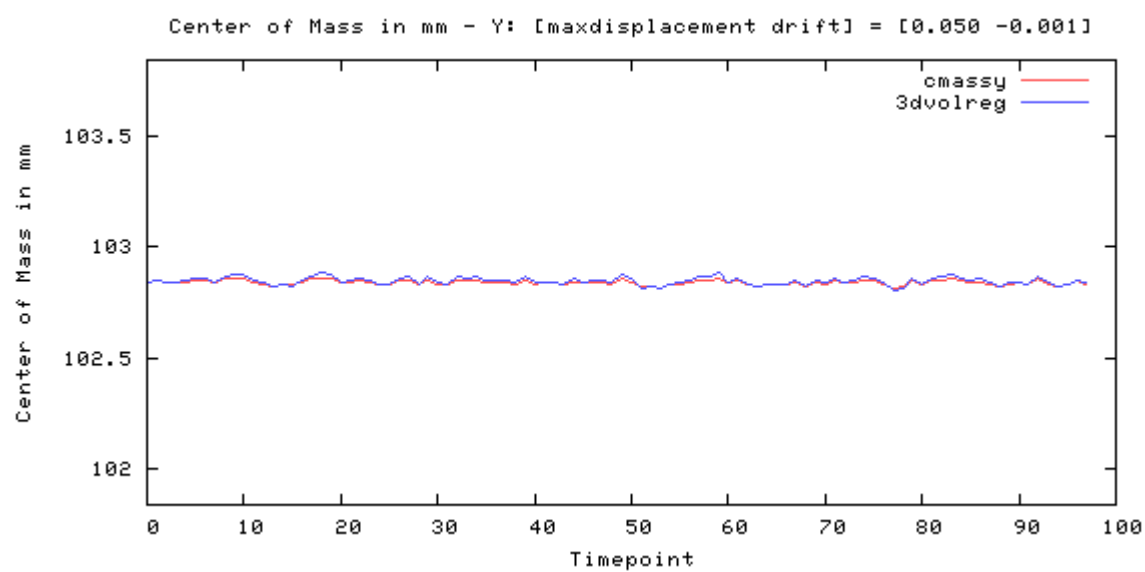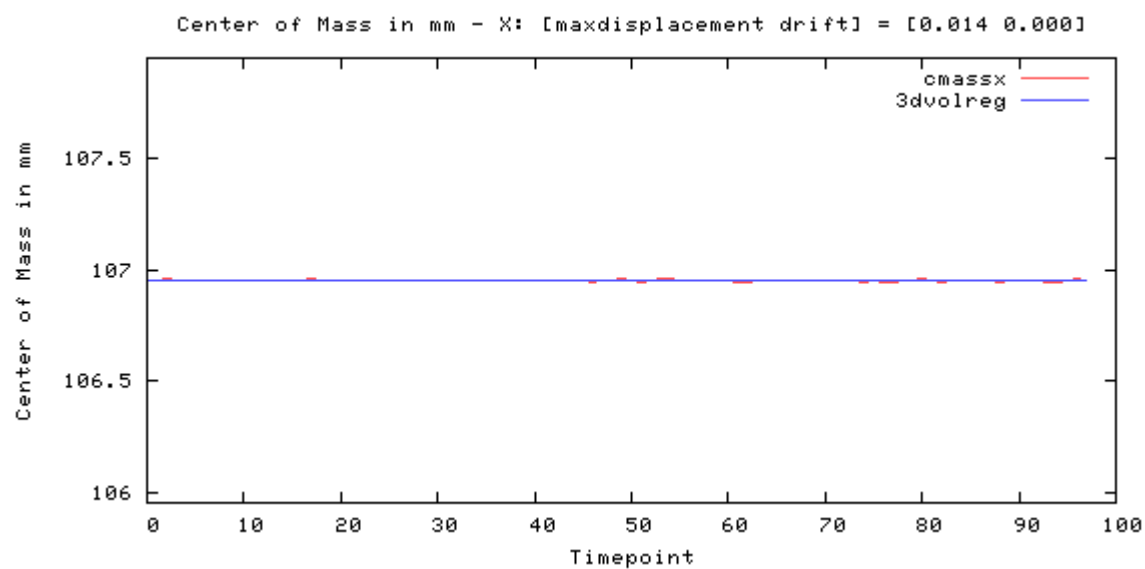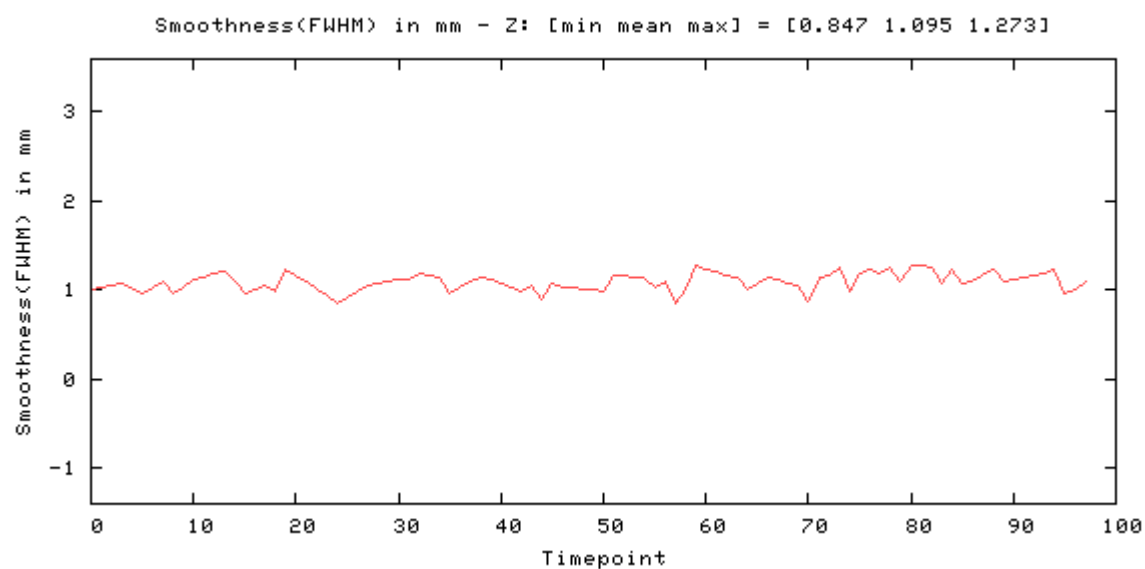


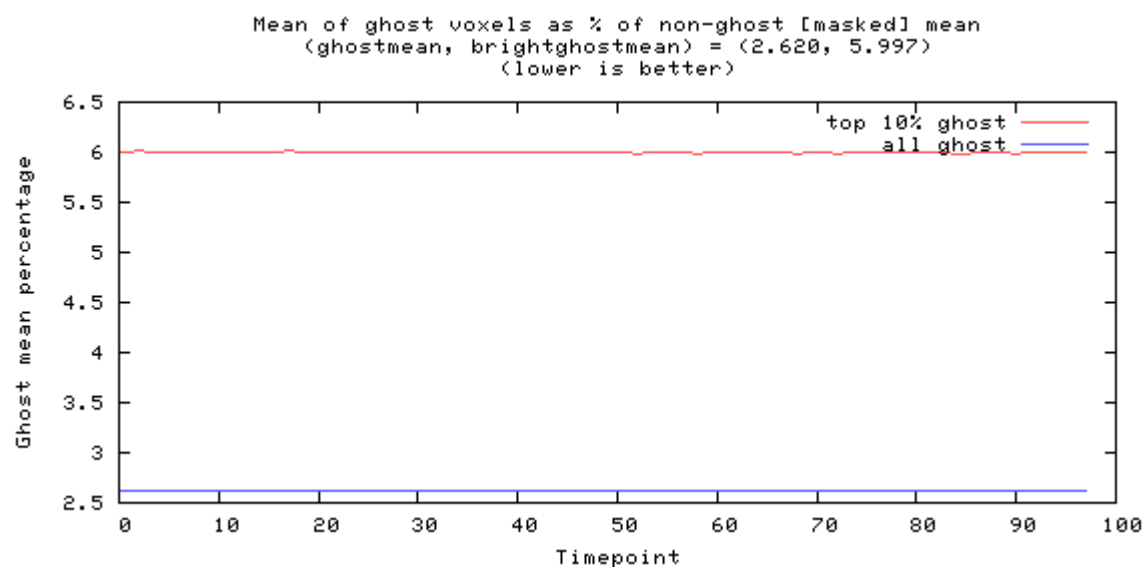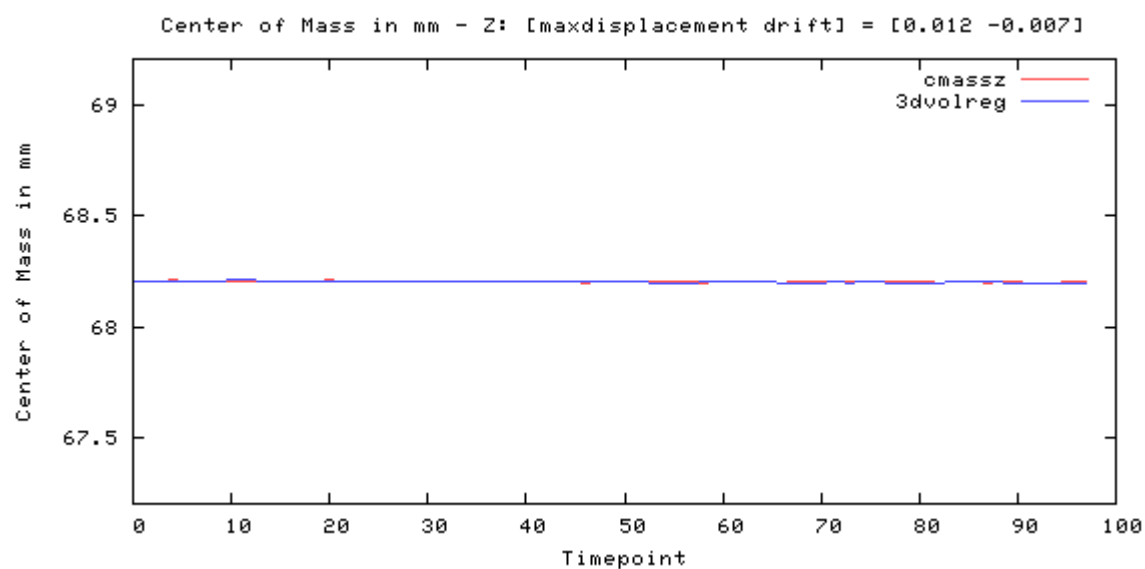

Odd-even difference image




image min: -7688, image max: 7184

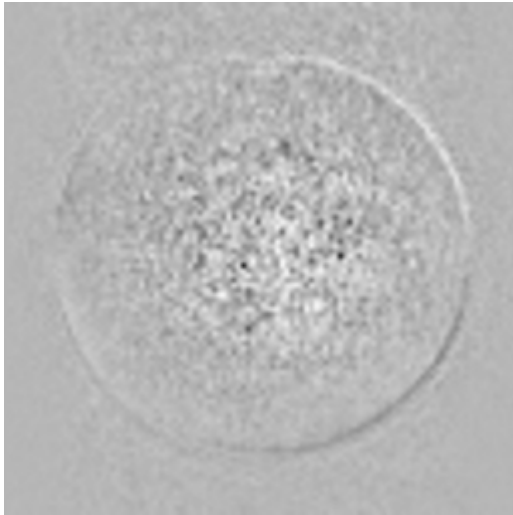

Mean image




image min: 0, image max: 44723.9

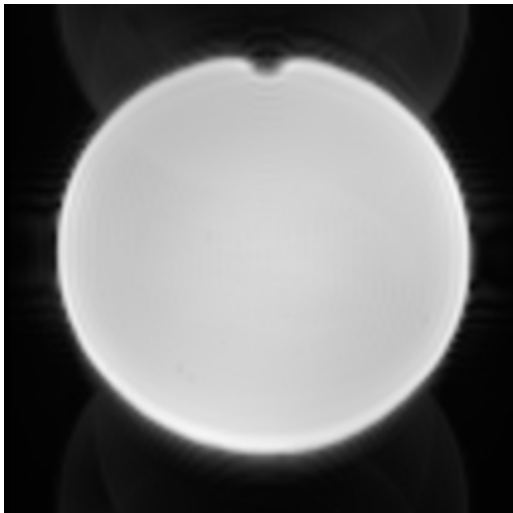

Standard Deviation image




image min: 0, image max: 363.57

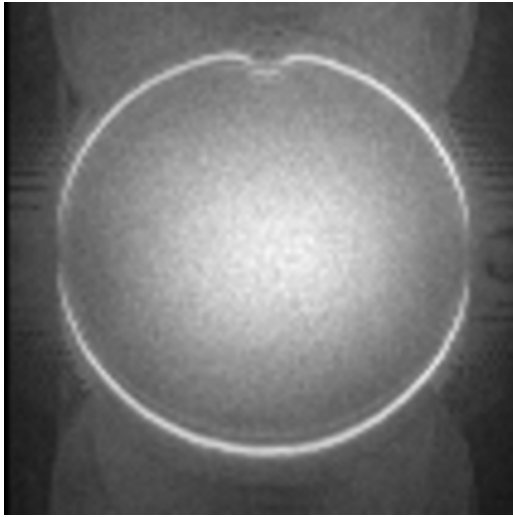

SFNR image



image min: 0, image max: 834.545

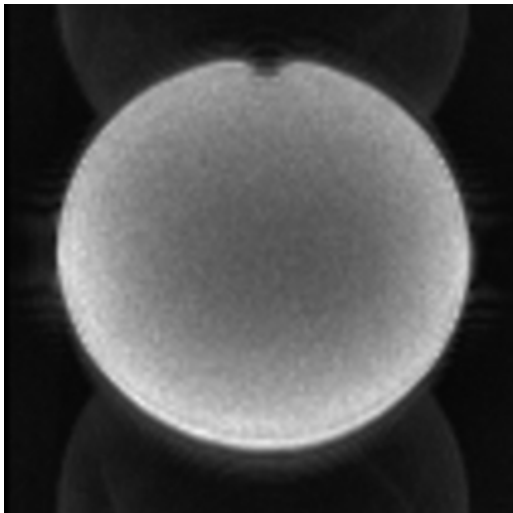

Acquisition parameters

dimensions	106x106x69x100
spacing	2mmx2mmx2mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	de5ccfd2-7e27-4f
studyid	Stability_20260930
seriesnumber	10
runnumber	1
magneticfield	3

description	ses-64ch_func-bold_acq-cmrrstd_task-tsnrxacpc
scandate	2026-09-30
scantime	13:16:47
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-64ch_func-bold_acq-cmrrstd_task-tsnrxacpc
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260930131645.535000
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	2
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	RESEARCH
acquisitionduration	217.892

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	100
imagecomments	Not for diagnostic use, Unaliased MB3/ PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	69
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00