






Phantom QA results for /fbirn_out/ses-64ch_func-bold_acq-cmrrstd_task-tsnr.xml, slice 34




`>bold_acq-cmrrstd_task-tsnr.xml` [percent fluct (trend removed), drift, driftfit]


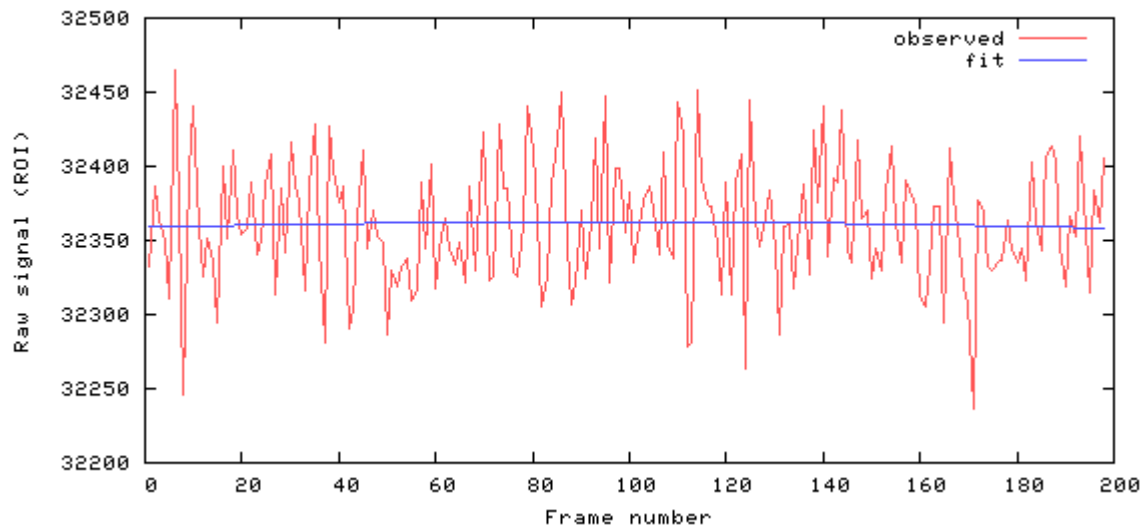


[mean, SNR, SFNR] = [32361.1 135.4 126.3]

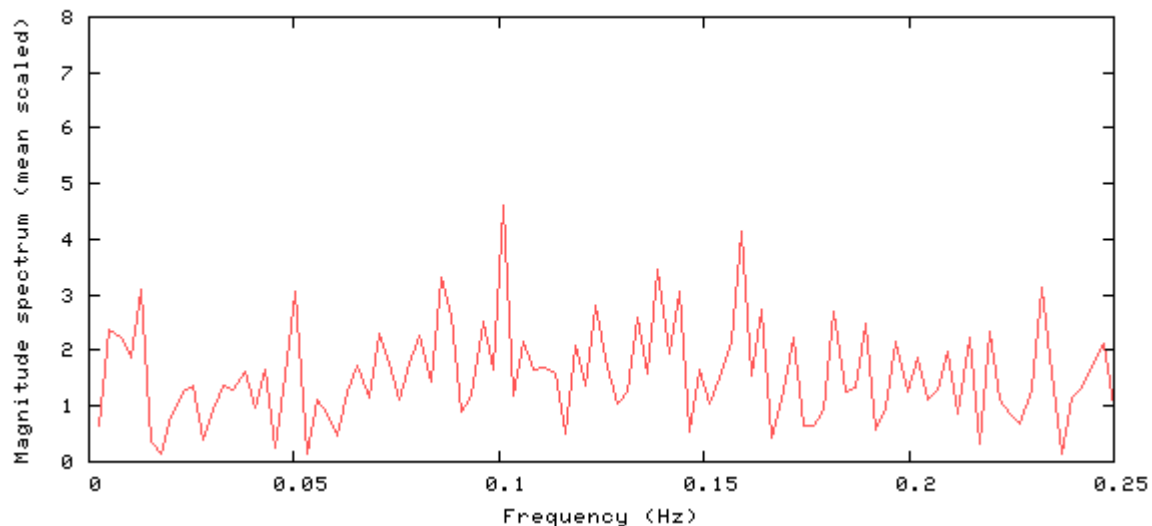

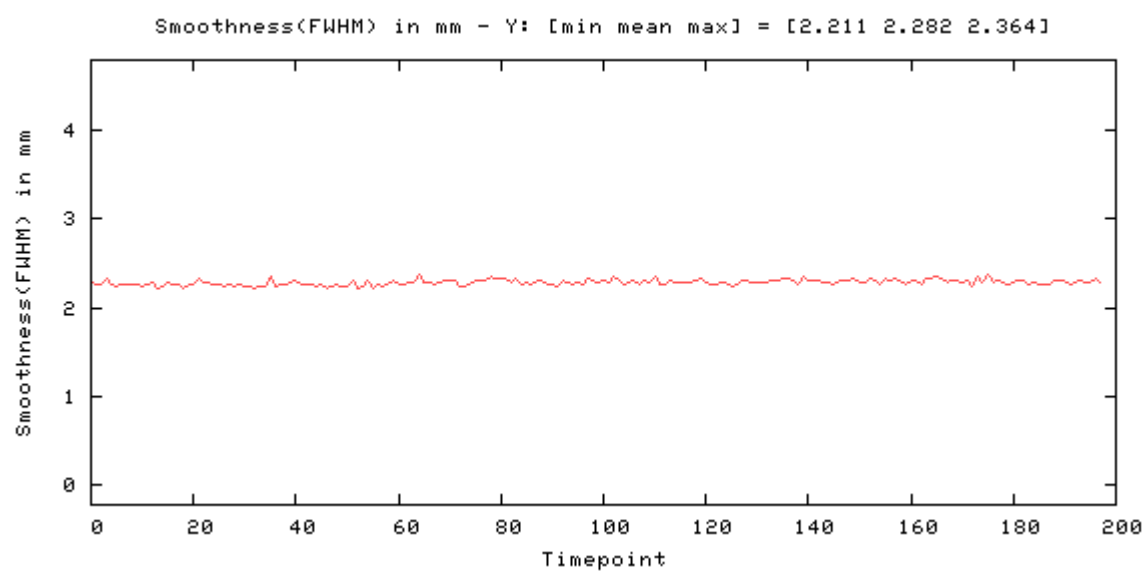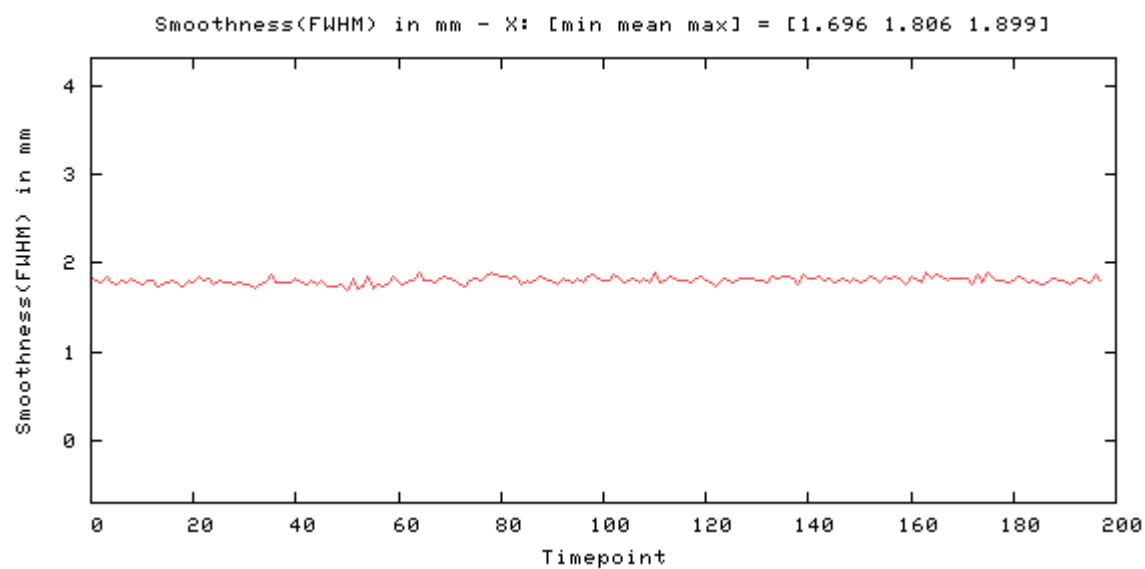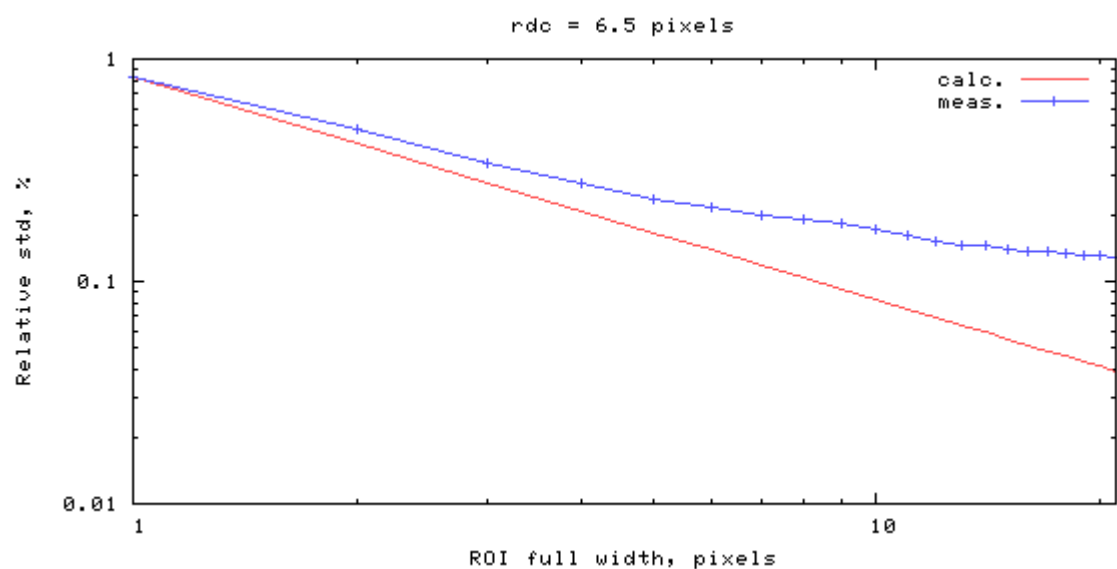

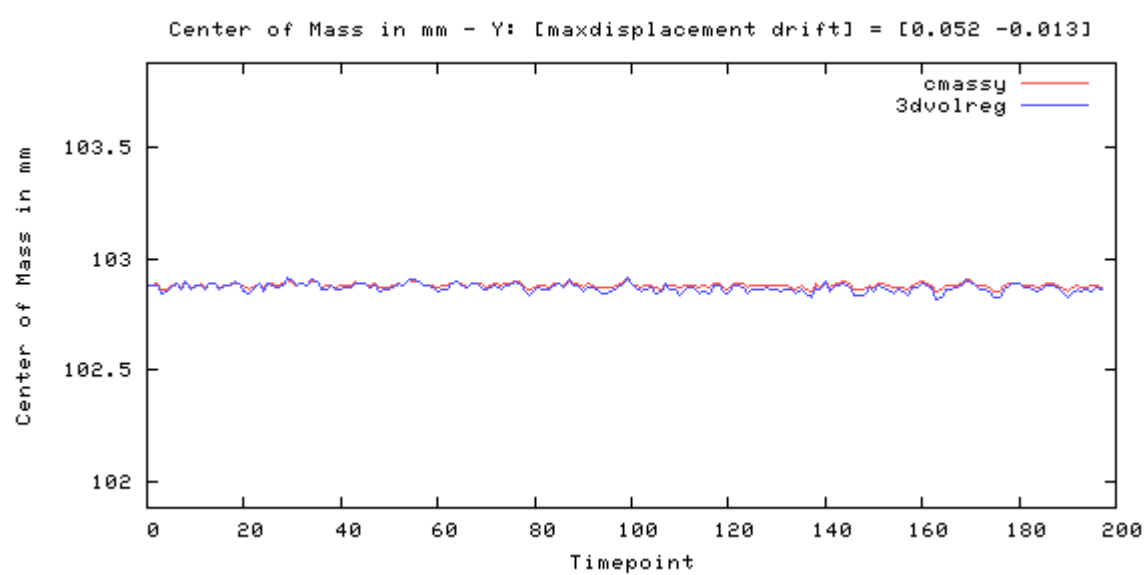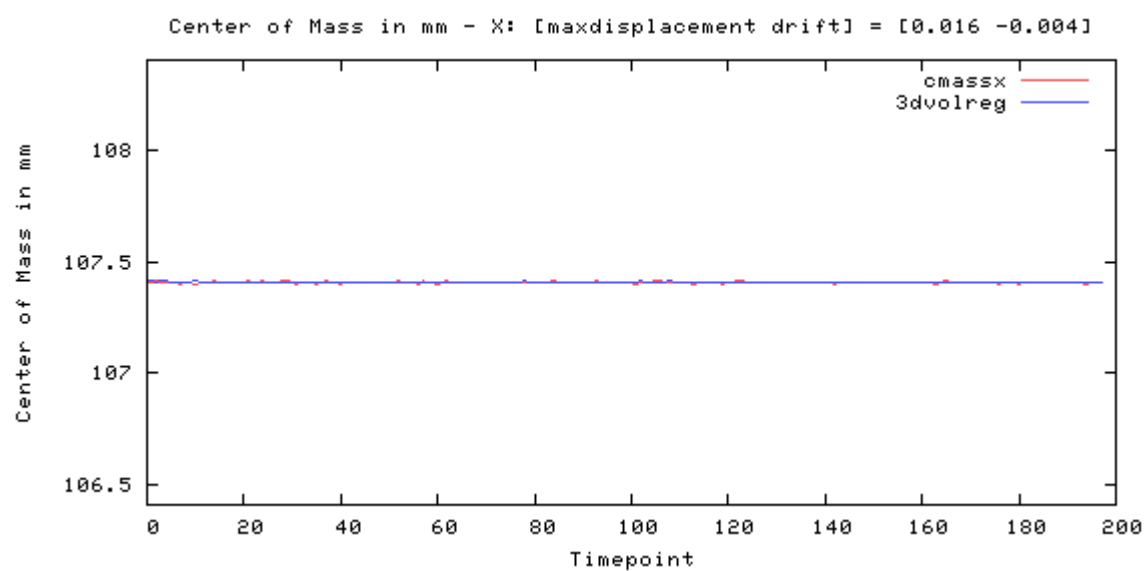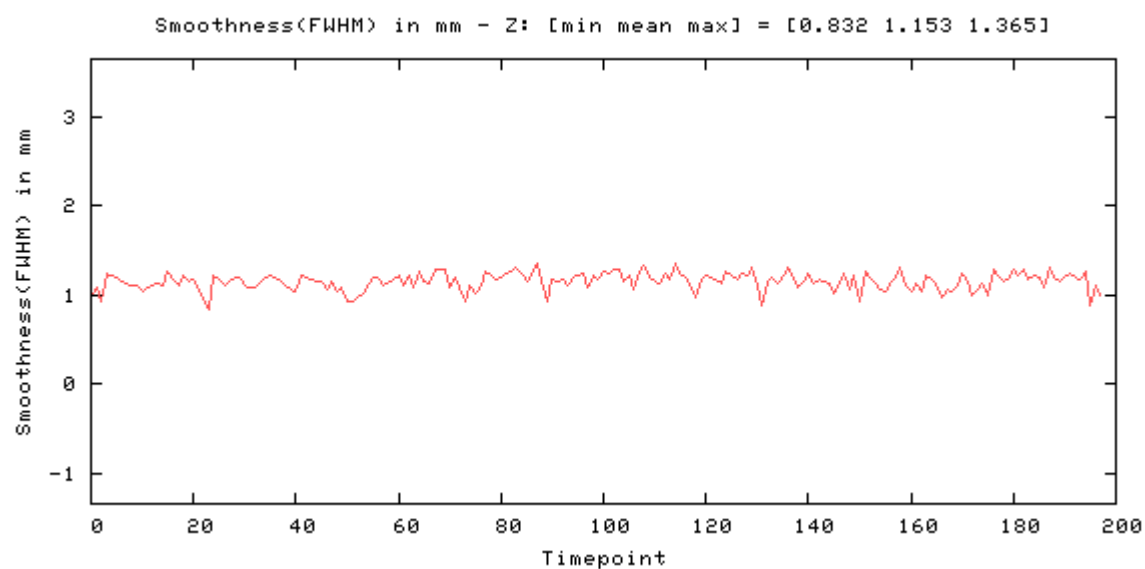


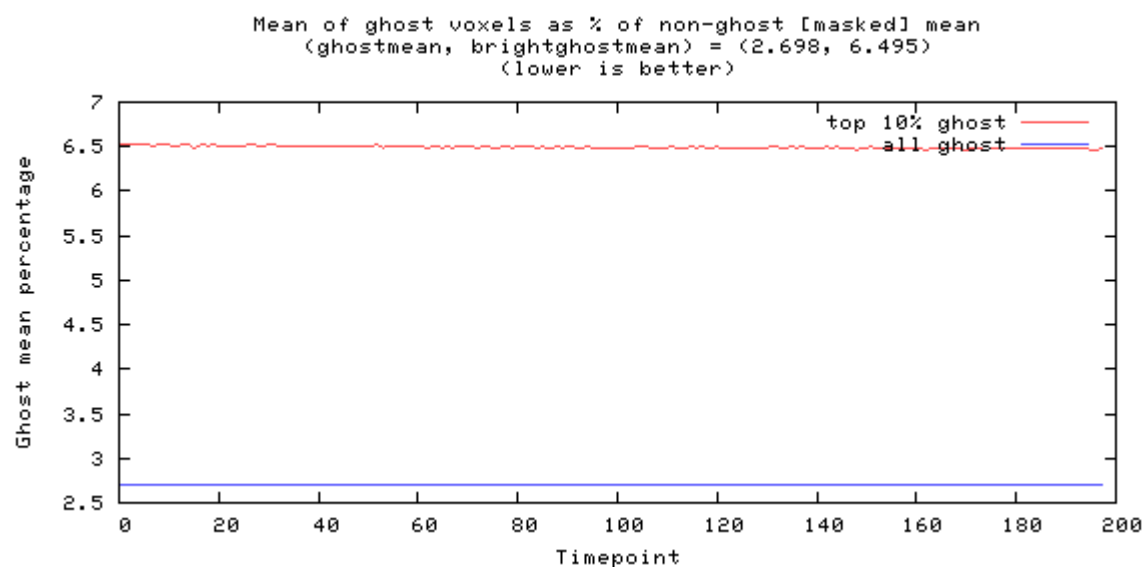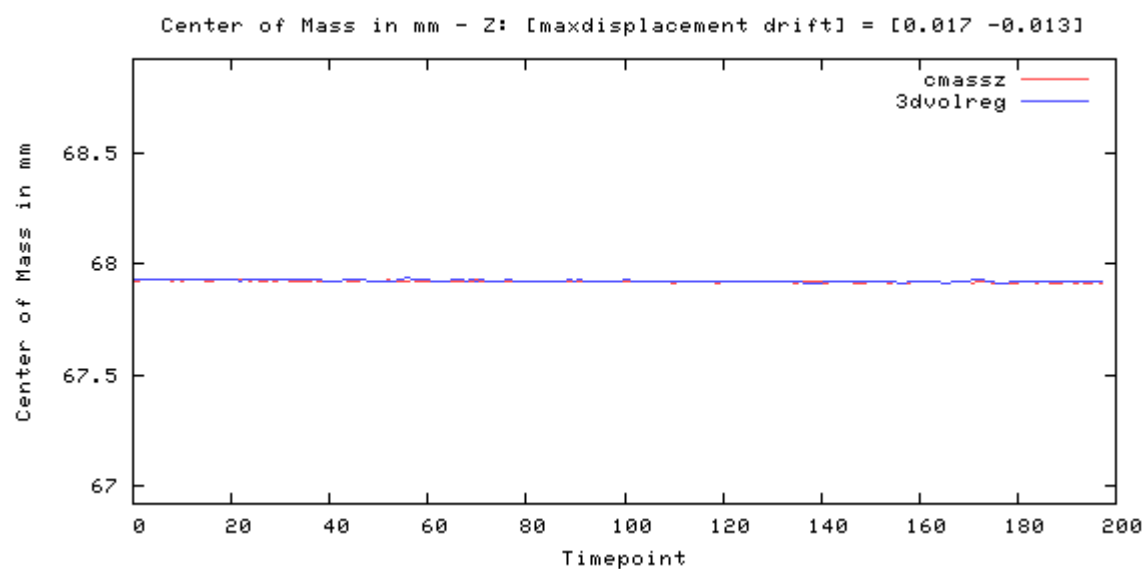

Odd-even difference image




image min: -10844, image max: 11768

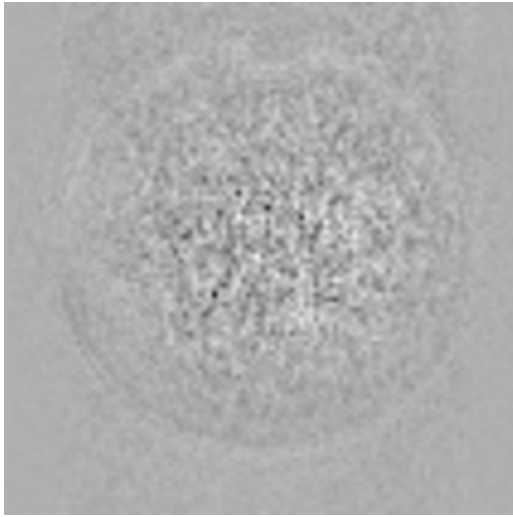

Mean image




image min: 0, image max: 44308.1

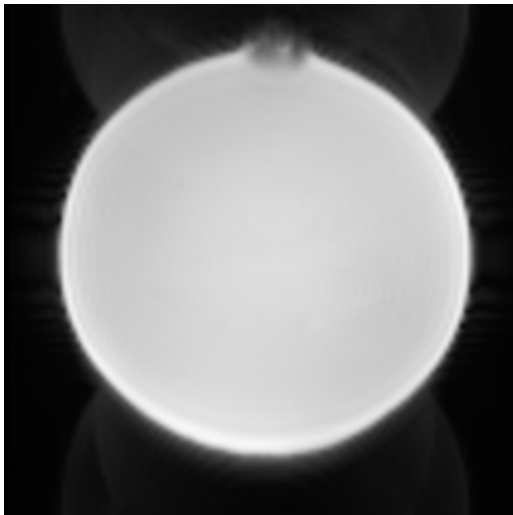

Standard Deviation image




image min: 0, image max: 318.076

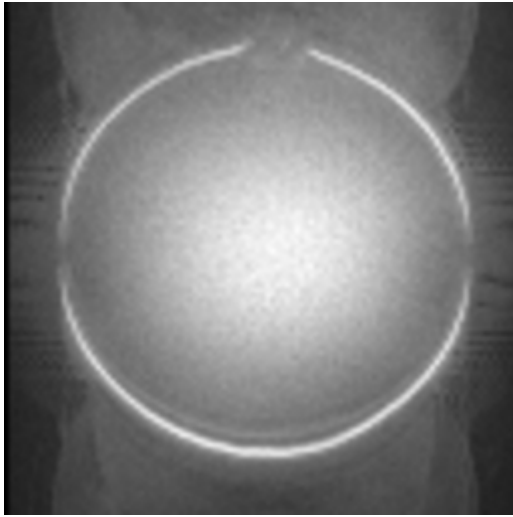

SFNR image



image min: 0, image max: 693.493

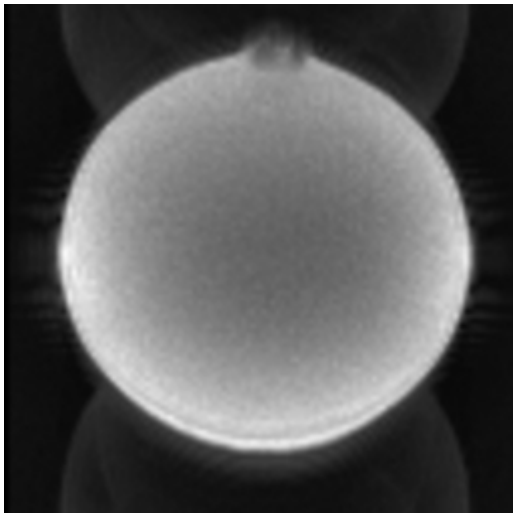

Acquisition parameters

dimensions	106x106x69x200
spacing	2mmx2mmx2mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	71482d6a-9d02-4c
studyid	Stability_20260909
seriesnumber	6
runnumber	1
magneticfield	3

description	ses-64ch_func-bold_acq-cmrrstd_task-tsnr
scandate	2026-09-09
scantime	09:38:44
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-64ch_func-bold_acq-cmrrstd_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260909093842.782500
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	2
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	RESEARCH
acquisitionduration	417.892
resonantnucleus	1H

applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
imagecomments	Not for diagnostic use, Unaliased MB3/ PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	69
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00