






Phantom QA results for / fbirn_out/ ses20ch_funcbold_acqcmrrstd_tasktsnr slice 34




.d acqcmrrstd tasktsnrxaopo 34.xml [percent fluct (trend removed), drift, driftfit


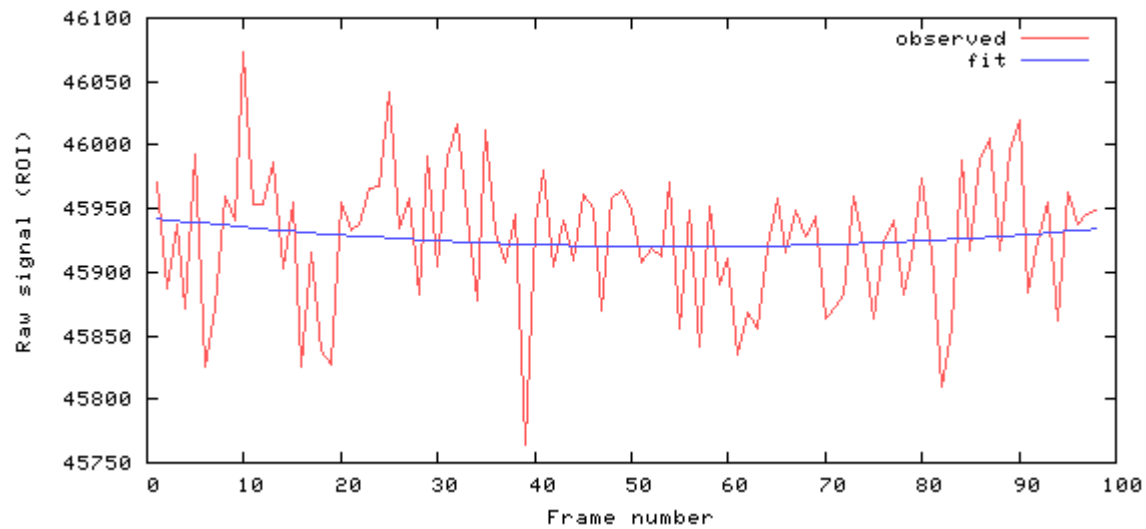


[mean, SNR, SFNR] = [45926.4 127.4 131.9]

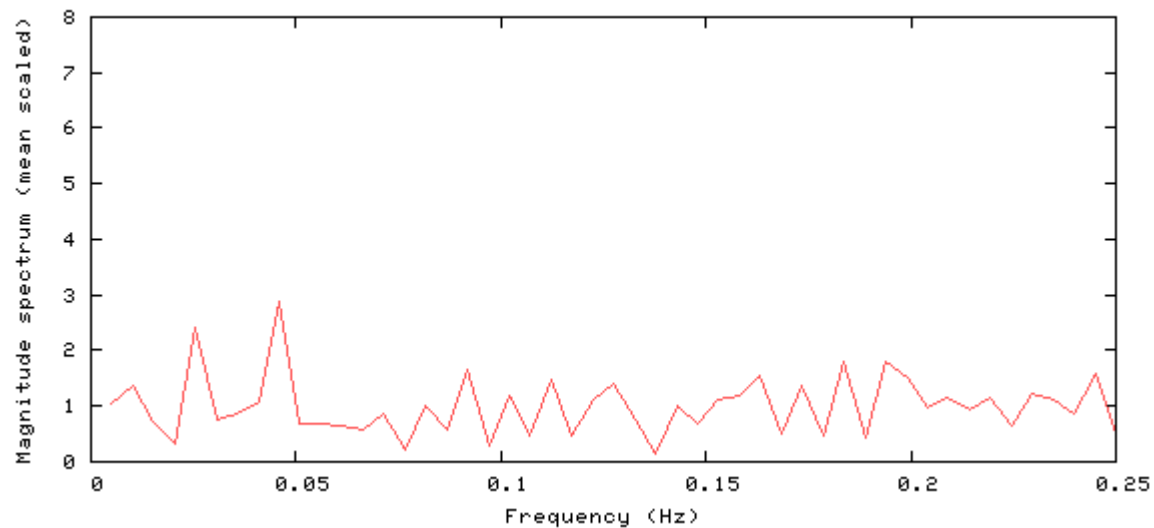

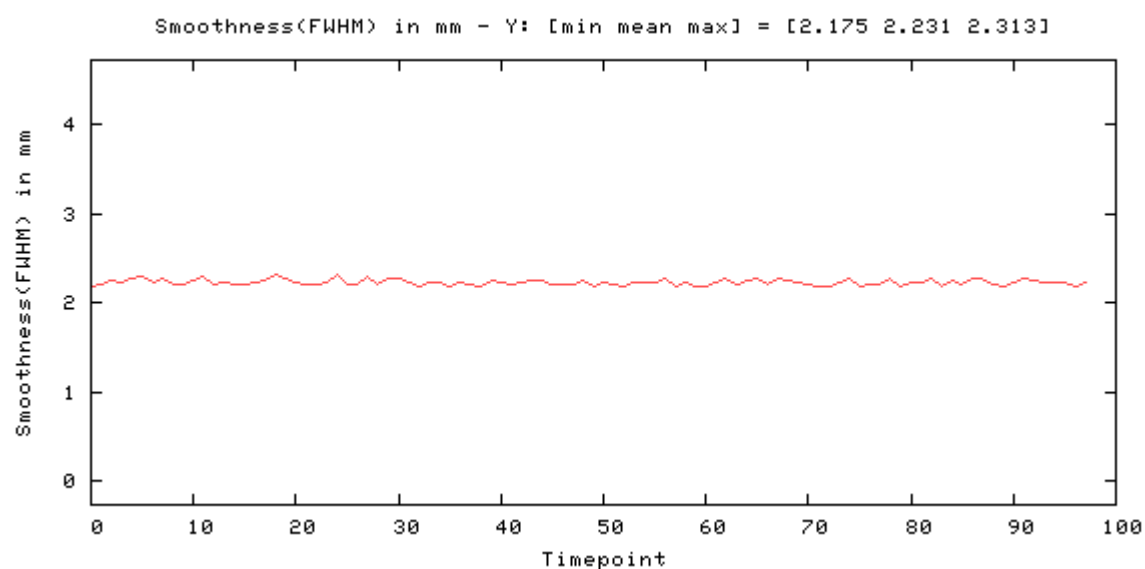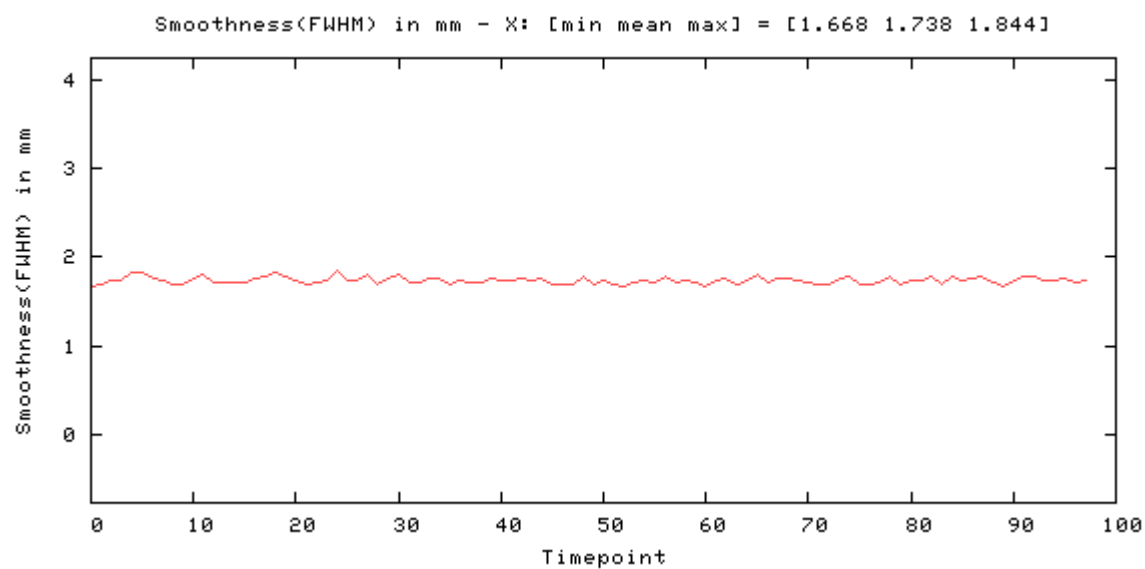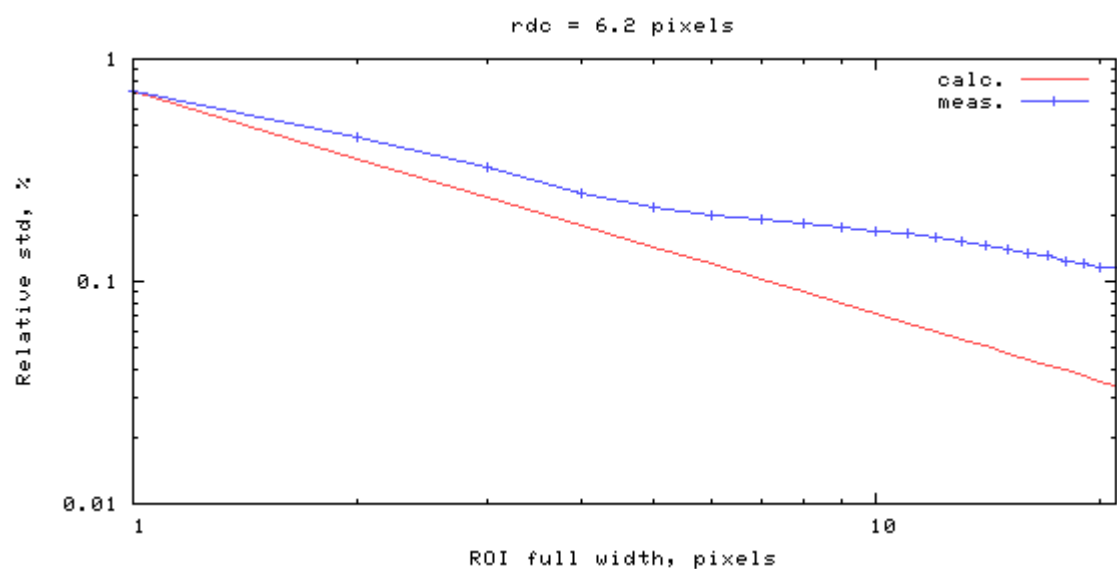

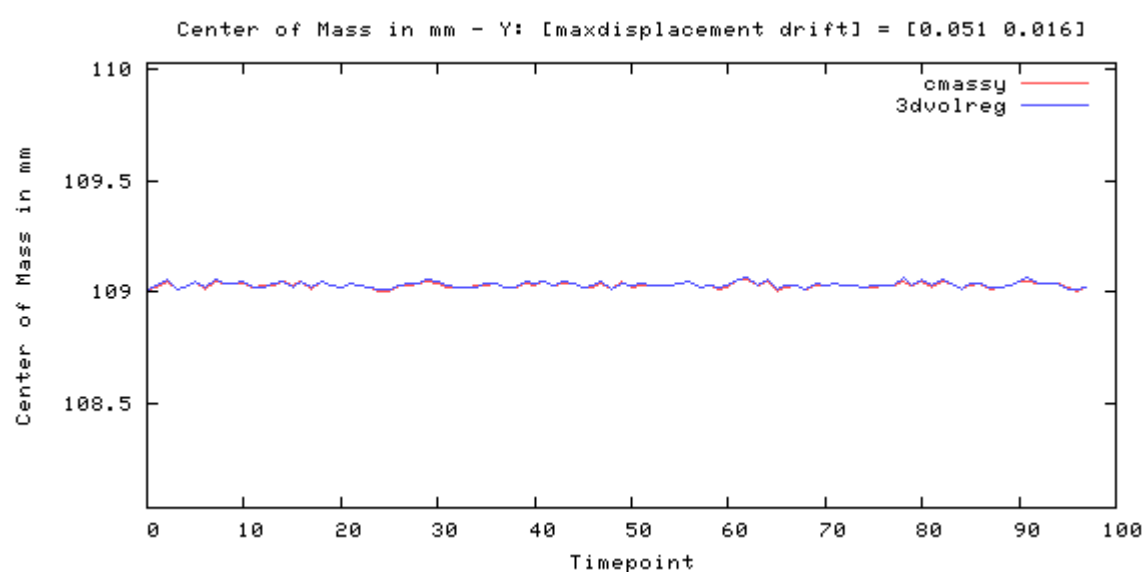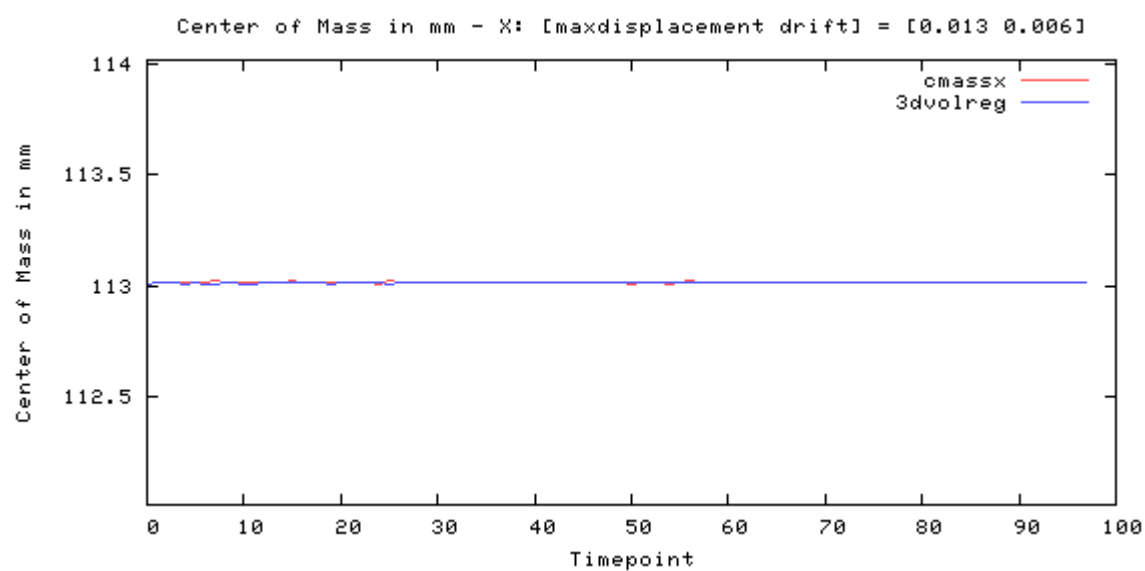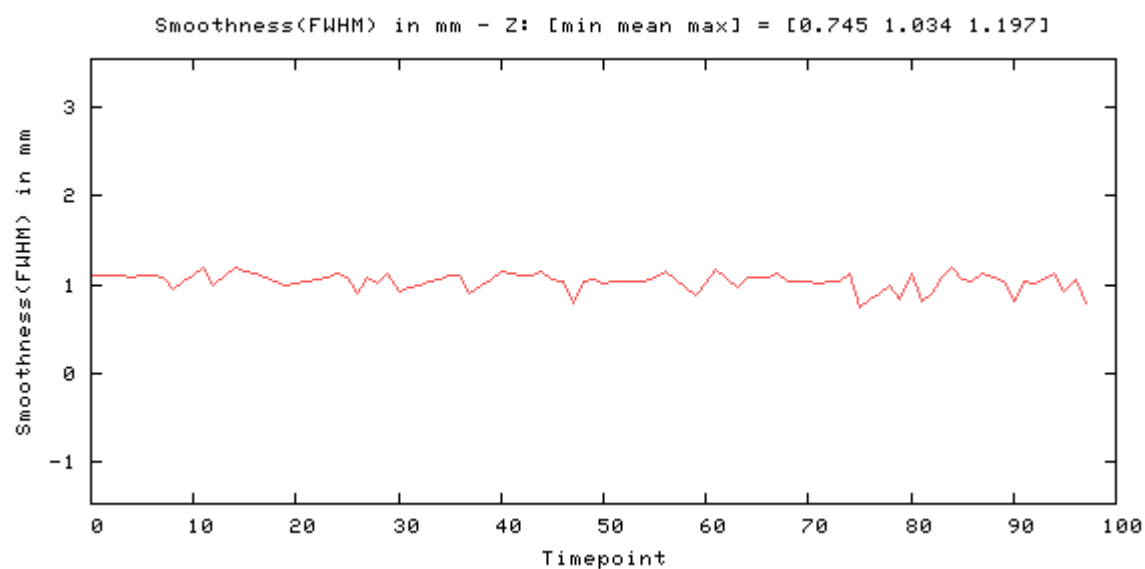

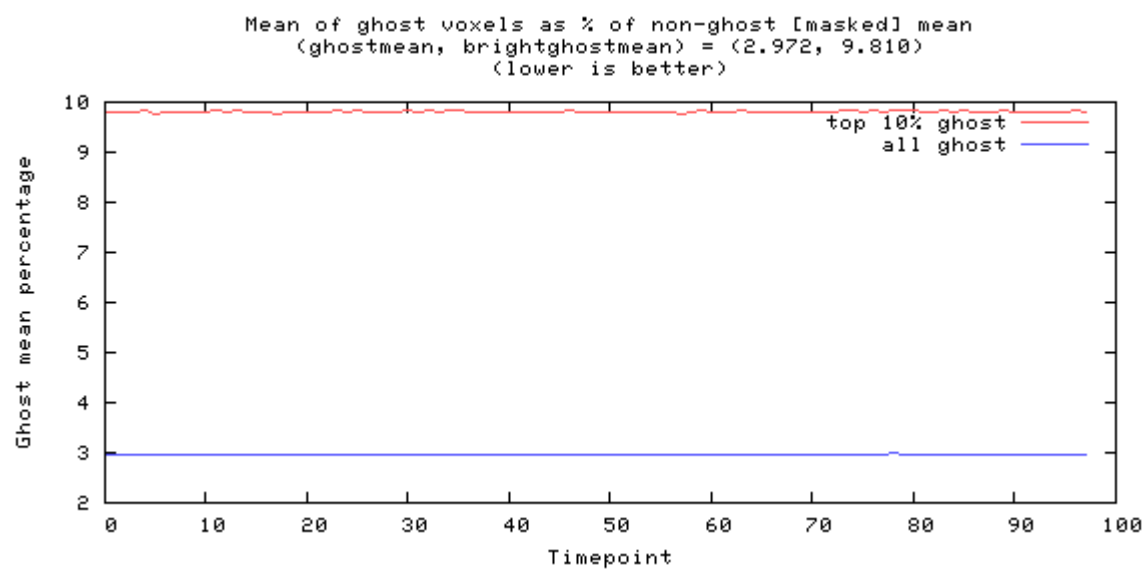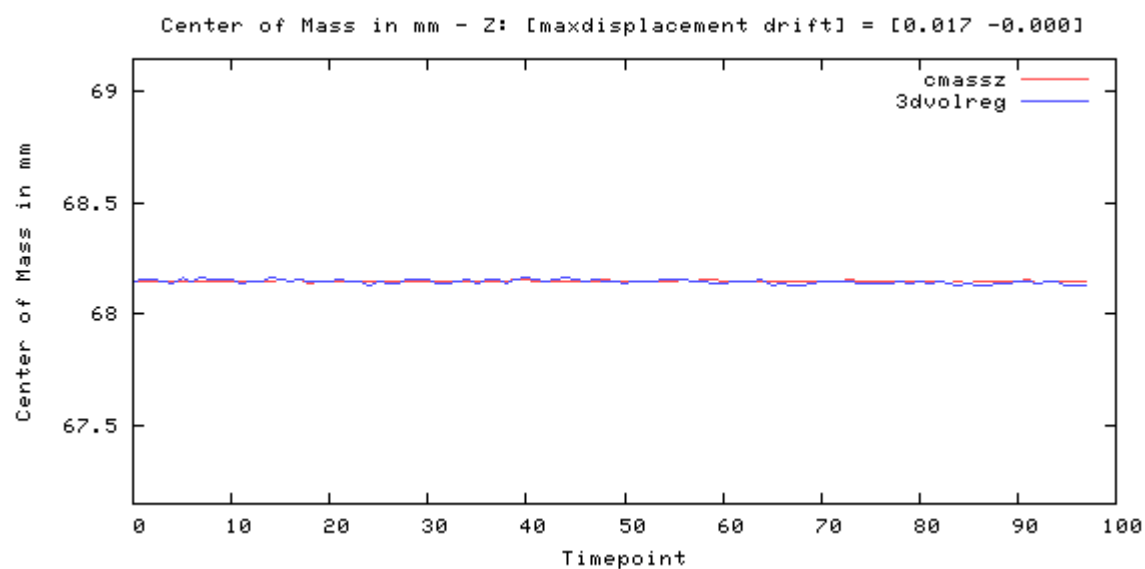


Odd-even difference image

-11807 14619




image min: -11807, image max: 14619

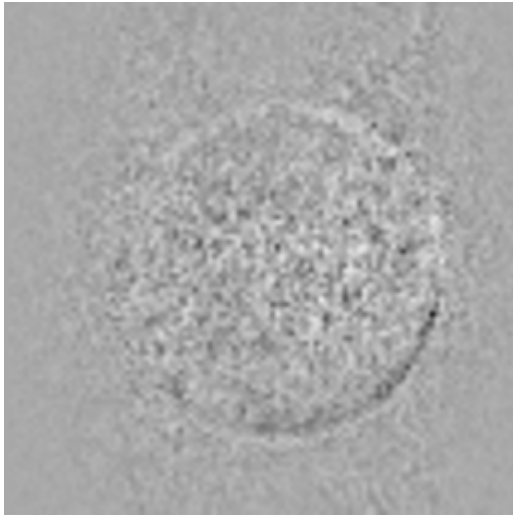

Mean image




image min: 0, image max: 52891.3

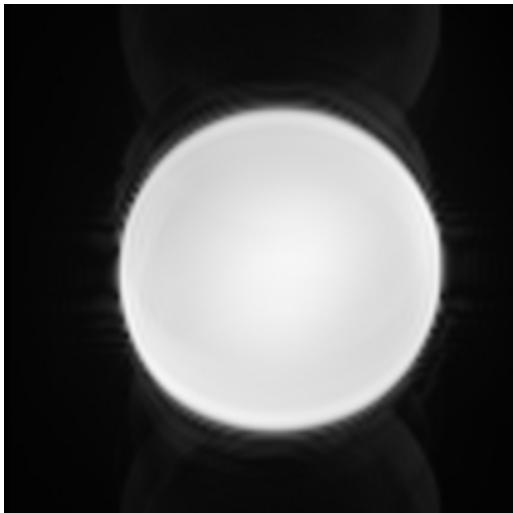

Standard Deviation image




image min: 0, image max: 526.029

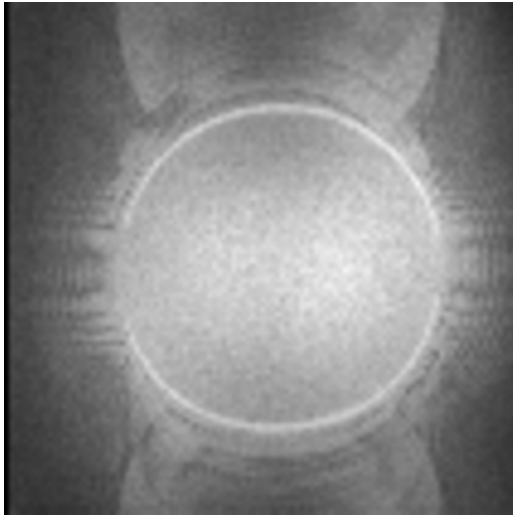

SFNR image



image min: 0, image max: 275.093

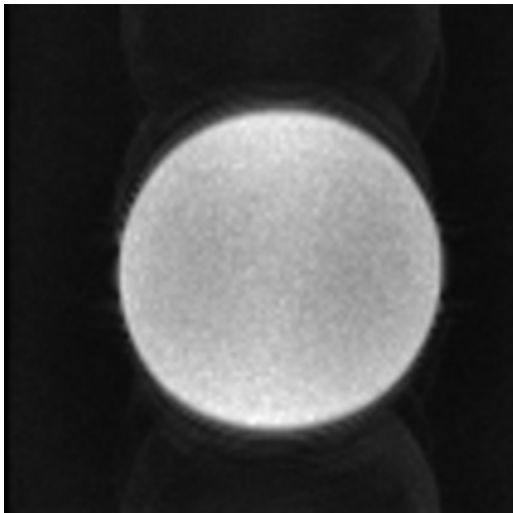

Acquisition parameters

dimensions	106x106x69x100
spacing	2mmx2mmx2mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	acdebd74-b79b-4f
studyid	Stability_20260615
seriesnumber	34
runnumber	1
magneticfield	3

description	ses-20ch_func-bold_acq-cmrrstd_task-tsnrxacpc
scandate	2026-06-15
scantime	10:49:33
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
receivecoilname	HeadNeck_20_TCS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-20ch_func-bold_acq-cmrrstd_task-tsnrxacpc
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260615104931.622500
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	2
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	RESEARCH
acquisitionduration	217.892

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	100
imagecomments	Not for diagnostic use, Unaliased MB3/ PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	69
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00