






Phantom QA results for /fbirn_out/ses32ch_funcbold_acqsiemensstd_taskk slice 17




.d acqsiemensstd tasktsnr 21.xml [percent fluct (trend removed), drift, driftfit]


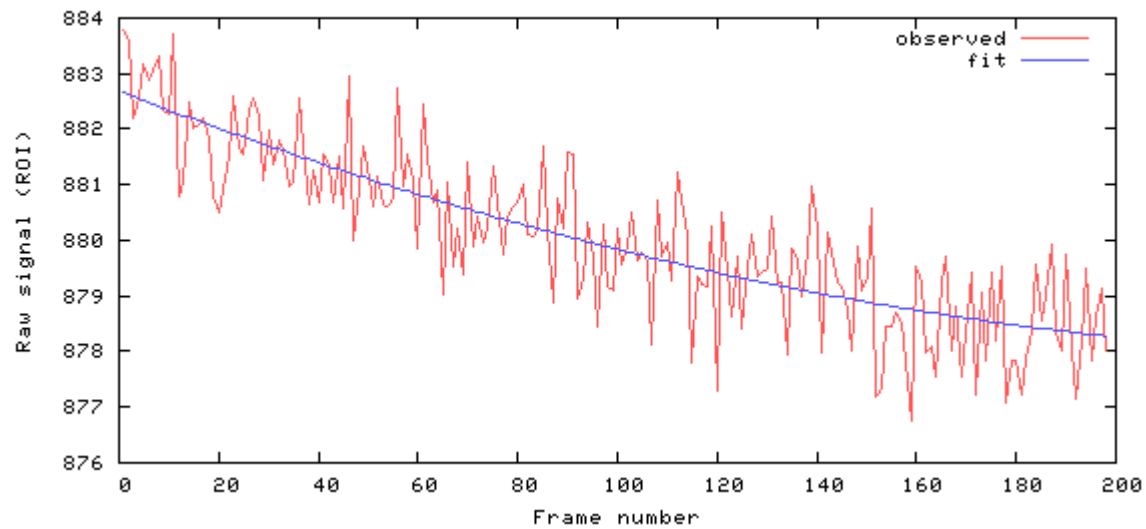


[mean, SNR, SFNR] = [880.1 113.3 125.1]

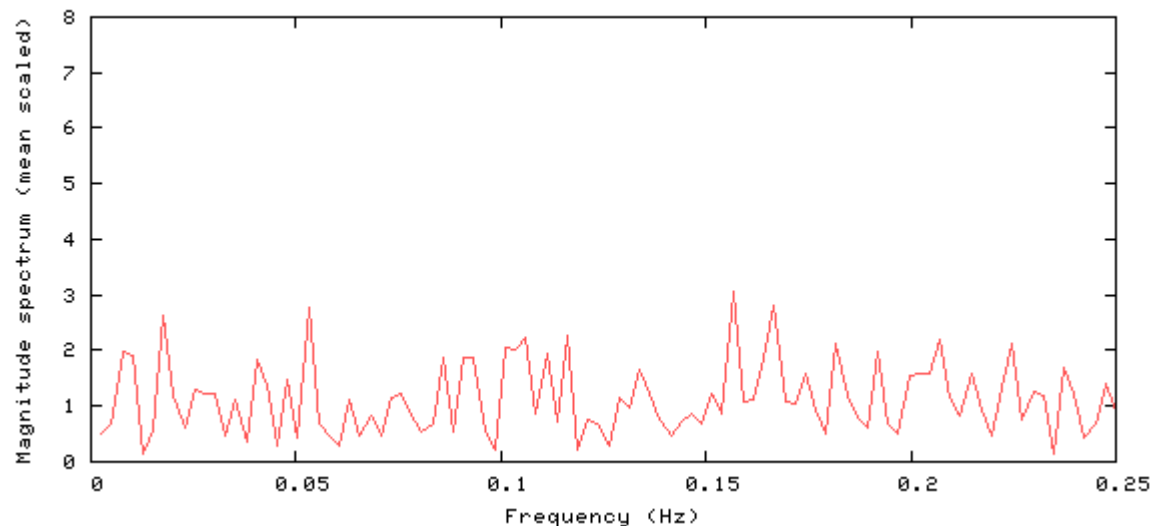

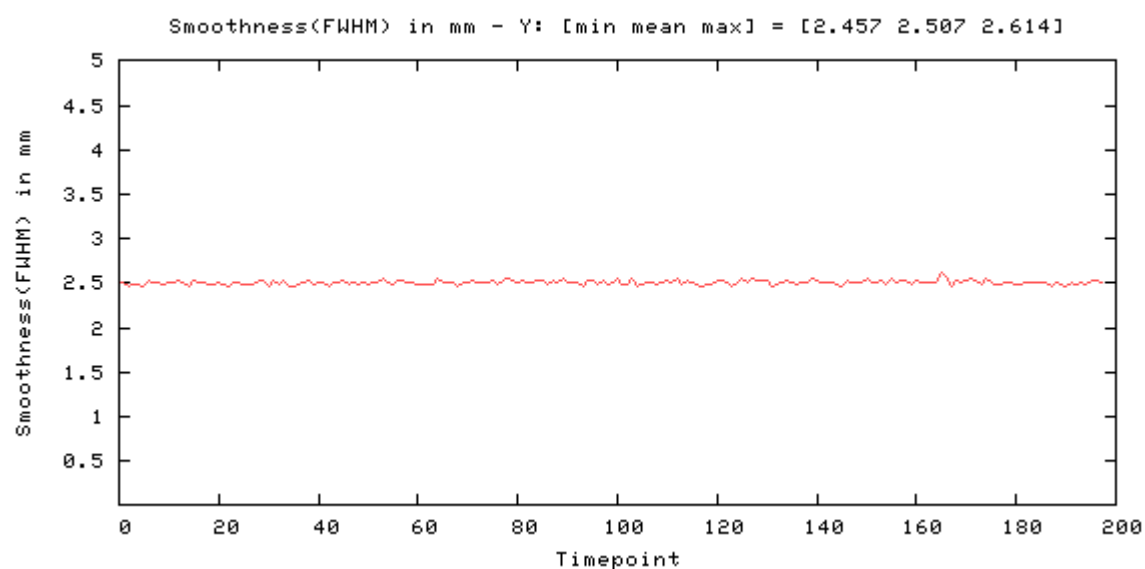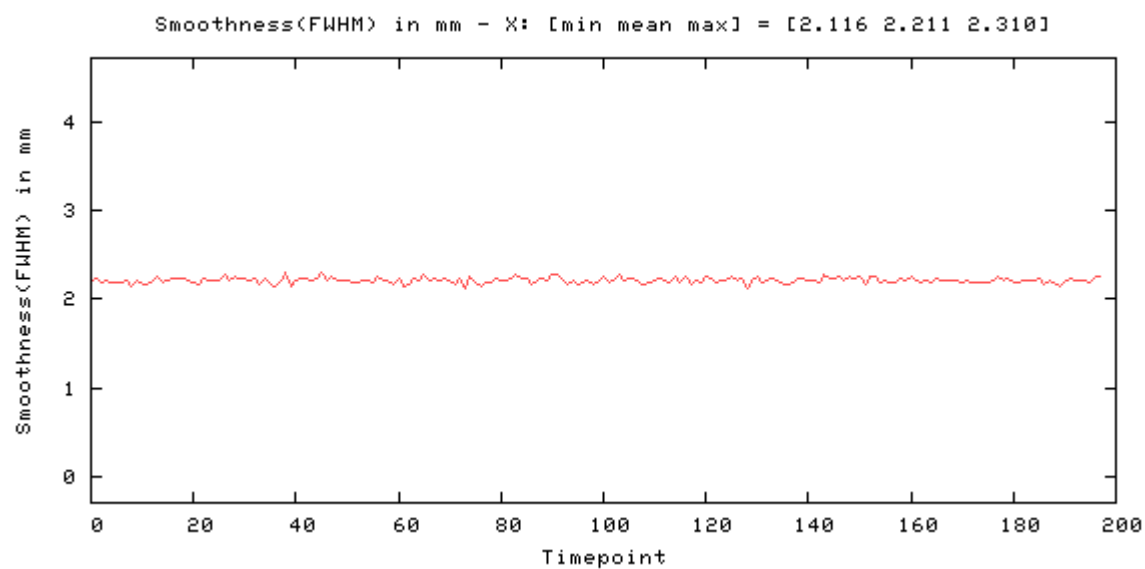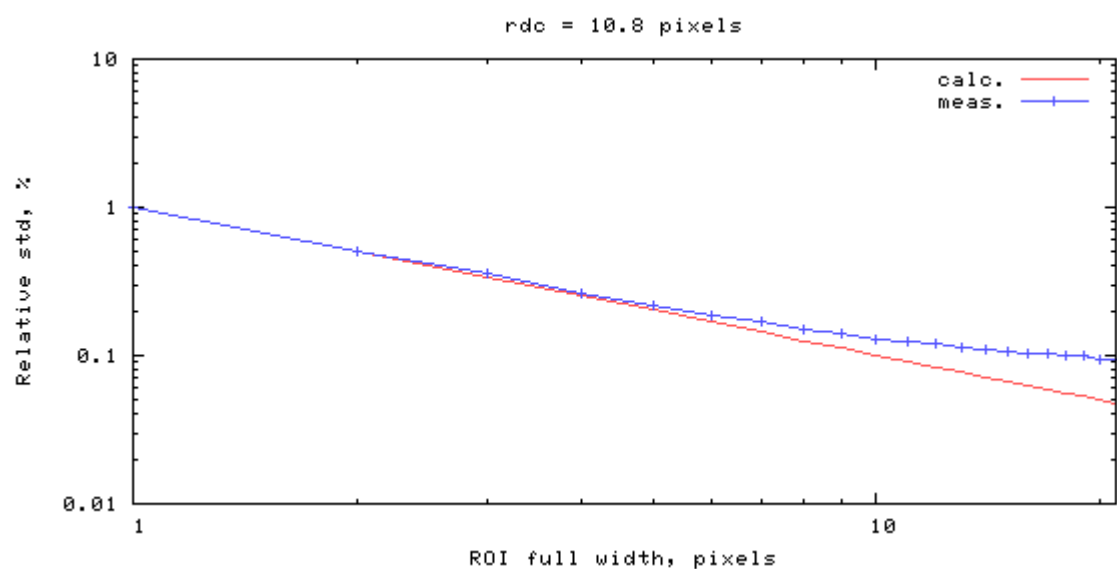

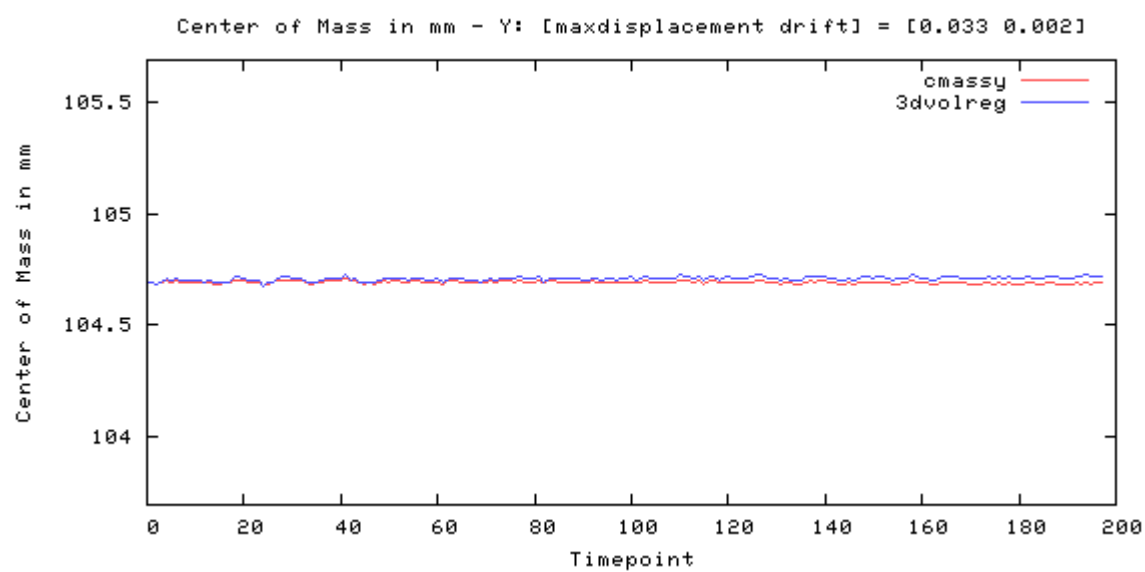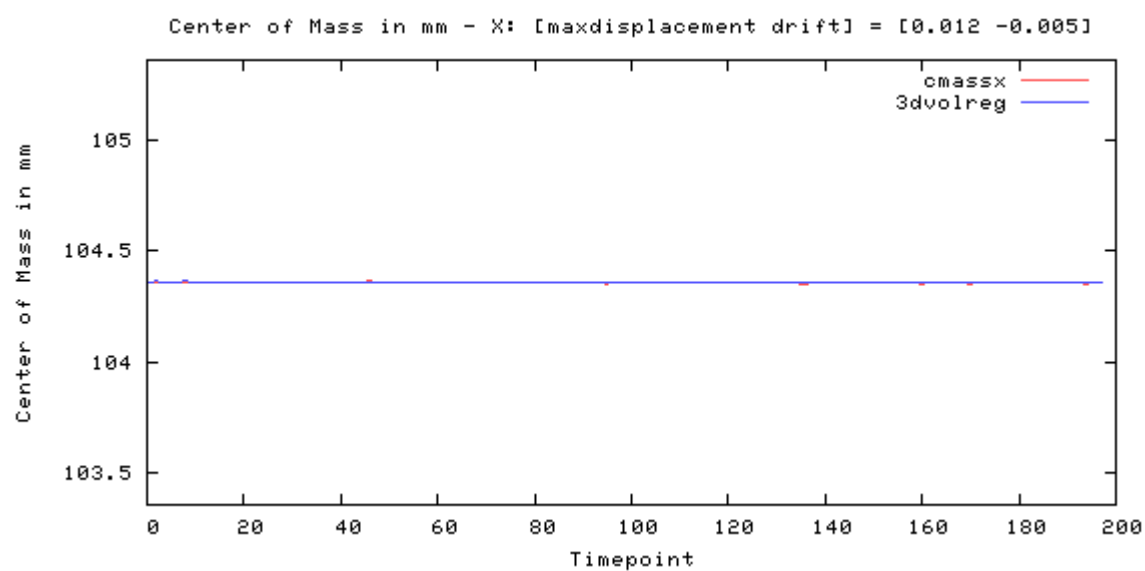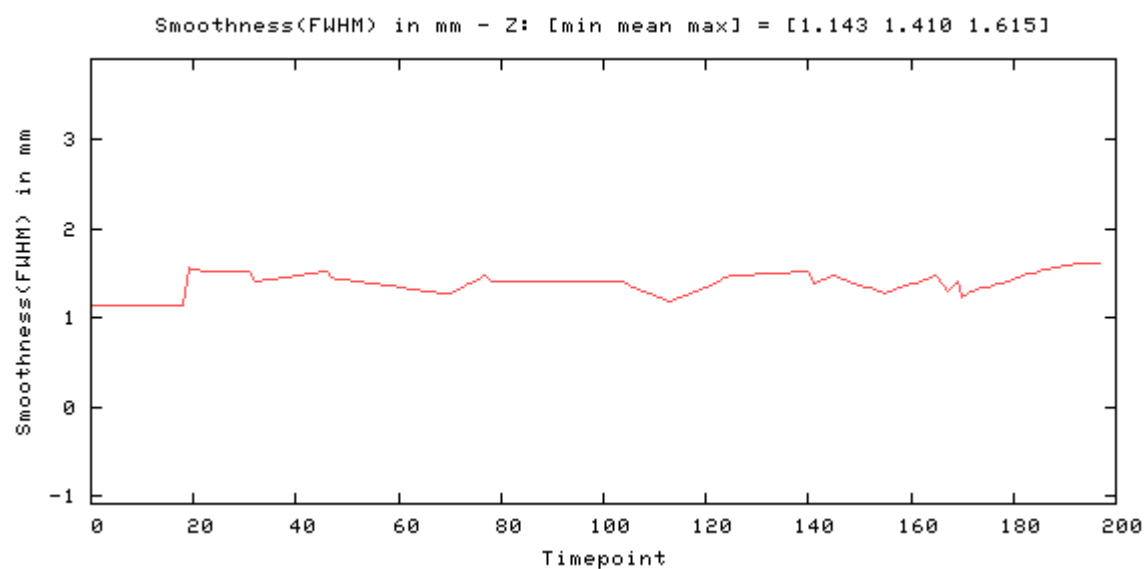


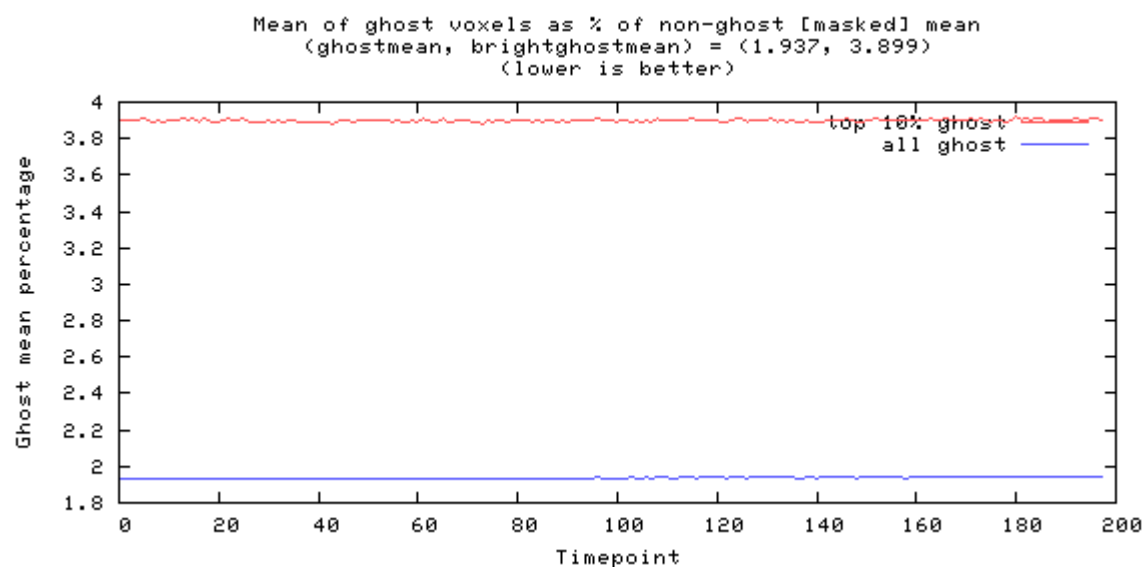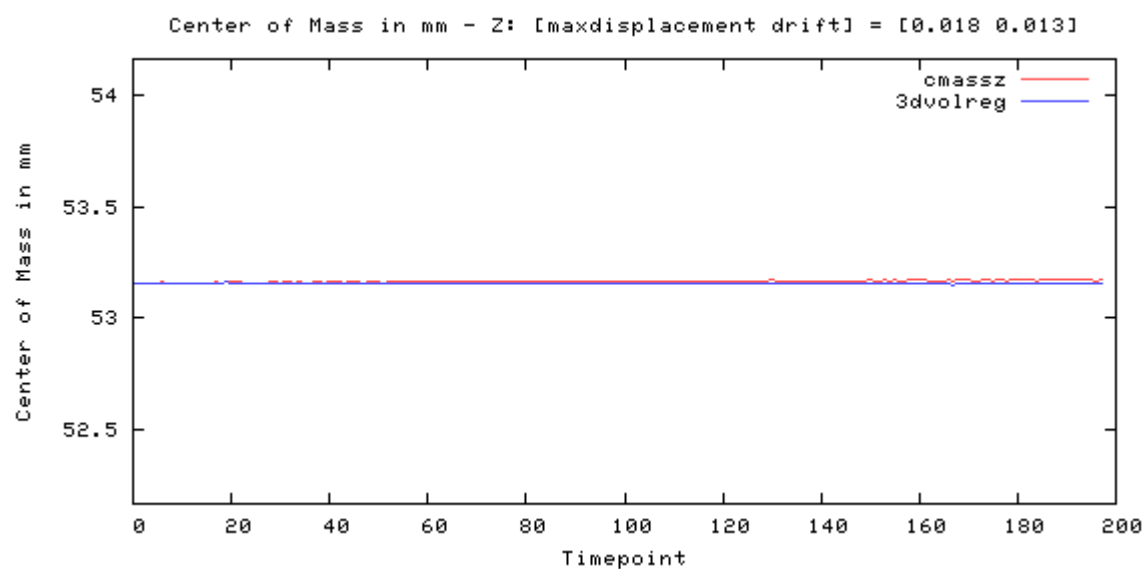

Odd-even difference image




image min: -400, image max: 298

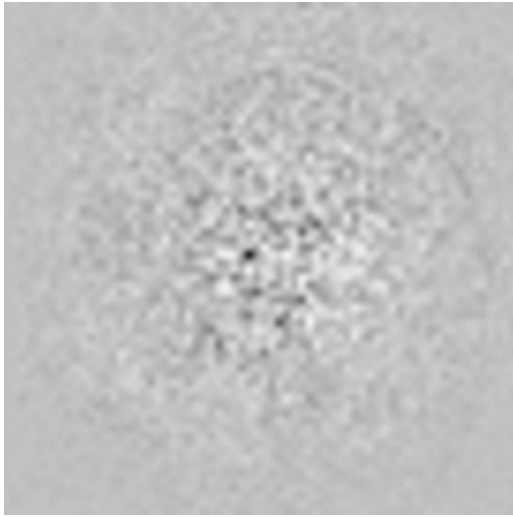

Mean image




image min: 0, image max: 1209.81

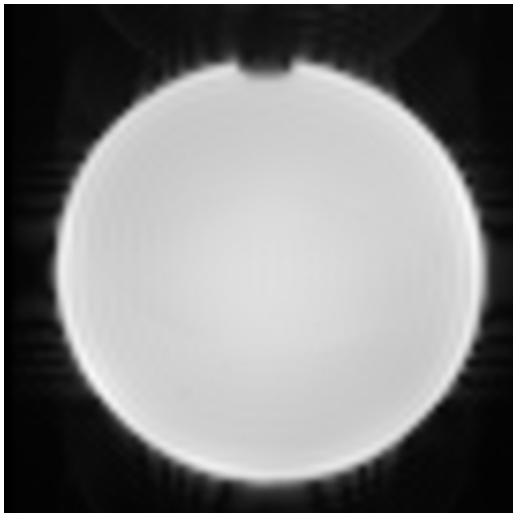

Standard Deviation image




image min: 0, image max: 9.57218

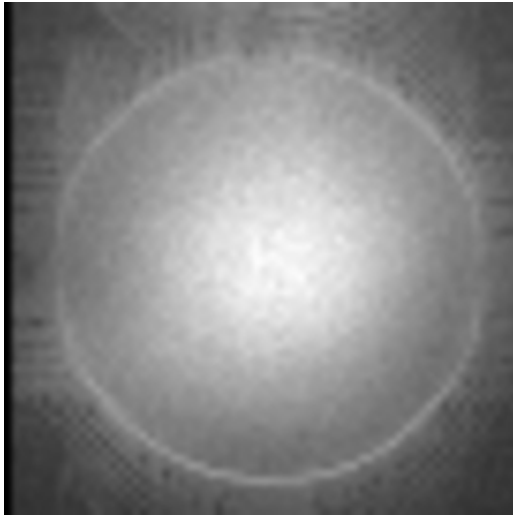

SFNR image



image min: 0, image max: 605.598

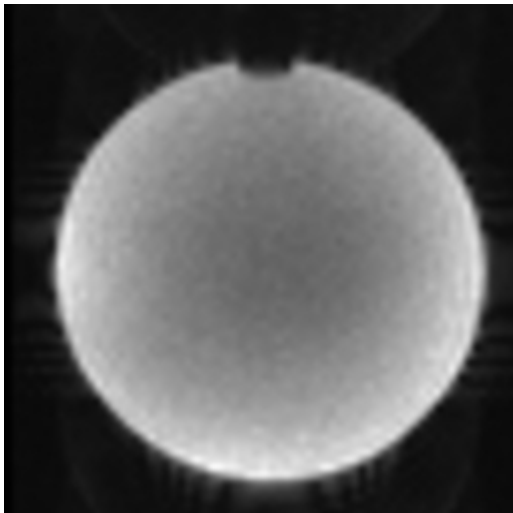

Acquisition parameters

dimensions	68x68x36x200
spacing	3mmx3mmx3mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	acdebd74-b79b-4f
studyid	Stability_20260615
seriesnumber	21
runnumber	1
magneticfield	3

description	ses-32ch_func-bold_acq-siemensstd_task-tsnr
scandate	2026-06-15
scantime	10:04:55
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	3
frequencydirection	1
receivecoilname	Head_32
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-32ch_func-bold_acq-siemensstd_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260615100453.237500
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	3
echotrainlength	33
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	PRODUCT
acquisitionduration	407.92500000000001

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	36
rows	68
columns	68
bitsallocated	16
bitsstored	12
highbit	11
lossyimagecompression	00