






Phantom QA results for /fbirn_out/ses-20ch_func-bold_acq-siemensstd_task-tsnr.xml, slice 17




old_acq-siemensstd_task-tsnr.xml [percent fluct (trend removed), drift, driftfit]


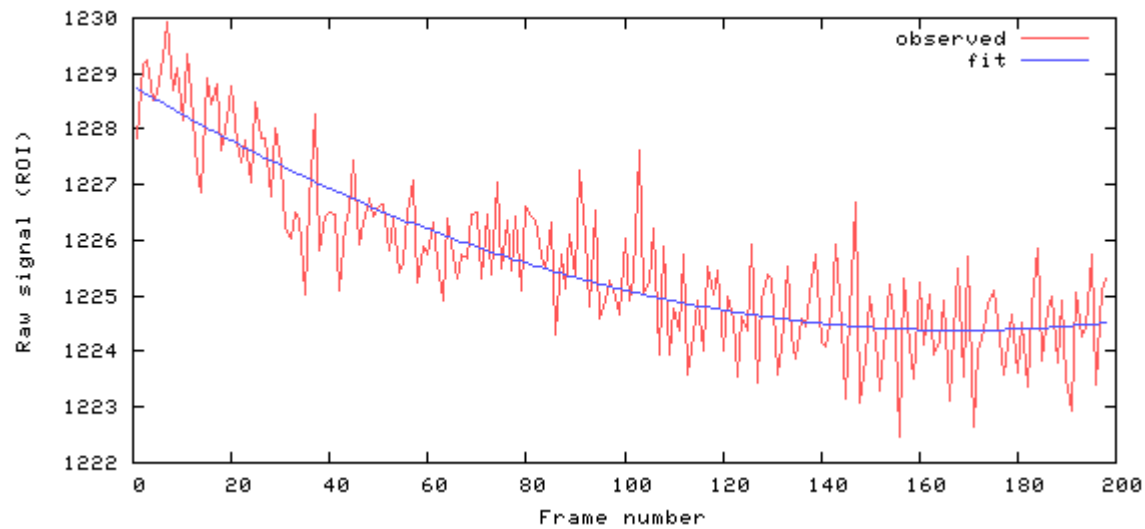


[mean, SNR, SFNR] = [1225.6 152.6 155.8]

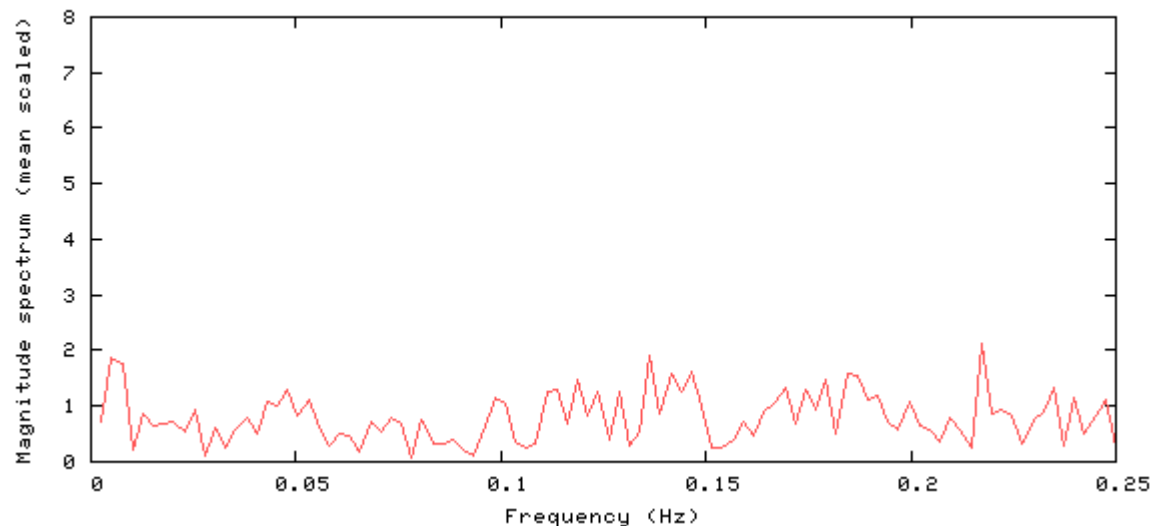

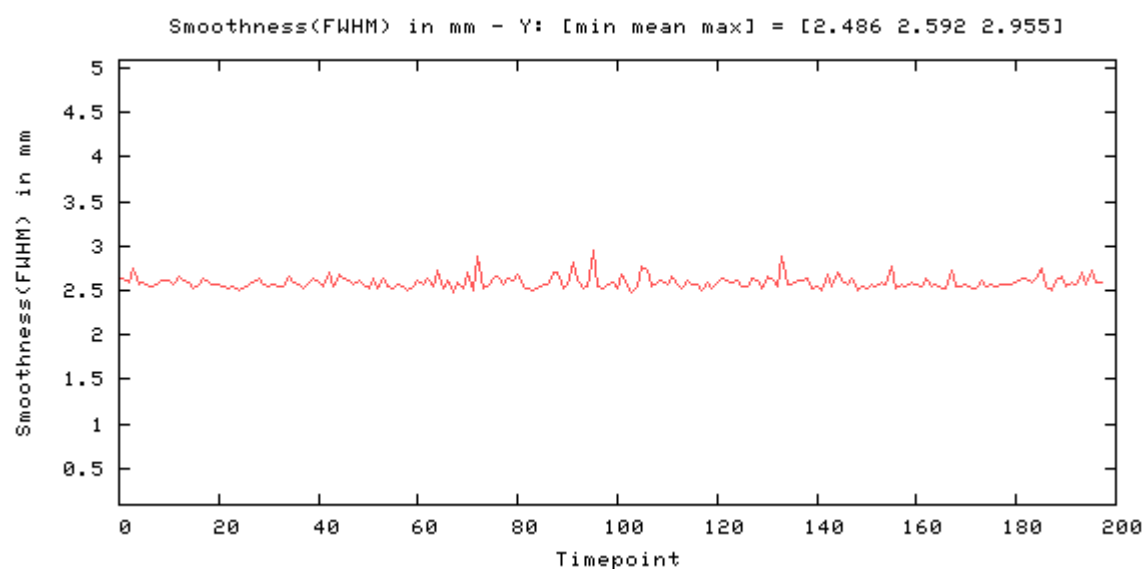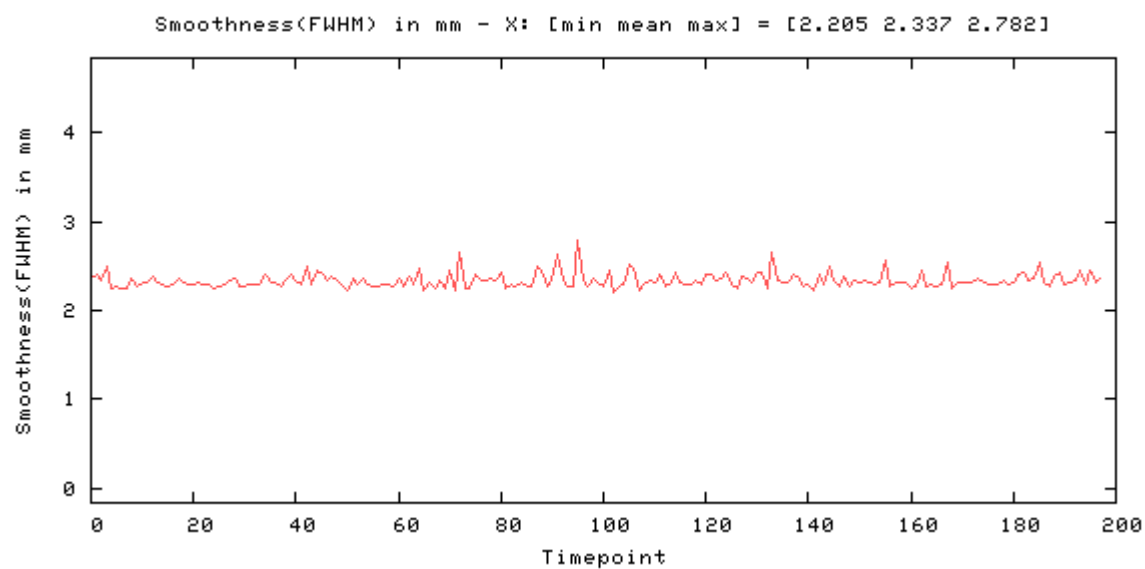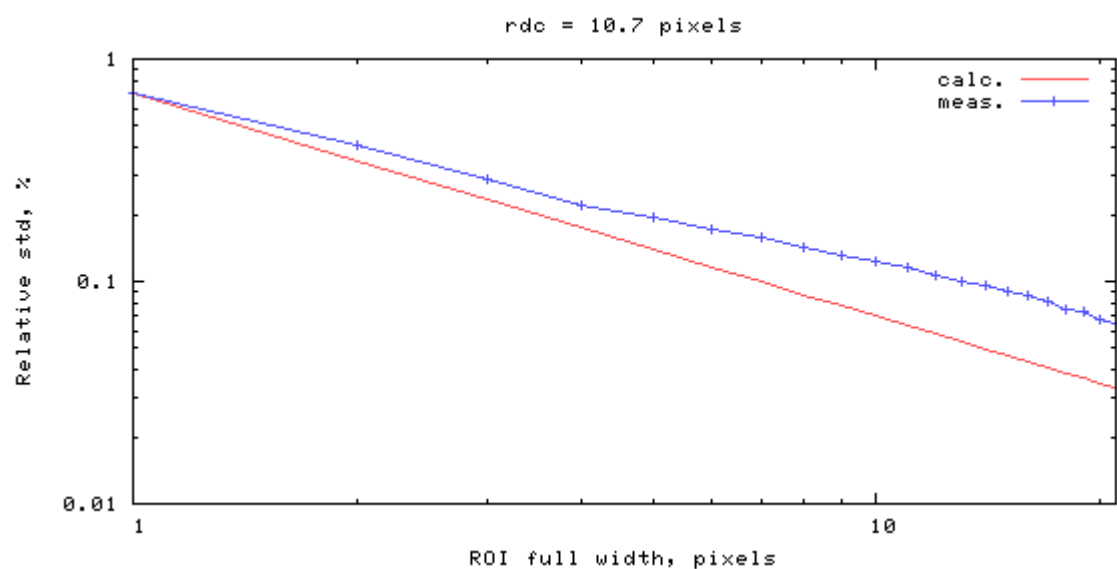

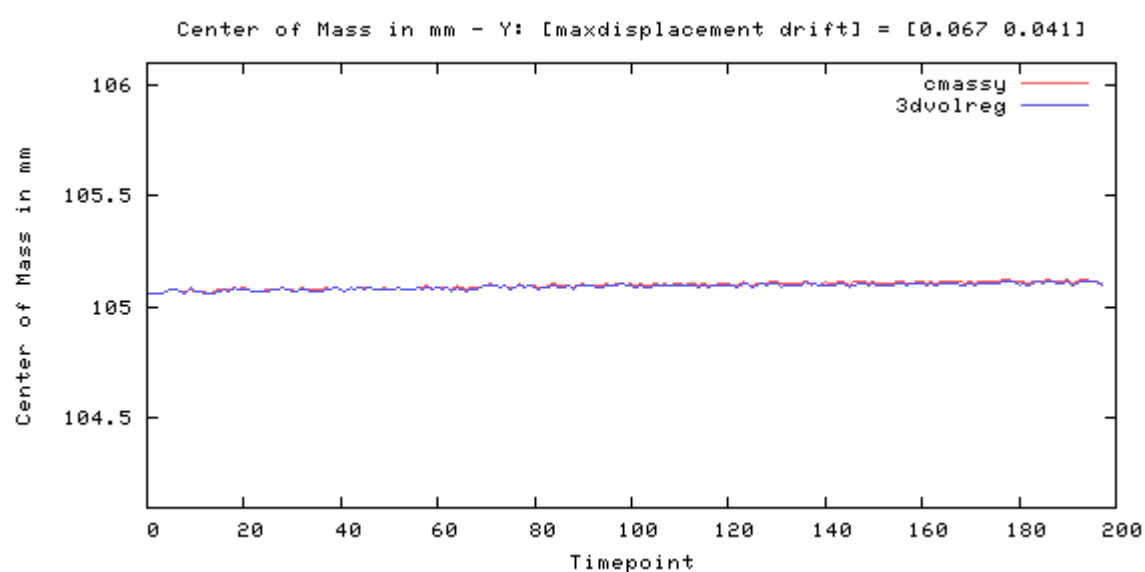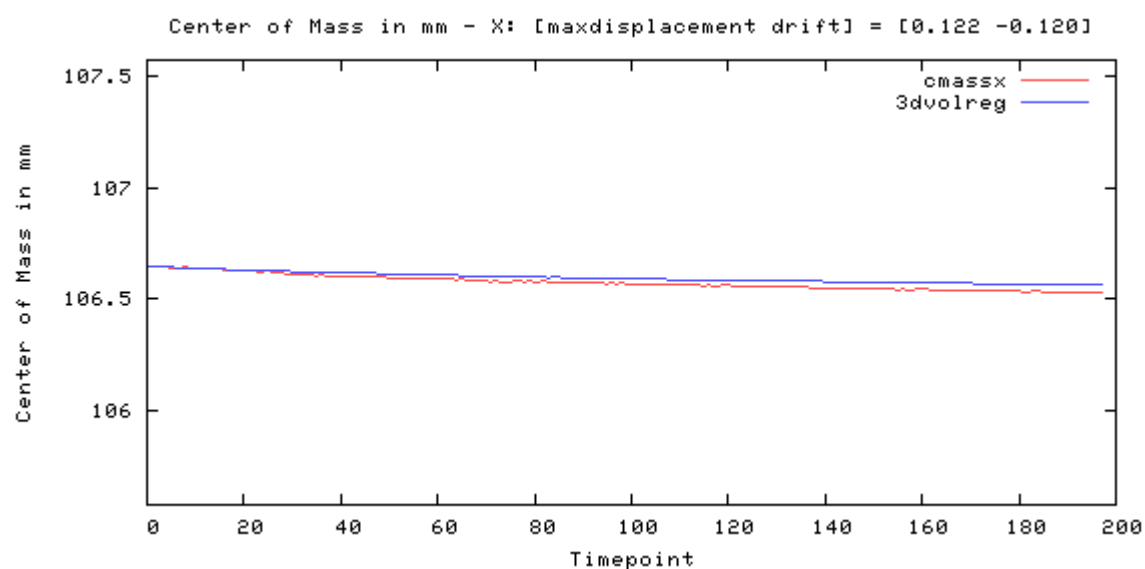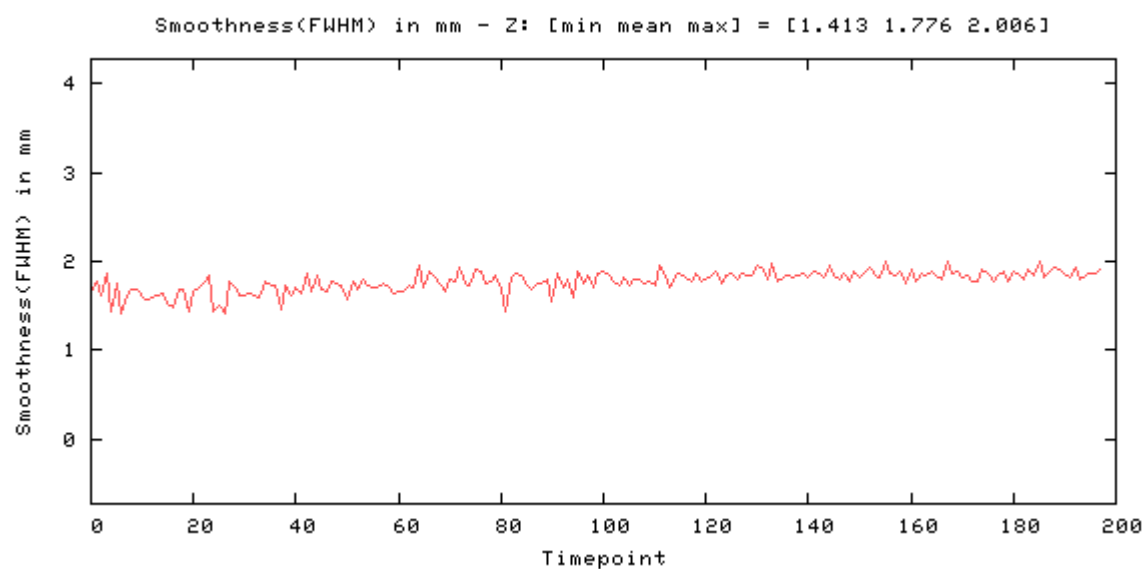


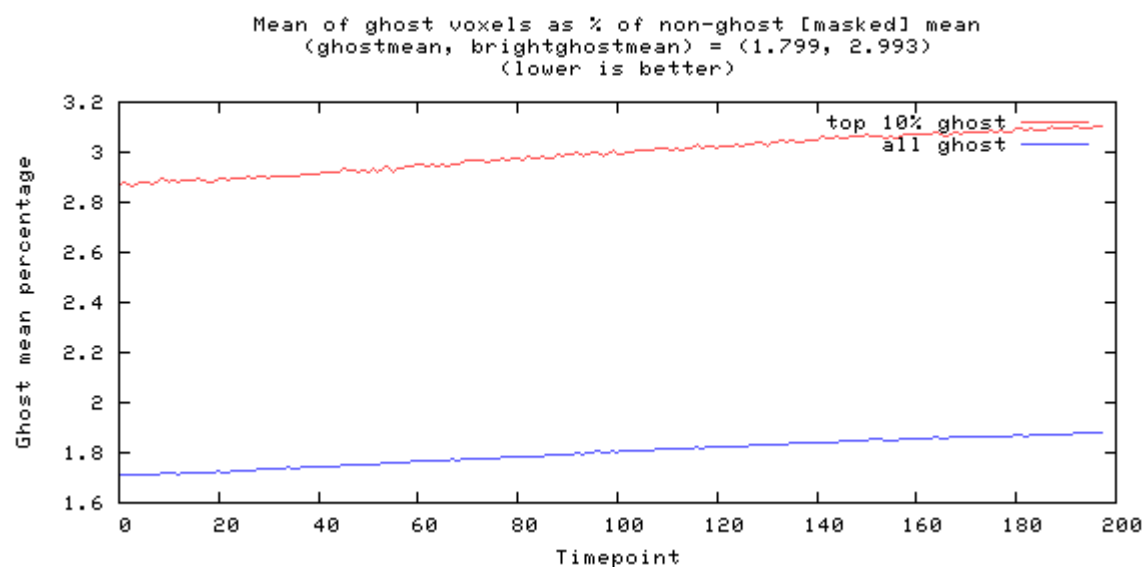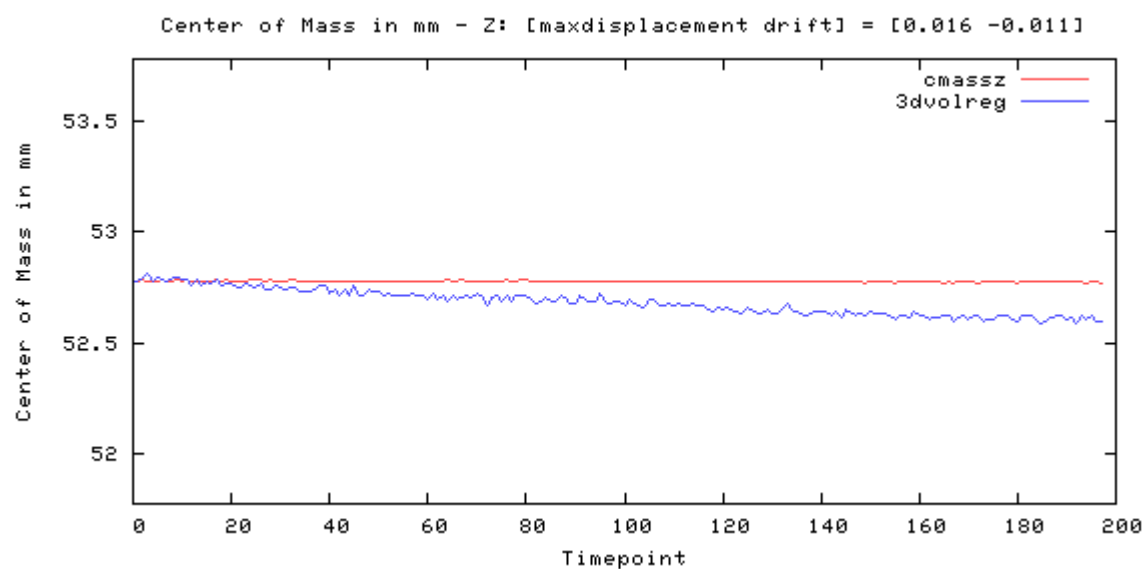

Odd-even difference image



image min: -313, image max: 383

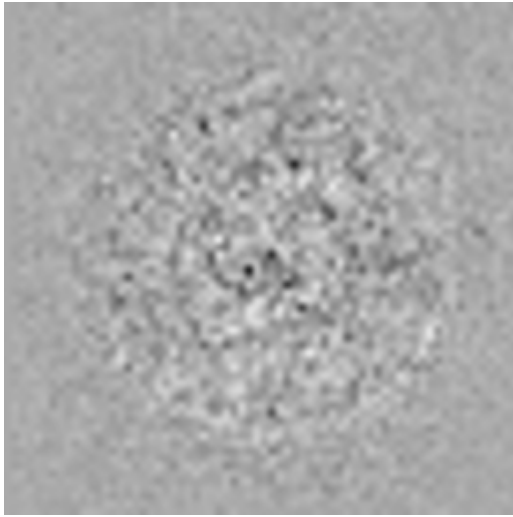

Mean image




image min: 0, image max: 1471.31



Standard Deviation image




image min: 0, image max: 10.1763

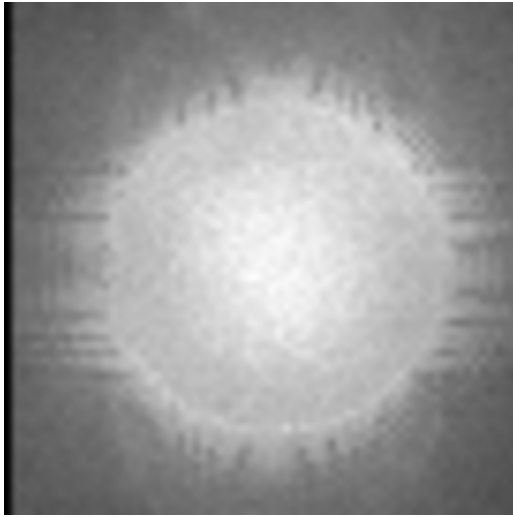

SFNR image



image min: 0, image max: 299.862

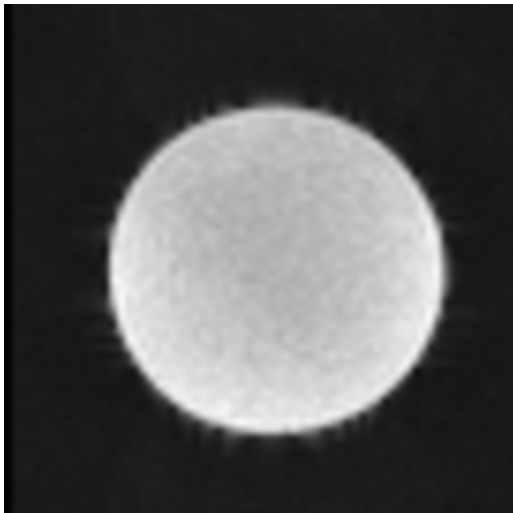

Acquisition parameters

dimensions	68x68x36x200
spacing	3mmx3mmx3mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	2743d5d9-e47d-45
studyid	Stability_20260831
seriesnumber	24
runnumber	1
magneticfield	3

description	ses-20ch_func-bold_acq-siemensstd_task-tsnr
scandate	2026-08-31
scantime	10:12:09
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	3
frequencydirection	1
receivecoilname	HeadNeck_20_TCS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-20ch_func-bold_acq-siemensstd_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260831101207.567500
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	3
echotrainlength	33
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	PRODUCT
acquisitionduration	407.92500000000001

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	36
rows	68
columns	68
bitsallocated	16
bitsstored	12
highbit	11
lossyimagecompression	00