






Phantom QA results for / fbirn_out/ ses64ch_funcbold_acqsiemensstd_task slice 17




```
ld acqsiemensstd tasktsnr 4.xml [percent fluct (trend removed), drift, driftfit] :
```


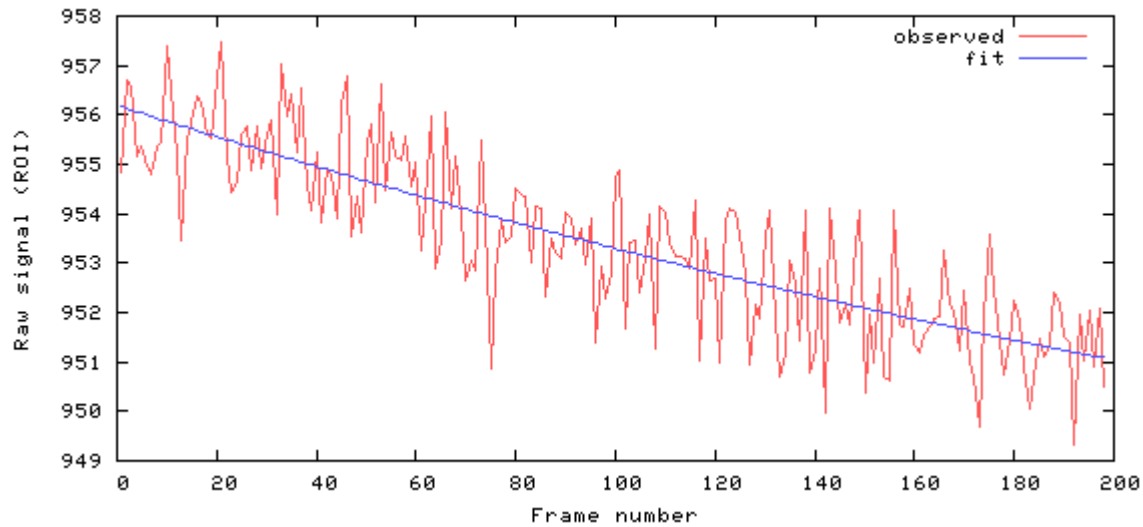


```
[mean, SNR, SFNR] = [953.4 111.2 111.9]
```

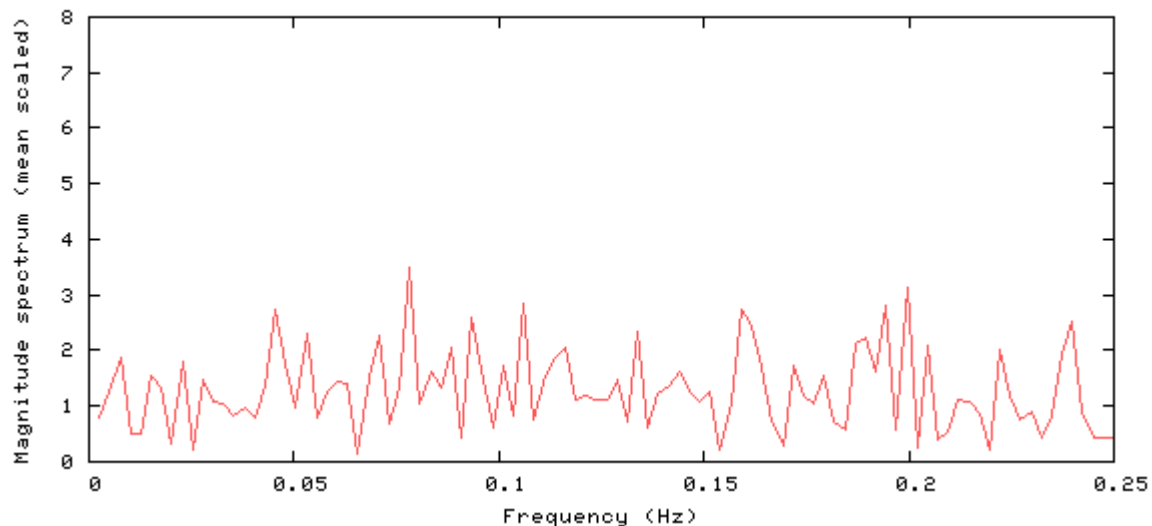

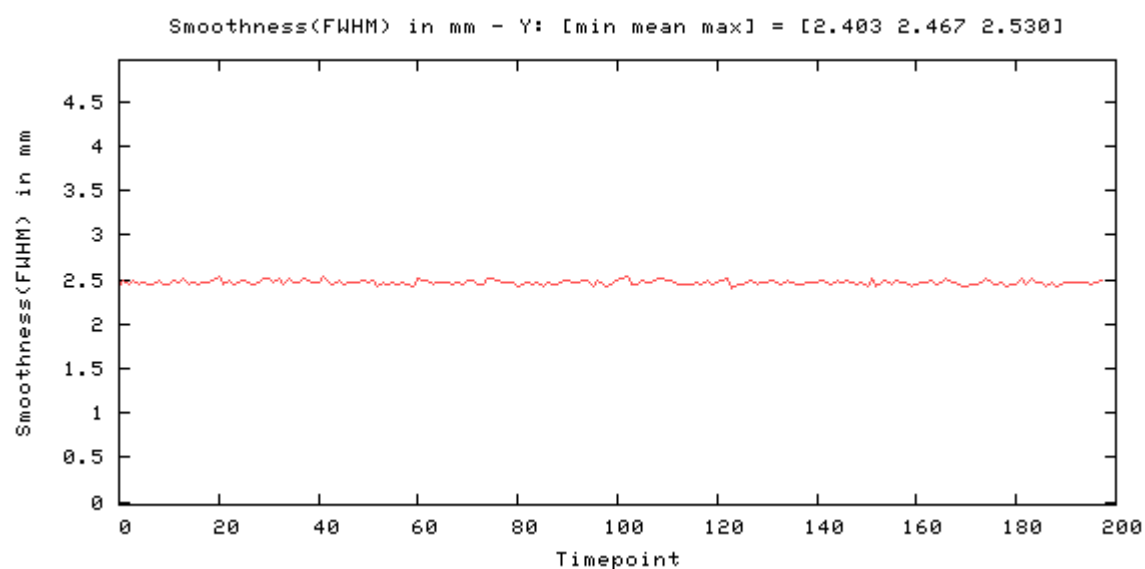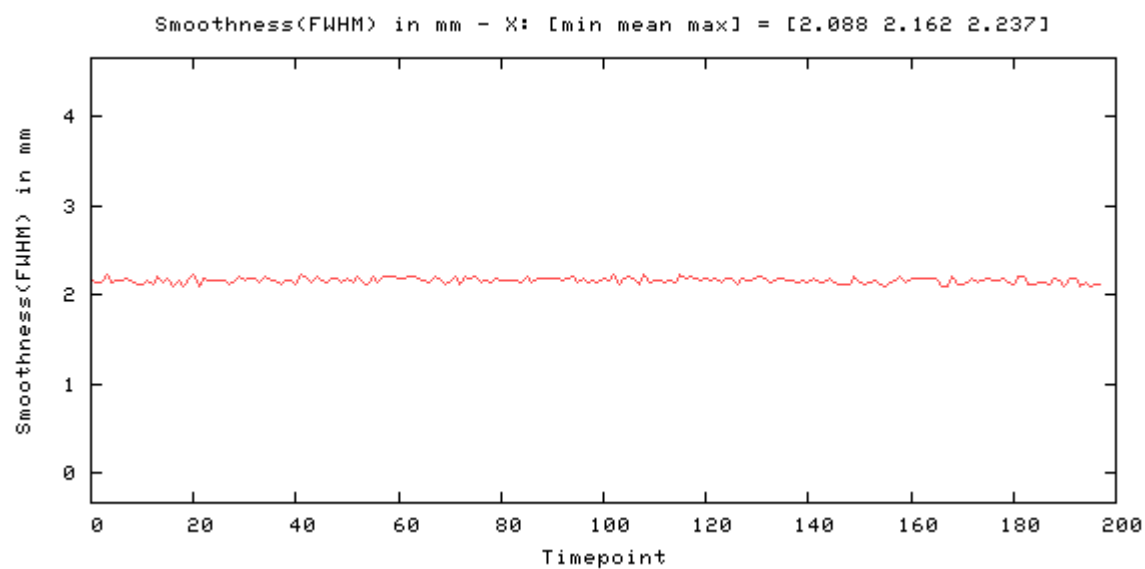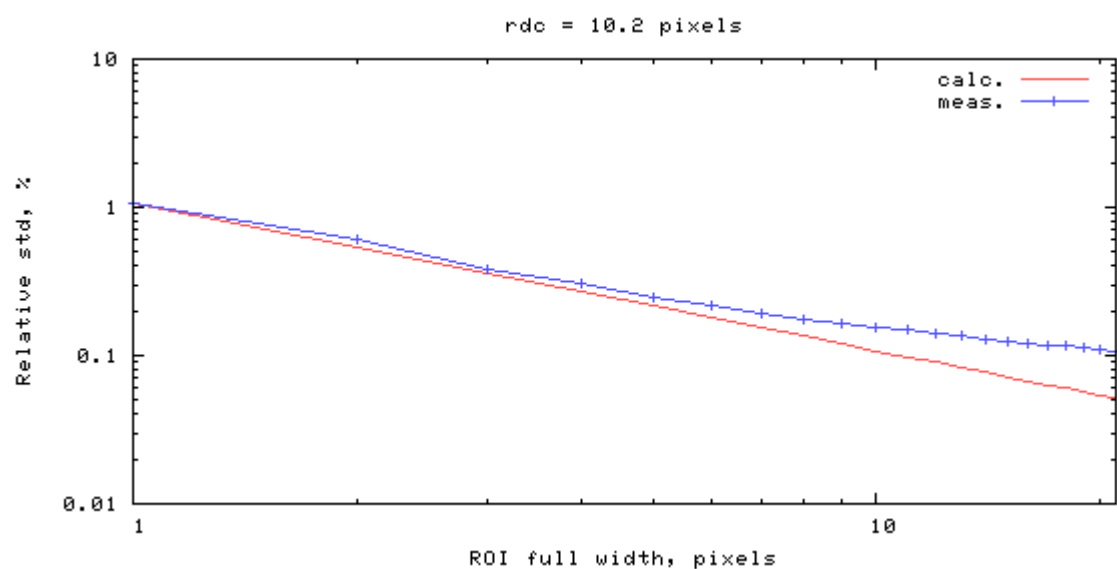

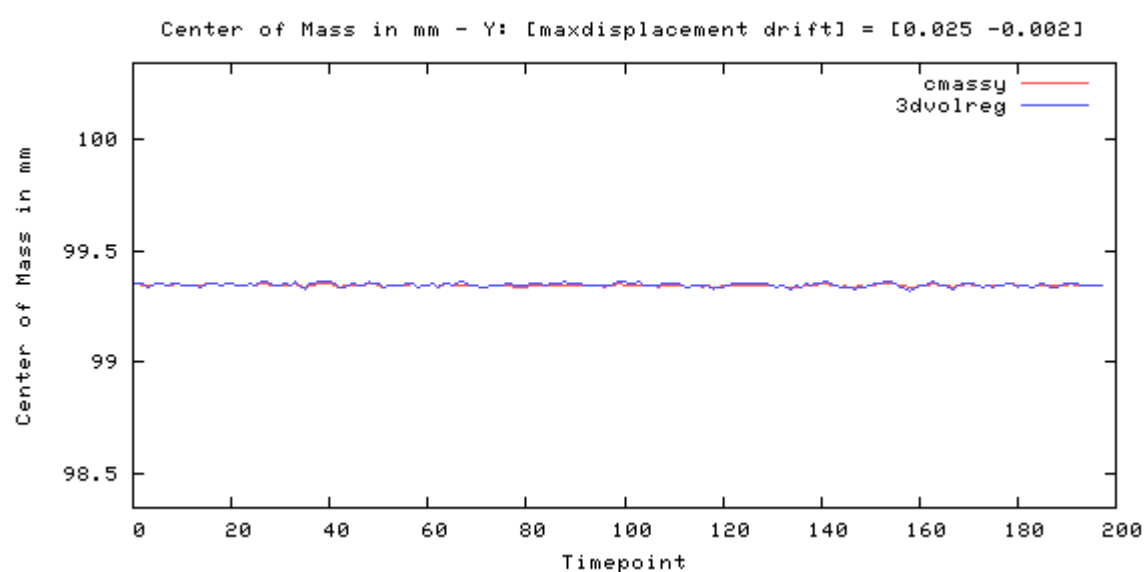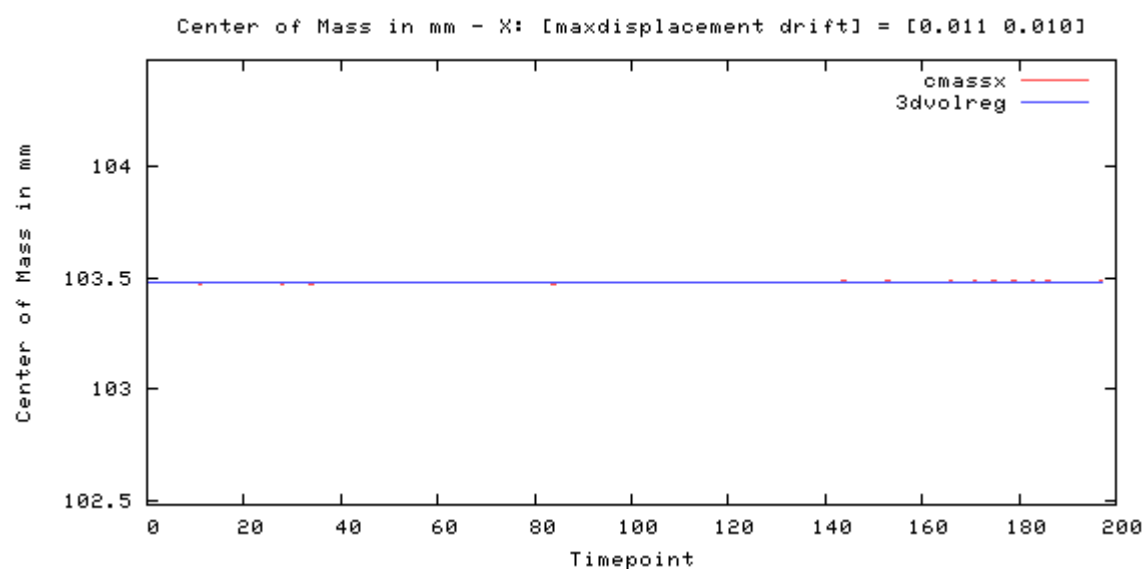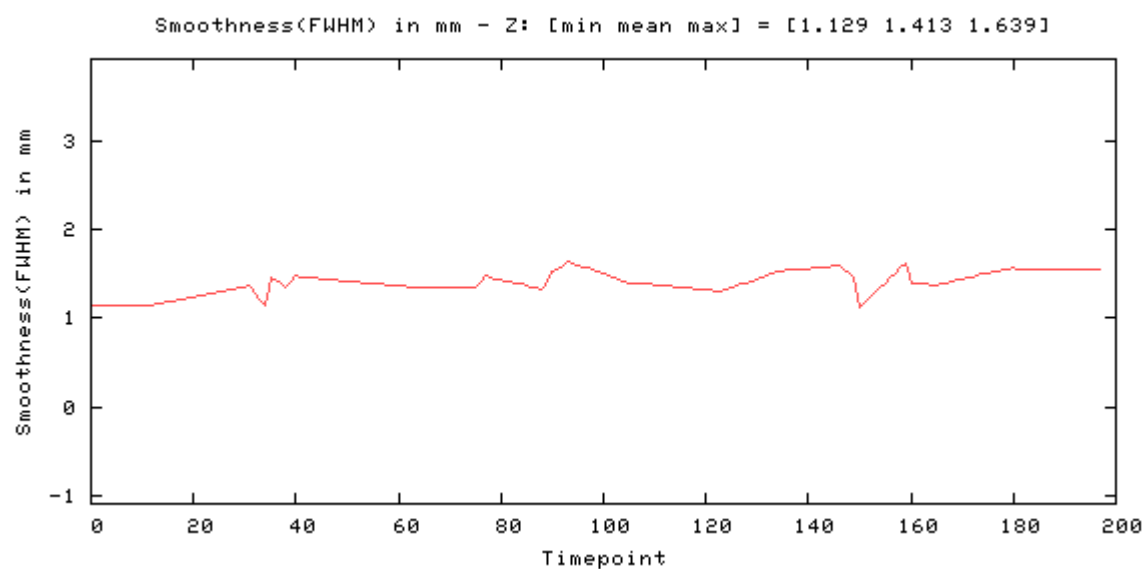


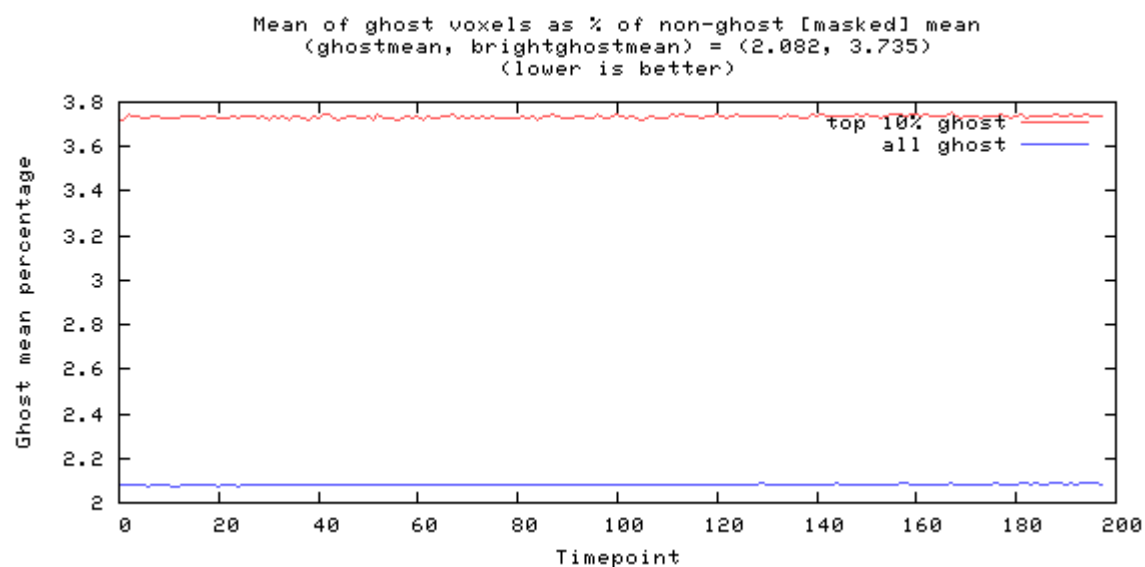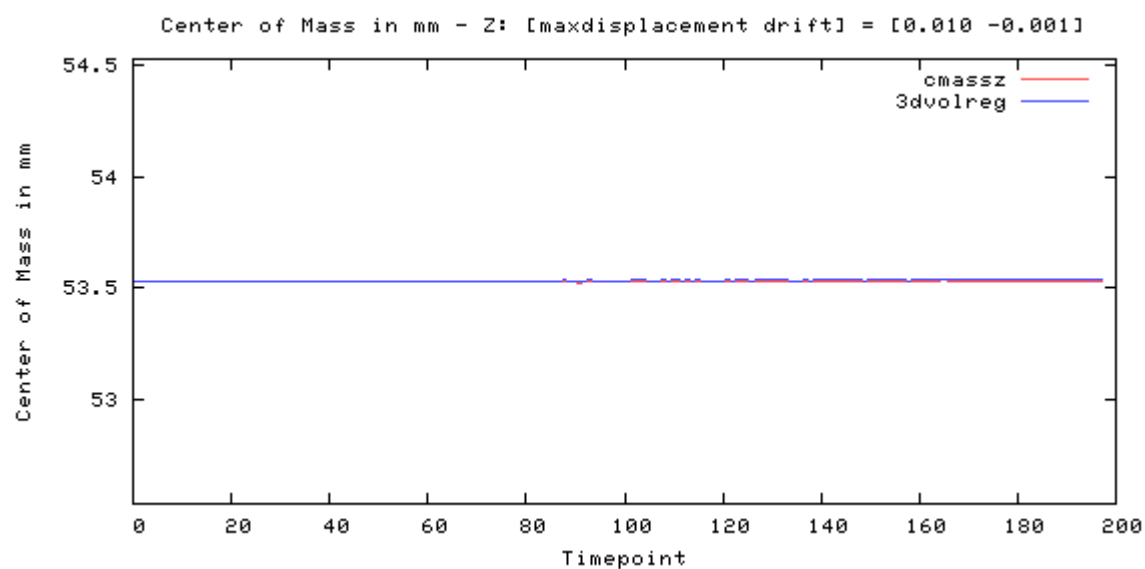

Odd-even difference image




image min: -303, image max: 358

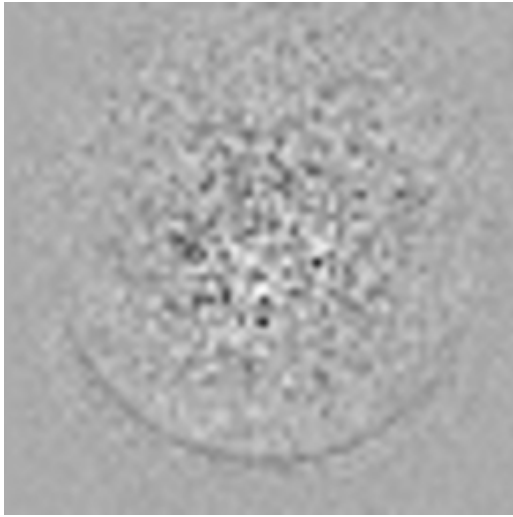

Mean image




image min: 0, image max: 1323.08

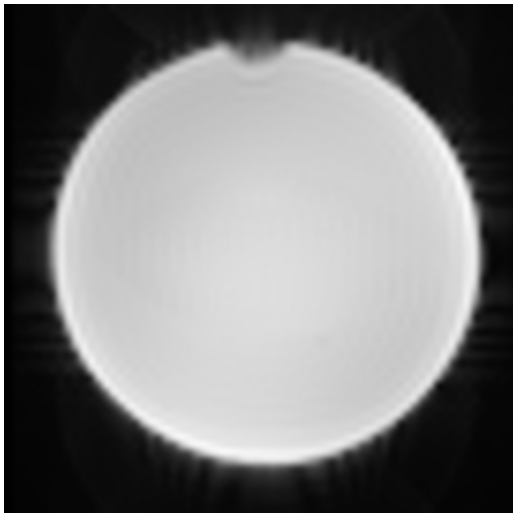

Standard Deviation image




image min: 0, image max: 12.3644

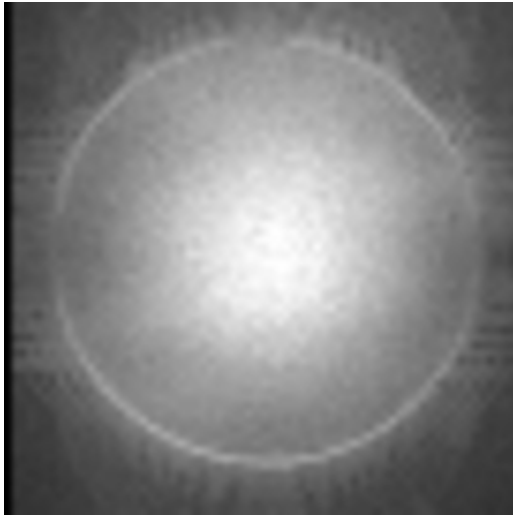

SFNR image



image min: 0, image max: 633.652

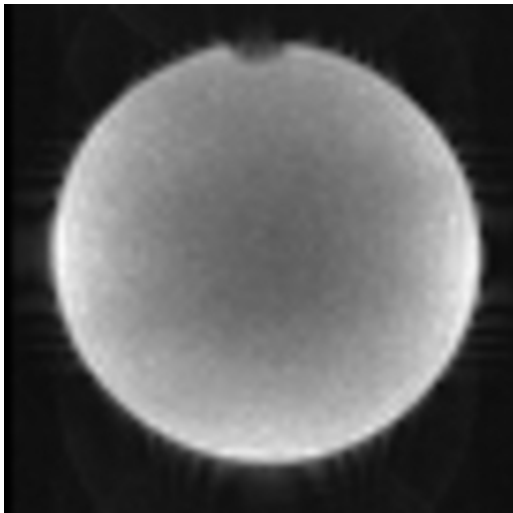

Acquisition parameters

dimensions	68x68x36x200
spacing	3mmx3mmx3mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	364e27d6-dd69-4a
studyid	Stability_20260427
seriesnumber	4
runnumber	1
magneticfield	3

description	ses-64ch_func-bold_acq-siemensstd_task-tsnr
scandate	2026-04-27
scantime	09:25:27
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	3
frequencydirection	1
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-64ch_func-bold_acq-siemensstd_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260427092525.547500
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	3
echotrainlength	33
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	PRODUCT
acquisitionduration	407.92500000000001

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	36
rows	68
columns	68
bitsallocated	16
bitsstored	12
highbit	11
lossyimagecompression	00