






Phantom QA results for / fbirn_out/ ses20ch_funcbold_acqcmrrstd_tasktsnr slice 34




bold_acqcmrrstd_tasktsnr_21.xml [percent fluct (trend removed), drift, driftfit] :


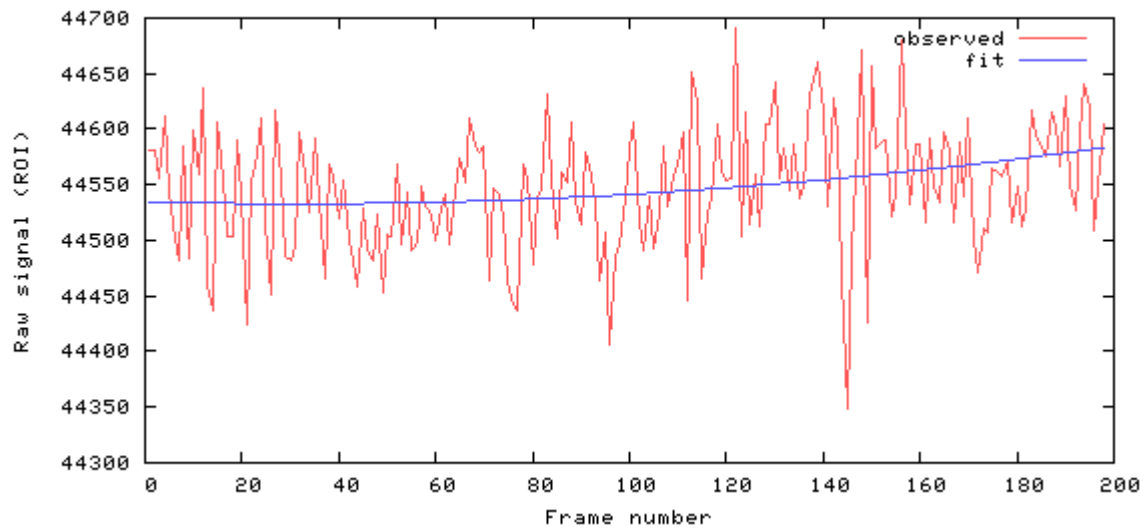


[mean, SNR, SFNR] = [44547.4 166.2 165.9]

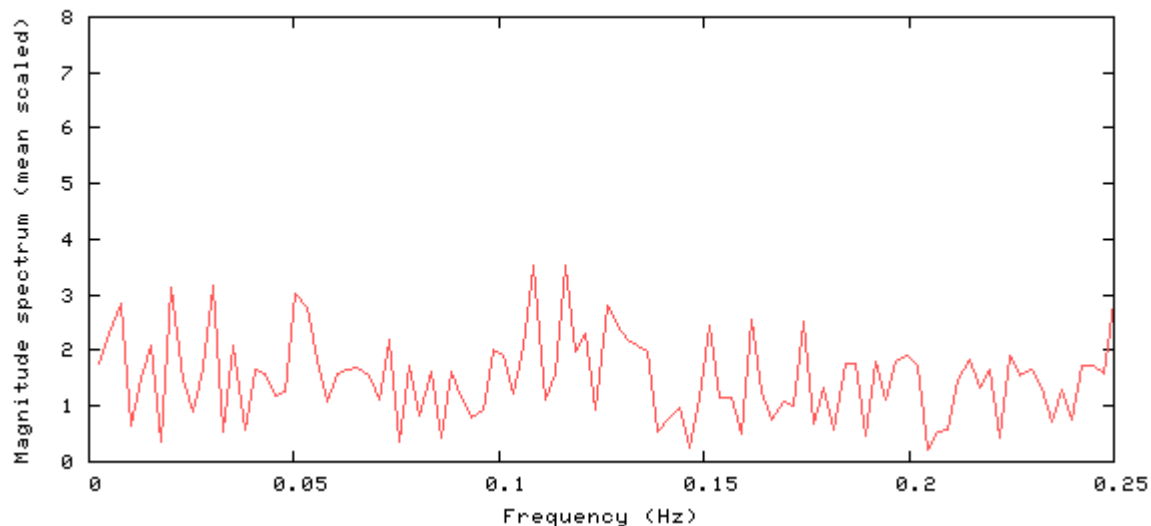

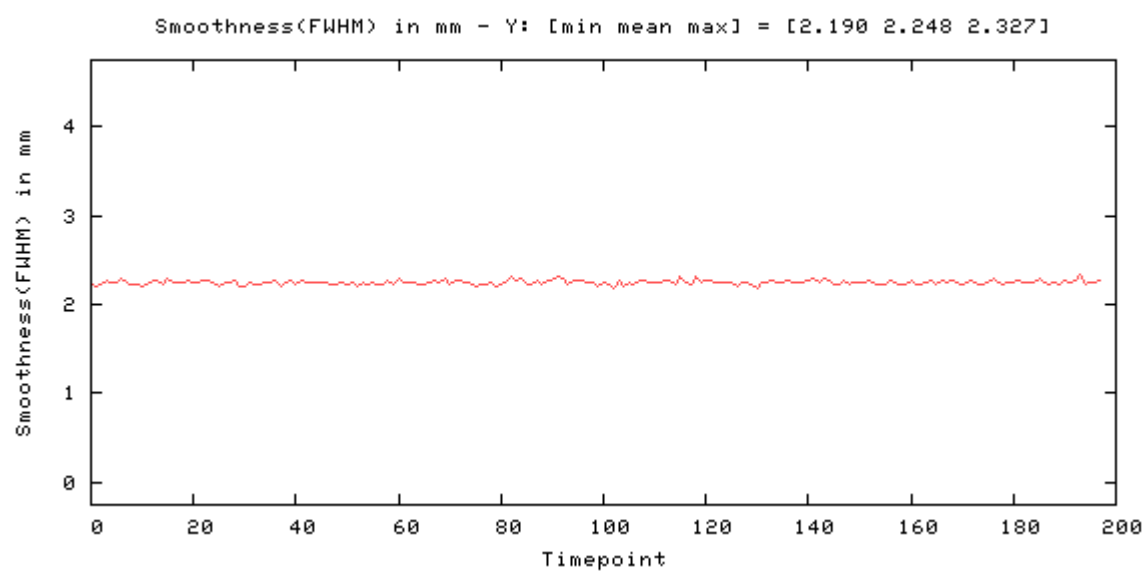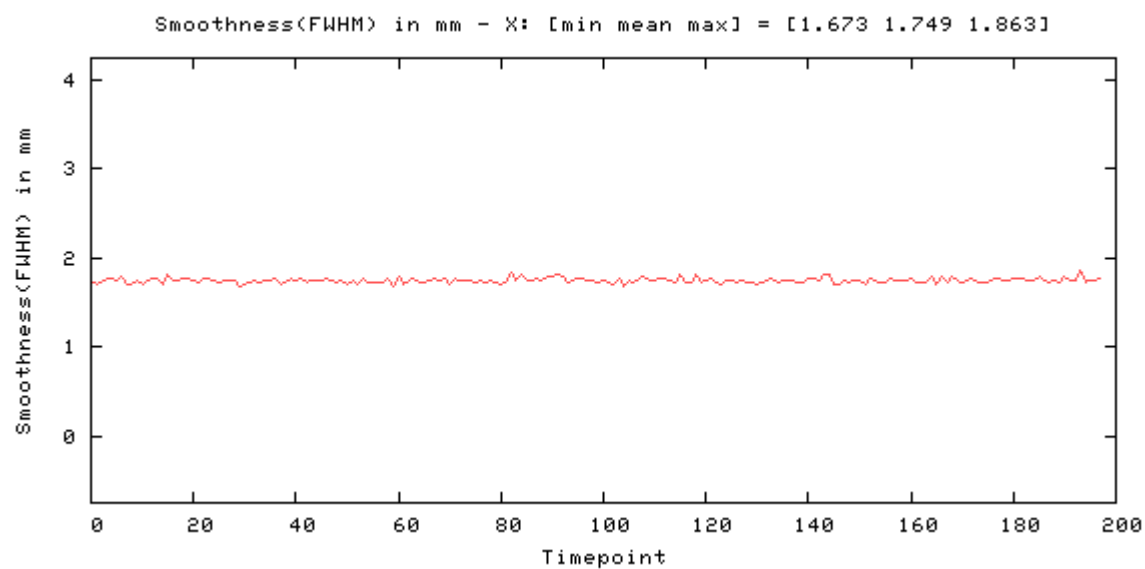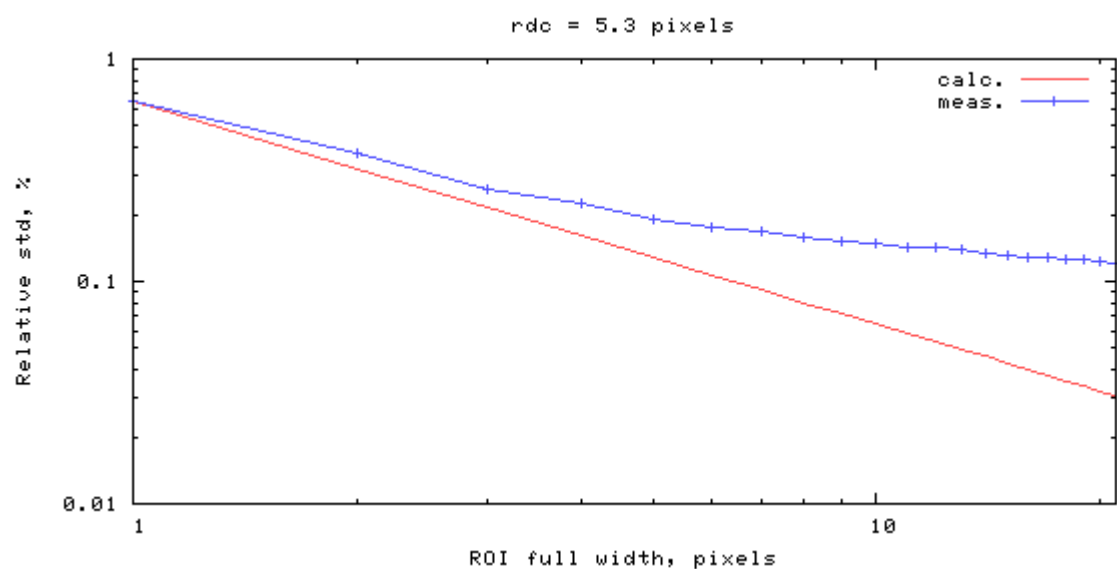

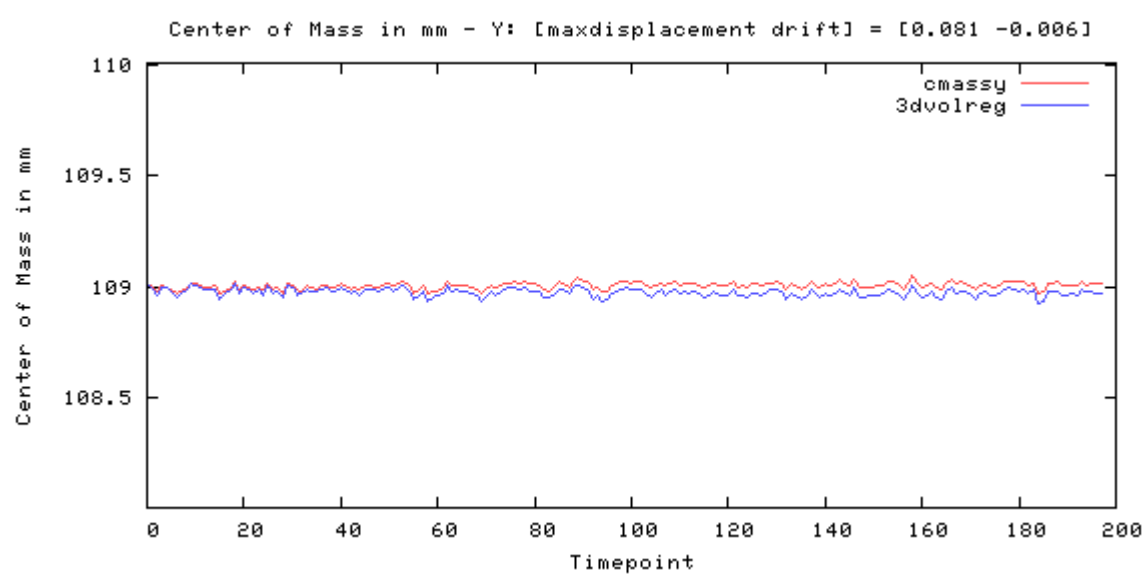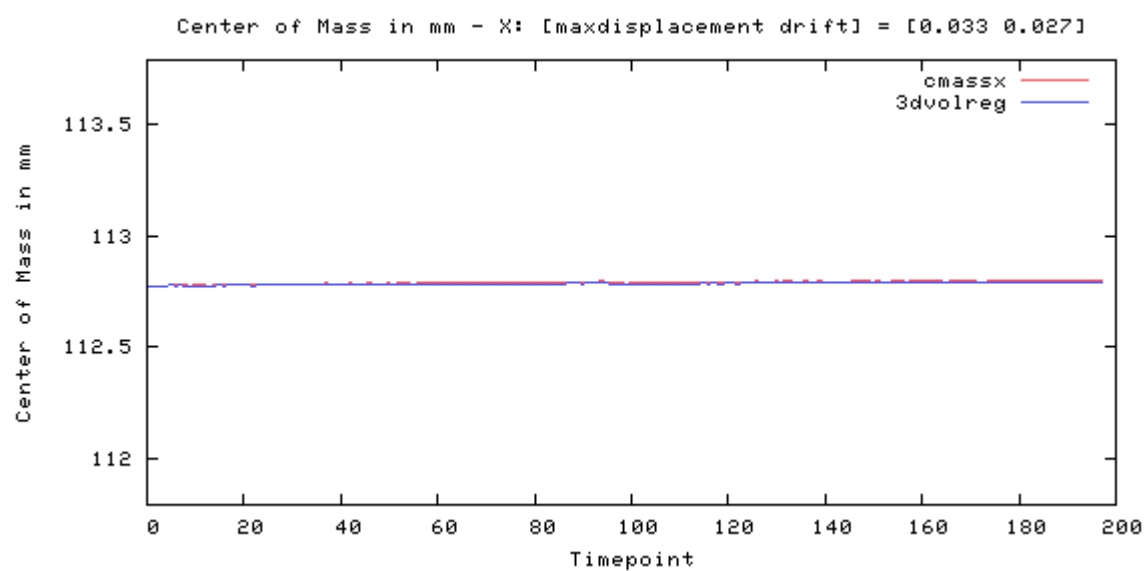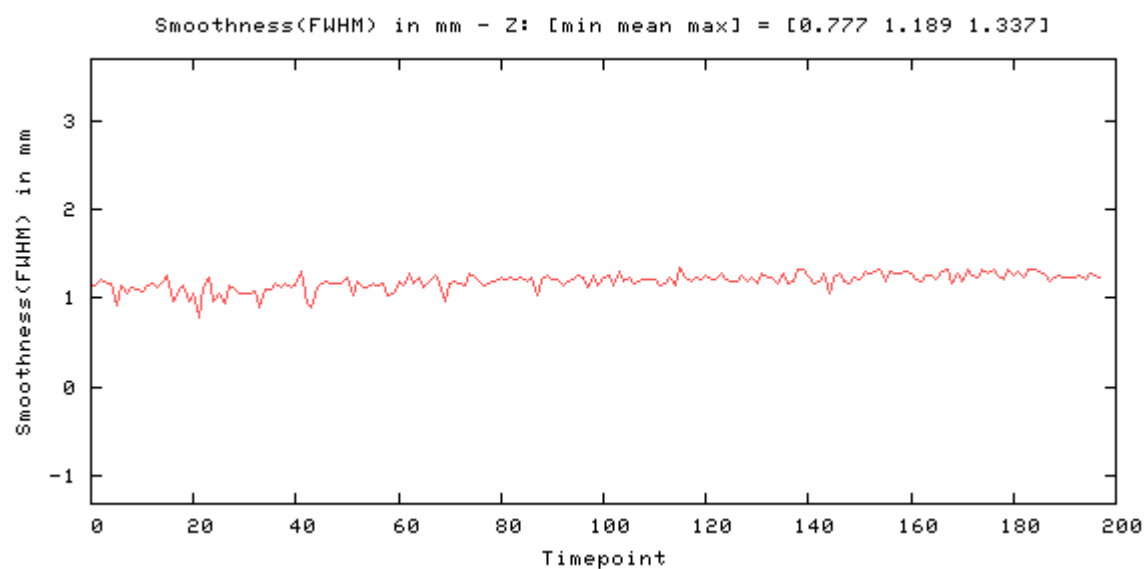

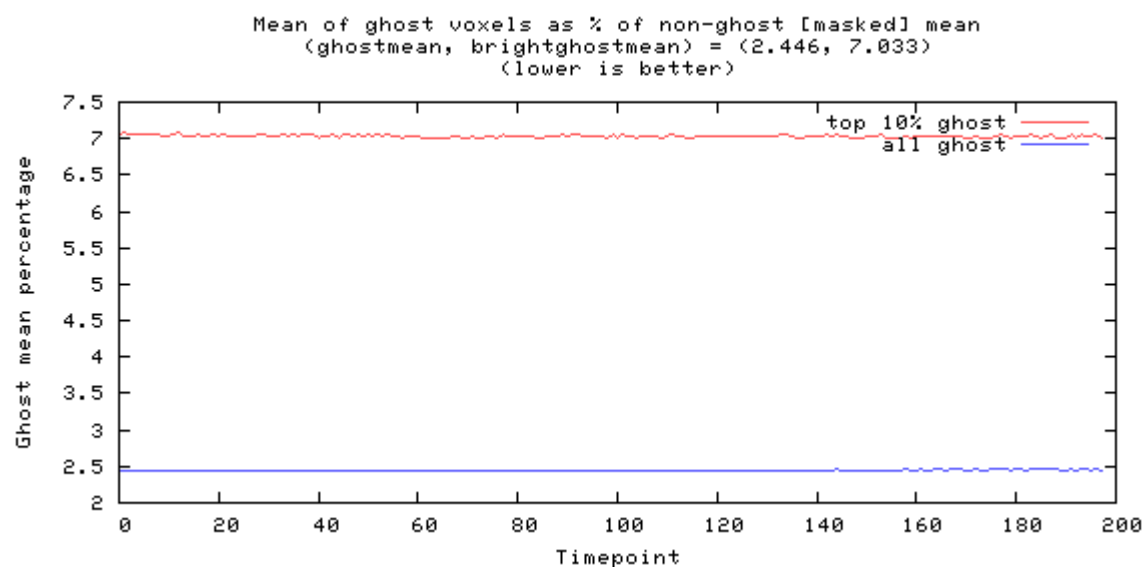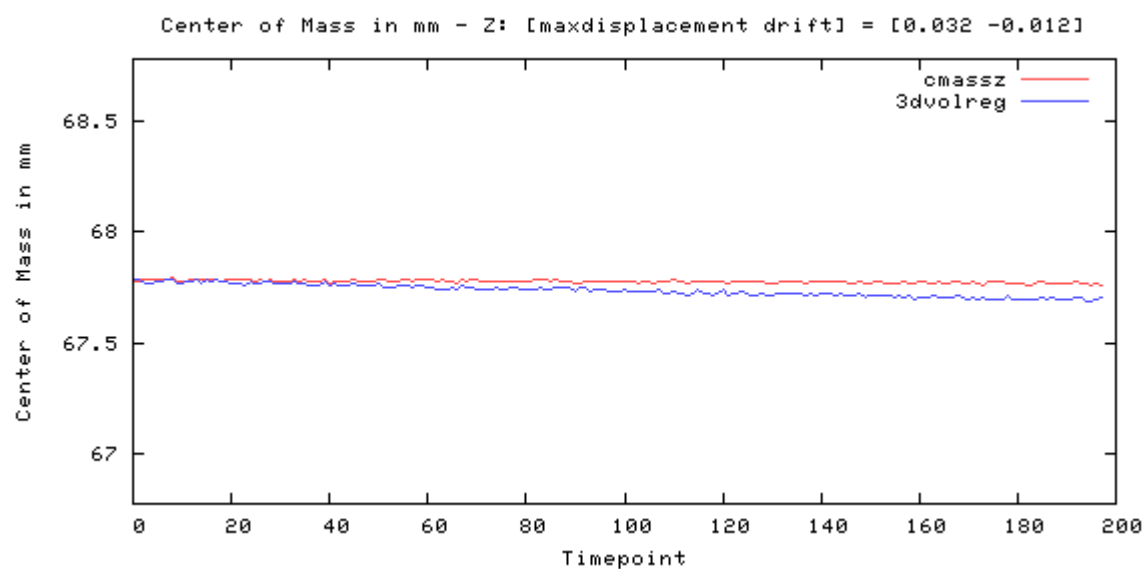


Odd-even difference image

-15039 11267




image min: -15039, image max: 11267

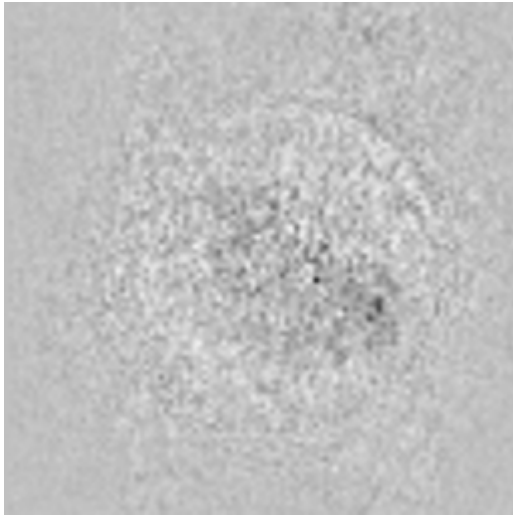

Mean image




image min: 0, image max: 52038.3

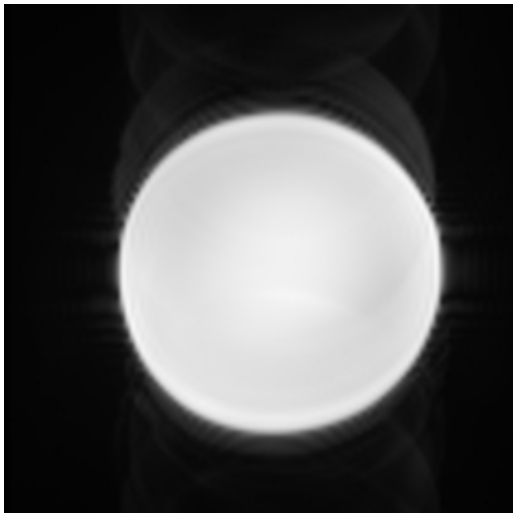

Standard Deviation image




image min: 0, image max: 449.493

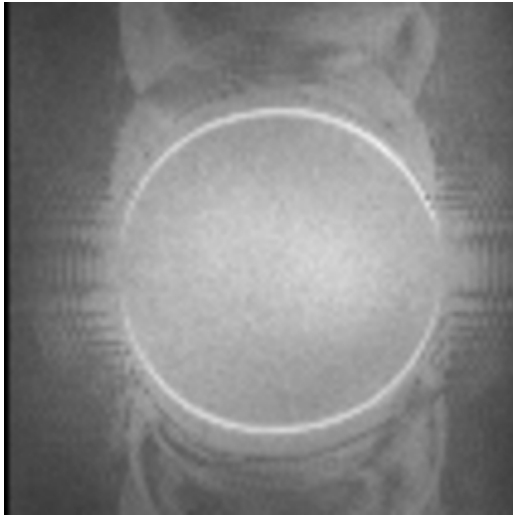

SFNR image



image min: 0, image max: 301.113

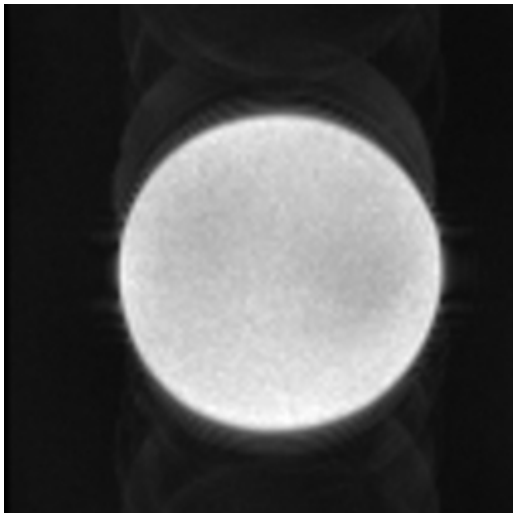

Acquisition parameters

dimensions	106x106x69x200
spacing	2mmx2mmx2mmx1[units?]
gap	0mmx0mmx0mmx0[units?]
scanner	AWP213018
psdname	
examnumber	68a520d1-425b-42
studyid	Stability_20250204
seriesnumber	21
runnumber	1-200
scanningsequence	GR/EP

sequencevariant	SK
magneticfield	3
description	ses-20ch_func-bold_acq-cmrrstd_task-tsnr
scandate	2025-02-04
scantime	12:24:25
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
acquisitionmatrix	106 106
bandwidth	192.284
receivecoilname	HeadNeck_20_TCS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
ti	
protocolname	ses-20ch_func-bold_acq-cmrrstd_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20250204122423.030000
volumetricproperties	DISTORTED
scanoptions	PFP/FS/PER
mracquisitiontype	2D
sequencename	epfid2d1_106
angioflag	N
numaverages	1
imagingfrequency	123.253102
imagednucleus	1H
spacingbetweenslices	2
numphaseencodingsteps	106
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61

triggertime	
beatrejectionflag	N
cardiacnumberofimages	1
inplane phase encoding direction	COL
variableflipangleflag	N
sar	0.01534216663899
contentqualification	RESEARCH
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
imagecomments	Not for diagnostic use, Unaliased MB3/PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00
echonumber	0
pixelbandwidth	1814