






Phantom QA results for / fbirn_out/ ses64ch_funcbold_acqsiemensbasic_ta slice 14




ld acqsiemensbasic tasktsnr 6.xml [percent fluct (trend removed), drift, driftfit:


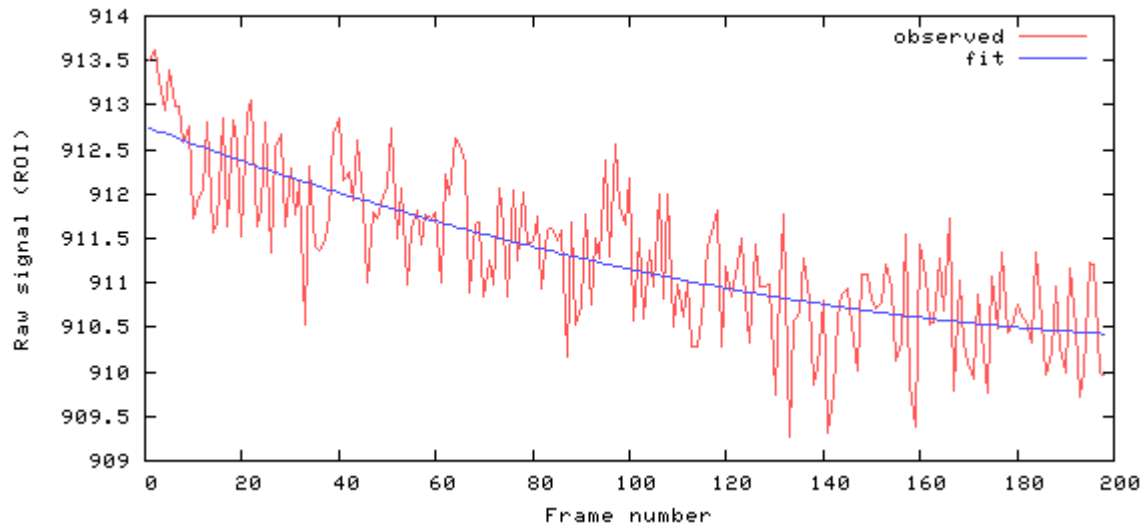


[mean, SNR, SFNR] = [911.3 159.3 169.2]

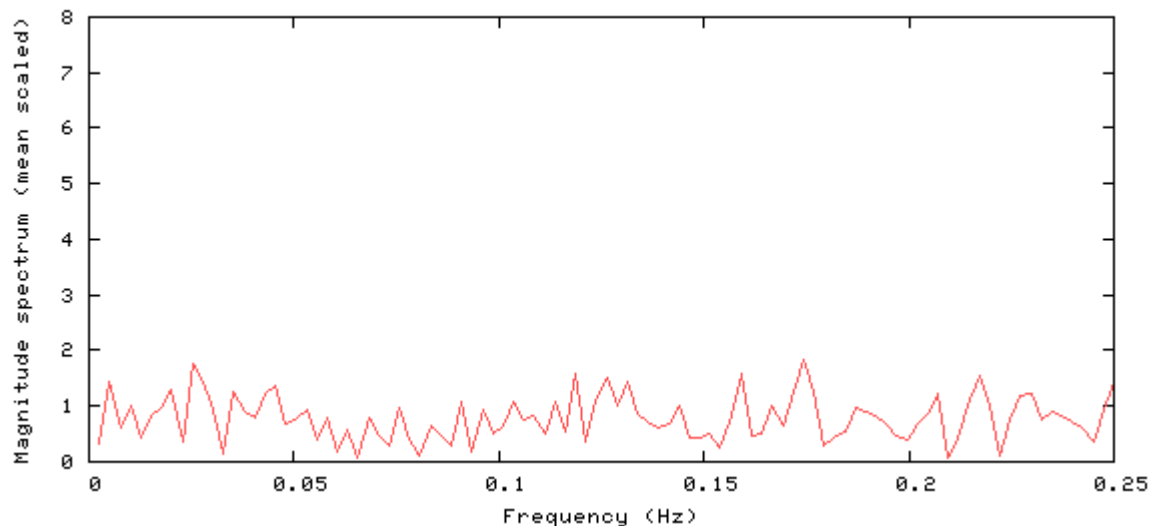

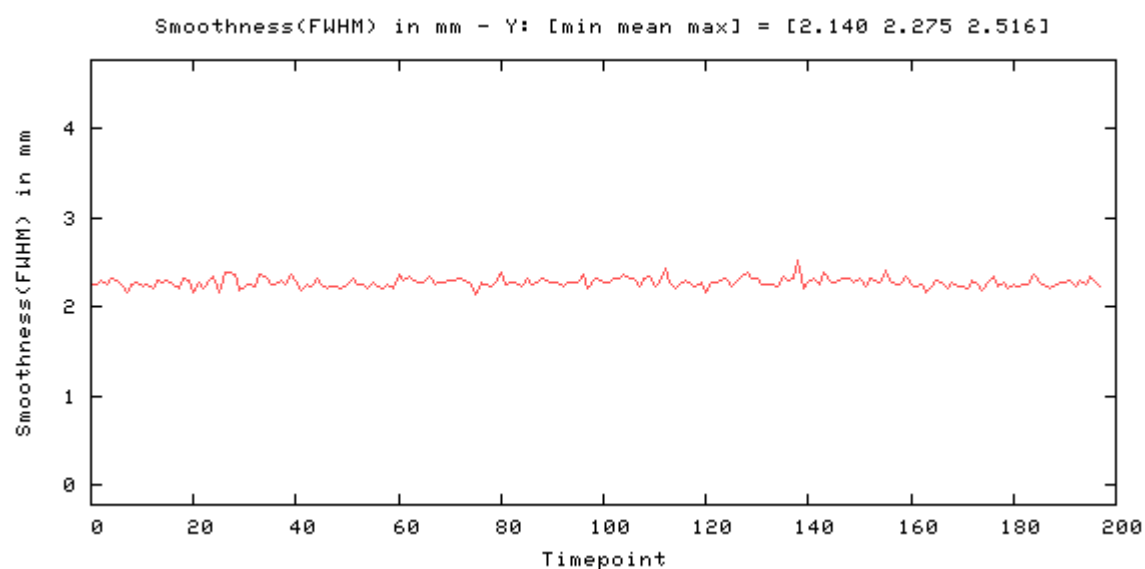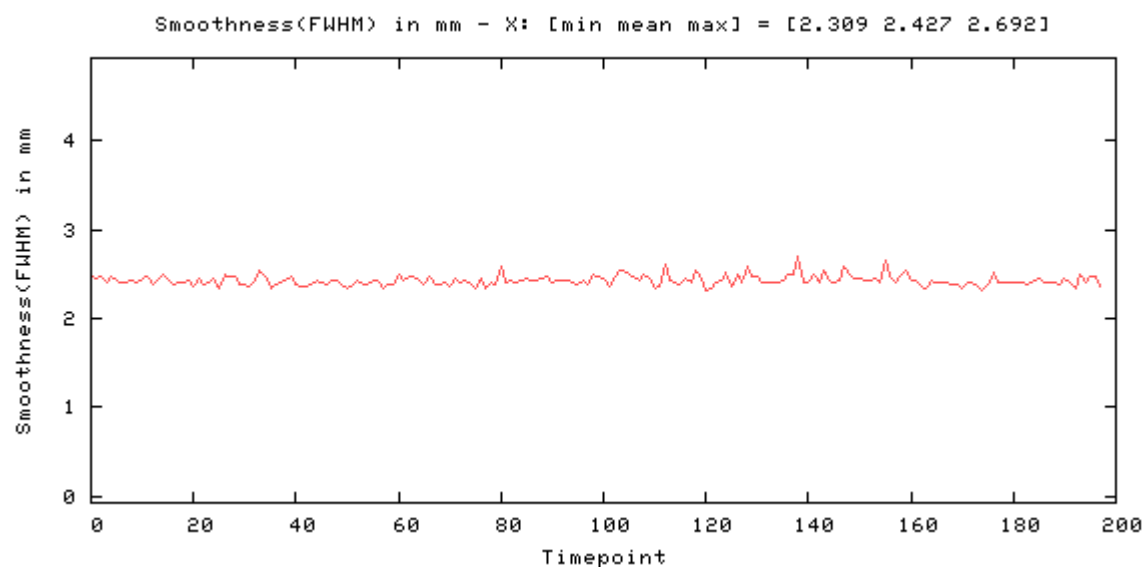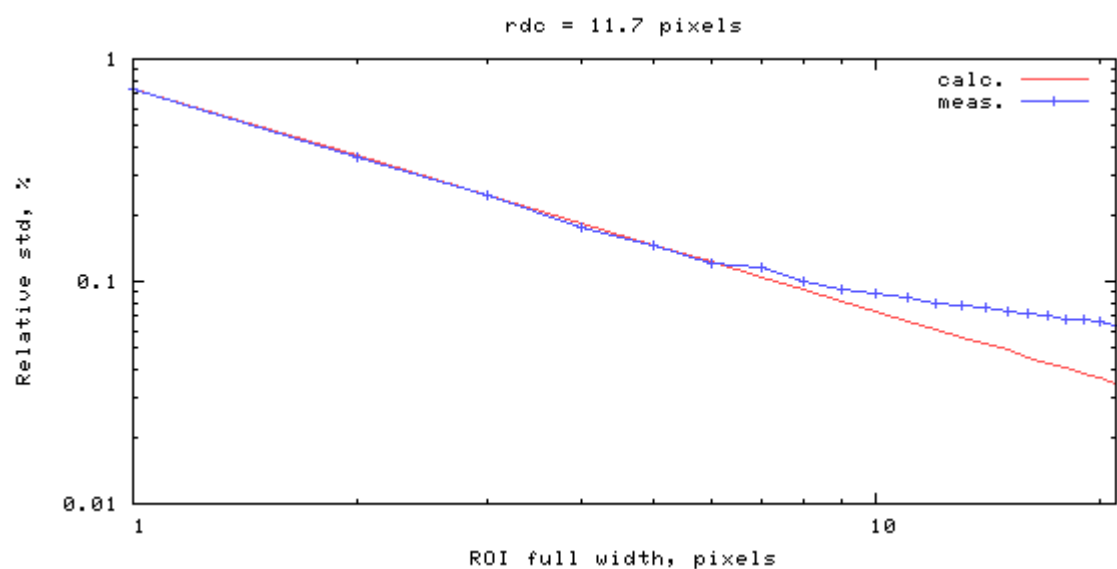

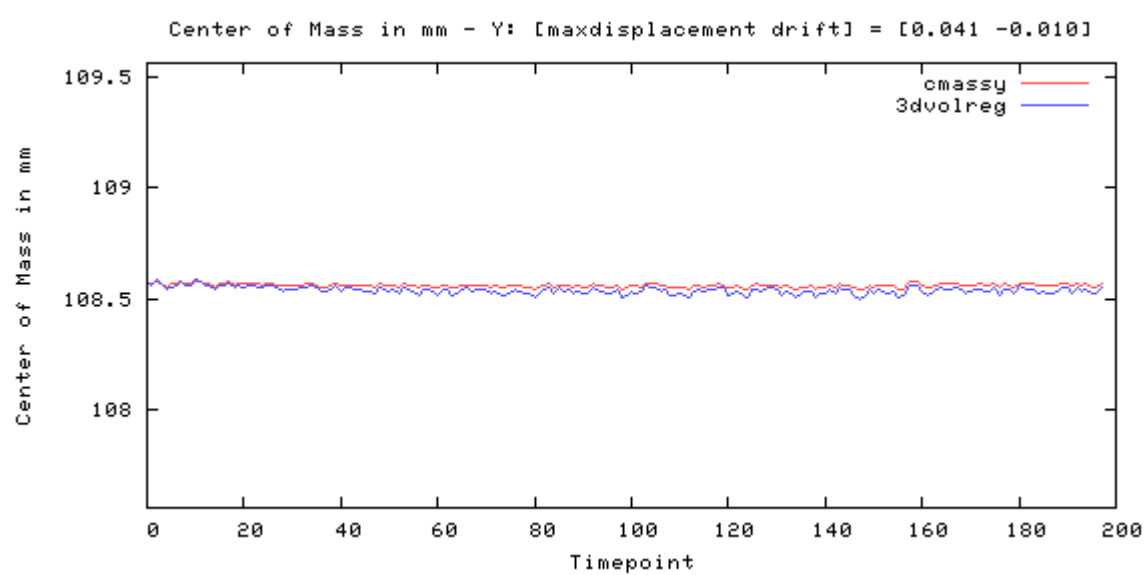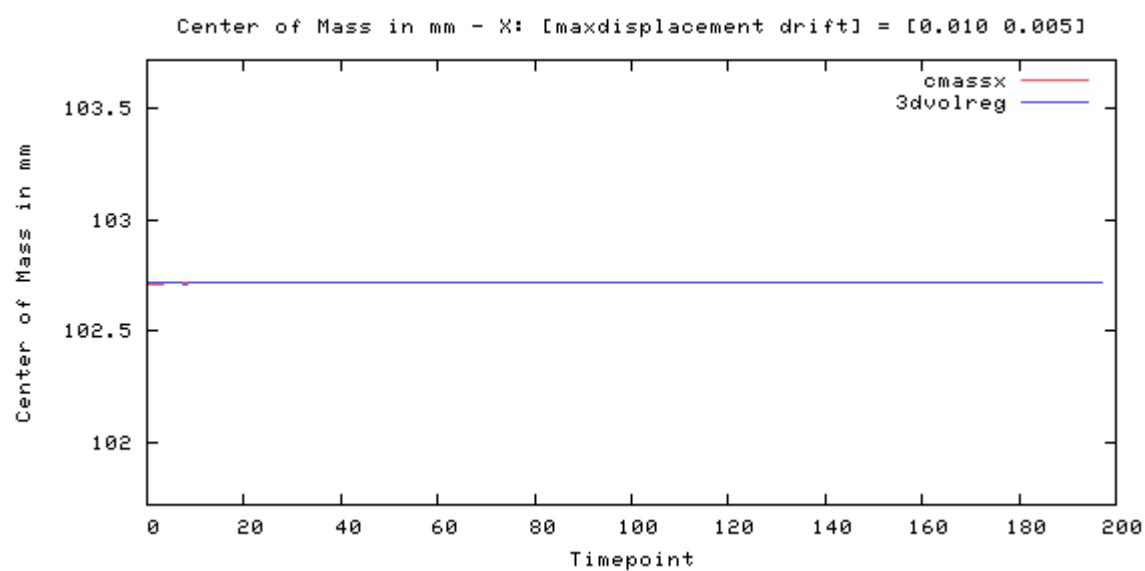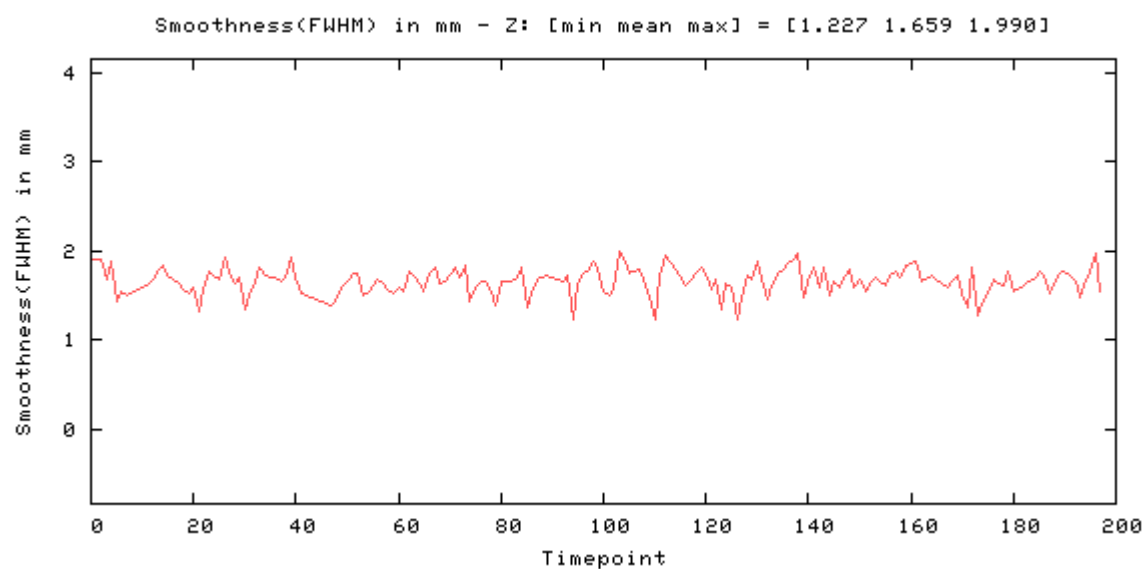


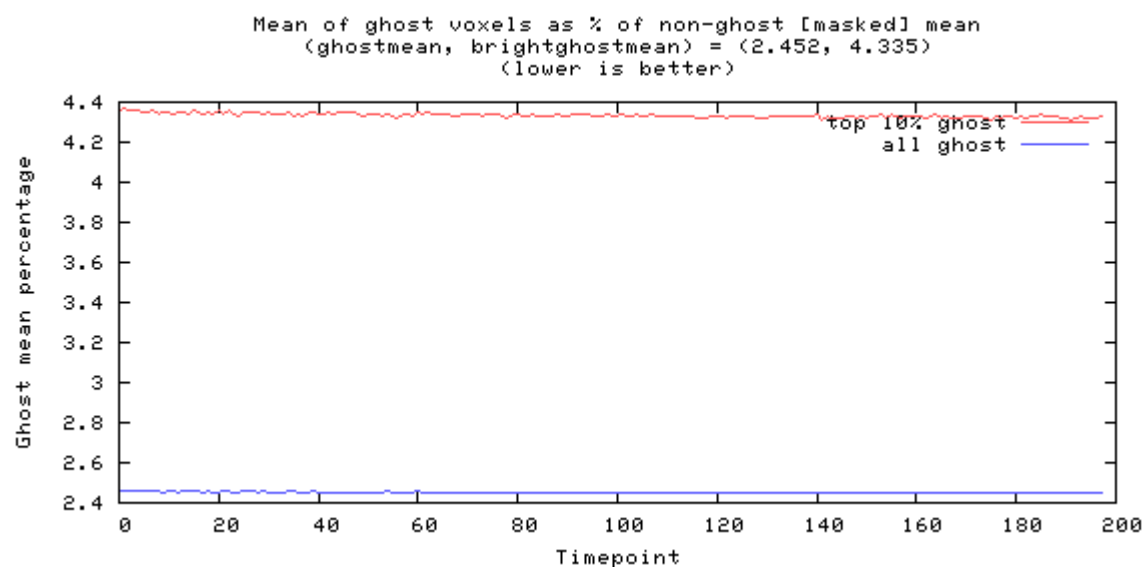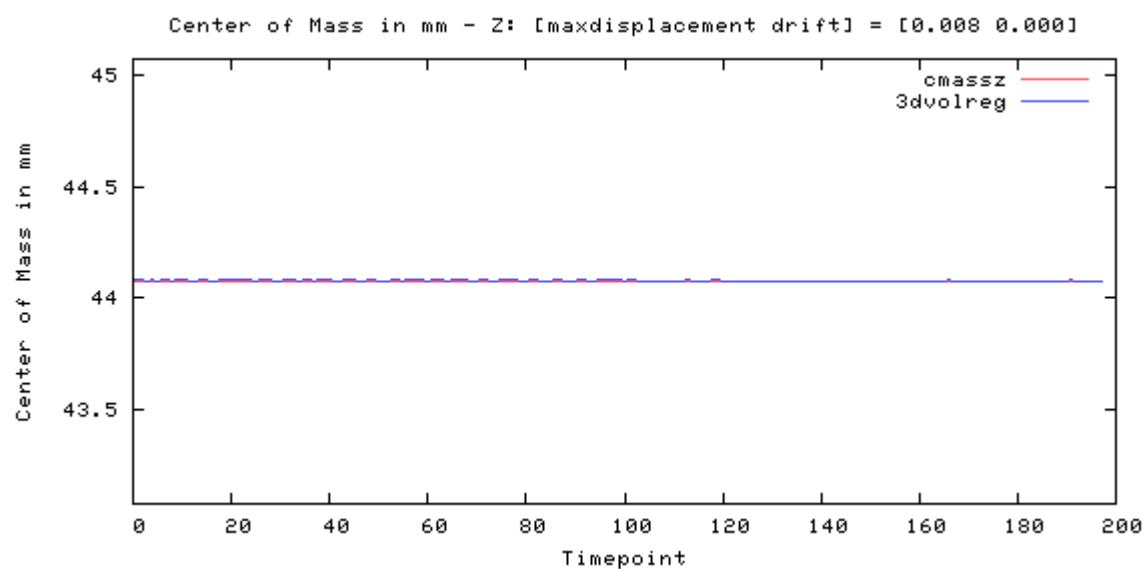

Odd-even difference image



image min: -235, image max: 301

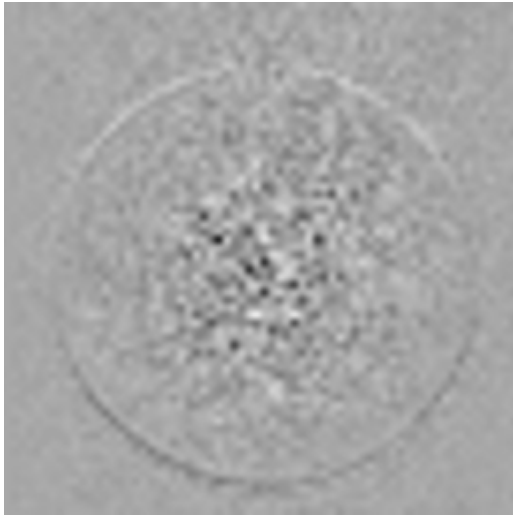

Mean image




image min: 0, image max: 1257.15



Standard Deviation image




image min: 0, image max: 7.99915

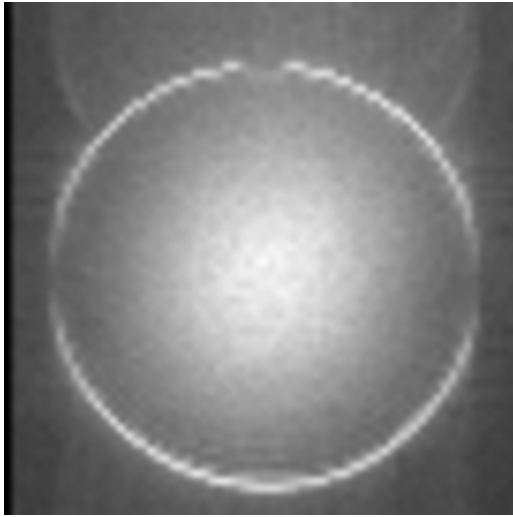

SFNR image



image min: 0, image max: 928.349

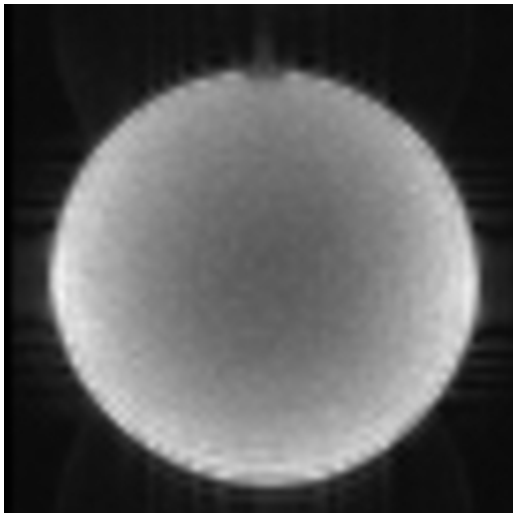

Acquisition parameters

dimensions	68x68x30x200
spacing	3mmx3mmx3mmx1[units?]
gap	0mmx0mmx0mmx0[units?]
scanner	AWP213018
psdname	
examnumber	2b45b240-1de4-44
studyid	Stability_20250318
seriesnumber	6
runnumber	1-200
scanningsequence	GR/EP

sequencevariant	SK
magneticfield	3
description	ses-64ch_func-bold_acq-siemensbasic_task-tsnr
scandate	2025-03-18
scantime	10:21:00
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	3
frequencydirection	1
acquisitionmatrix	68 68
bandwidth	166.668
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
ti	
protocolname	ses-64ch_func-bold_acq-siemensbasic_task-tsnr
sliceorder	2,4,6,8,10,12,14,16,18,20,22,24,26,28,30,1,3,5,7,9,11,
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,6
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20250318102059.350000
volumetricproperties	DISTORTED
scanoptions	FS/PER
mracquisitiontype	2D
sequencename	*epfid2d1_68
angioflag	N
numaverages	1
imagingfrequency	123.253268
imagednucleus	1H
spacingbetweenslices	3
numphaseencodingsteps	68
echotrainlength	68
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61

triggertime	
beatrejectionflag	N
cardiacnumberofimages	1
inplane phase encoding direction	COL
variableflipangleflag	N
sar	0.06788646892361
contentqualification	PRODUCT
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
samplesperpixel	1
photometricinterpretation	MONOCHROME2
rows	68
columns	68
bitsallocated	16
bitsstored	12
highbit	11
lossyimagecompression	00
echonumber	0
pixelbandwidth	2451