






# Phantom QA results for / fbirn\_out/ ses20ch\_funcbold\_acqsiemensstd\_task slice 17




1d acqsiemensstd tasktsnr 24.xml [percent fluct (trend removed), drift, driftfit]


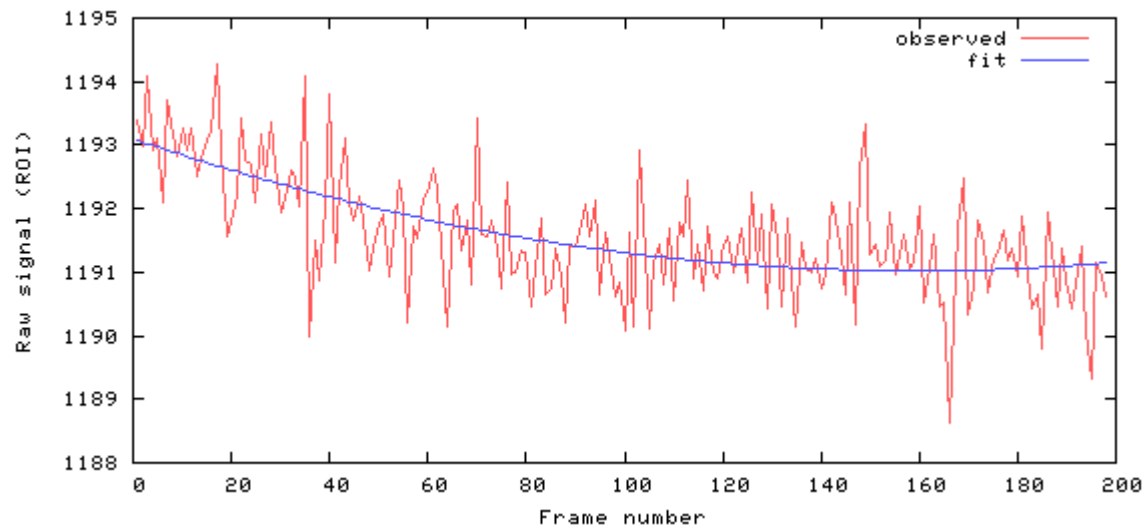


[mean, SNR, SFNR] = [1191.6 148.4 157.5]

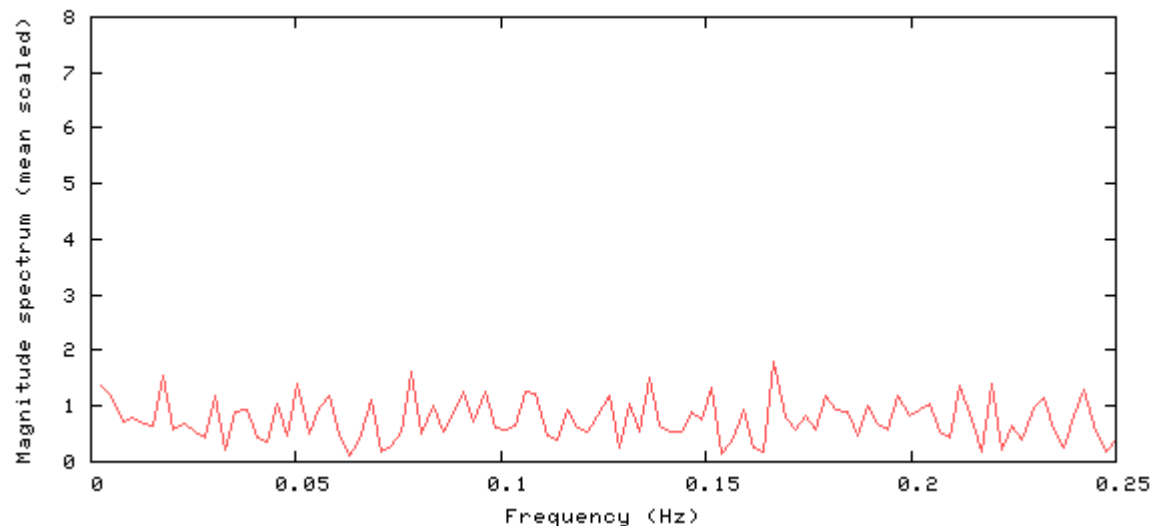

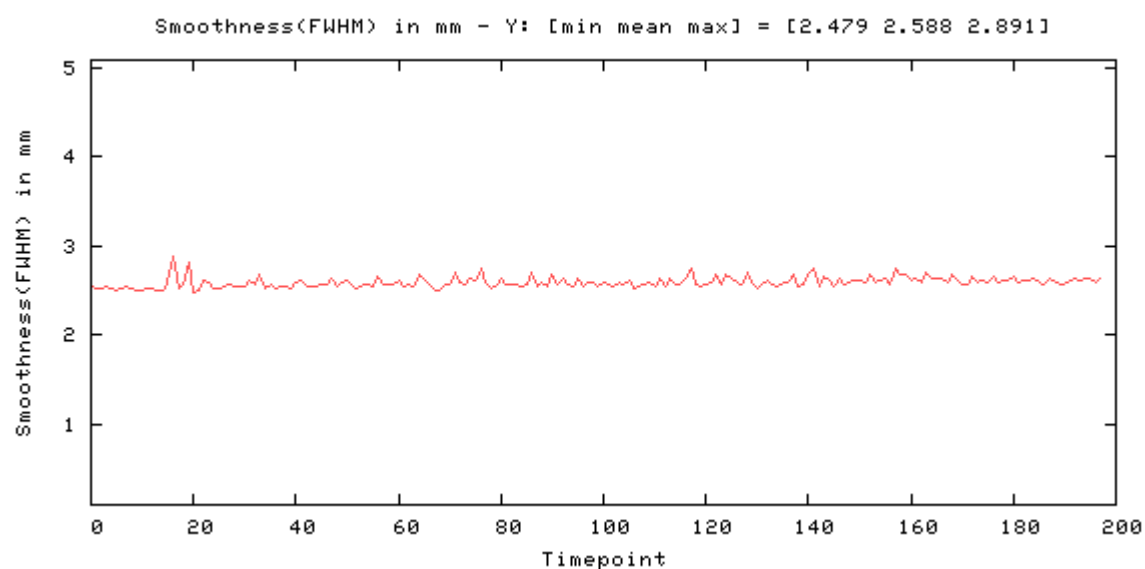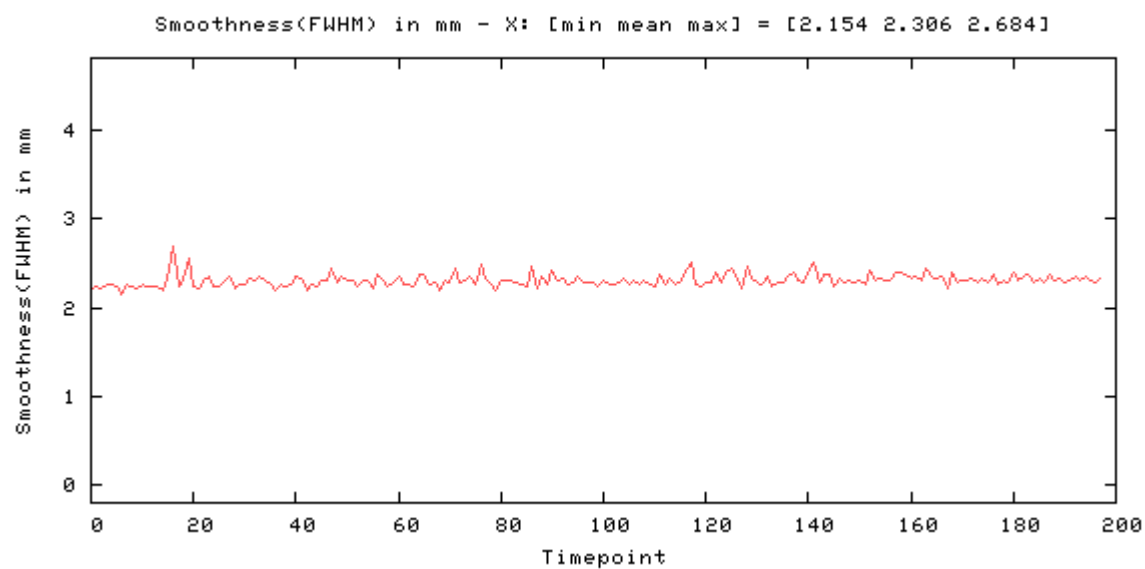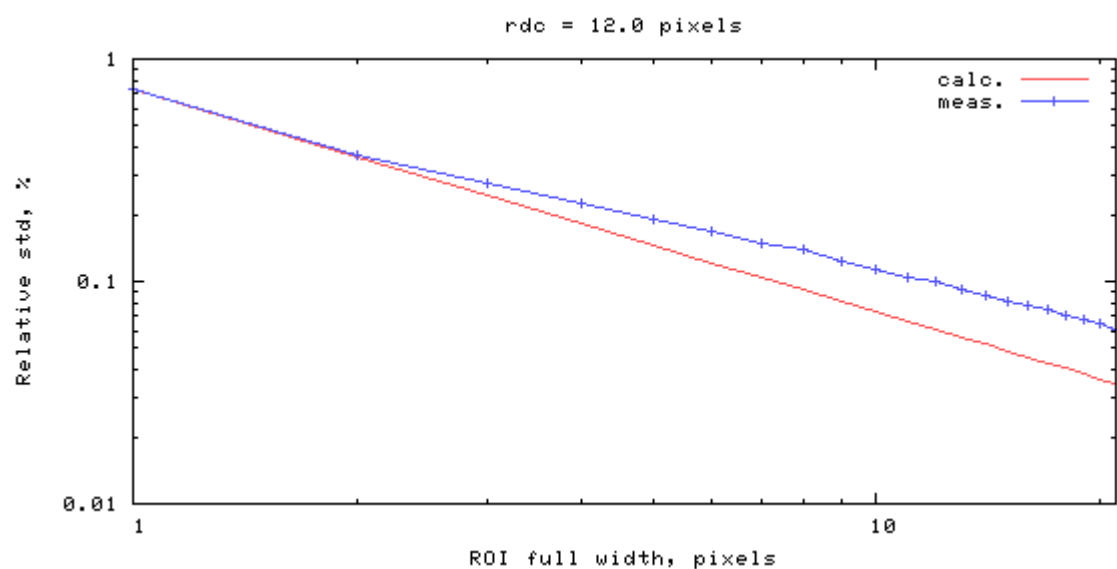

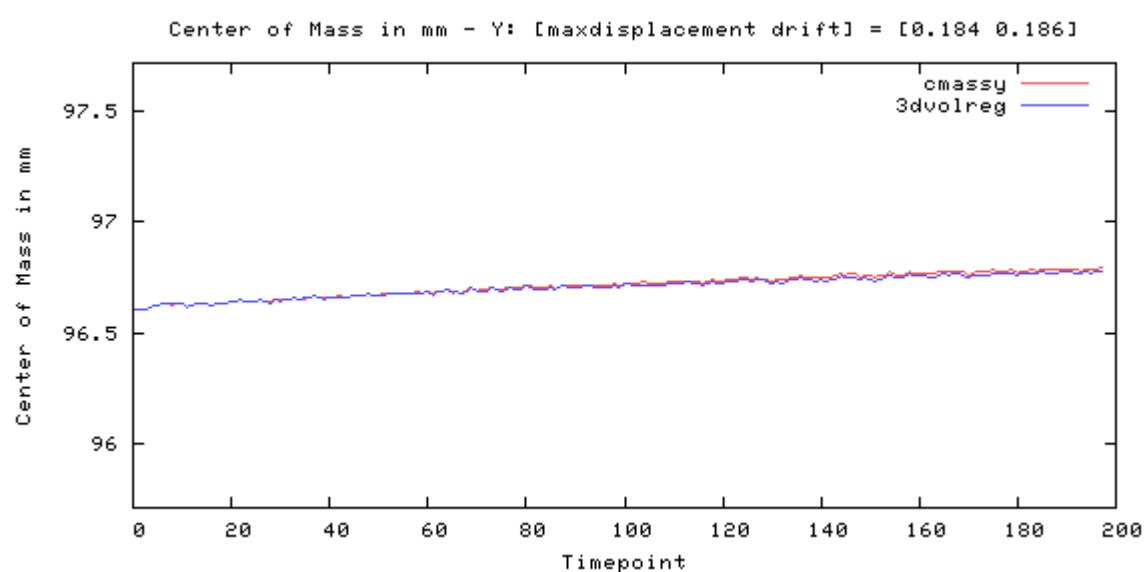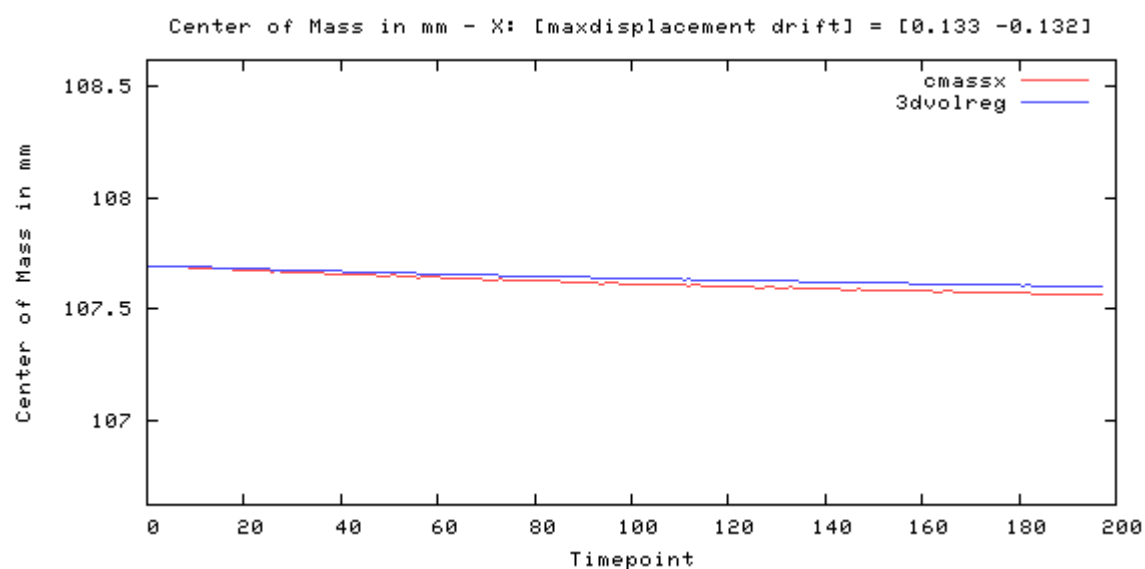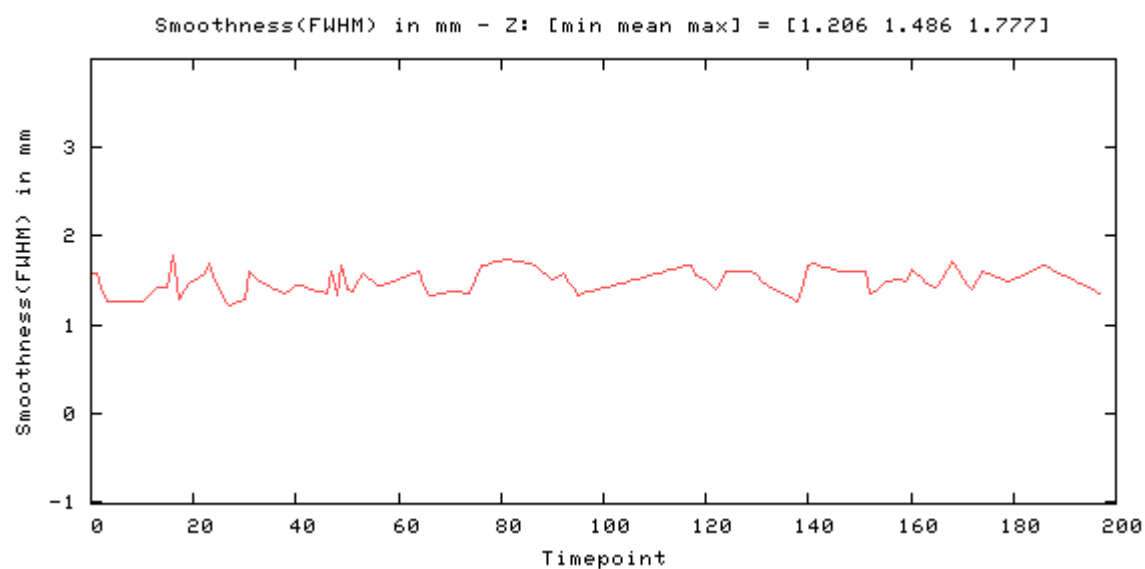

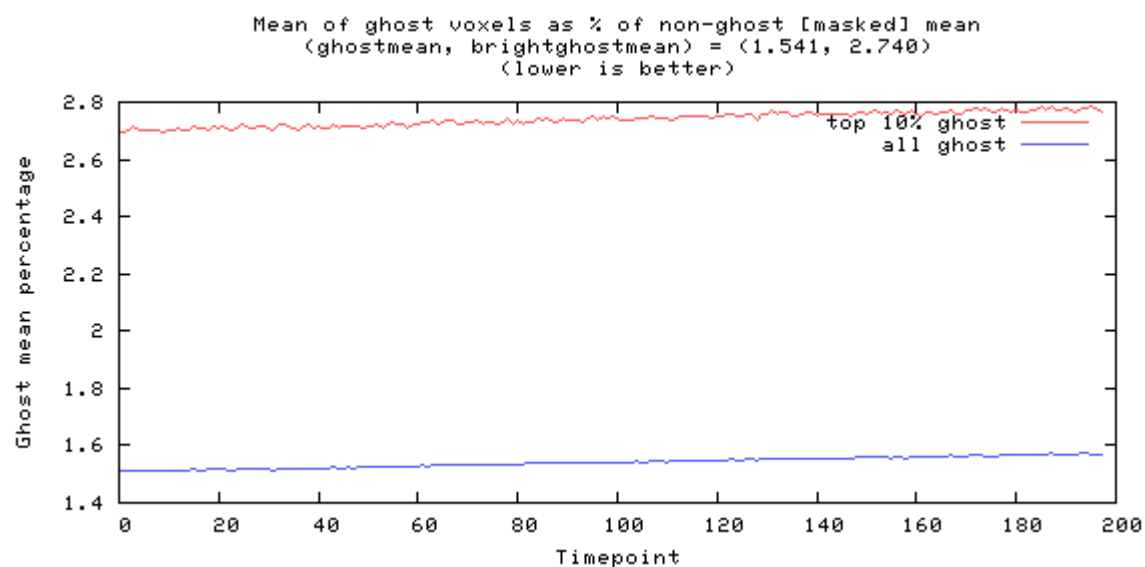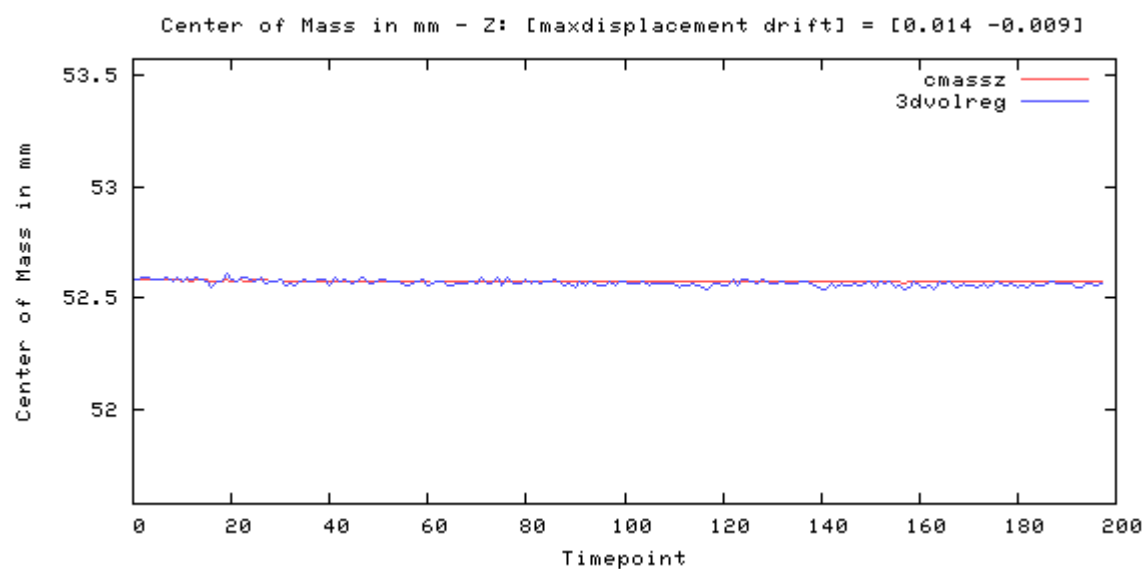


---

## Odd-even difference image



image min: -346, image max: 421

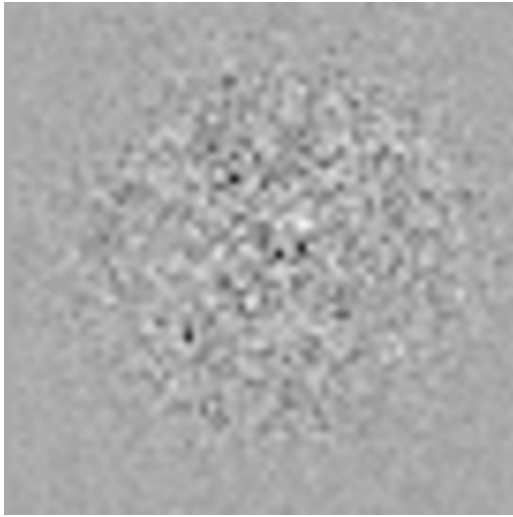

---

### Mean image



image min: 0, image max: 1465.12




---

### Standard Deviation image



image min: 0, image max: 9.65176

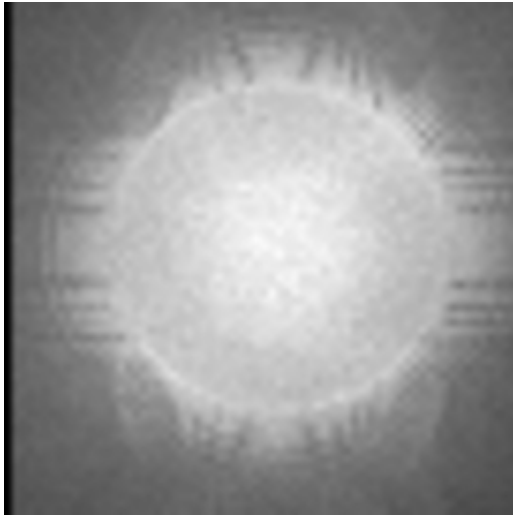


---

### SFNR image



image min: 0, image max: 282.214

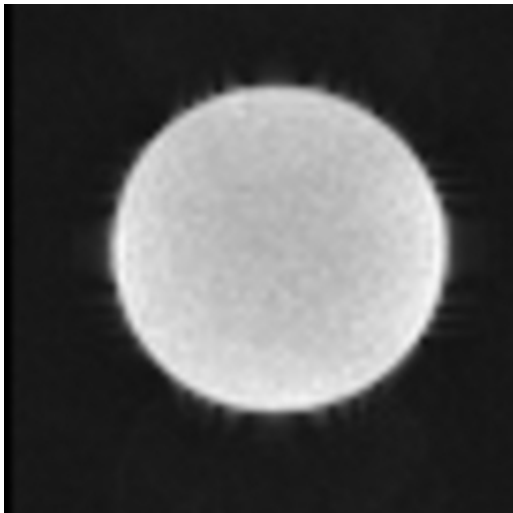

---

### Acquisition parameters

|                  |                       |
|------------------|-----------------------|
| dimensions       | 68x68x36x200          |
| spacing          | 3mmx3mmx3mmx1[units?] |
| gap              | 0mmx0mmx0mmx0[units?] |
| scanner          | AWP213018             |
| psdname          |                       |
| examnumber       | 8895a78e-5707-41      |
| studyid          | Stability_07072025    |
| seriesnumber     | 24                    |
| runnumber        | 1-200                 |
| scanningsequence | GR/EP                 |

|                             |                                                       |
|-----------------------------|-------------------------------------------------------|
| sequencevariant             | SK                                                    |
| magneticfield               | 3                                                     |
| description                 | ses-20ch_func-bold_acq-siemensstd_task-tsnr           |
| scandate                    | 2025-07-07                                            |
| scantime                    | 14:13:13                                              |
| tr                          | 2000                                                  |
| te                          | 30                                                    |
| operator                    |                                                       |
| flipangle                   | 79                                                    |
| prescribedslicespacing      | 3                                                     |
| frequencydirection          | 1                                                     |
| acquisitionmatrix           | 68 68                                                 |
| bandwidth                   | 166.668                                               |
| receivecoilname             | HeadNeck_20_TCS                                       |
| transmitcoilname            | Body                                                  |
| institution                 | THE UNIVERSITY OF PENNSYLVANIA                        |
| examdescription             | ACTIVE KIRWAN                                         |
| scanner                     | AWP213018                                             |
| scannermanufacturer         | Siemens Healthineers                                  |
| scannermodelname            | MAGNETOM Cima.X                                       |
| ti                          |                                                       |
| protocolname                | ses-20ch_func-bold_acq-siemensstd_task-tsnr           |
| sliceorder                  | 2,4,6,8,10,12,14,16,18,20,22,24,26,28,30,32,34,36,1,3 |
| institutionaddress          | Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,6 |
| institutionaldepartmentname | MindCORE Neuroimaging Suite 140                       |
| acquisitiondatetime         | 20250707141312.962500                                 |
| volumetricproperties        | DISTORTED                                             |
| scanoptions                 | FS/PER                                                |
| mracquisitiontype           | 2D                                                    |
| sequencename                | *epfid2d1_68                                          |
| angioflag                   | N                                                     |
| numaverages                 | 1                                                     |
| imagingfrequency            | 123.252981                                            |
| imagednucleus               | 1H                                                    |
| spacingbetweenslices        | 3                                                     |
| numphaseencodingsteps       | 68                                                    |
| echoTRAINlength             | 33                                                    |
| percentsampling             | 100                                                   |
| percentphasefieldofview     | 100                                                   |
| scannerserialnumber         | 213018                                                |
| softwareversions            | syngo MR XA61                                         |

|                                  |                  |
|----------------------------------|------------------|
| triggertime                      |                  |
| beatrejectionflag                | N                |
| cardiacnumberofimages            | 1                |
| inplane phase encoding direction | COL              |
| variableflipangleflag            | N                |
| sar                              | 0.11859437080037 |
| contentqualification             | PRODUCT          |
| acquisitionnumber                | 1                |
| instancenumber                   | 1                |
| numtemporalpositions             | 200              |
| samplesperpixel                  | 1                |
| photometricinterpretation        | MONOCHROME2      |
| rows                             | 68               |
| columns                          | 68               |
| bitsallocated                    | 16               |
| bitsstored                       | 12               |
| highbit                          | 11               |
| lossyimagecompression            | 00               |
| echonumber                       | 0                |
| pixelbandwidth                   | 2451             |