






# Phantom QA results for / fbirn\_out/ ses64ch\_funcbold\_acqcmrrstd\_tasktsnr slice 34




```
bold acqcmrrstd tasktsnr 6.xml [percent fluct (trend removed), drift, driftfit] =
```


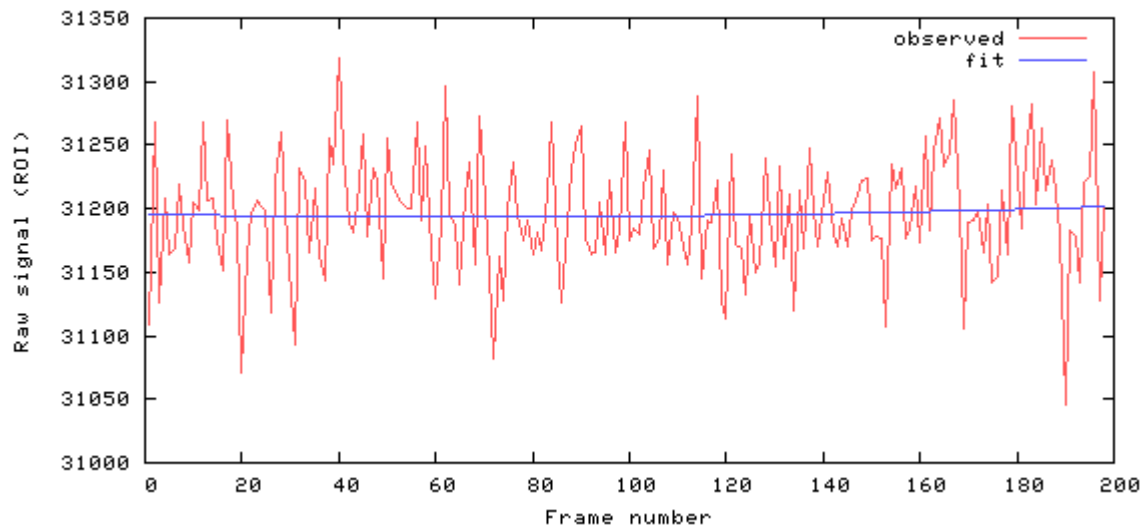


```
[mean, SNR, SFNR] = [31195.6 129.3 122.7]
```

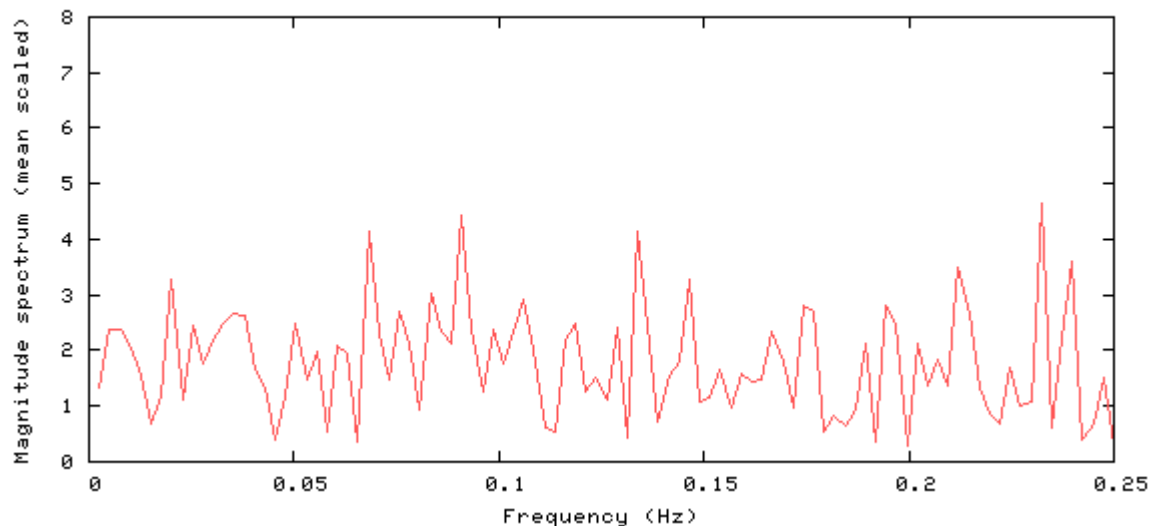

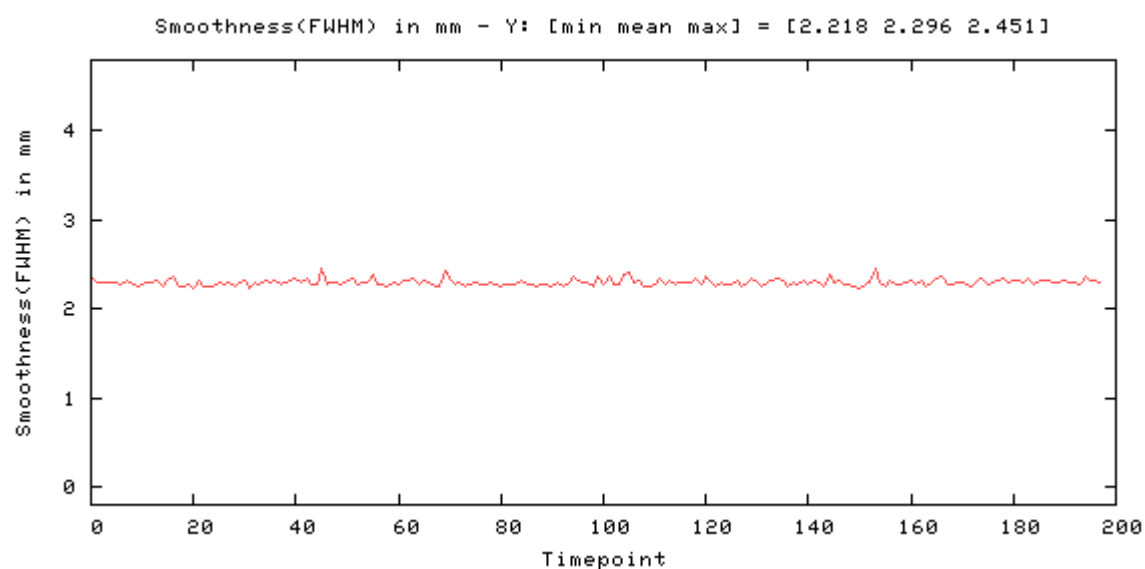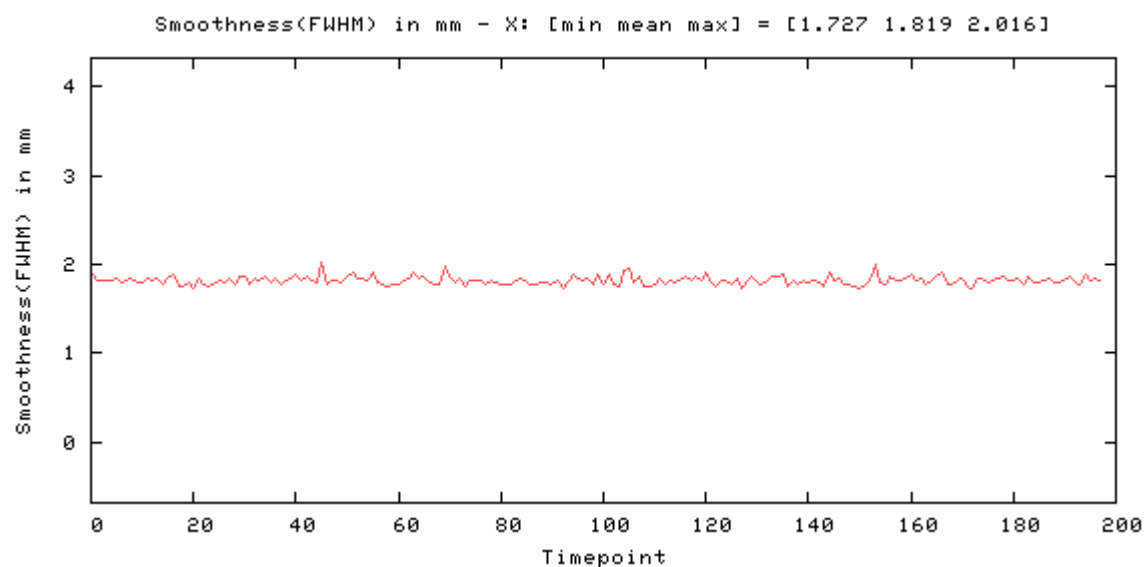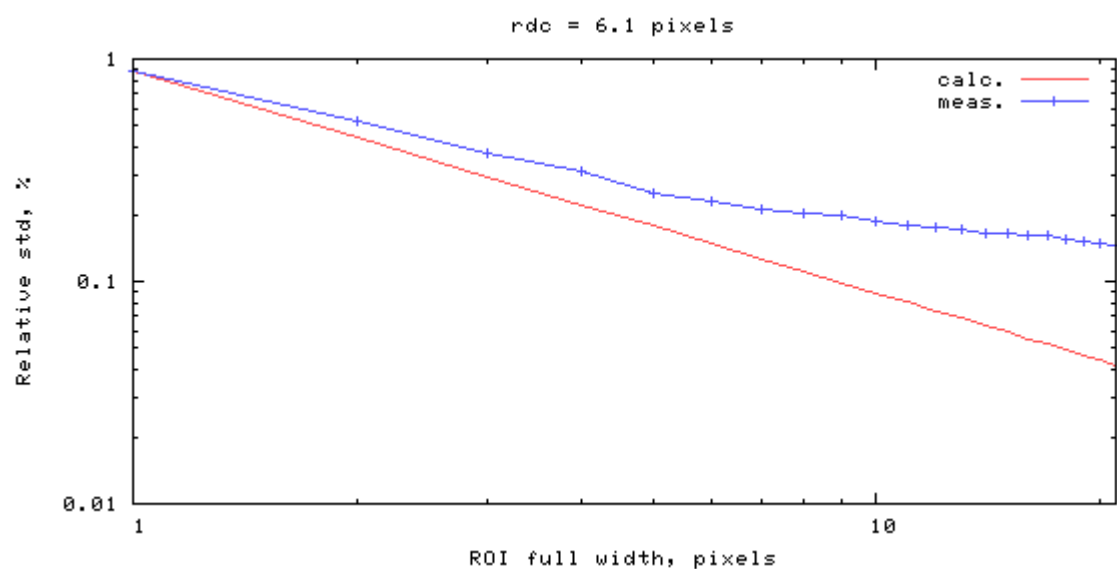

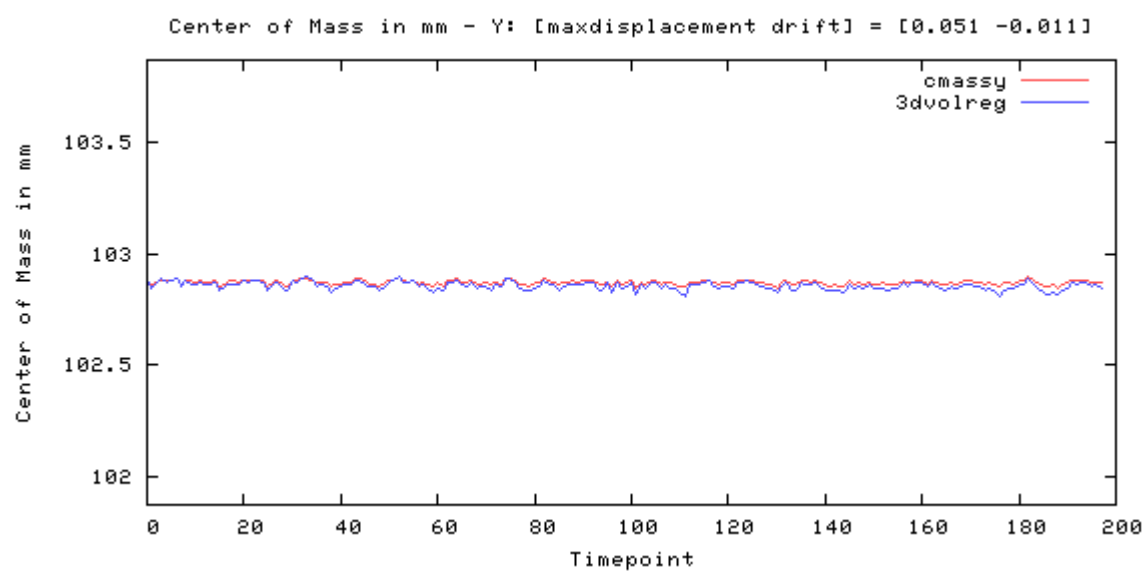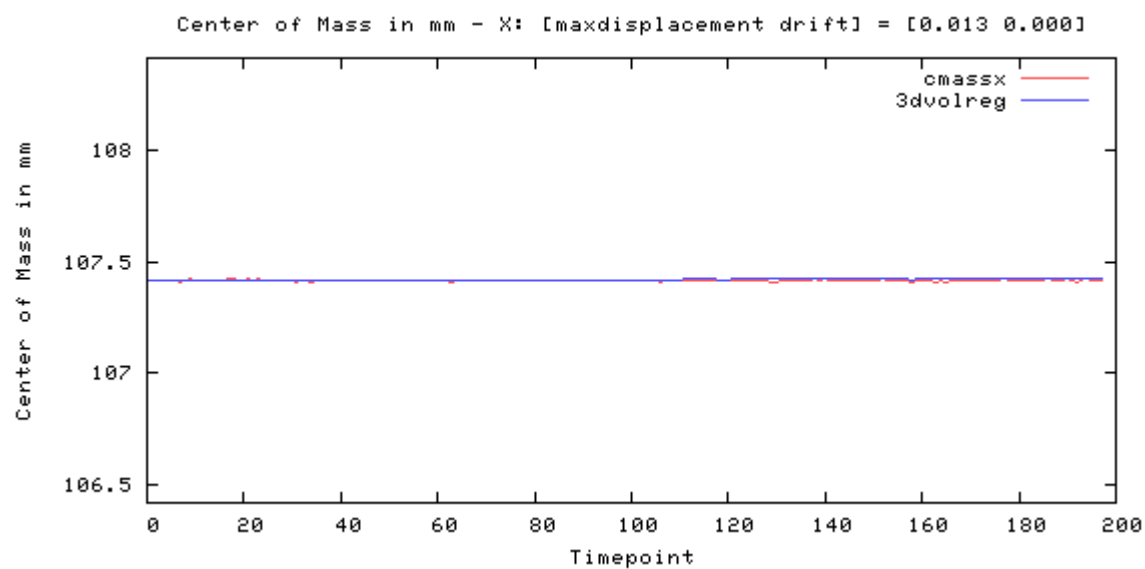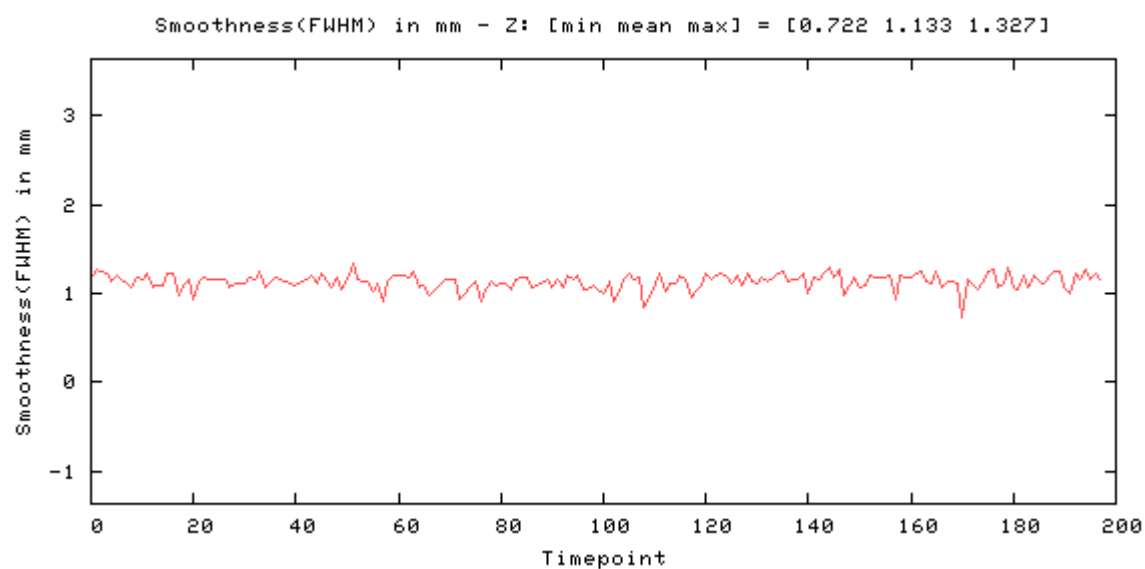

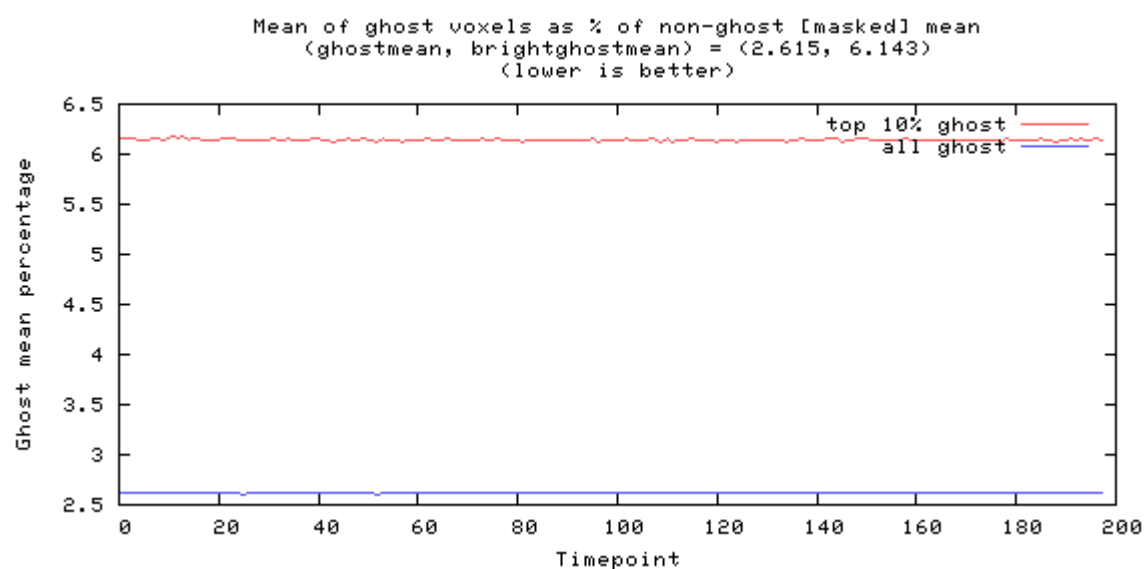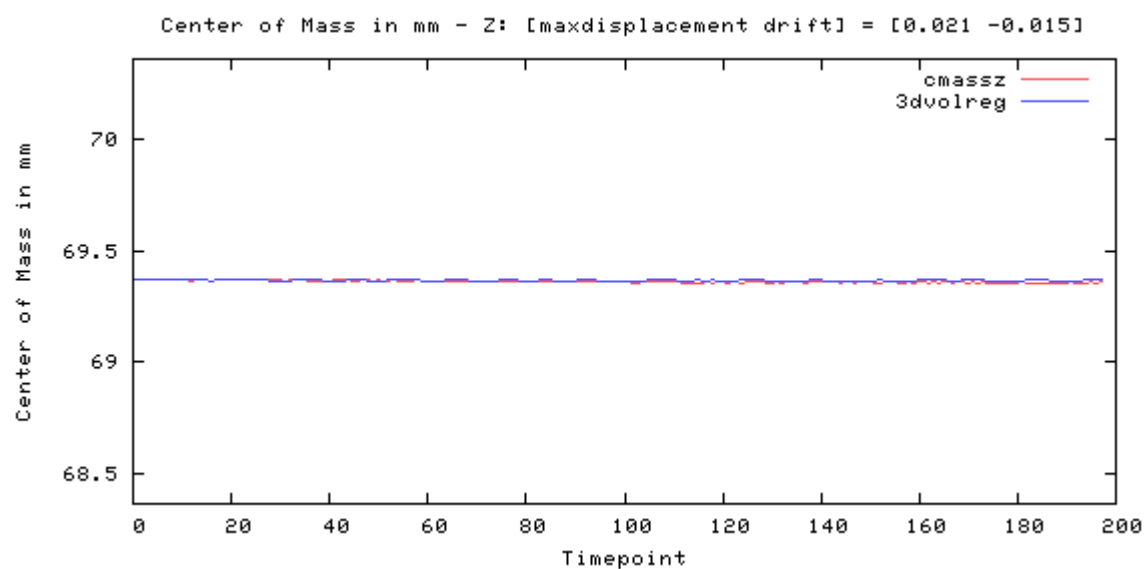


---

## Odd-even difference image

-10577 11046



image min: -10577, image max: 11046

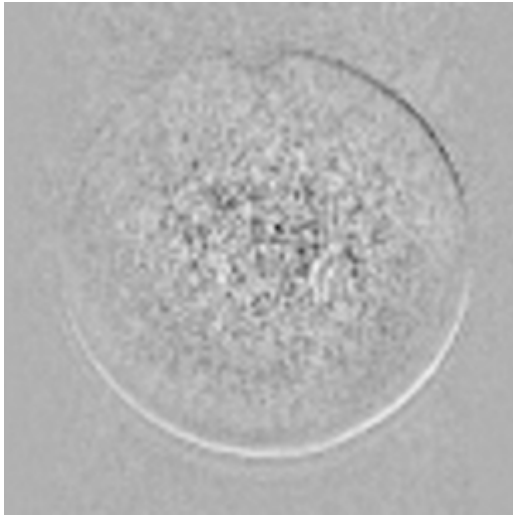


---

### Mean image



image min: 0, image max: 43161.2

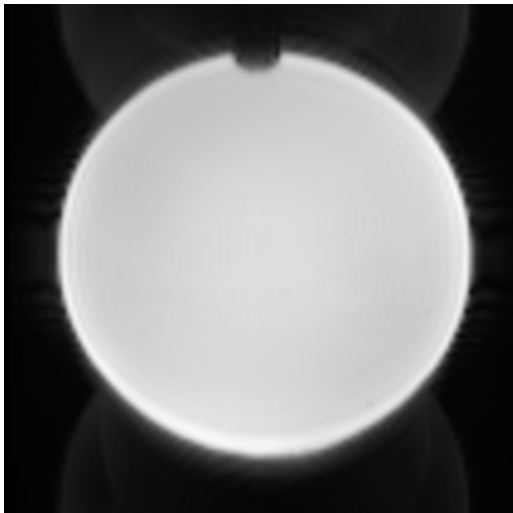


---

### Standard Deviation image



image min: 0, image max: 316.697

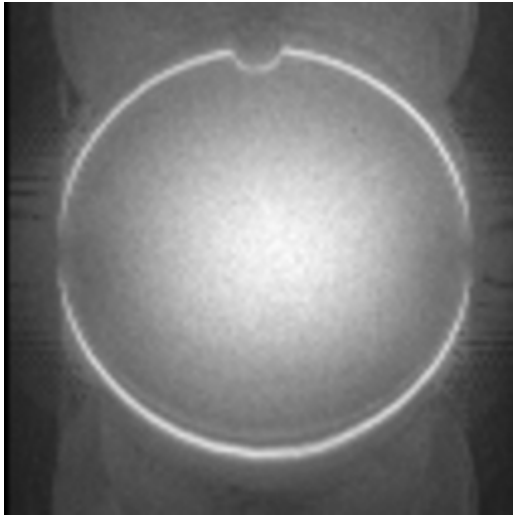


---

### SFNR image



image min: 0, image max: 731.512

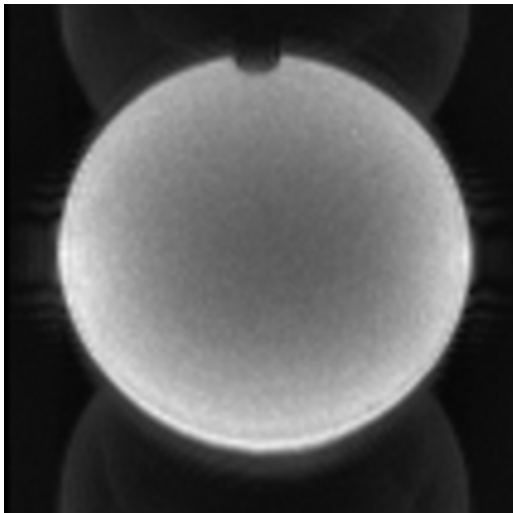

---

### Acquisition parameters

|                  |                       |
|------------------|-----------------------|
| dimensions       | 106x106x69x200        |
| spacing          | 2mmx2mmx2mmx1[units?] |
| gap              | 0mmx0mmx0mmx0[units?] |
| scanner          | AWP213018             |
| psdname          |                       |
| examnumber       | 9bf9b951-052d-4f      |
| studyid          | Stability_20250908    |
| seriesnumber     | 6                     |
| runnumber        | 1-200                 |
| scanningsequence | GR/EP                 |

|                             |                                                                               |
|-----------------------------|-------------------------------------------------------------------------------|
| sequencevariant             | SK                                                                            |
| magneticfield               | 3                                                                             |
| description                 | ses-64ch_func-bold_acq-cmrrstd_task-tsnr                                      |
| scandate                    | 2025-09-08                                                                    |
| scantime                    | 10:20:16                                                                      |
| tr                          | 2000                                                                          |
| te                          | 30                                                                            |
| operator                    |                                                                               |
| flipangle                   | 79                                                                            |
| prescribedslicespacing      | 2                                                                             |
| frequencydirection          | 1                                                                             |
| acquisitionmatrix           | 106 106                                                                       |
| bandwidth                   | 192.284                                                                       |
| receivecoilname             | HeadNeck_64_CS                                                                |
| transmitcoilname            | Body                                                                          |
| institution                 | THE UNIVERSITY OF PENNSYLVANIA                                                |
| examdescription             | ACTIVE KIRWAN                                                                 |
| scanner                     | AWP213018                                                                     |
| scannermanufacturer         | Siemens Healthineers                                                          |
| scannermodelname            | MAGNETOM Cima.X                                                               |
| ti                          |                                                                               |
| protocolname                | ses-64ch_func-bold_acq-cmrrstd_task-tsnr                                      |
| institutionaddress          | Grays Ferry Ave. BLDG 250 STE 140<br>3401,Philadelphia,69 - KEYSTONE,US,19146 |
| institutionaldepartmentname | MindCORE Neuroimaging Suite 140                                               |
| acquisitiondatetime         | 20250908102014.412500                                                         |
| volumetricproperties        | DISTORTED                                                                     |
| scanoptions                 | PFP/FS/PER                                                                    |
| mracquisitiontype           | 2D                                                                            |
| sequencename                | epfid2d1_106                                                                  |
| angioflag                   | N                                                                             |
| numaverages                 | 1                                                                             |
| imagingfrequency            | 123.253209                                                                    |
| imagednucleus               | 1H                                                                            |
| spacingbetweenslices        | 2                                                                             |
| numphaseencodingsteps       | 106                                                                           |
| echoTRAINlength             | 80                                                                            |
| percentsampling             | 100                                                                           |
| percentphasefieldofview     | 100                                                                           |
| scannerserialnumber         | 213018                                                                        |
| softwareversions            | syngo MR XA61                                                                 |

|                                  |                                              |
|----------------------------------|----------------------------------------------|
| triggertime                      |                                              |
| beatrejectionflag                | N                                            |
| cardiacnumberofimages            | 1                                            |
| inplane phase encoding direction | COL                                          |
| variableflipangleflag            | N                                            |
| sar                              | 0.01725800797120                             |
| contentqualification             | RESEARCH                                     |
| acquisitionnumber                | 1                                            |
| instancenumber                   | 1                                            |
| numtemporalpositions             | 200                                          |
| imagecomments                    | Not for diagnostic use, Unaliased MB3/PE3/LB |
| samplesperpixel                  | 1                                            |
| photometricinterpretation        | MONOCHROME2                                  |
| rows                             | 106                                          |
| columns                          | 106                                          |
| bitsallocated                    | 16                                           |
| bitsstored                       | 16                                           |
| highbit                          | 15                                           |
| lossyimagecompression            | 00                                           |
| echonumber                       | 0                                            |
| pixelbandwidth                   | 1814                                         |