






Phantom QA results for /fbirn_out/ses64ch_funcbold_acqcmrrstd_taskwarmslice 34




old acqcmrrstd taskwarmup 2.xml [percent fluct (trend removed), drift, driftfit]


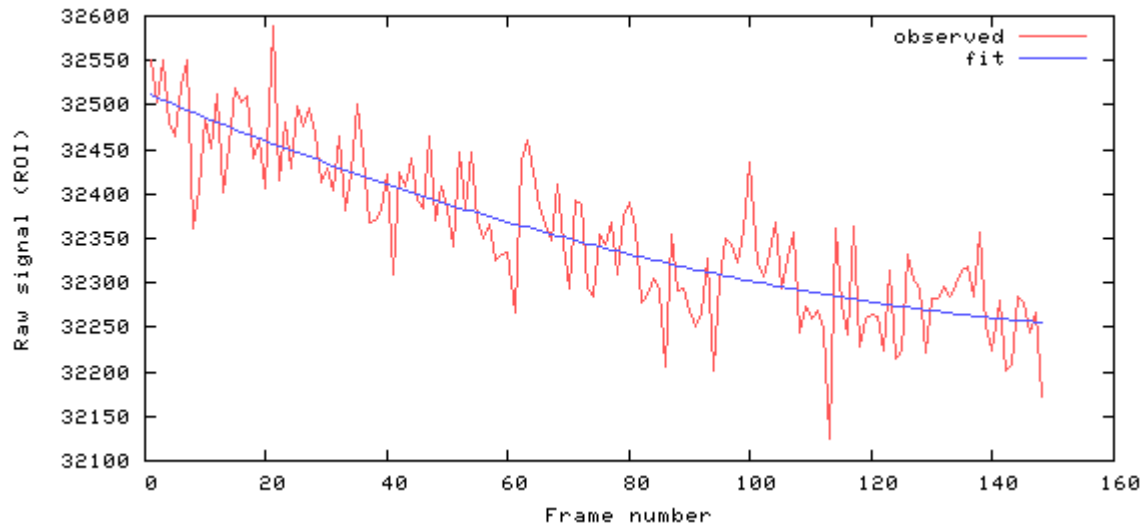


[mean, SNR, SFNR] = [32355.7 125.6 125.2]

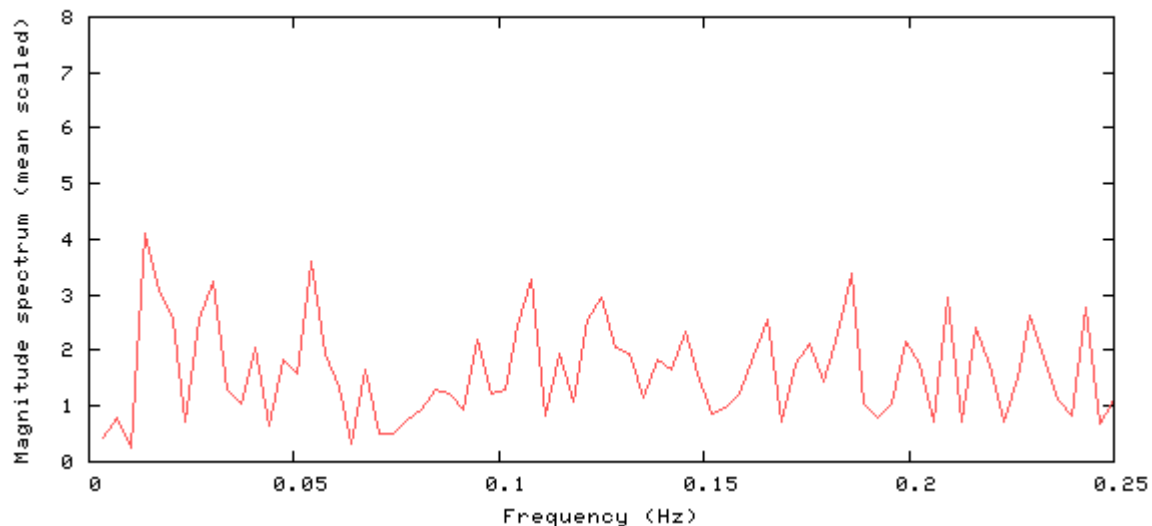

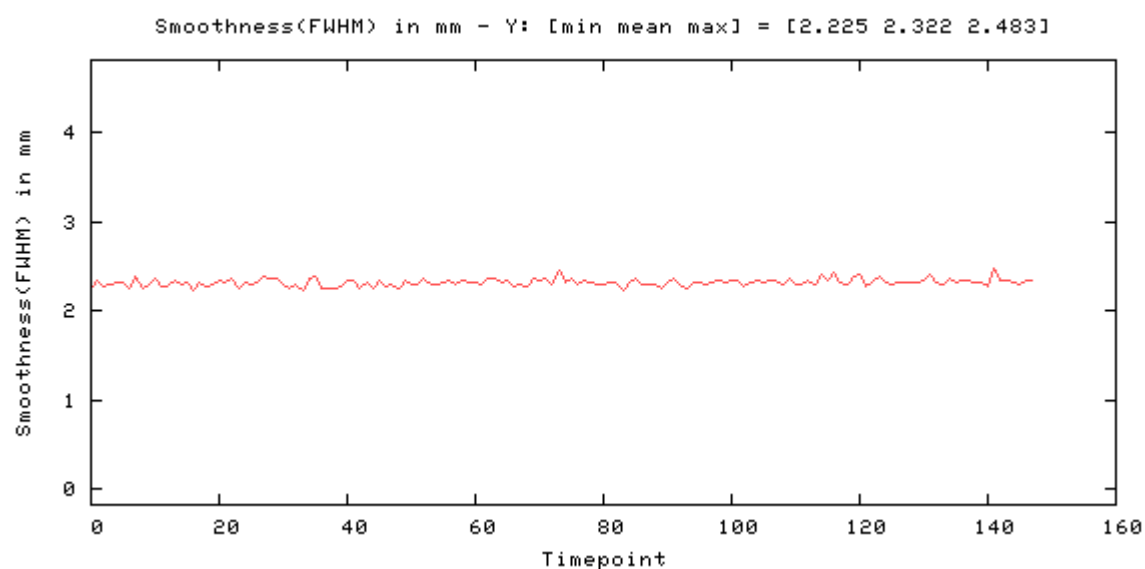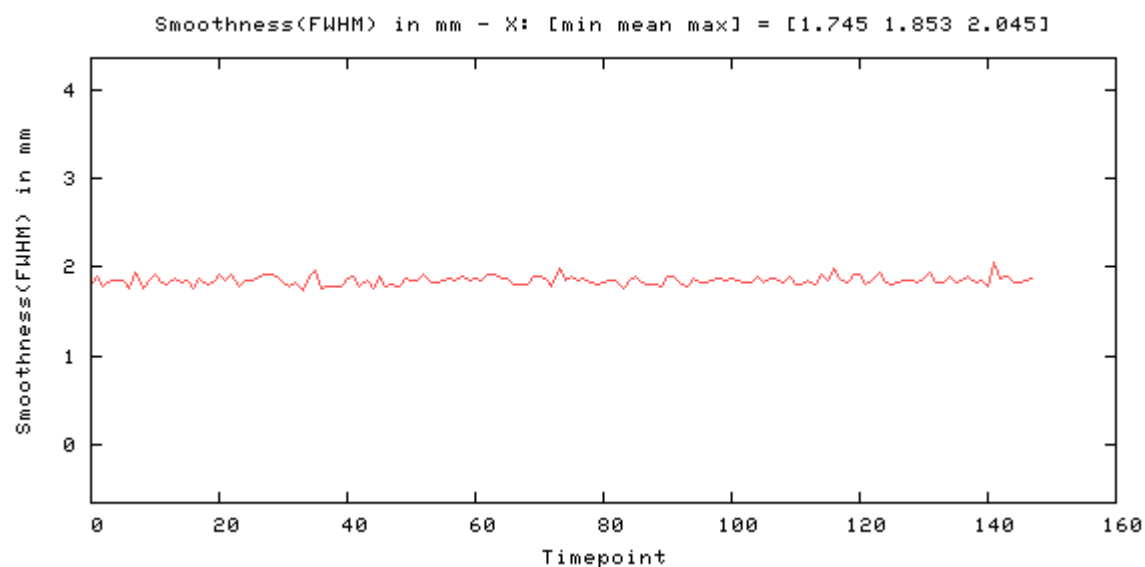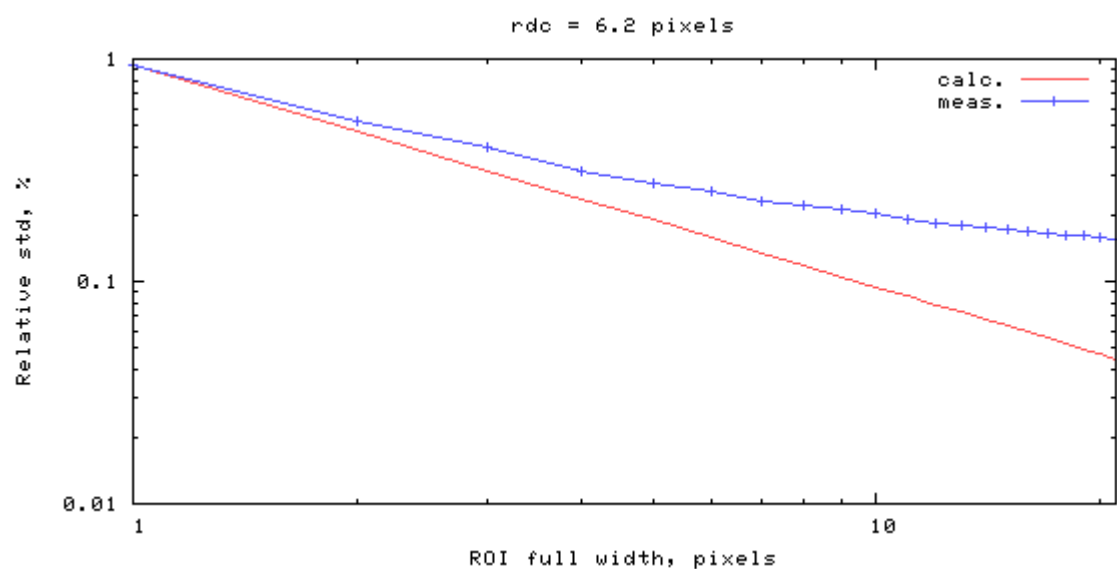

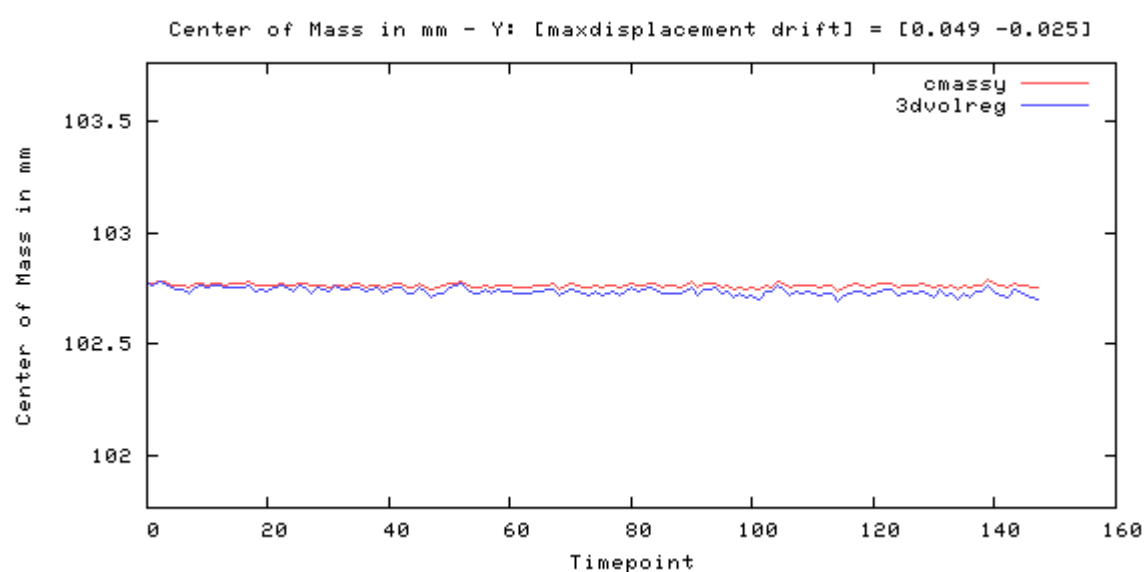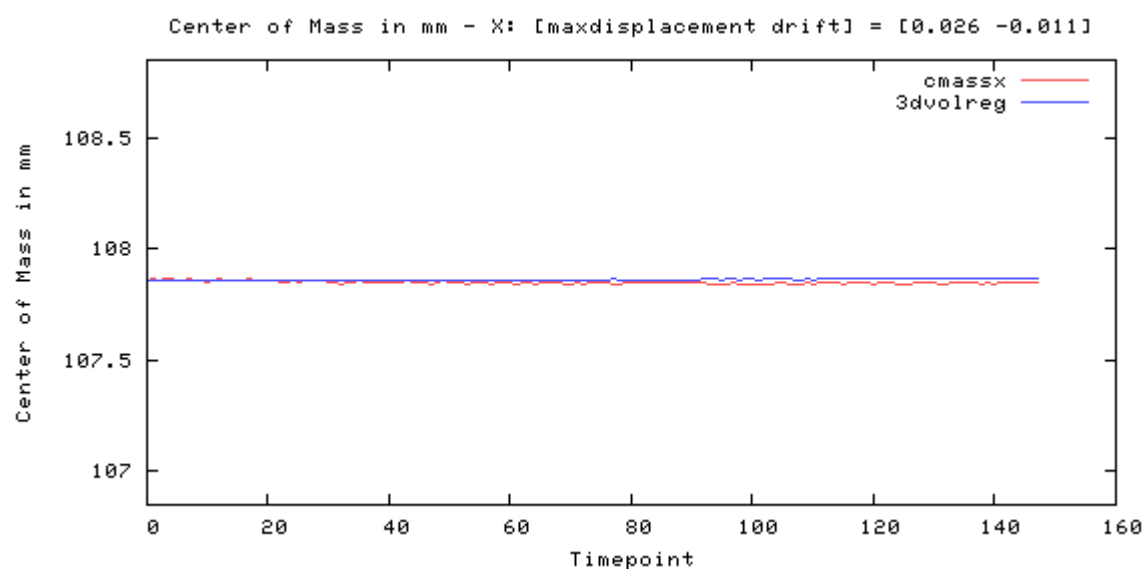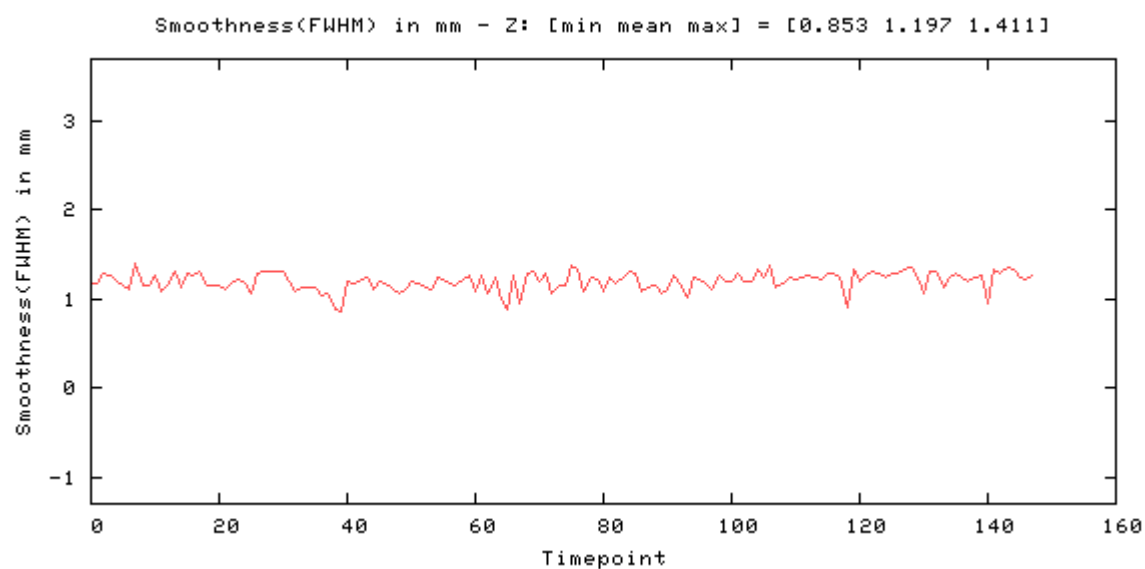

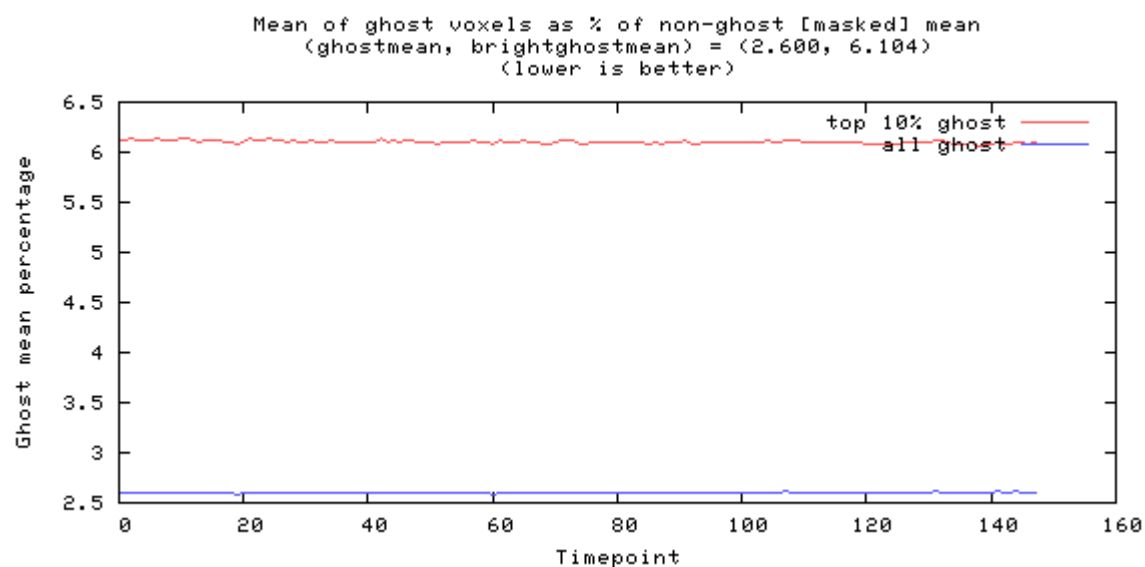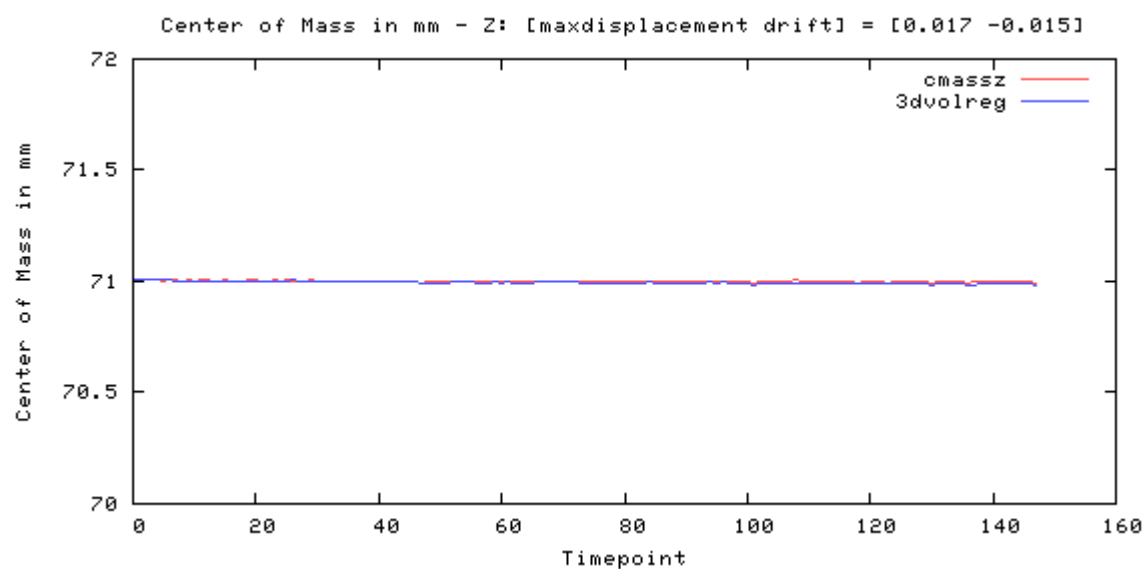


Odd-even difference image

-9441 12164




image min: -9441, image max: 12164

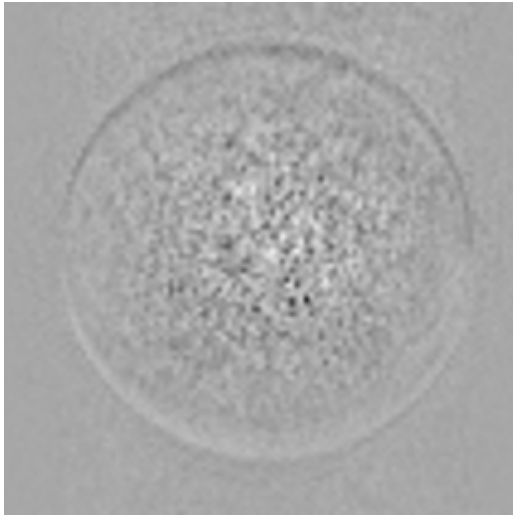

Mean image




image min: 0, image max: 43594.1

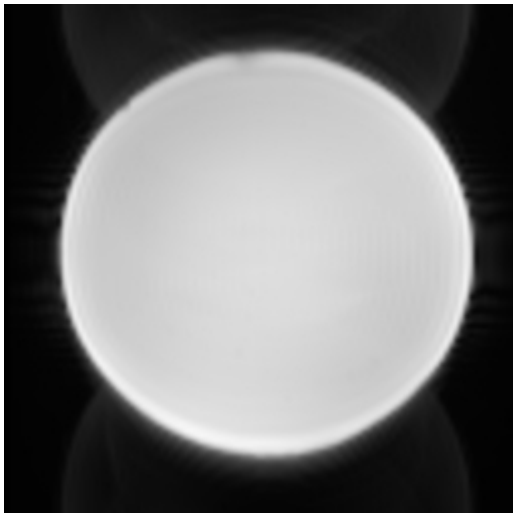

Standard Deviation image




image min: 0, image max: 335.966

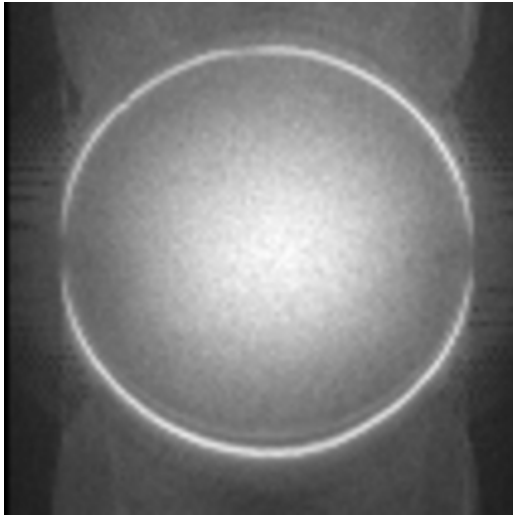

SFNR image



image min: 0, image max: 761.93

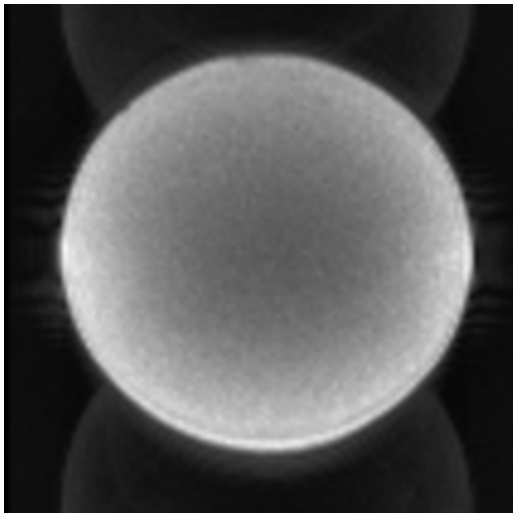

Acquisition parameters

dimensions	106x106x69x150
spacing	2mmx2mmx2mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	9b632d65-dffd-49
studyid	Stability_20250930
seriesnumber	2
runnumber	1
magneticfield	3

description	ses-64ch_func-bold_acq-cmrrstd_task-warmup
scandate	2025-09-30
scantime	09:53:05
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-64ch_func-bold_acq-cmrrstd_task-warmup
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20250930095303.935000
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	2
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	RESEARCH
acquisitionduration	317.89499999999998

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	150
imagecomments	Not for diagnostic use, Unaliased MB3/ PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	69
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00