






Phantom QA results for / fbirn_out/ ses64ch_funcbold_acqsiemensbasic_ta slice 14




1d acqsiemensbasic tasktsnr 9.xml [percent fluct (trend removed), drift, driftfit:


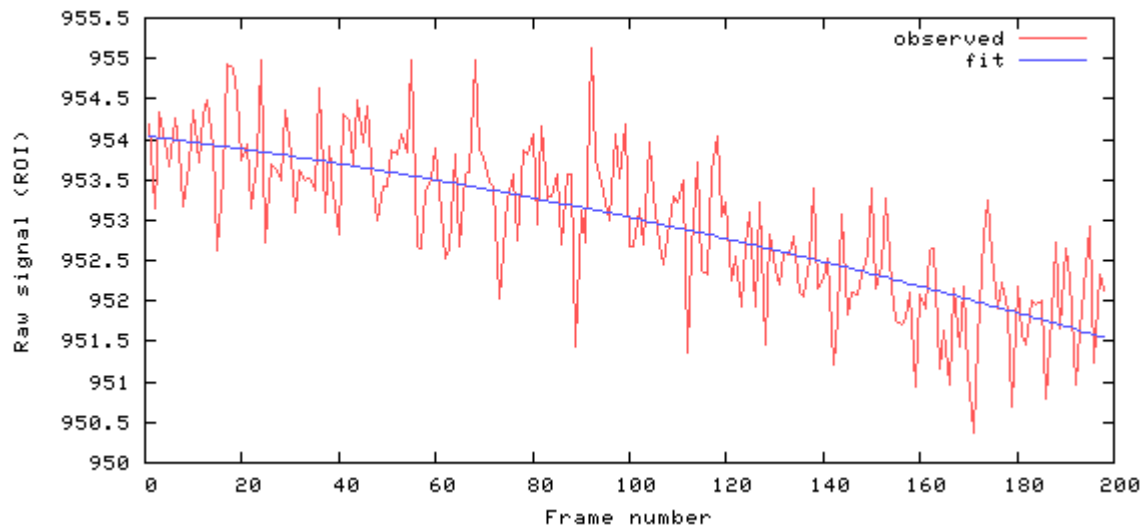


[mean, SNR, SFNR] = [953.0 162.5 171.2]

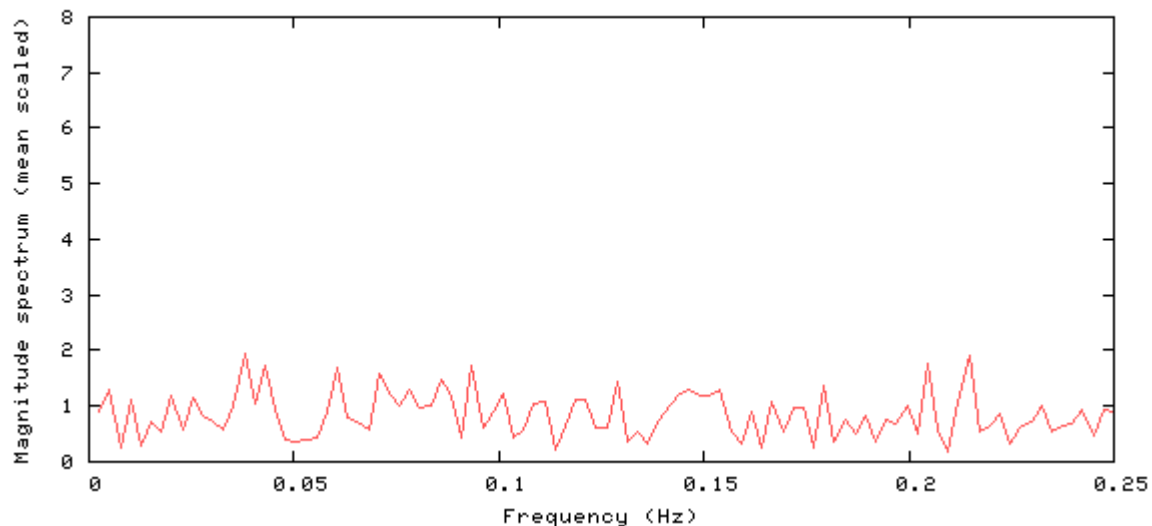

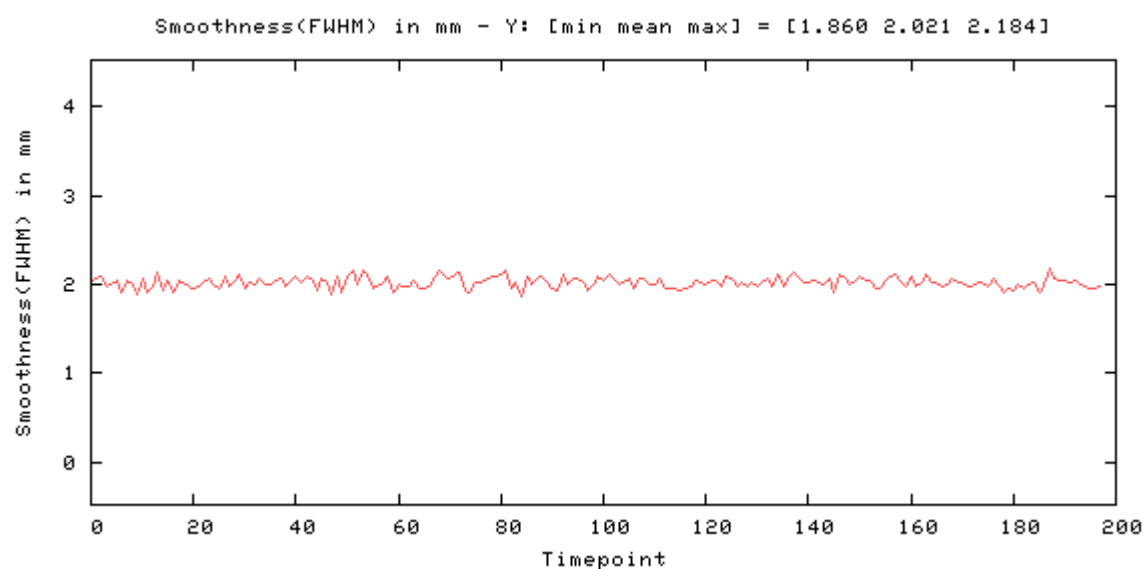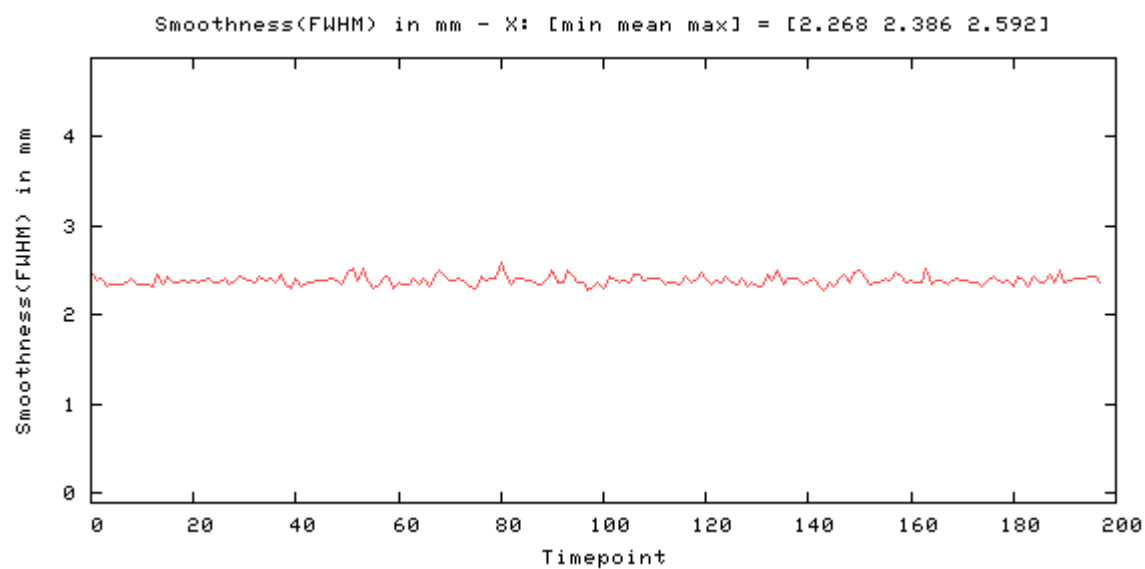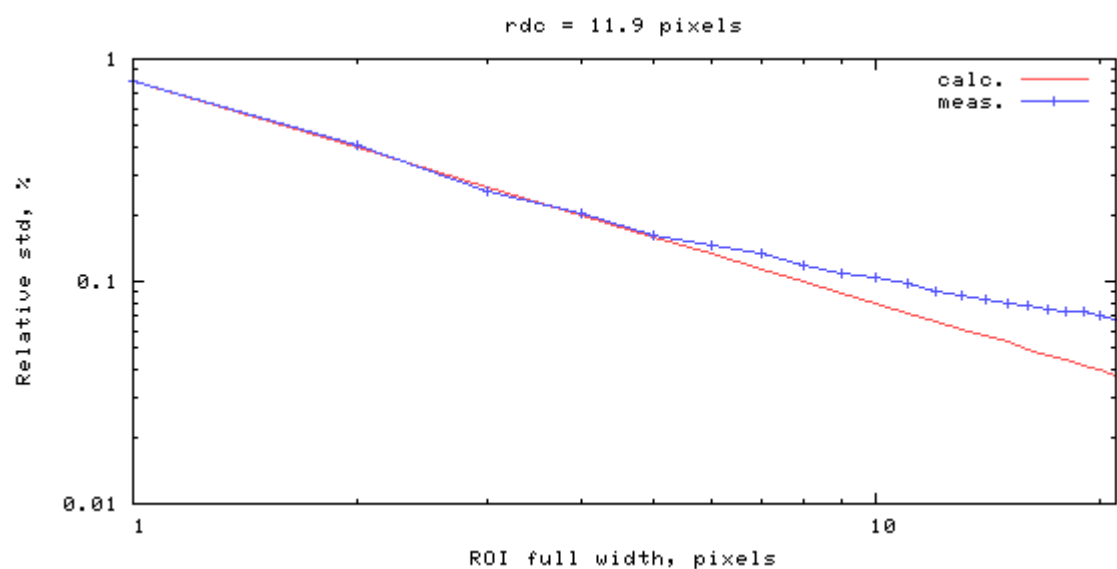

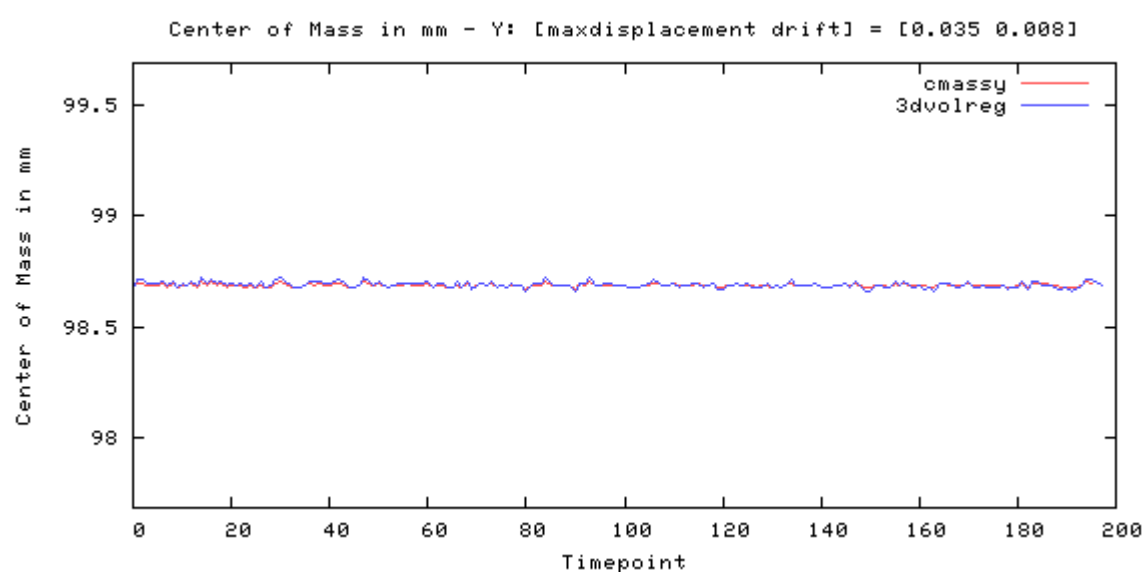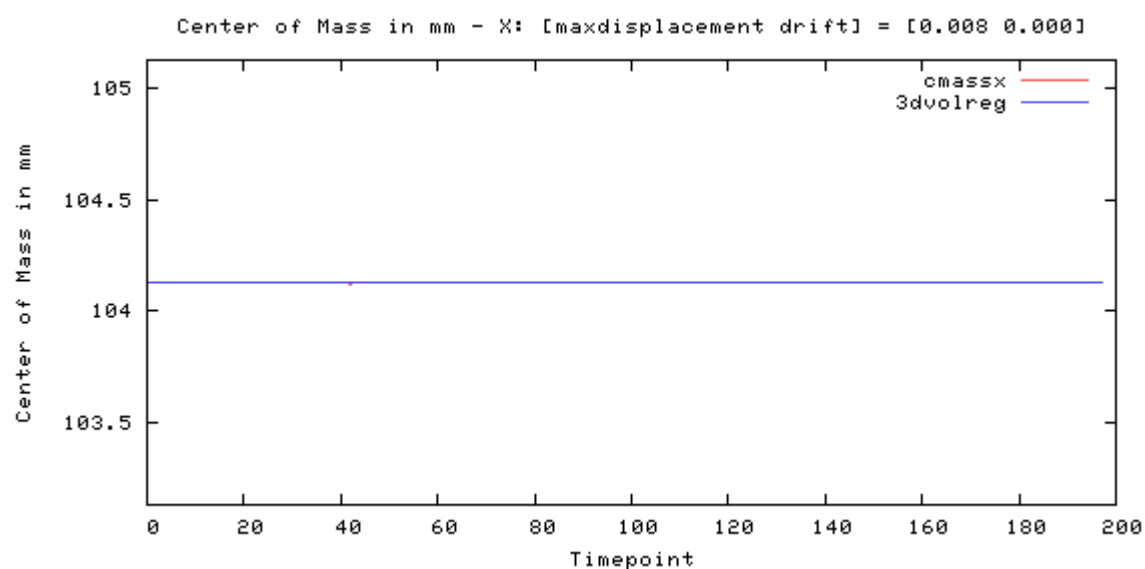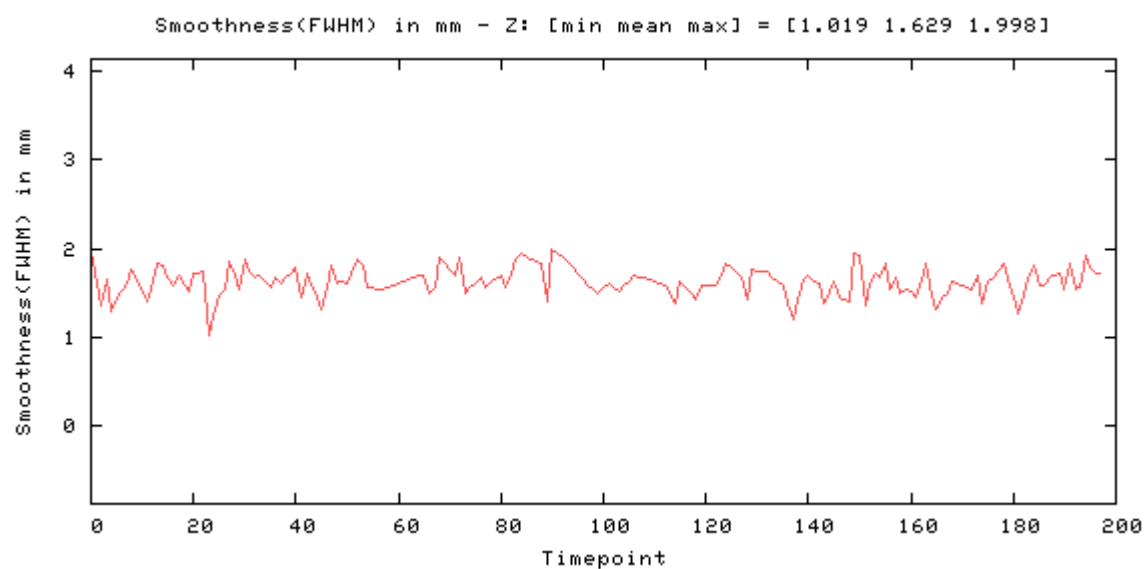


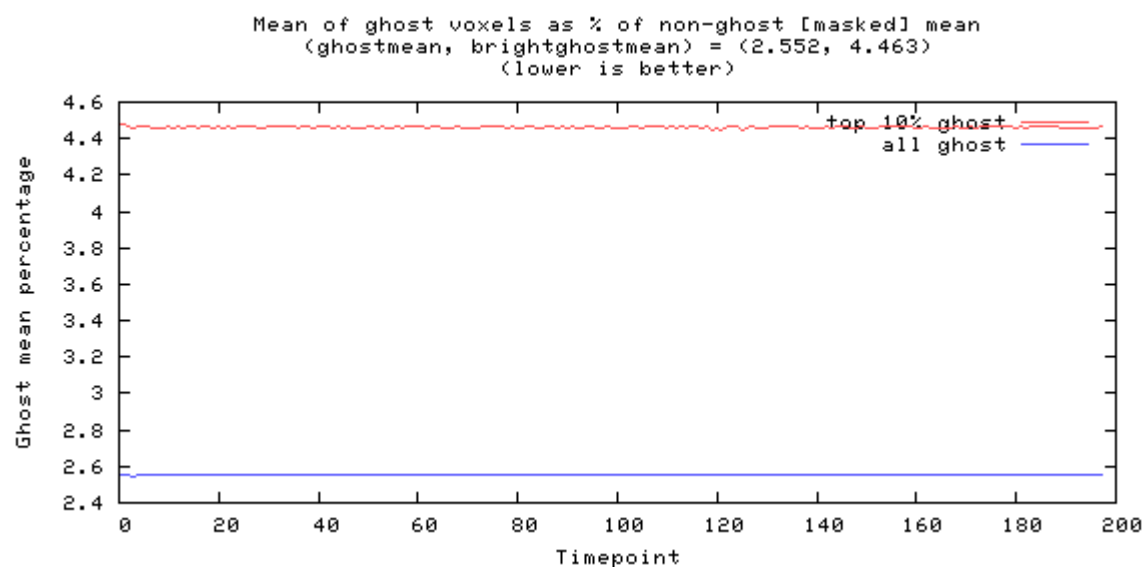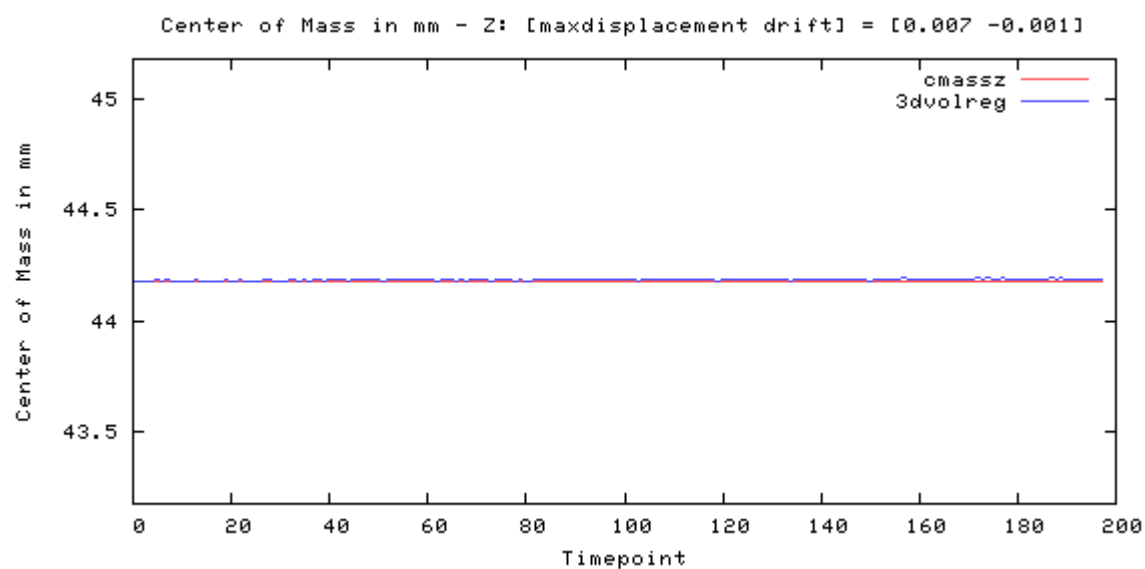

Odd-even difference image



image min: -270, image max: 223

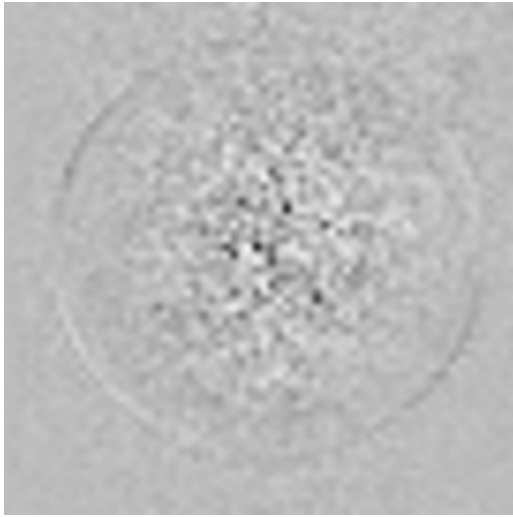

Mean image




image min: 0, image max: 1318.83



Standard Deviation image




image min: 0, image max: 7.86308

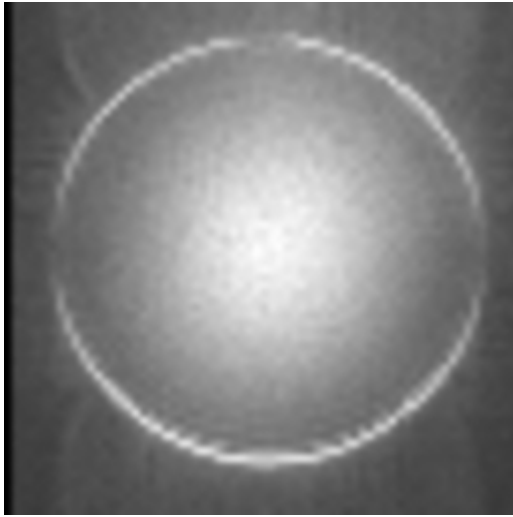

SFNR image



image min: 0, image max: 987.659

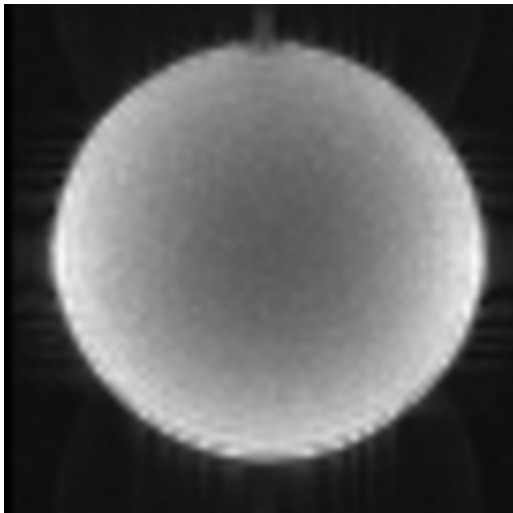

Acquisition parameters

dimensions	68x68x30x200
spacing	3mmx3mmx3mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	8a2add62-cd07-43
studyid	Stability_20251117
seriesnumber	9
runnumber	1
magneticfield	3

description	ses-64ch_func-bold_acq-siemensbasic_task-tsnr
scandate	2025-11-17
scantime	09:46:40
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	3
frequencydirection	1
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-64ch_func-bold_acq-siemensbasic_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20251117094638.422500
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	3
echotrainlength	68
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	PRODUCT
acquisitionduration	405.91500000000002

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	30
rows	68
columns	68
bitsallocated	16
bitsstored	12
highbit	11
lossyimagecompression	00