






Phantom QA results for / fbirn_out/ ses32ch_funcbold_acqcmrrstd_tasktsnr slice 34




bold_acqcmrrstd_tasktsnr_19.xml [percent fluct (trend removed), drift, driftfit] :


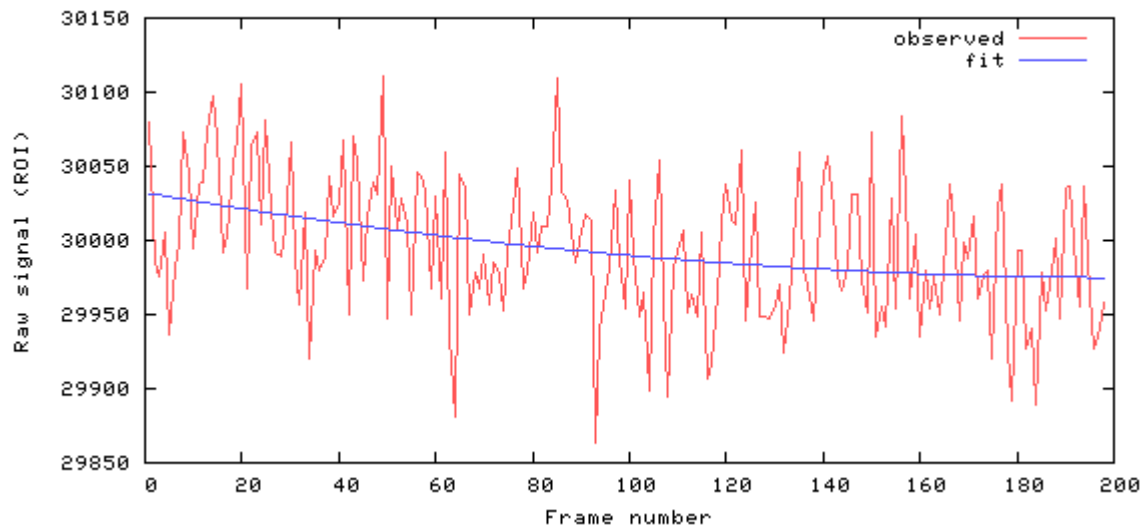


[mean, SNR, SFNR] = [29994.5 137.4 137.1]

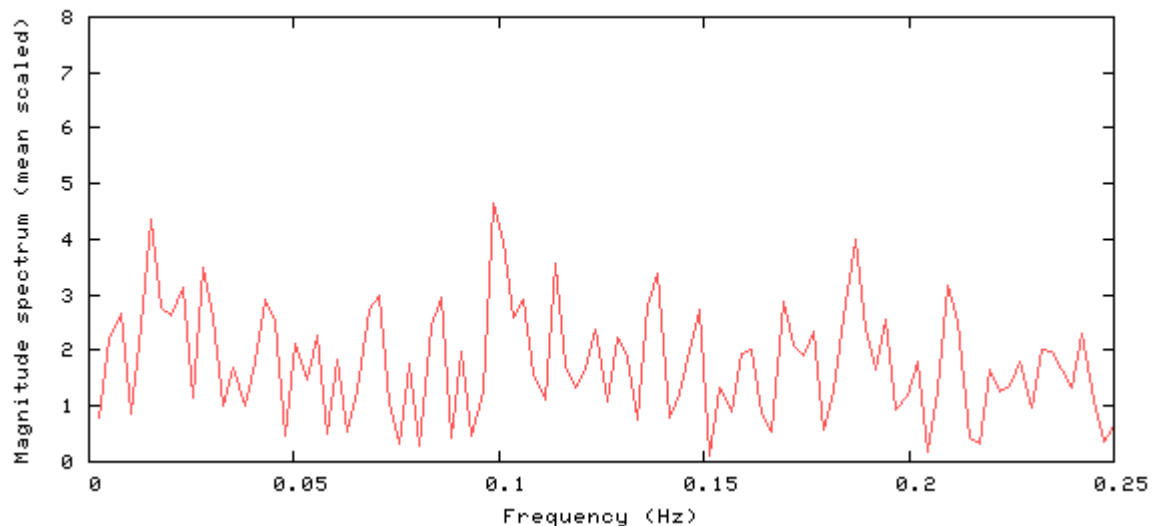

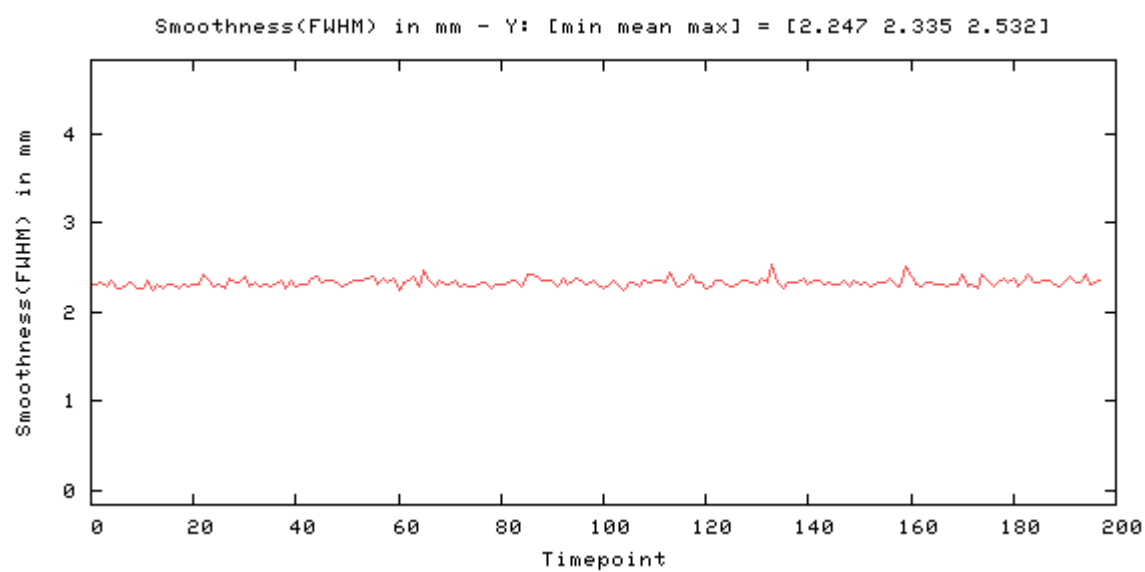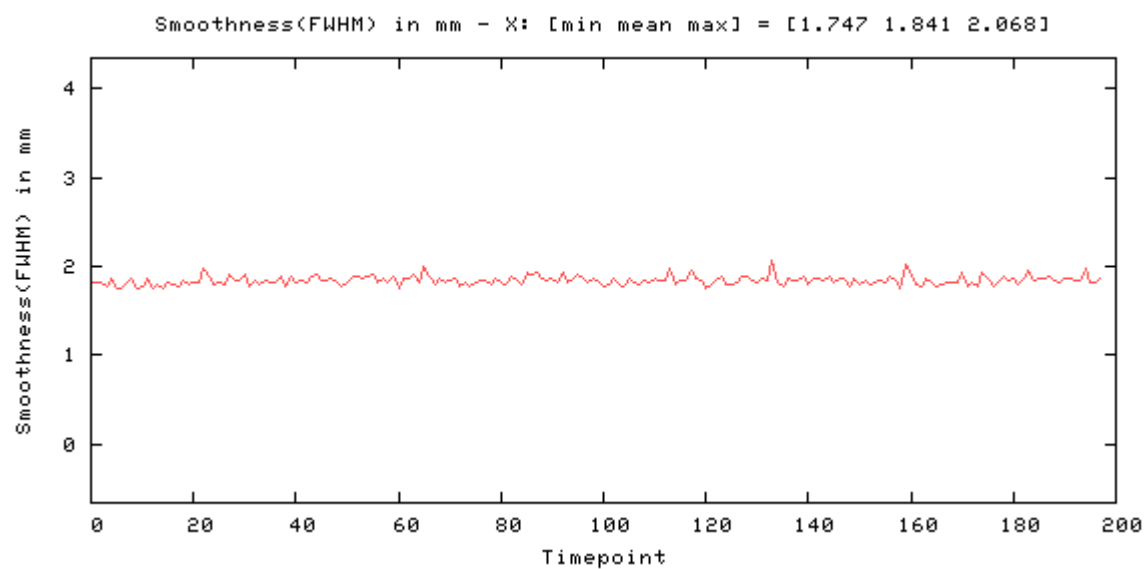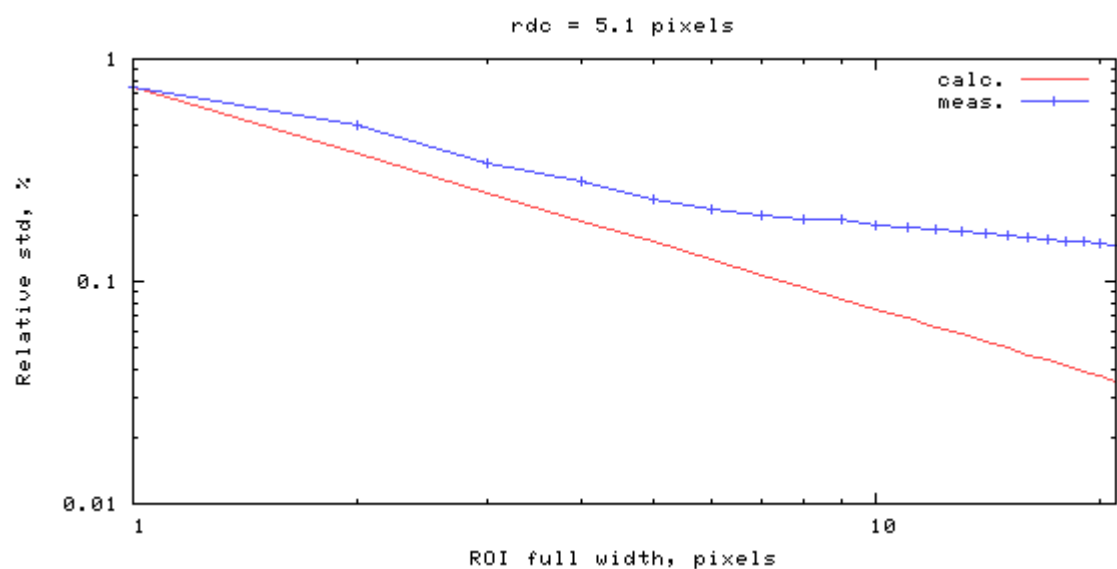

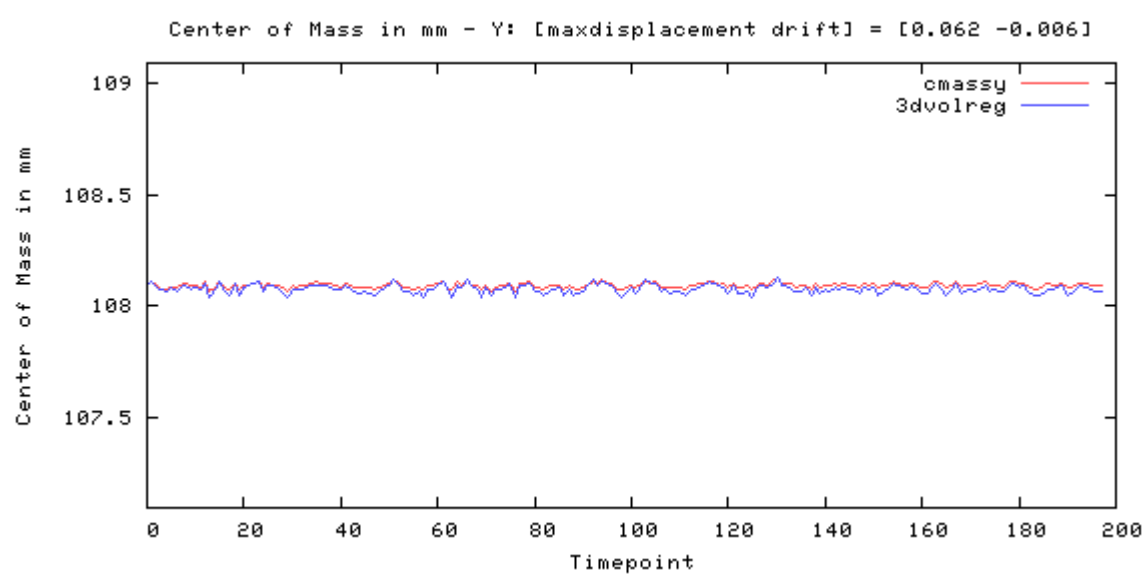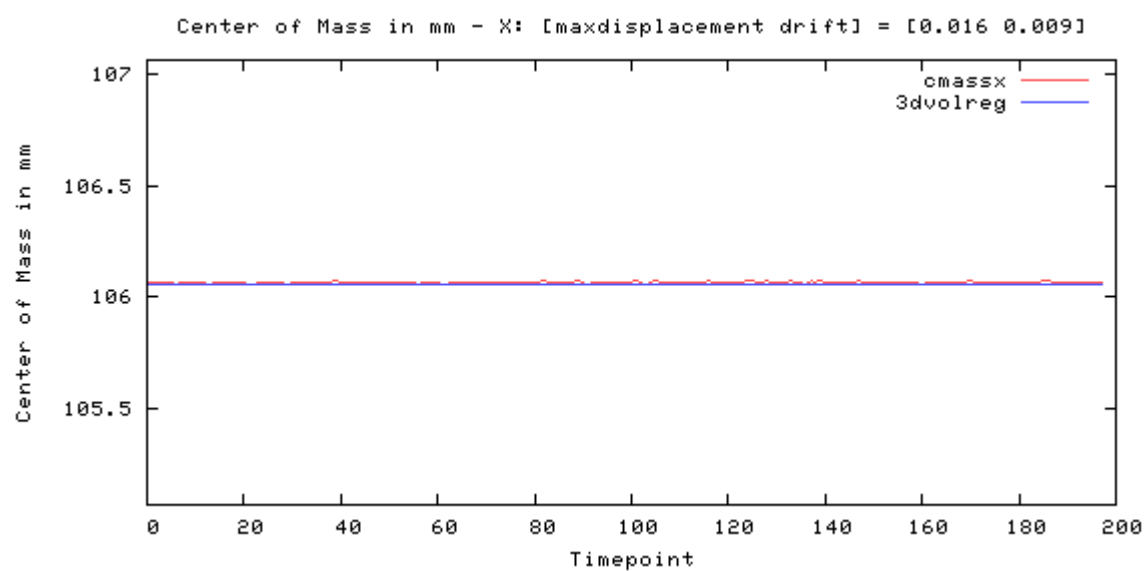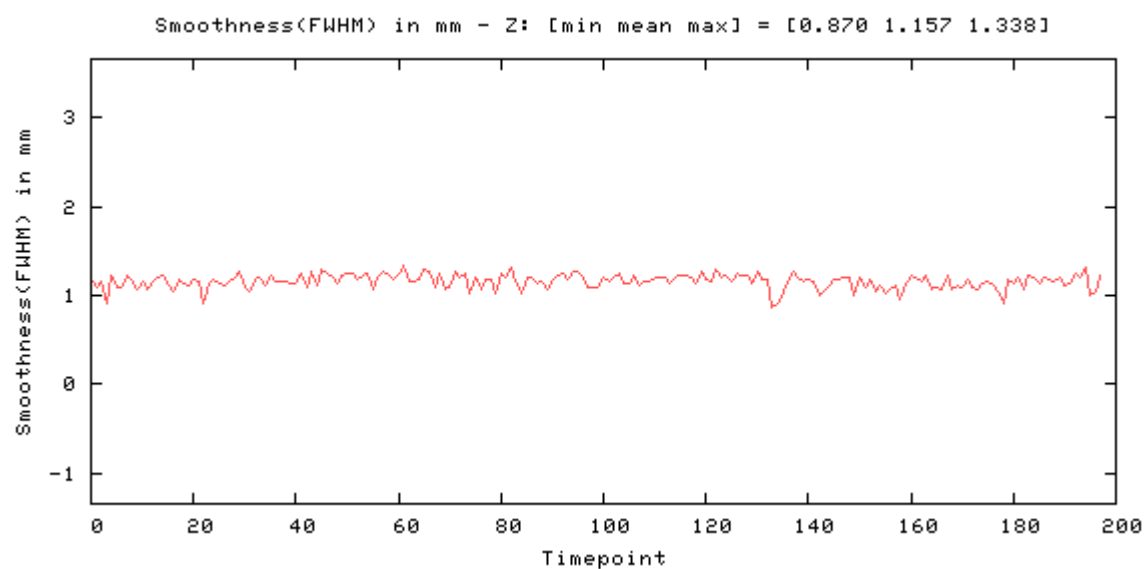


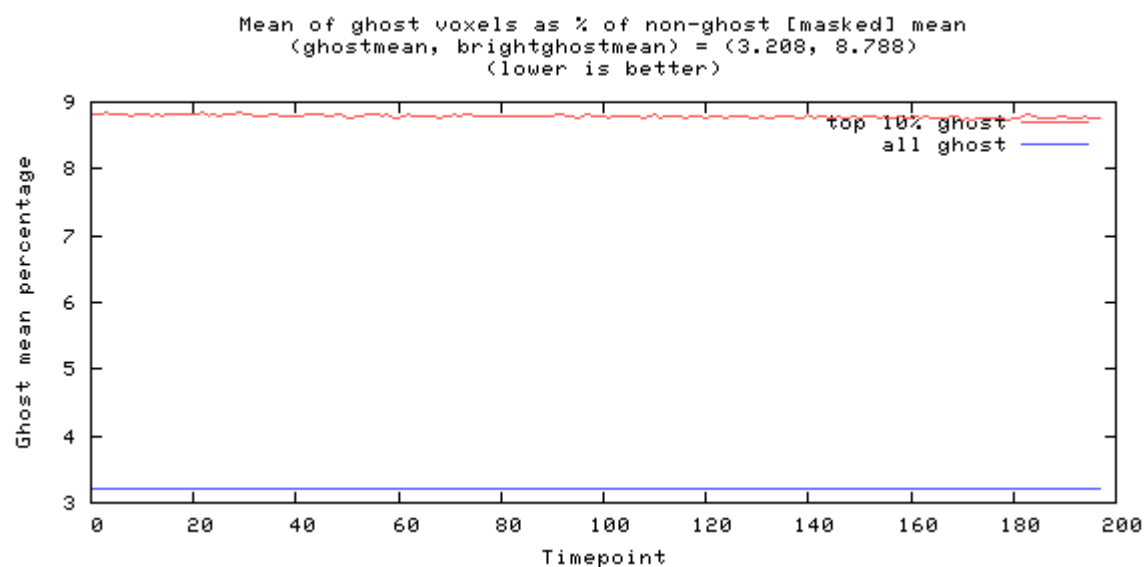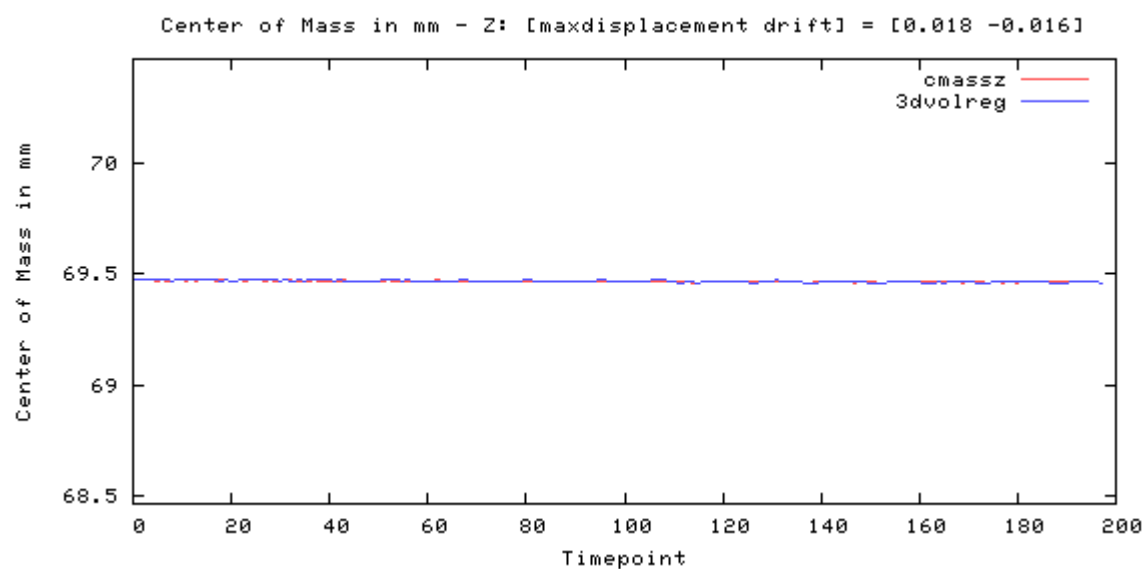

Odd-even difference image




image min: -11623, image max: 9068

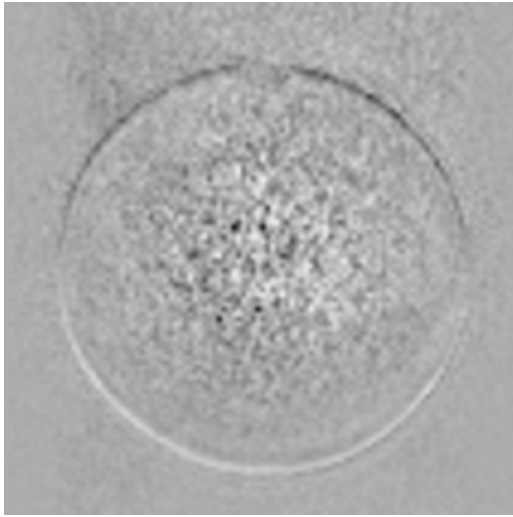

Mean image




image min: 0, image max: 40307.1

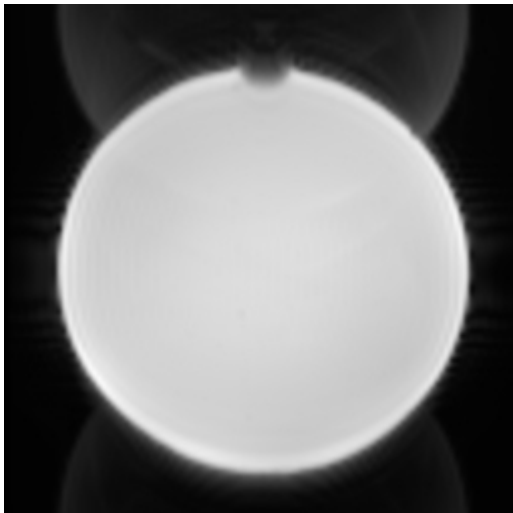

Standard Deviation image




image min: 0, image max: 286.814

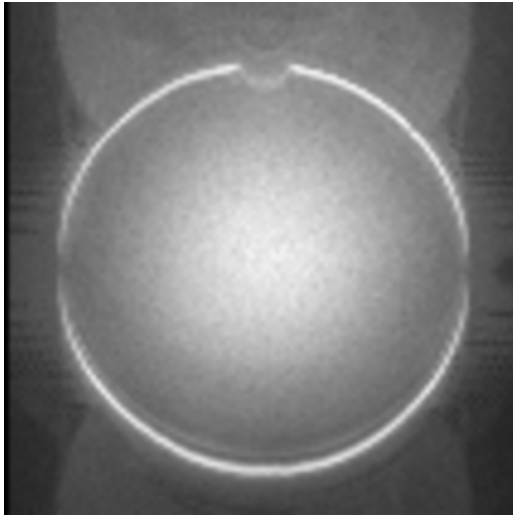

SFNR image



image min: 0, image max: 757.452

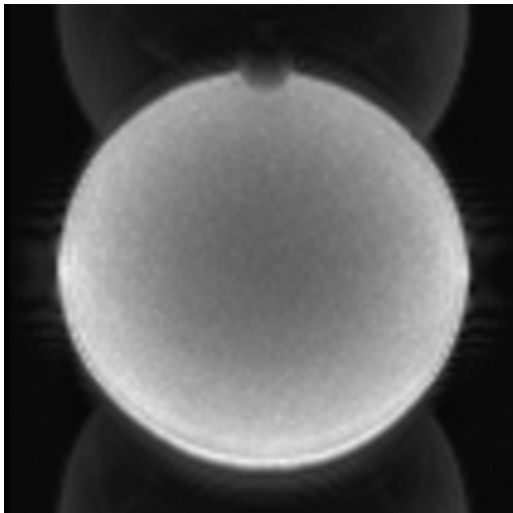

Acquisition parameters

dimensions	106x106x69x200
spacing	2mmx2mmx2mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	998d4285-9c7f-4d
studyid	Stability_20251201
seriesnumber	19
runnumber	1
magneticfield	3

description	ses-32ch_func-bold_acq-cmrrstd_task-tsnr
scandate	2025-12-01
scantime	11:00:21
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
receivecoilname	Head_32
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-32ch_func-bold_acq-cmrrstd_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20251201110019.220000
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	2
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	RESEARCH
acquisitionduration	417.89499999999998
resonantnucleus	1H

applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
imagecomments	Not for diagnostic use, Unaliased MB3/ PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	69
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00