






Phantom QA results for / fbirn_out/ ses20ch_funcbold_acqcmrrstd_tasktsnr slice 34




bold acqcmrrstd tasktsnr 26.xml [percent fluct (trend removed), drift, driftfit] :


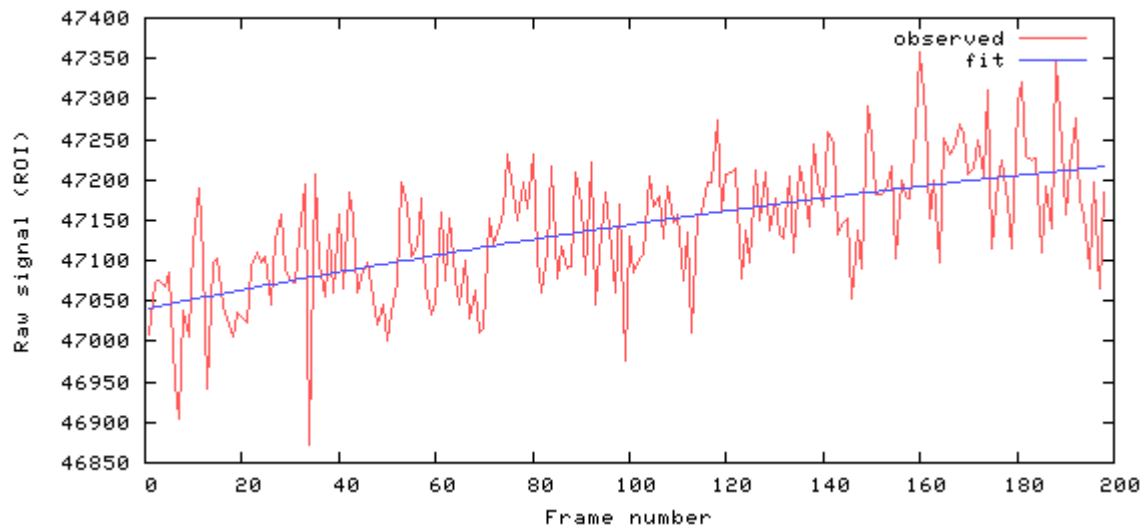


[mean, SNR, SFNR] = [47138.9 167.6 168.9]

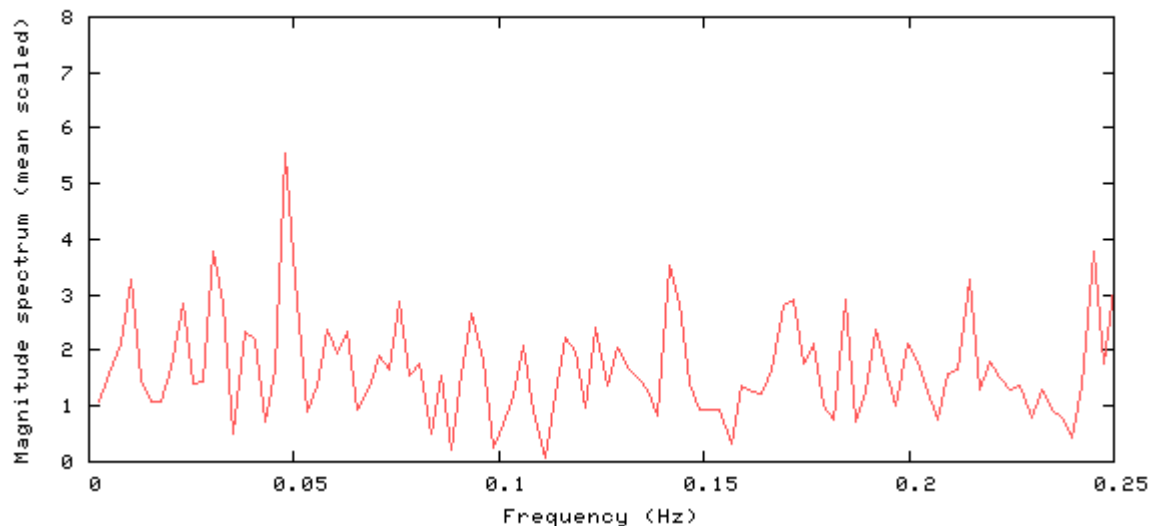

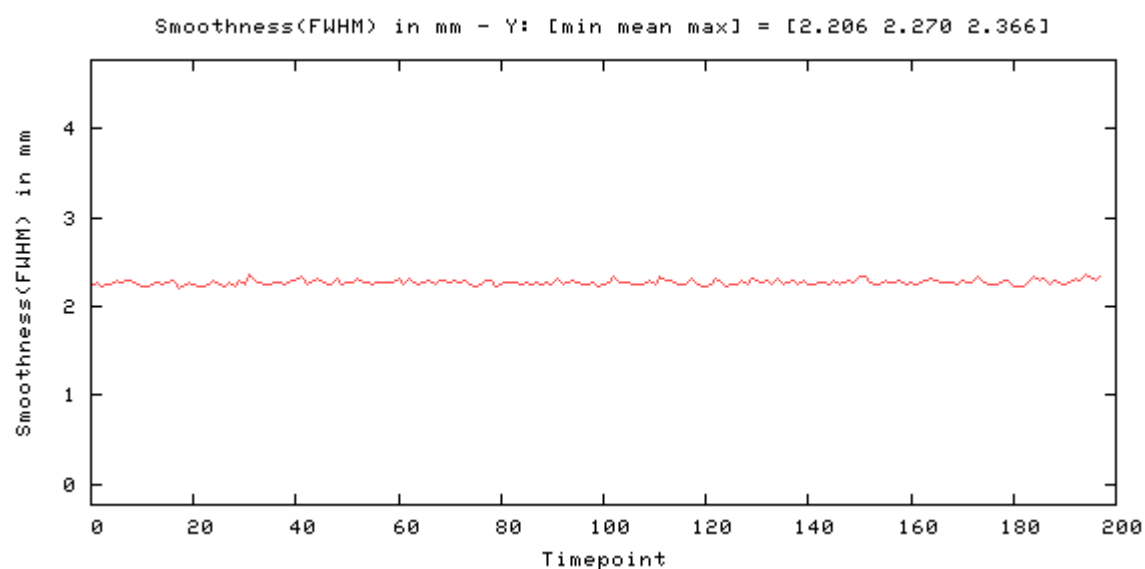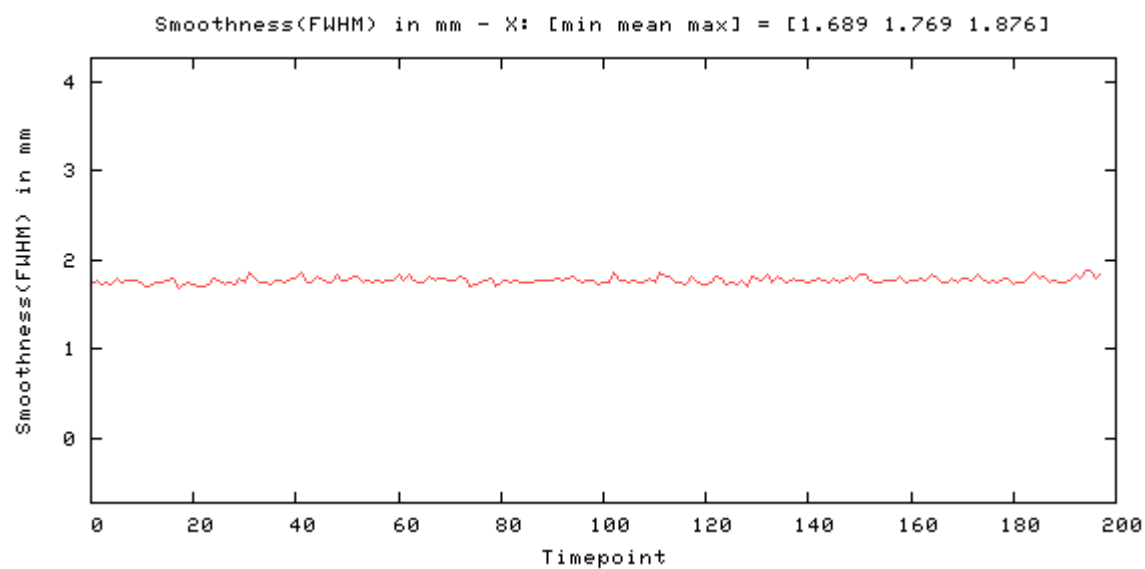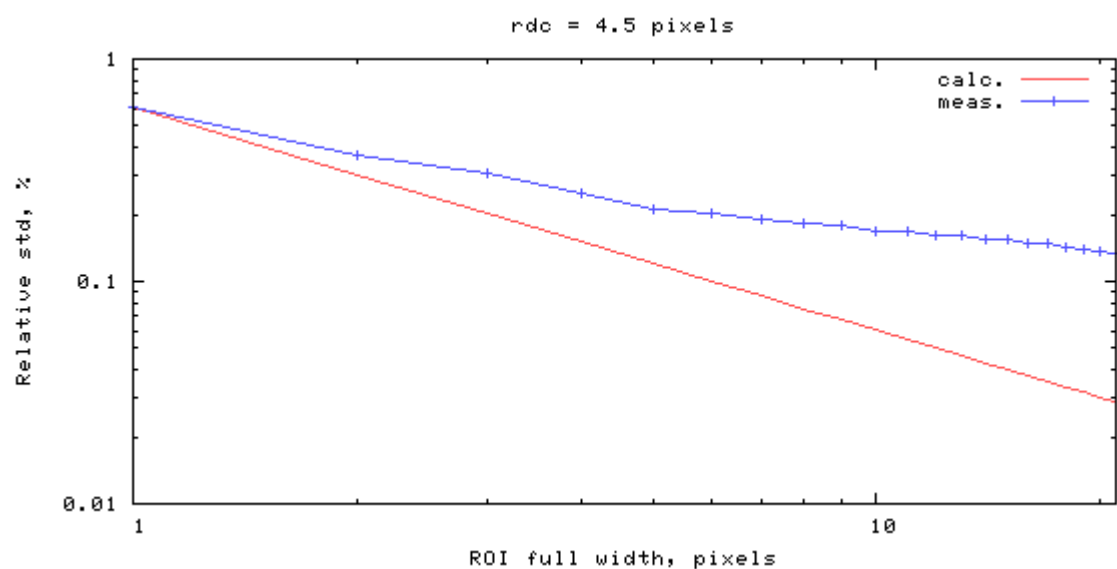

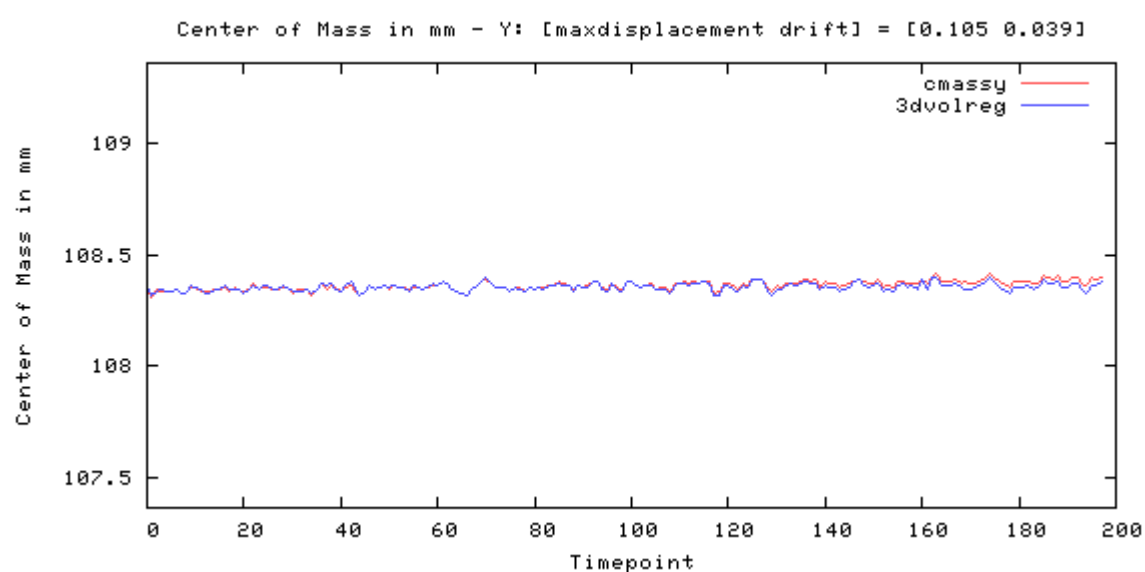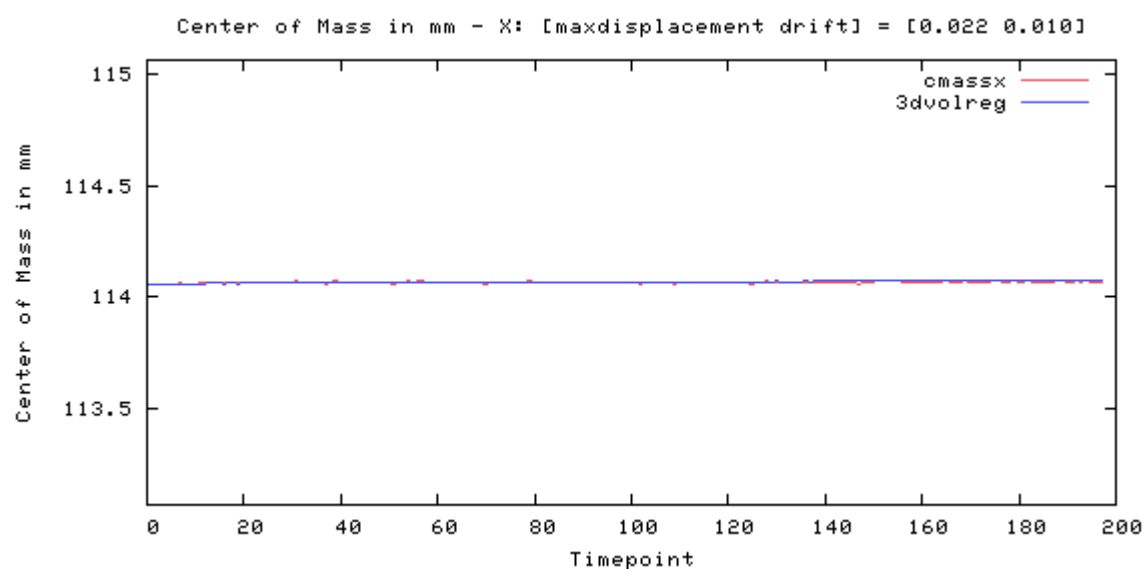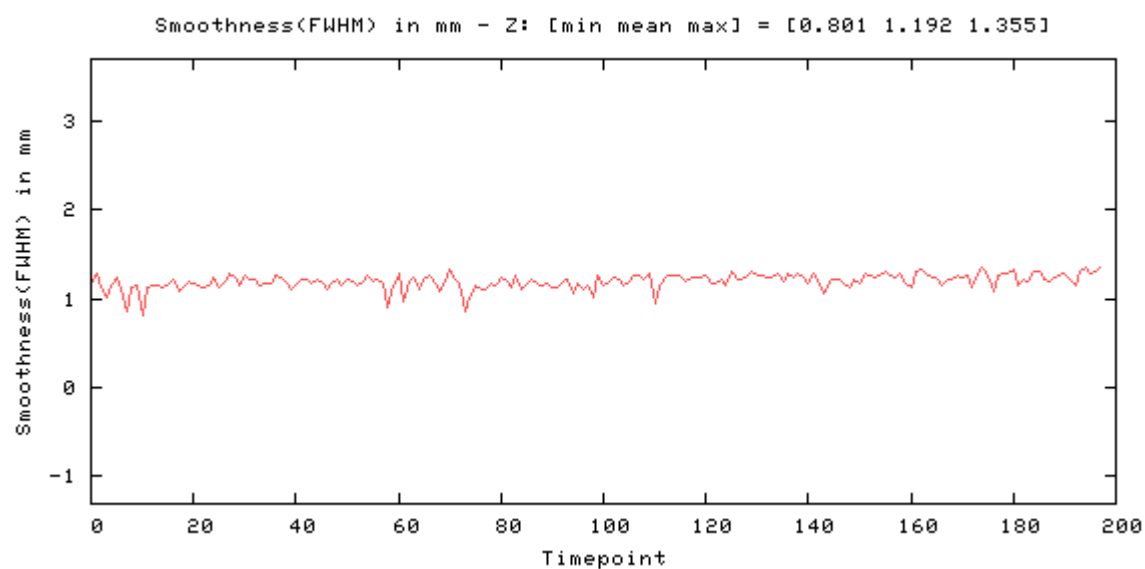

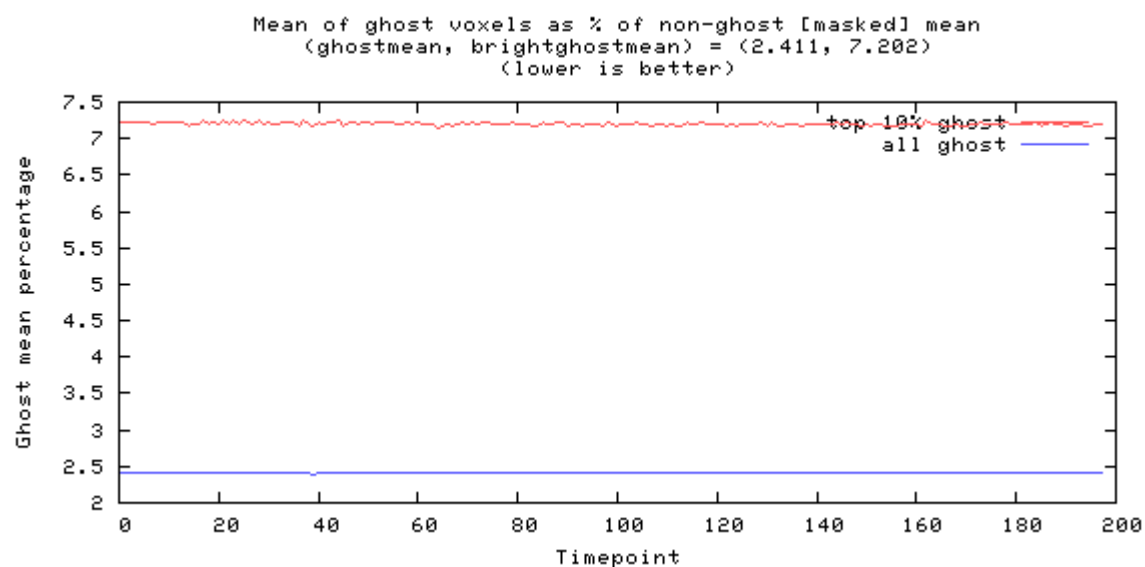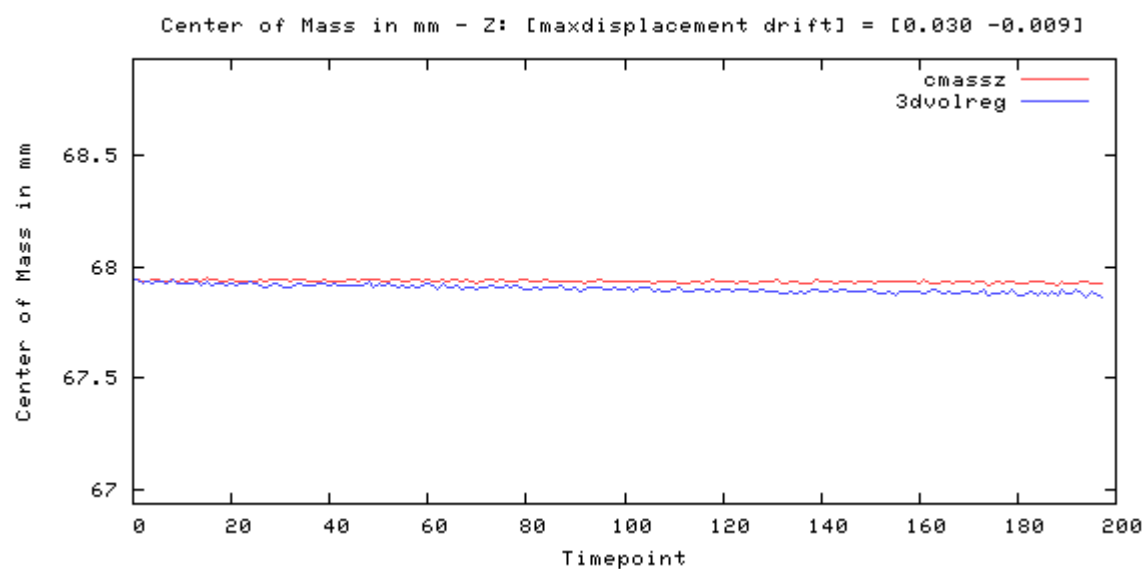


Odd-even difference image

-14538 11106



image min: -14538, image max: 11106

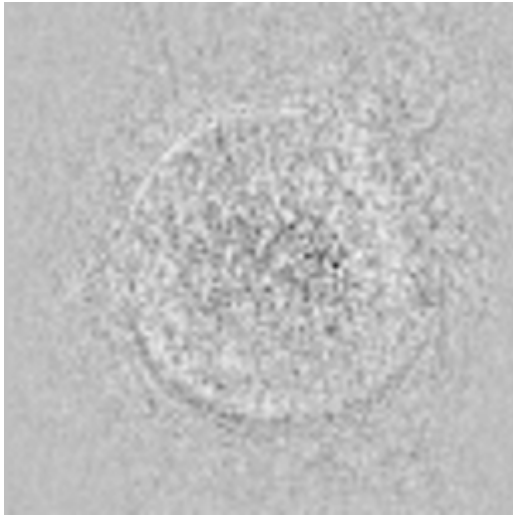

Mean image




image min: 0, image max: 54912.4



Standard Deviation image




image min: 0, image max: 406.878

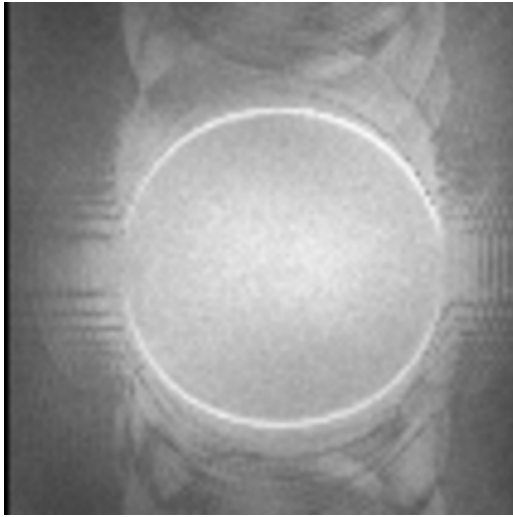

SFNR image



image min: 0, image max: 329.649

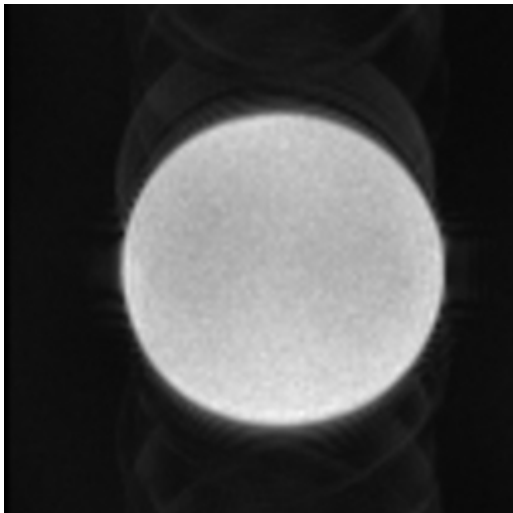

Acquisition parameters

dimensions	106x106x69x200
spacing	2mmx2mmx2mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	e90588c0-a506-4a
studyid	Stability_20251208
seriesnumber	26
runnumber	1
magneticfield	3

description	ses-20ch_func-bold_acq-cmrrstd_task-tsnr
scandate	2025-12-08
scantime	10:55:38
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
receivecoilname	HeadNeck_20_TCS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-20ch_func-bold_acq-cmrrstd_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20251208105536.505000
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	2
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	RESEARCH
acquisitionduration	417.892
resonantnucleus	1H

applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
imagecomments	Not for diagnostic use, Unaliased MB3/ PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	69
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00