






Phantom QA results for /fbirn_out/ses32ch_funcbold_acqsiemensstd_task slice 17




.d acqsiemensstd tasktsnr 21.xml [percent fluct (trend removed), drift, driftfit]


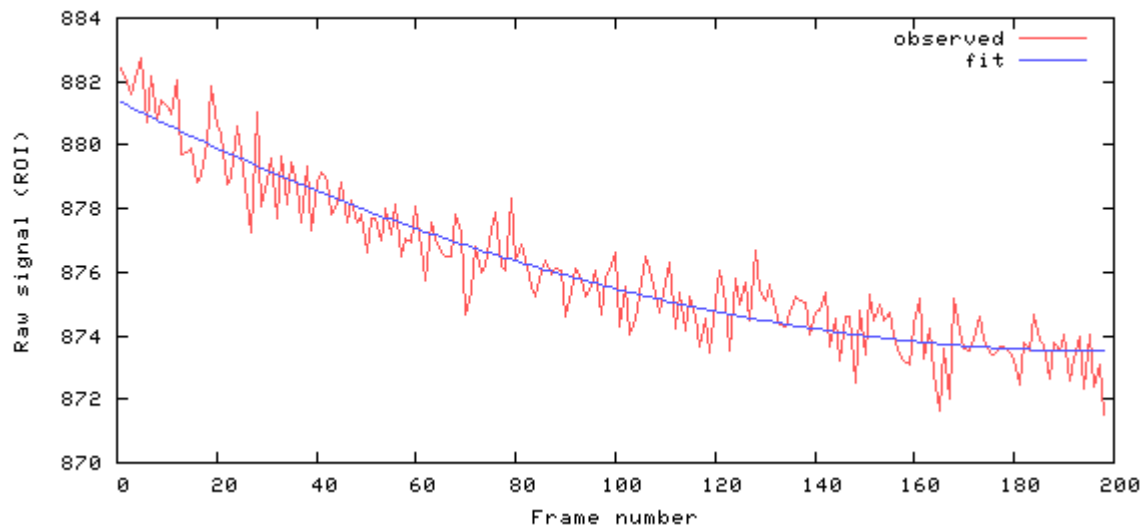


[mean, SNR, SFNR] = [876.2 124.3 123.6]

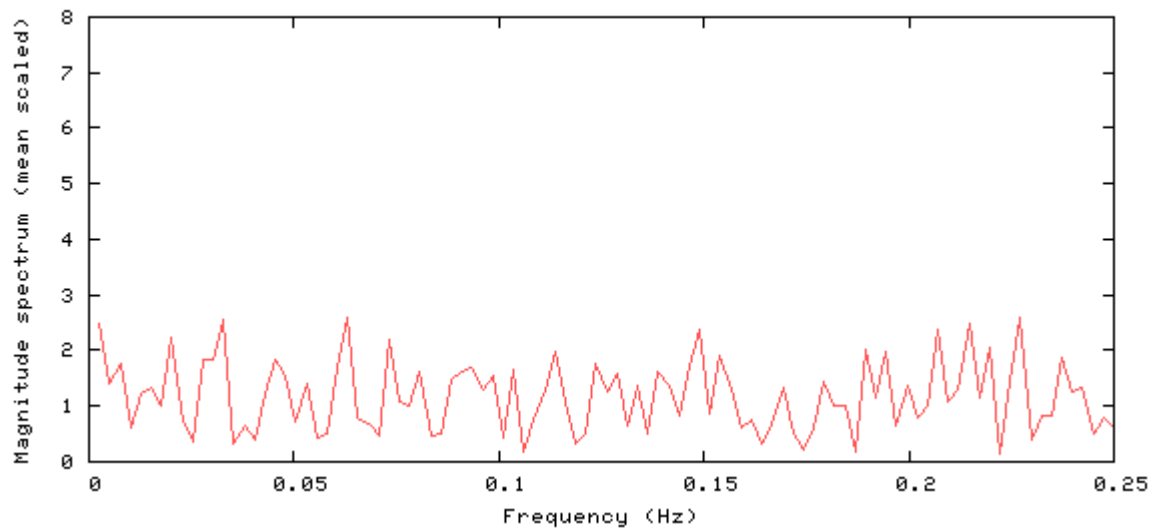

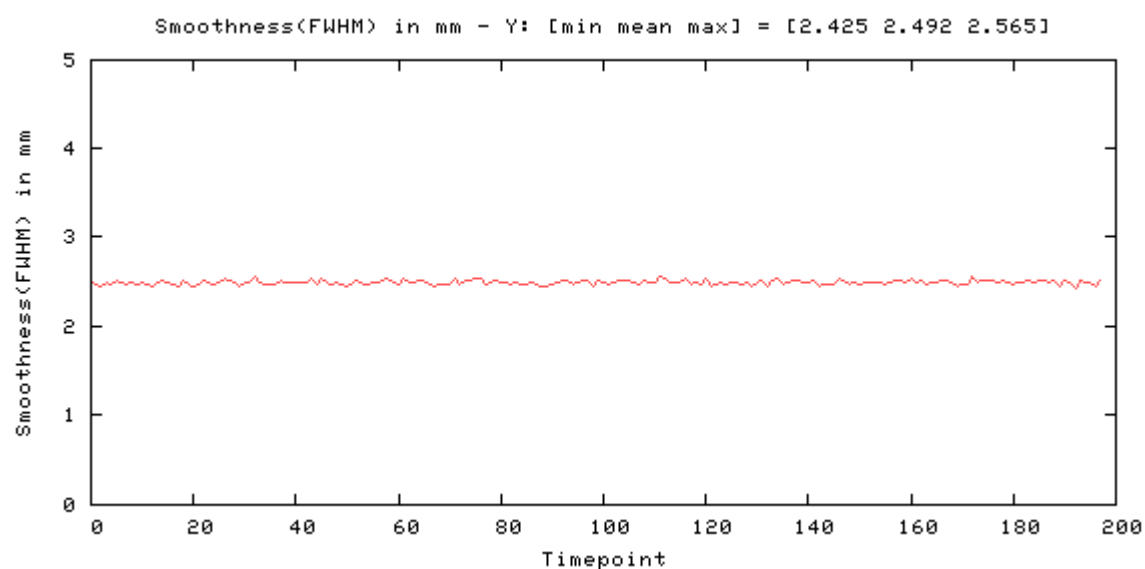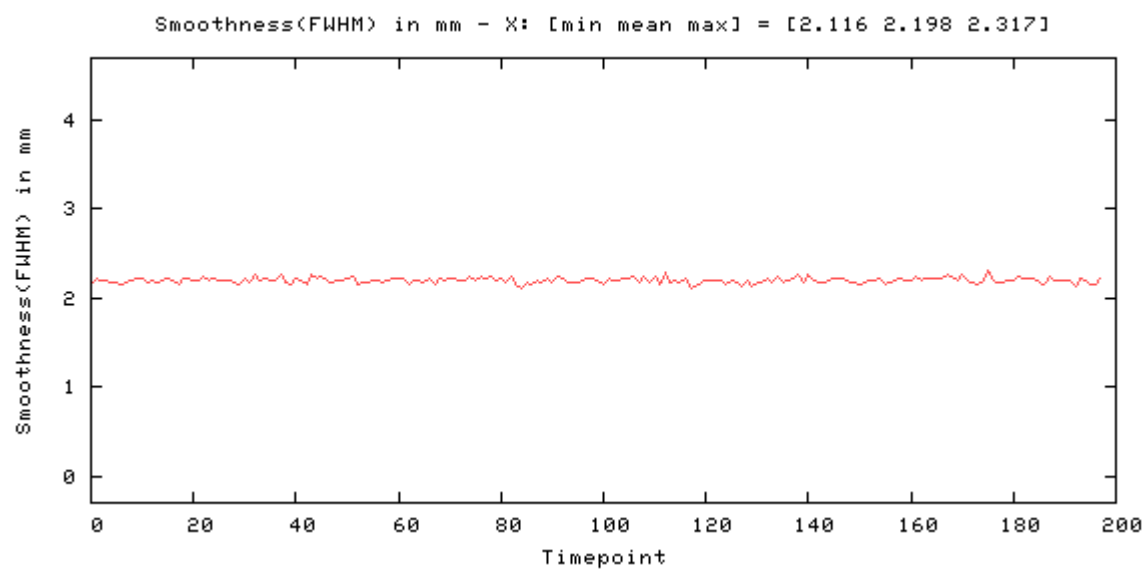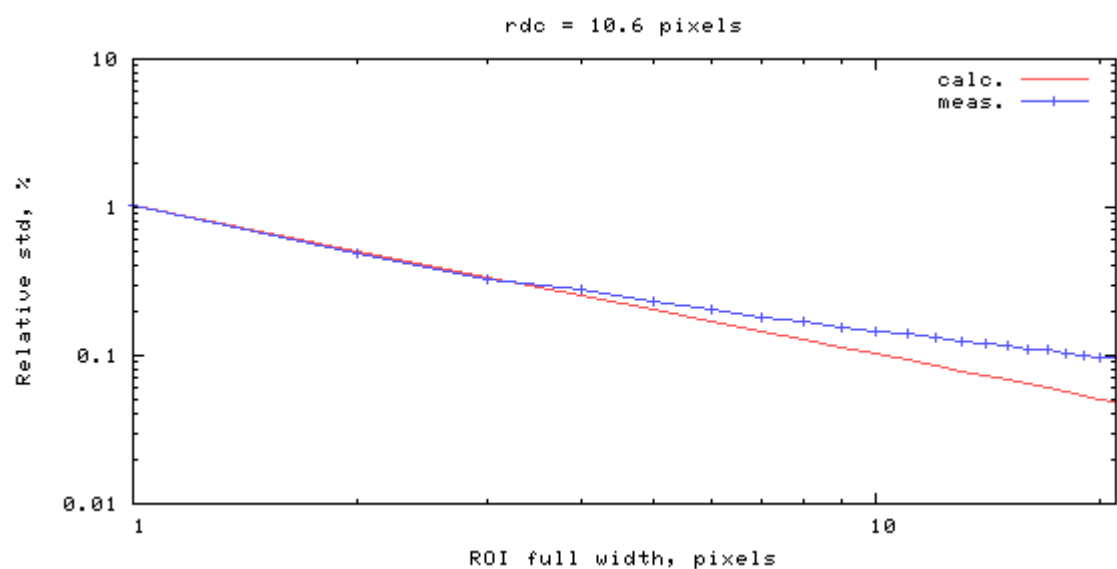

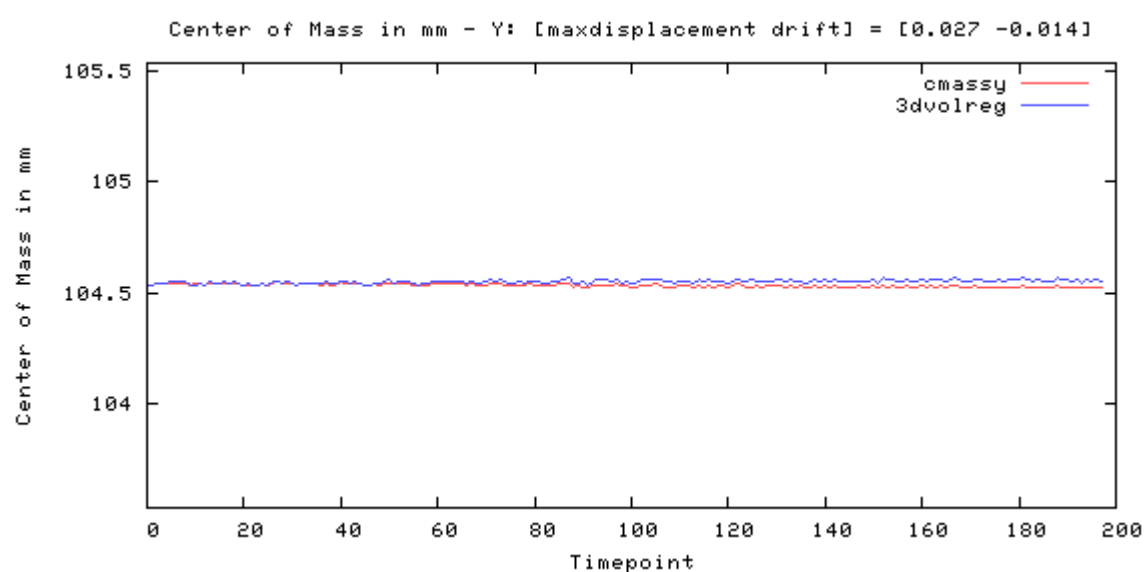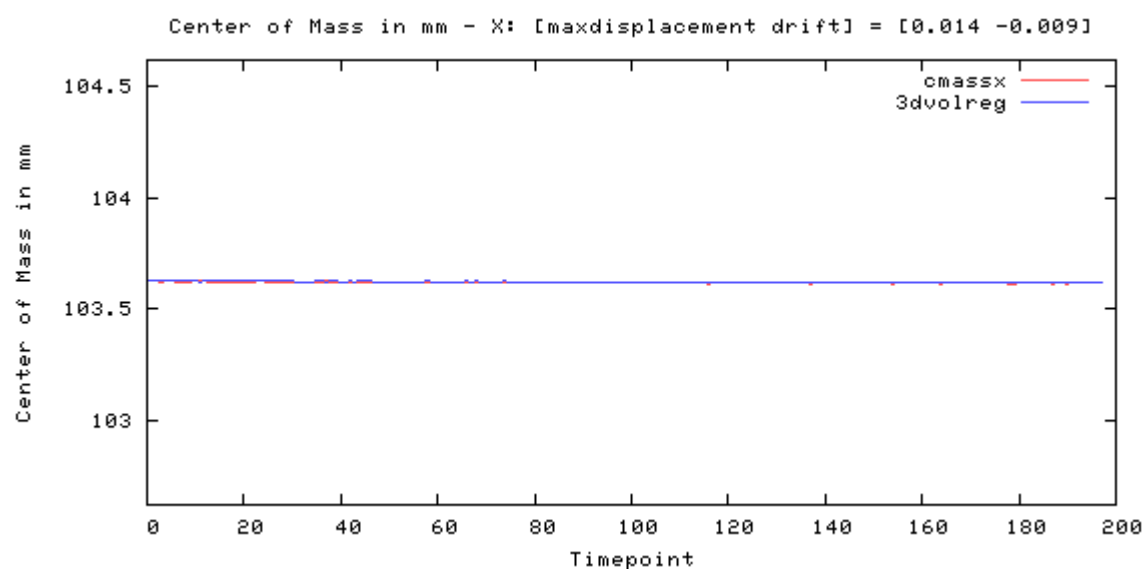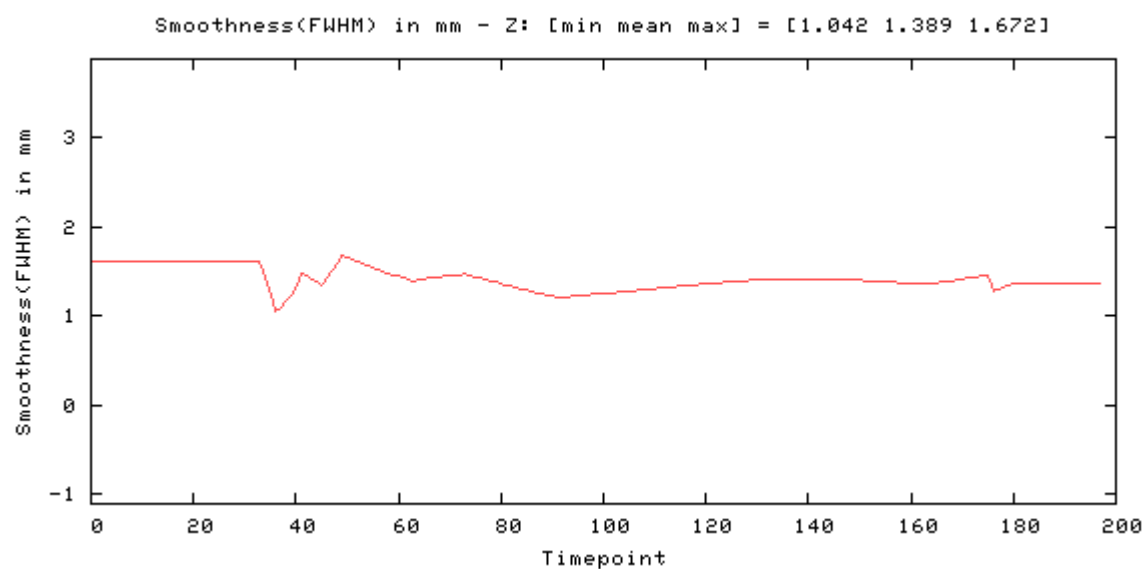


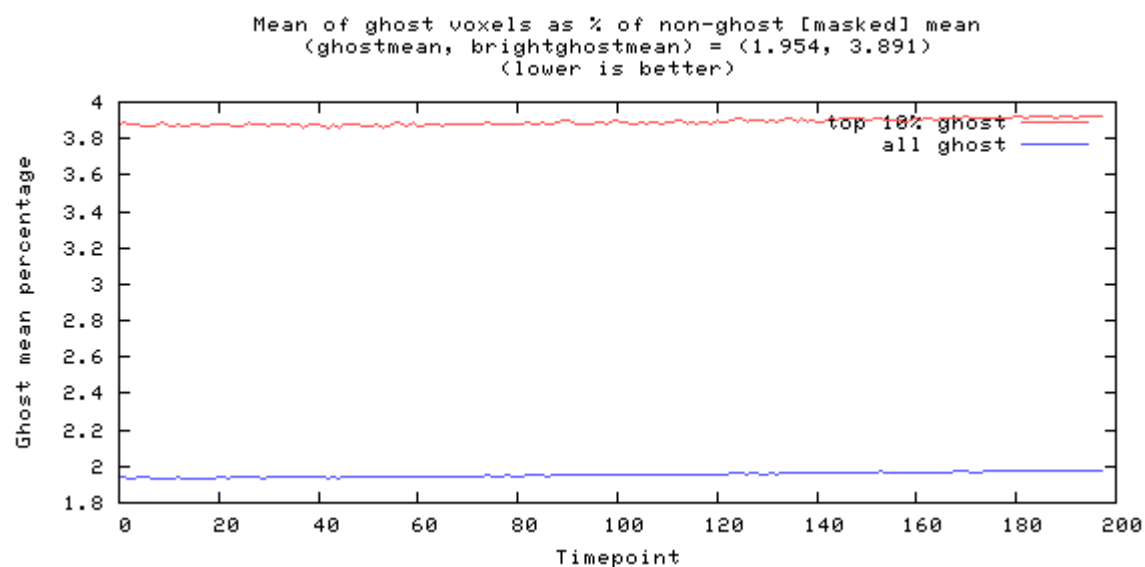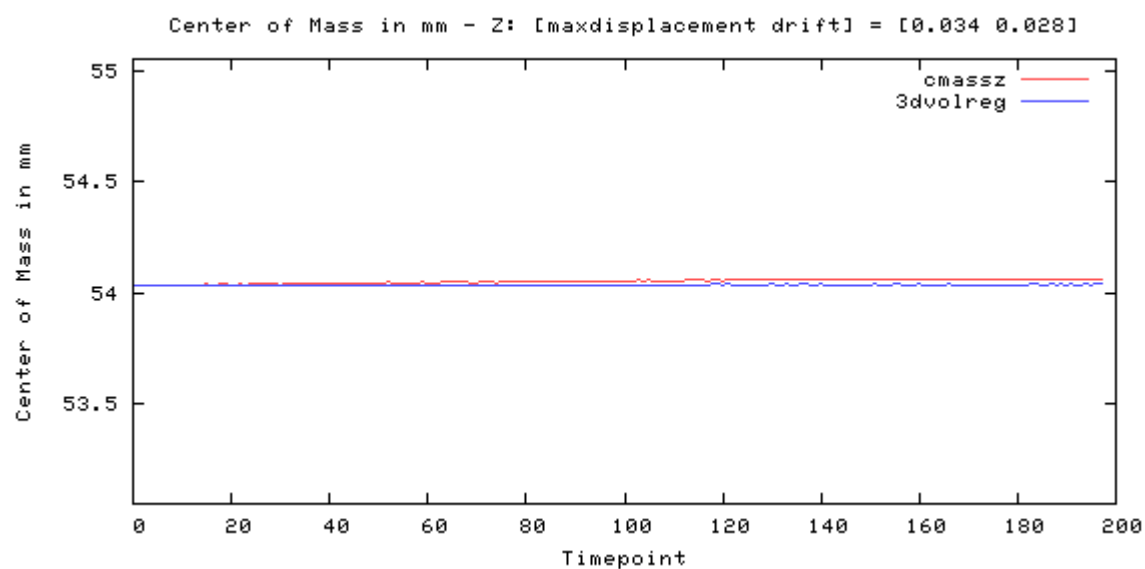

Odd-even difference image



image min: -297, image max: 290

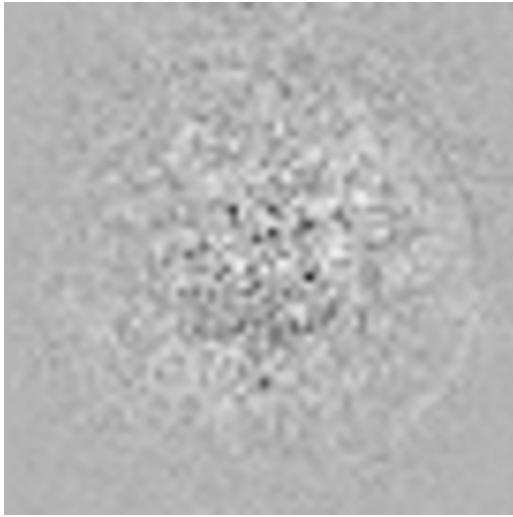

Mean image




image min: 0, image max: 1202.15



Standard Deviation image




image min: 0, image max: 9.28731

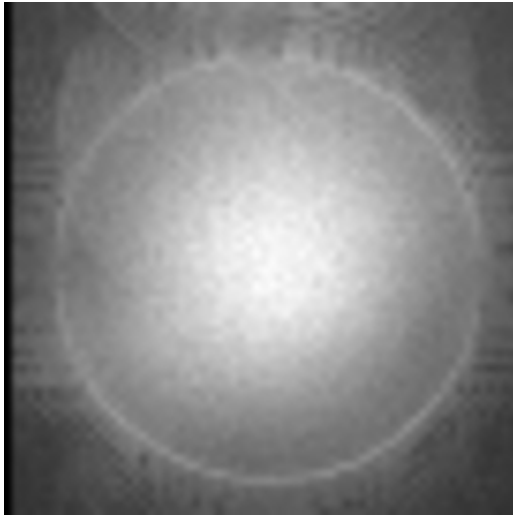

SFNR image



image min: 0, image max: 594.087

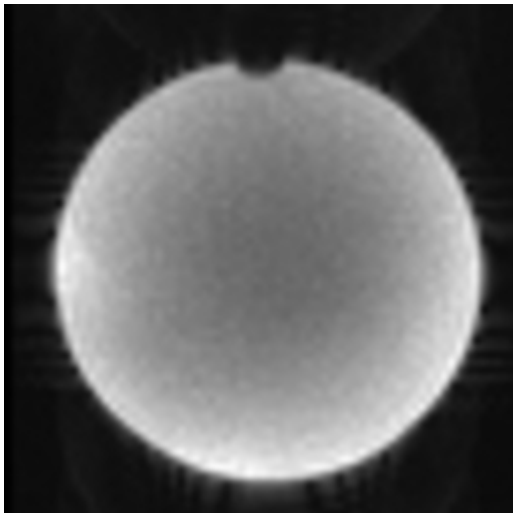

Acquisition parameters

dimensions	68x68x36x200
spacing	3mmx3mmx3mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	0adb9016-f407-46
studyid	2025.12.15-10:13:10-STD-1.3.12.2.1107.5.99.3
seriesnumber	21
runnumber	1

magneticfield	3
description	ses-32ch_func-bold_acq-siemensstd_task-tsnr
scandate	2025-12-15
scantime	11:07:20
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	3
frequencydirection	1
receivecoilname	Head_32
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-32ch_func-bold_acq-siemensstd_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20251215110718.405000
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	3
echotrainlength	33
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	PRODUCT

acquisitionduration	407.92500000000001
resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	36
rows	68
columns	68
bitsallocated	16
bitsstored	12
highbit	11
lossyimagecompression	00