






Phantom QA results for / fbirn_out/ ses20ch_funcbold_acqcmrrstd_tasktsnr slice 34




bold_acqcmrrstd_tasktsnr_28.xml [percent fluct (trend removed), drift, driftfit] :


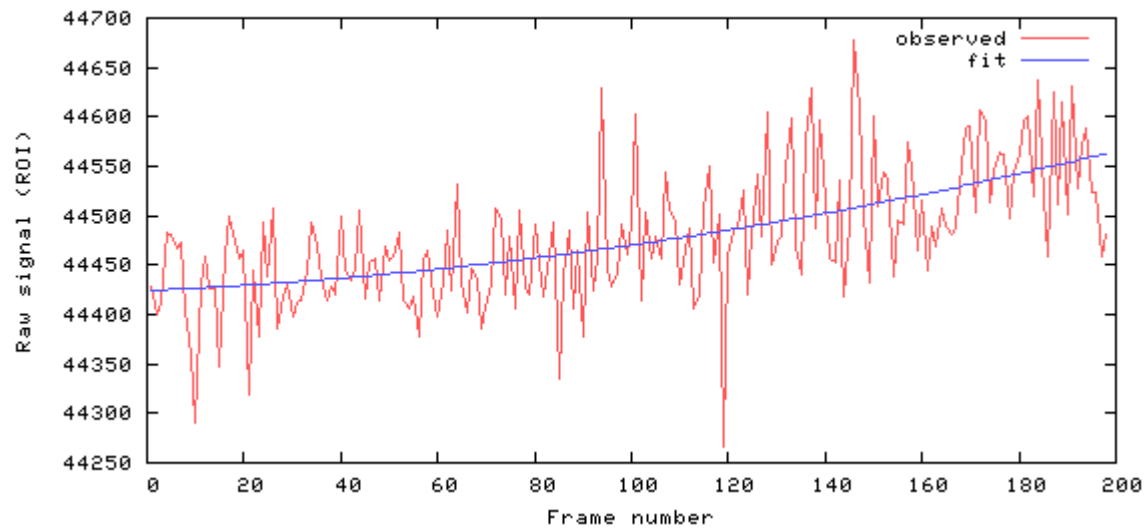


[mean, SNR, SFNR] = [44478.0 163.2 165.3]

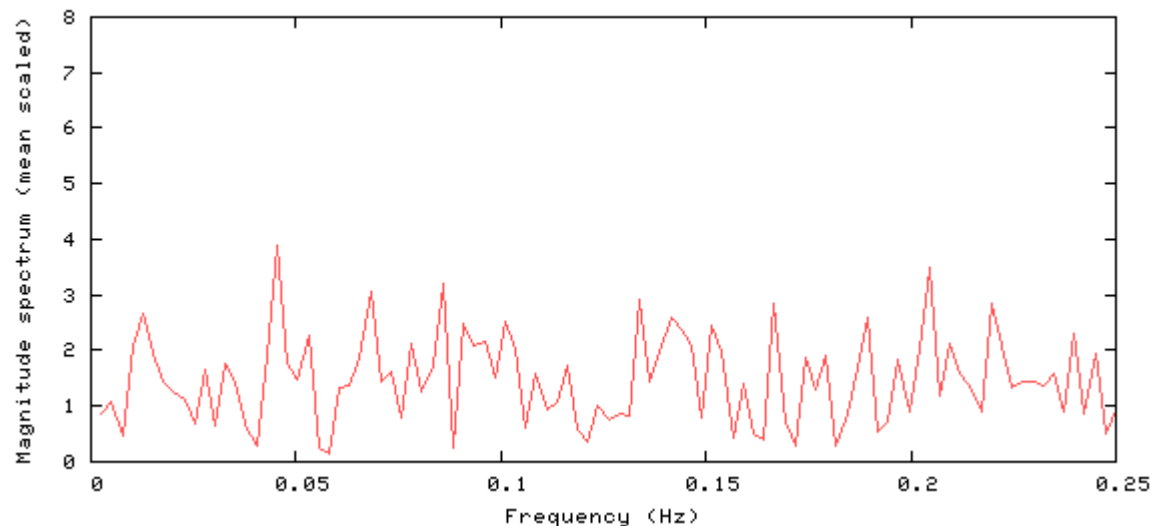

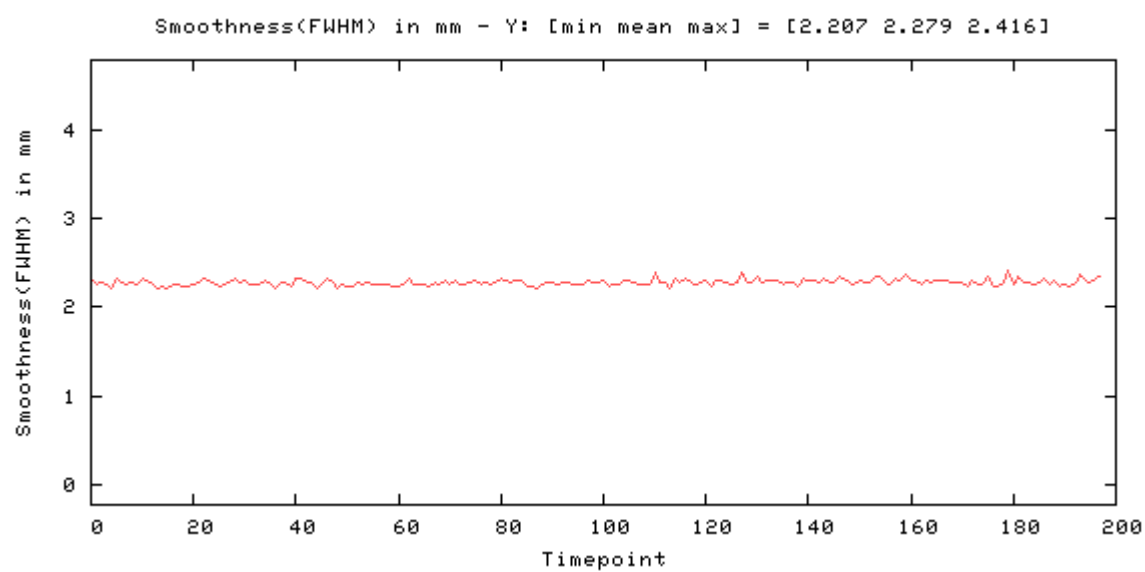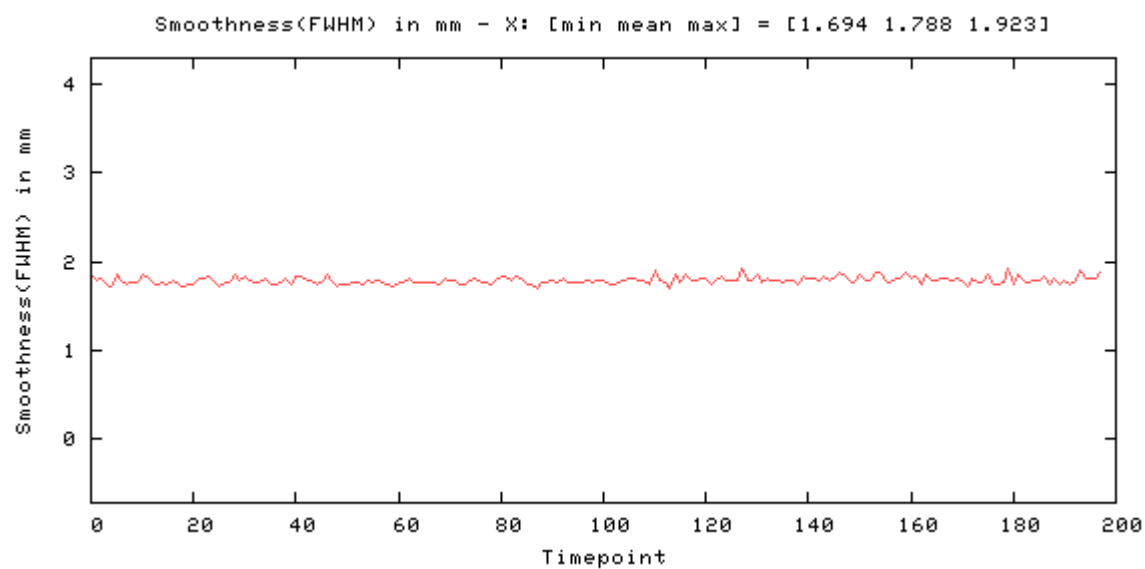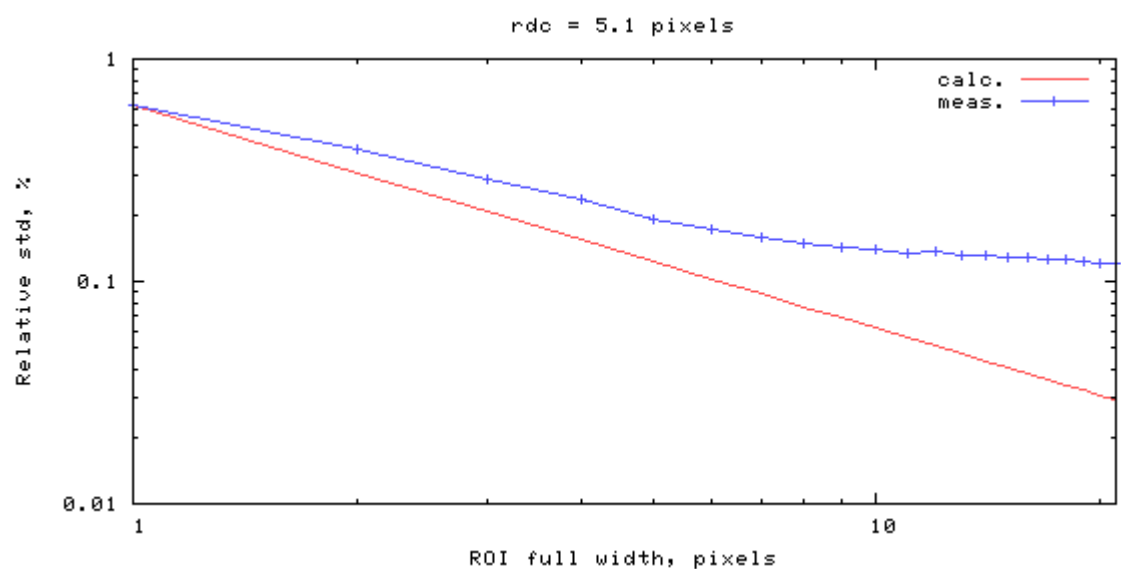

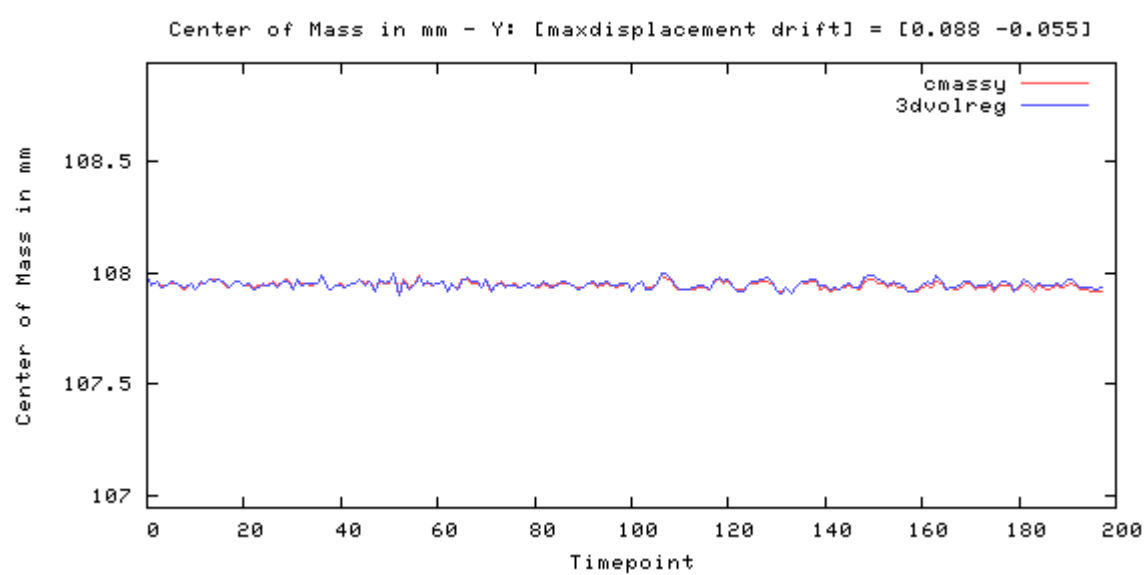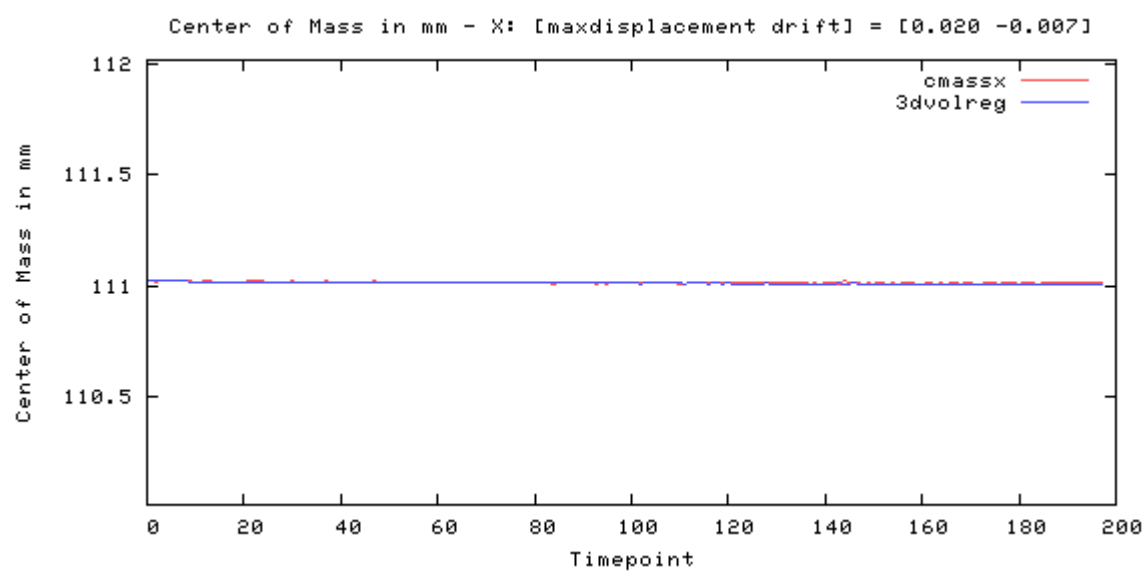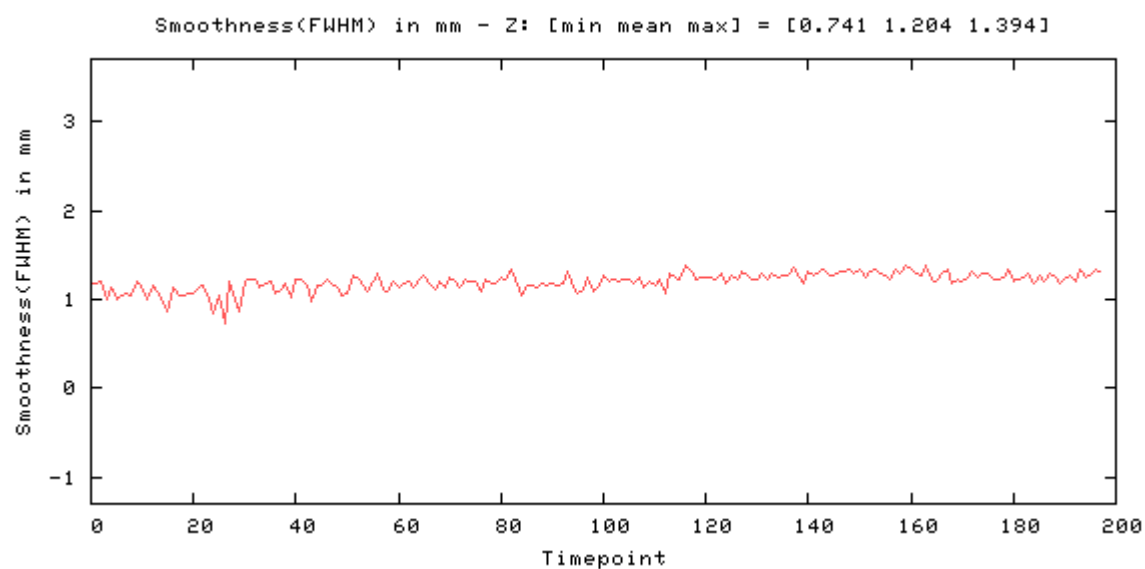

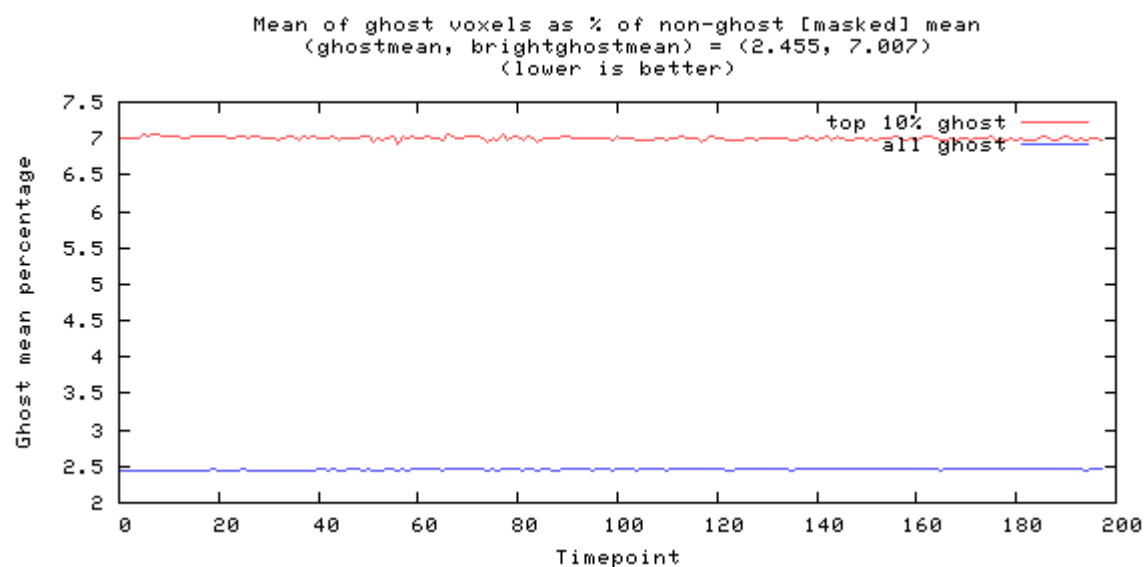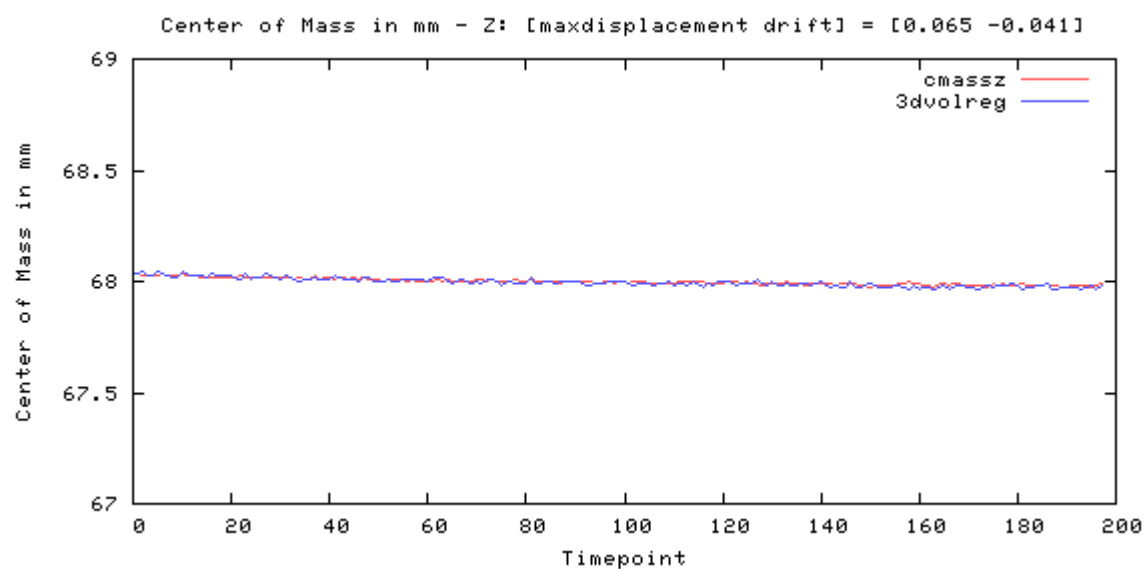


Odd-even difference image

-12342 13840



image min: -12342, image max: 13840

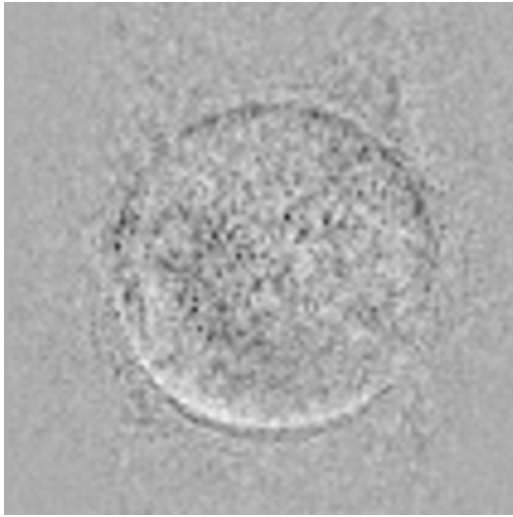

Mean image




image min: 0, image max: 54372.4



Standard Deviation image




image min: 0, image max: 416.698

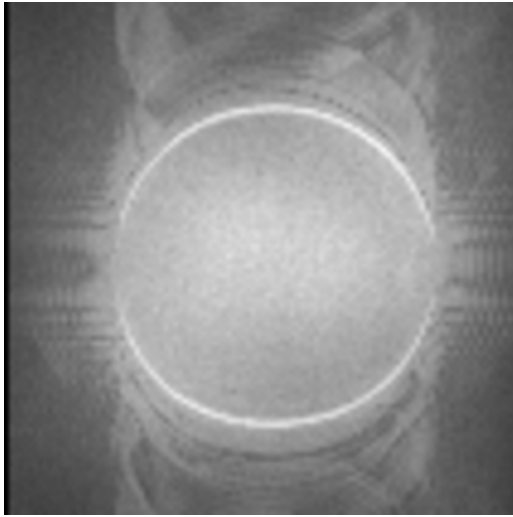

SFNR image



image min: 0, image max: 318.137

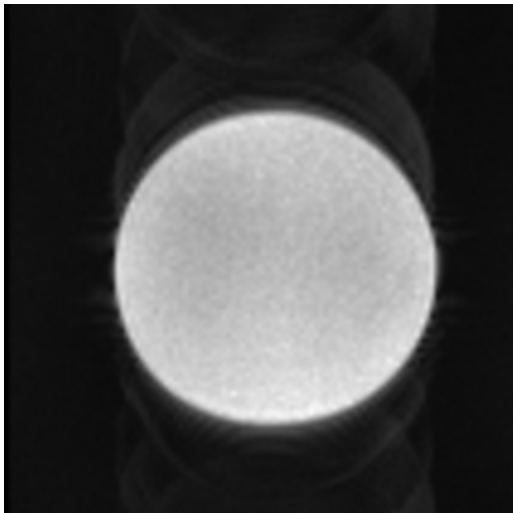

Acquisition parameters

dimensions	106x106x69x200
spacing	2mmx2mmx2mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	b62e6846-9e2a-46
studyid	Stability_20260224
seriesnumber	28
runnumber	1
magneticfield	3

description	ses-20ch_func-bold_acq-cmrrstd_task-tsnr
scandate	2026-02-24
scantime	10:46:11
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
receivecoilname	HeadNeck_20_TCS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-20ch_func-bold_acq-cmrrstd_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260224104609.047500
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	2
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	RESEARCH
acquisitionduration	417.892
resonantnucleus	1H

applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
imagecomments	Not for diagnostic use, Unaliased MB3/ PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	69
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00