






Phantom QA results for / fbirn_out/ ses64ch_funcbold_acqsiemensstd_task slice 17




```
ld acqsiemensstd tasktsnr 4.xml [percent fluct (trend removed), drift, driftfit] :
```


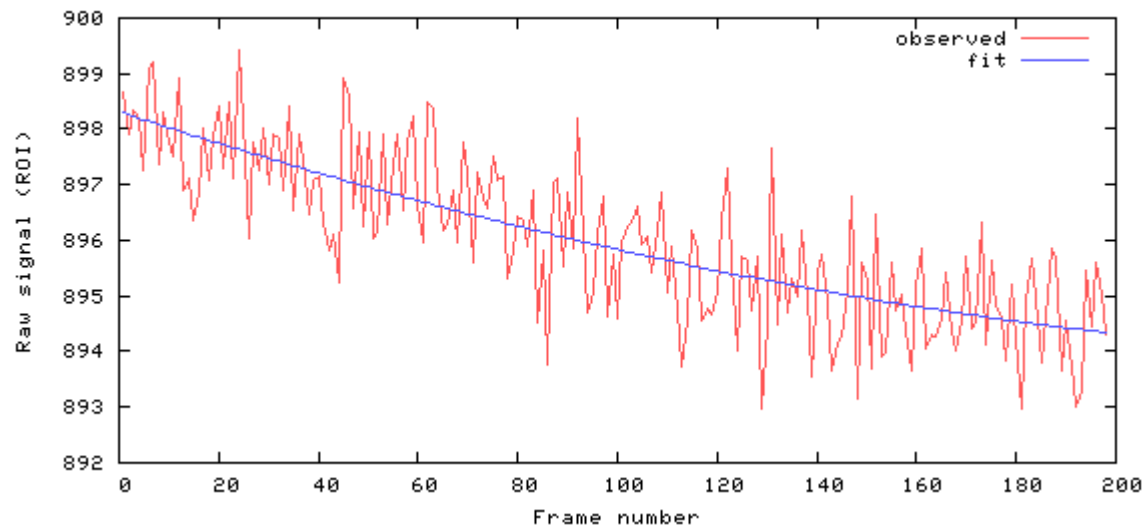


```
[mean, SNR, SFNR] = [896.0 107.2 109.1]
```

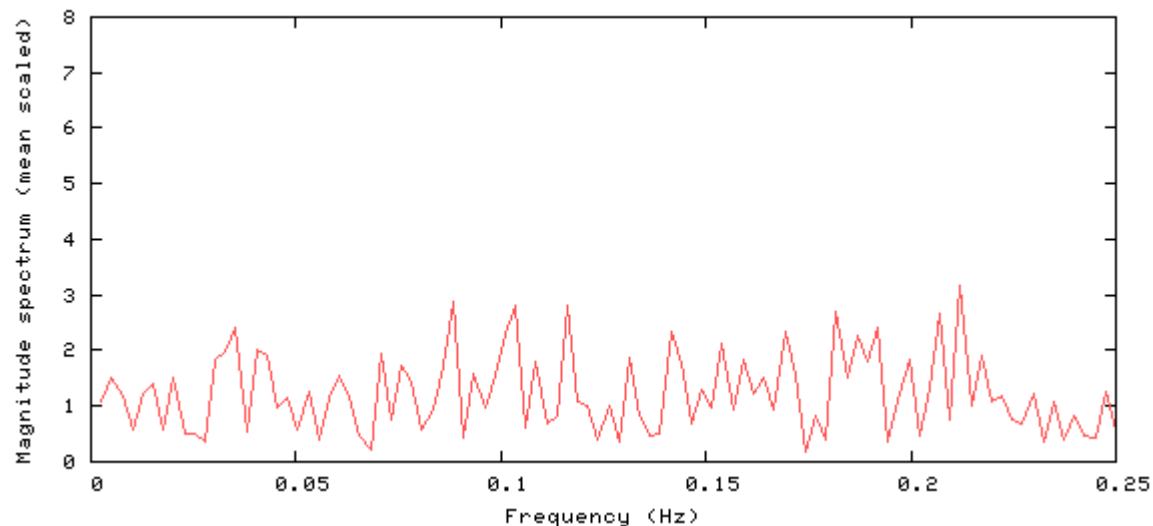

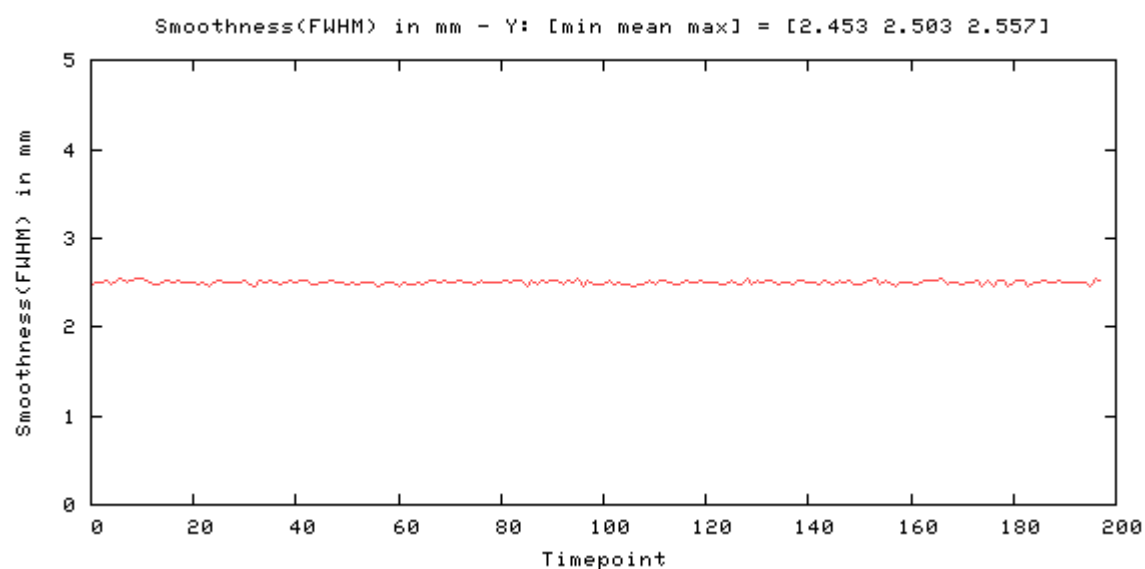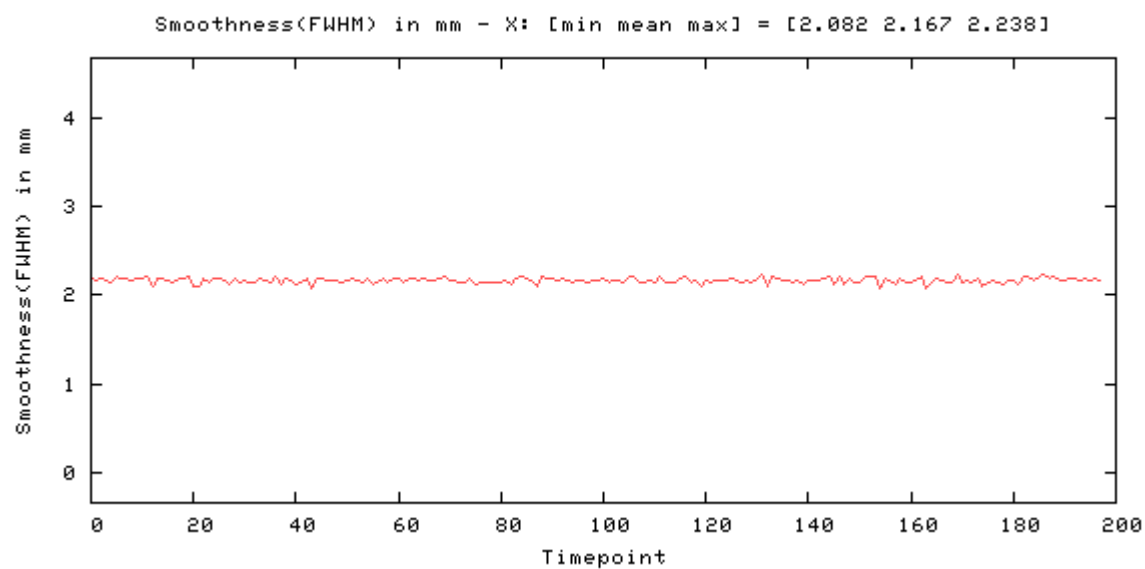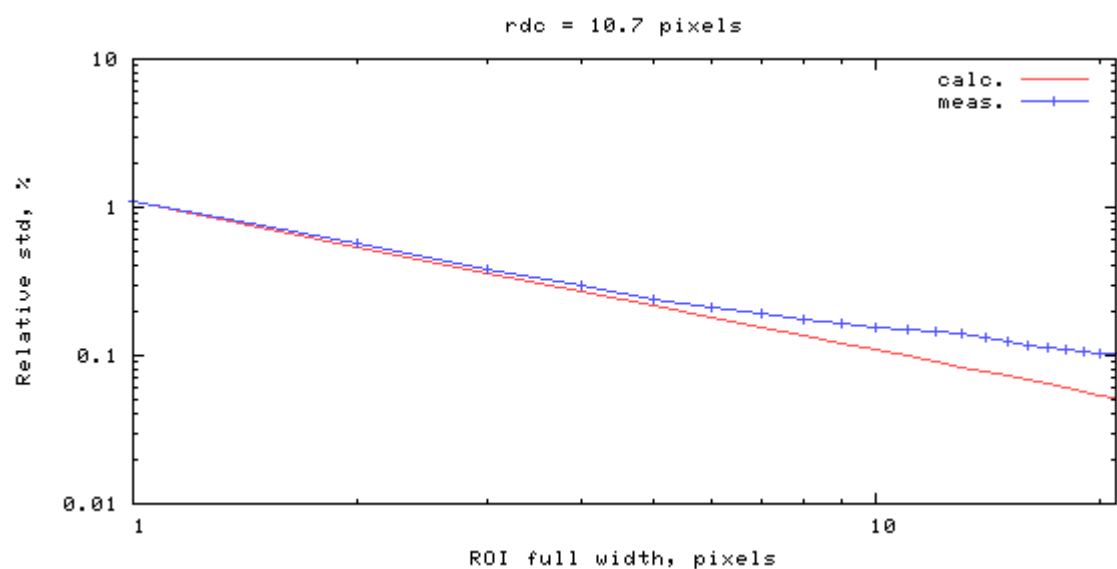

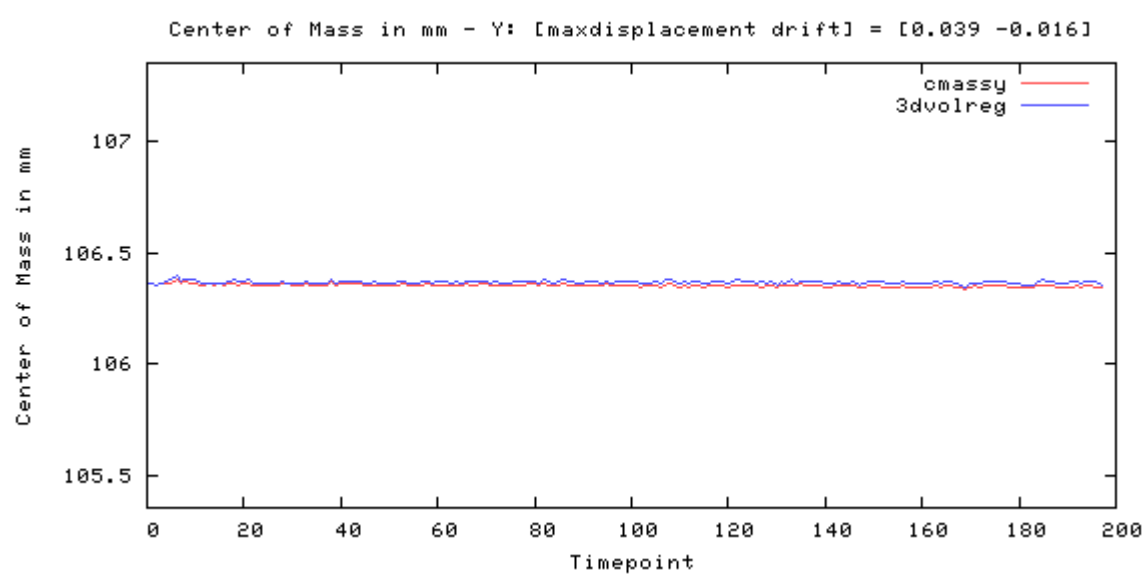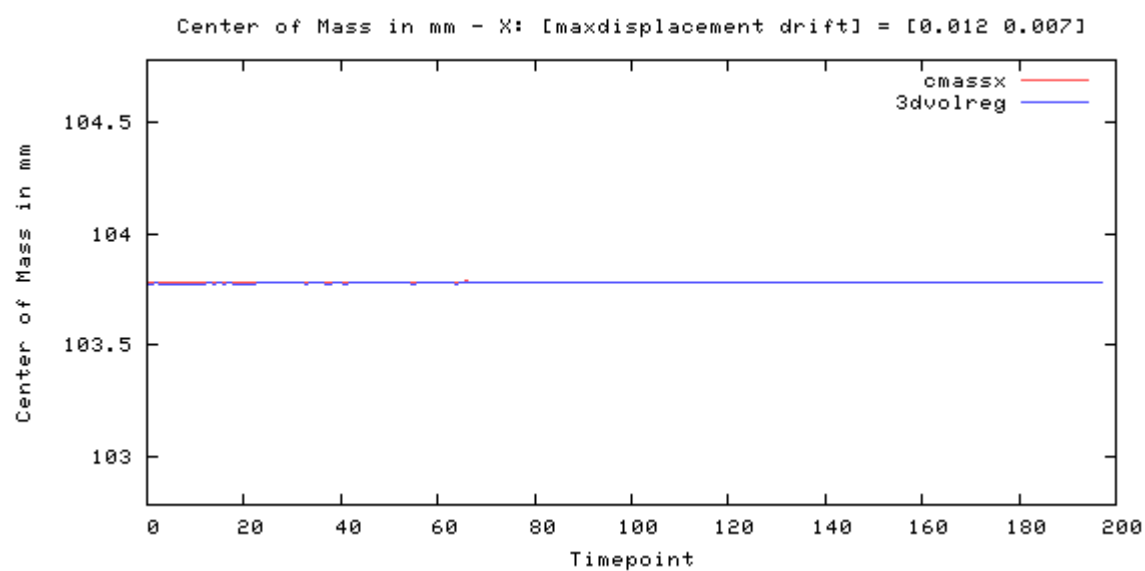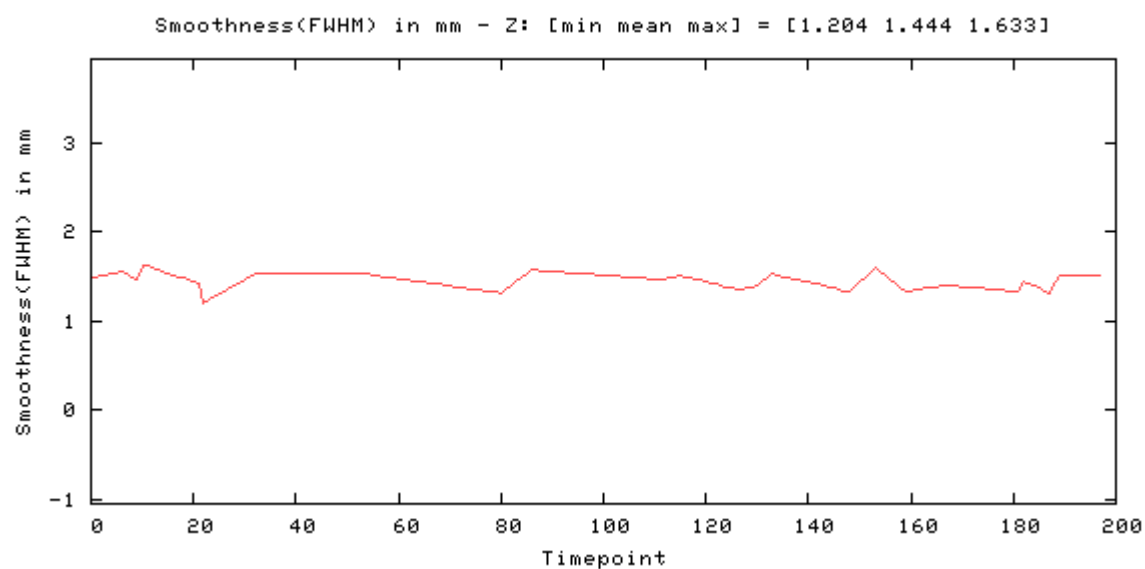


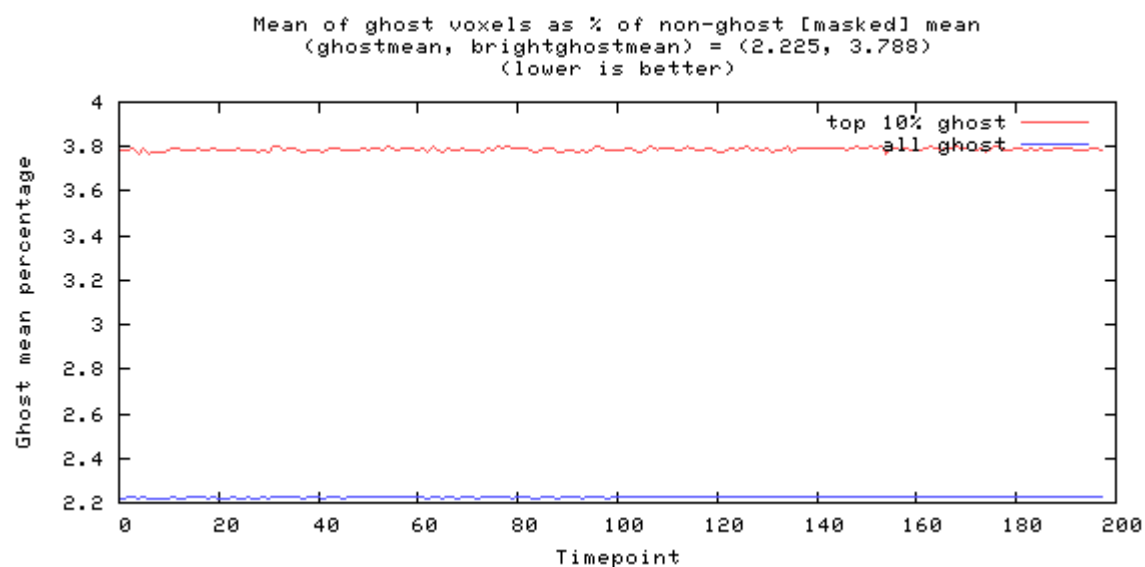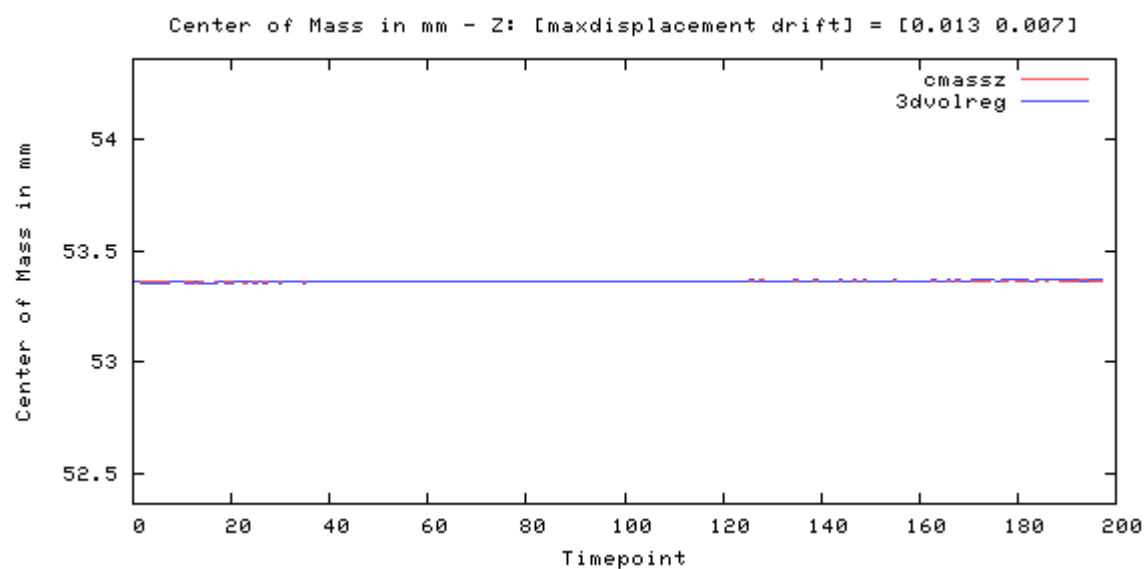

Odd-even difference image



image min: -340, image max: 334

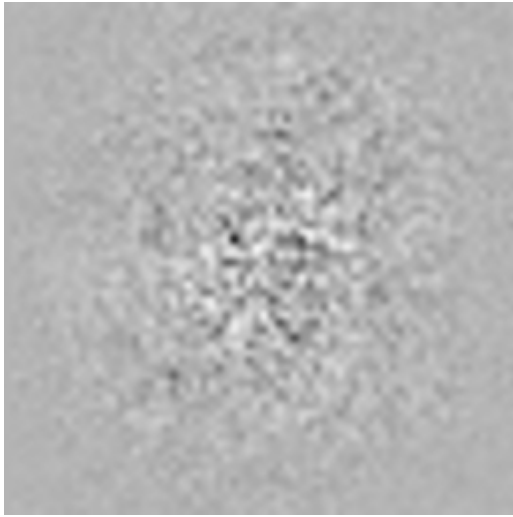

Mean image




image min: 0, image max: 1213.67



Standard Deviation image




image min: 0, image max: 11.8407

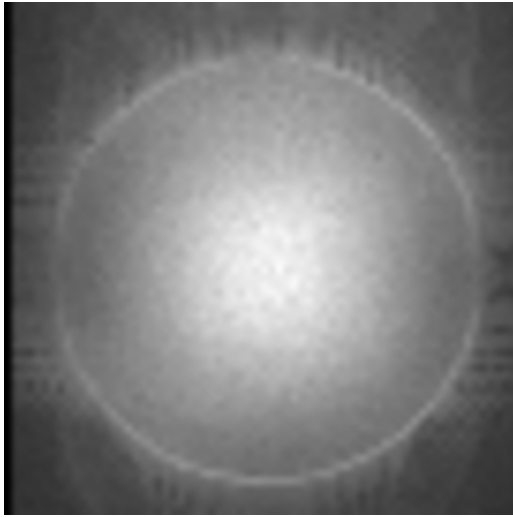

SFNR image



image min: 0, image max: 583.156

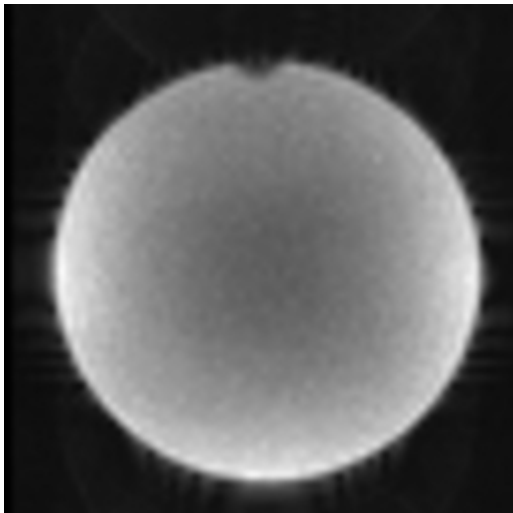

Acquisition parameters

dimensions	68x68x36x200
spacing	3mmx3mmx3mmx1[units?]
gap	0mmx0mmx0mmx0[units?]
scanner	AWP213018
psdname	
examnumber	8fe769f4-63e9-43
studyid	Stability_20250304
seriesnumber	4
runnumber	1-200
scanningsequence	GR/EP

sequencevariant	SK
magneticfield	3
description	ses-64ch_func-bold_acq-siemensstd_task-tsnr
scandate	2025-03-04
scantime	09:58:51
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	3
frequencydirection	1
acquisitionmatrix	68 68
bandwidth	166.668
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
ti	
protocolname	ses-64ch_func-bold_acq-siemensstd_task-tsnr
sliceorder	2,4,6,8,10,12,14,16,18,20,22,24,26,28,30,32,34,36,1,3
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,6
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20250304095850.445000
volumetricproperties	DISTORTED
scanoptions	FS/PER
mracquisitiontype	2D
sequencename	*epfid2d1_68
angioflag	N
numaverages	1
imagingfrequency	123.253209
imagednucleus	1H
spacingbetweenslices	3
numphaseencodingsteps	68
echoTRAINlength	33
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61

triggertime	
beatrejectionflag	N
cardiacnumberofimages	1
inplane phase encoding direction	COL
variableflipangleflag	N
sar	0.08150404198022
contentqualification	PRODUCT
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
samplesperpixel	1
photometricinterpretation	MONOCHROME2
rows	68
columns	68
bitsallocated	16
bitsstored	12
highbit	11
lossyimagecompression	00
echonumber	0
pixelbandwidth	2451