






# Phantom QA results for / fbirn\_out/ ses64ch\_funcbold\_acqsiemensbasic\_ta slice 14




1d acqsiemensbasic tasktsnr 8.xml [percent fluct (trend removed), drift, driftfit:


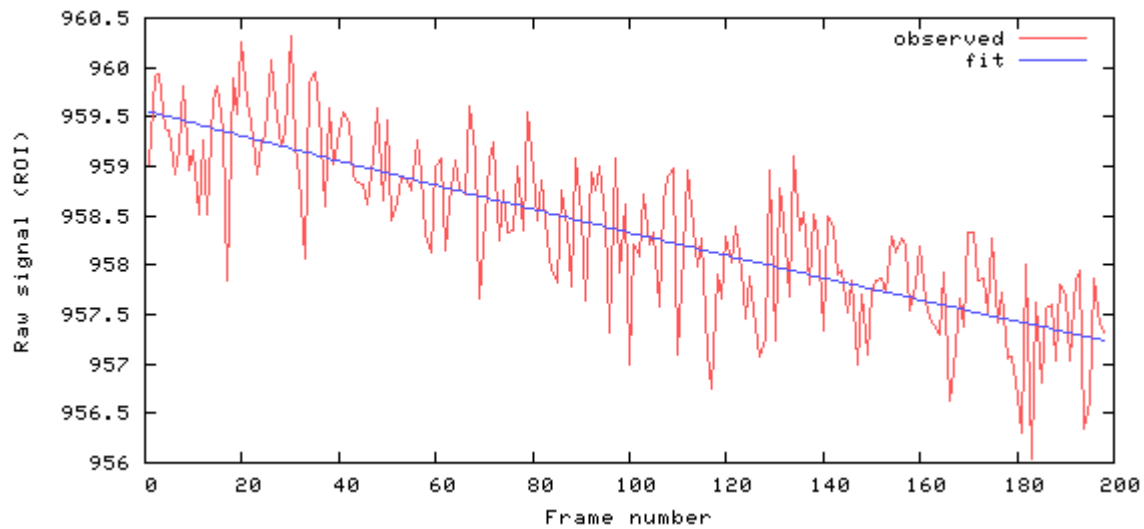


[mean, SNR, SFNR] = [958.4 164.4 169.4]

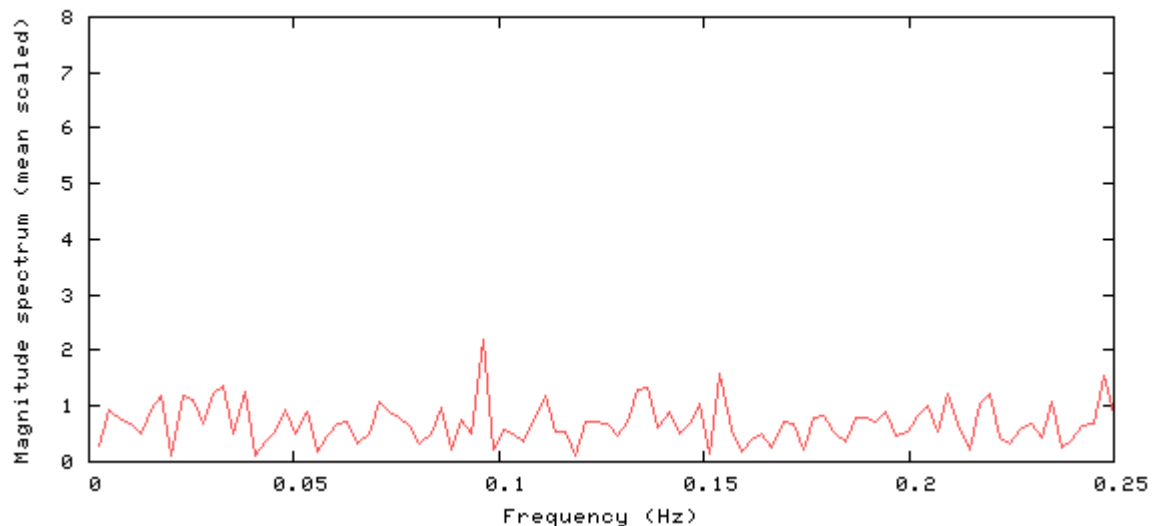

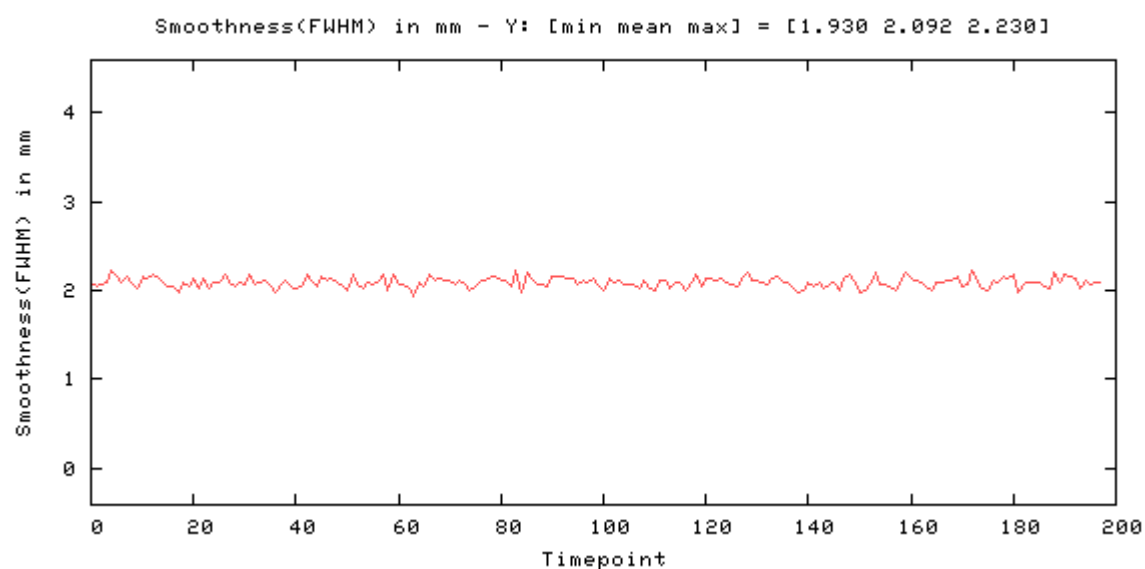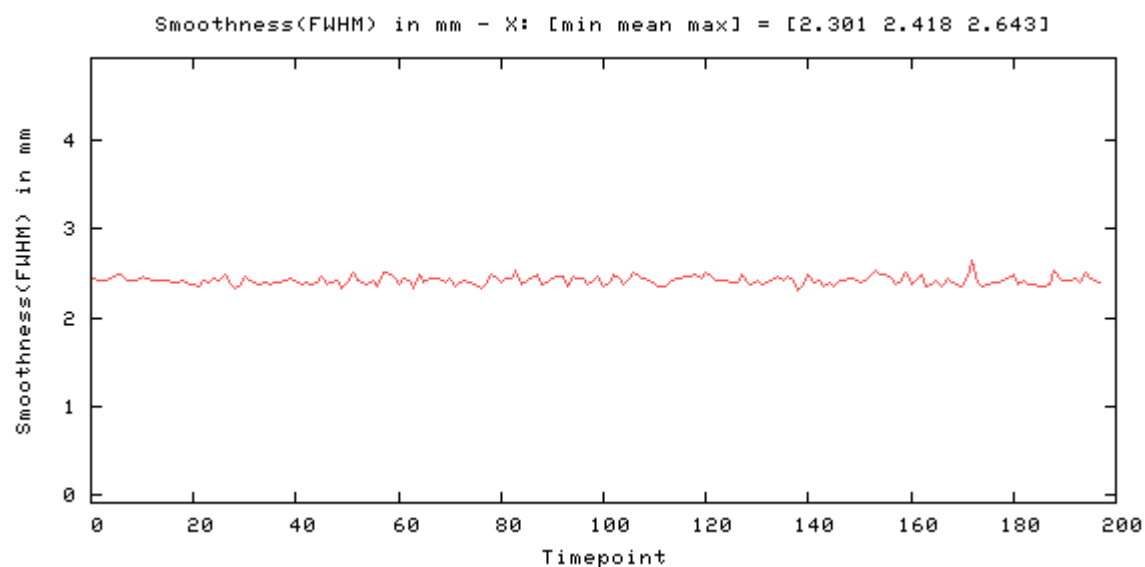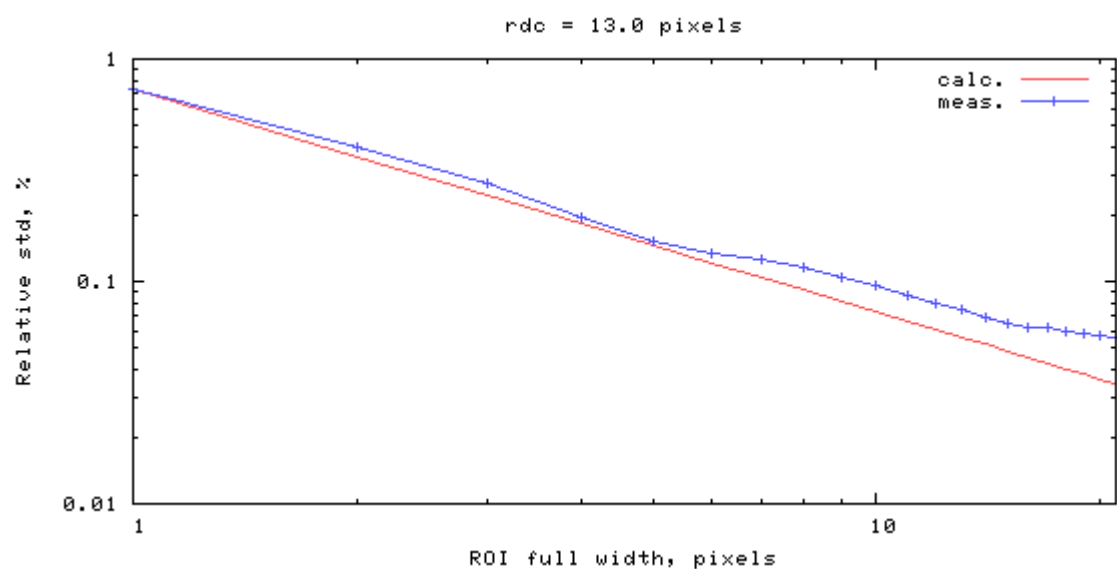

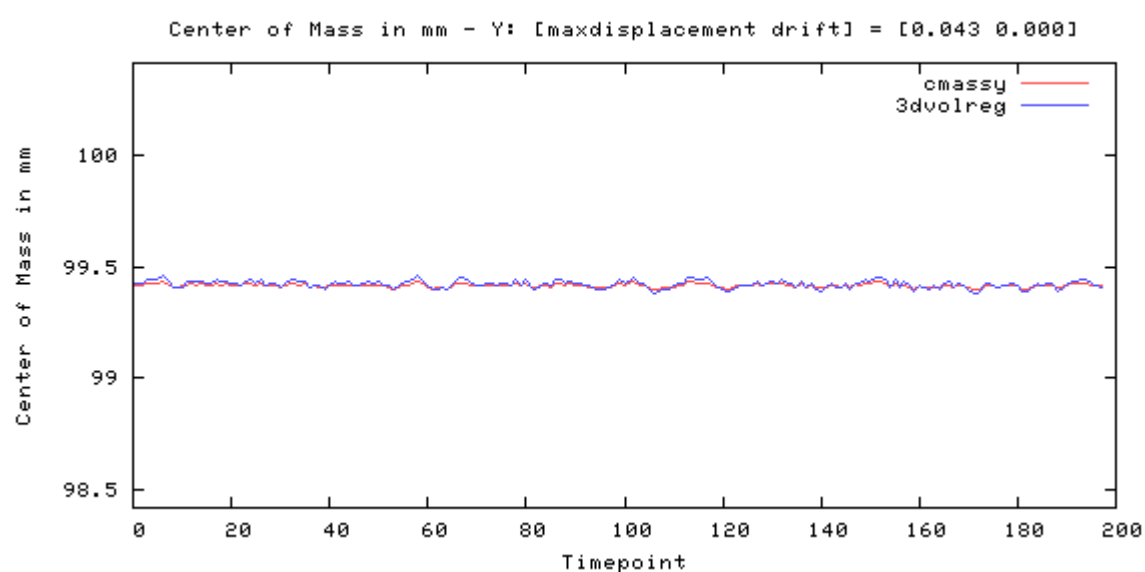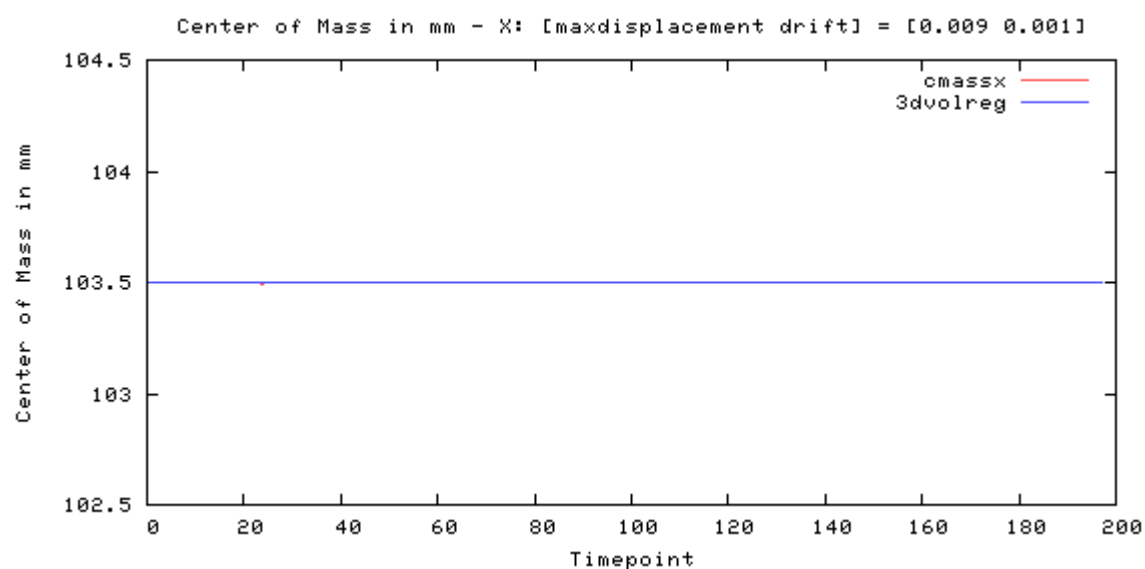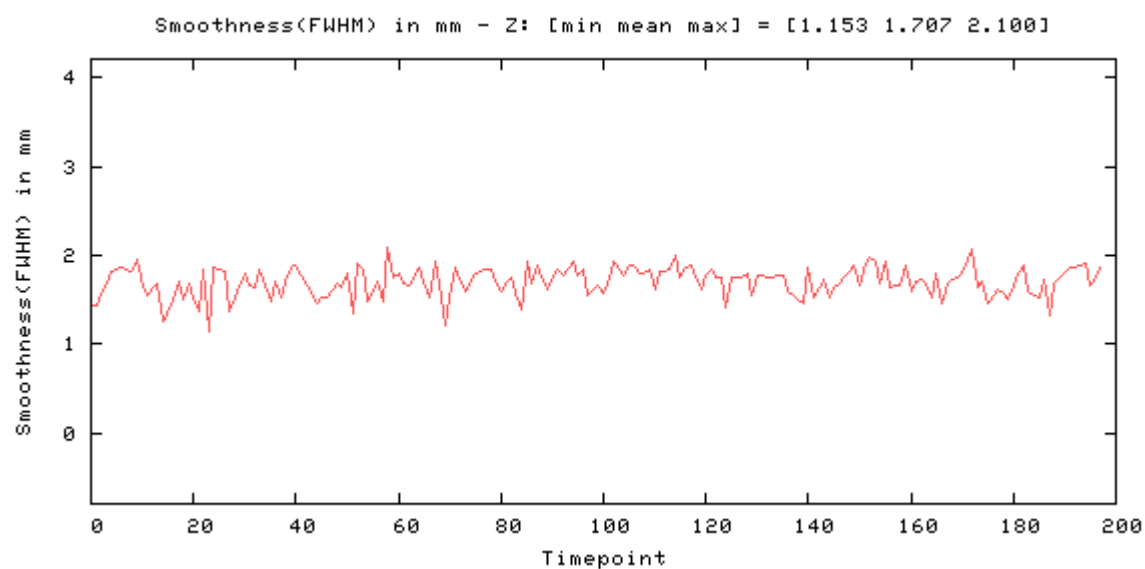

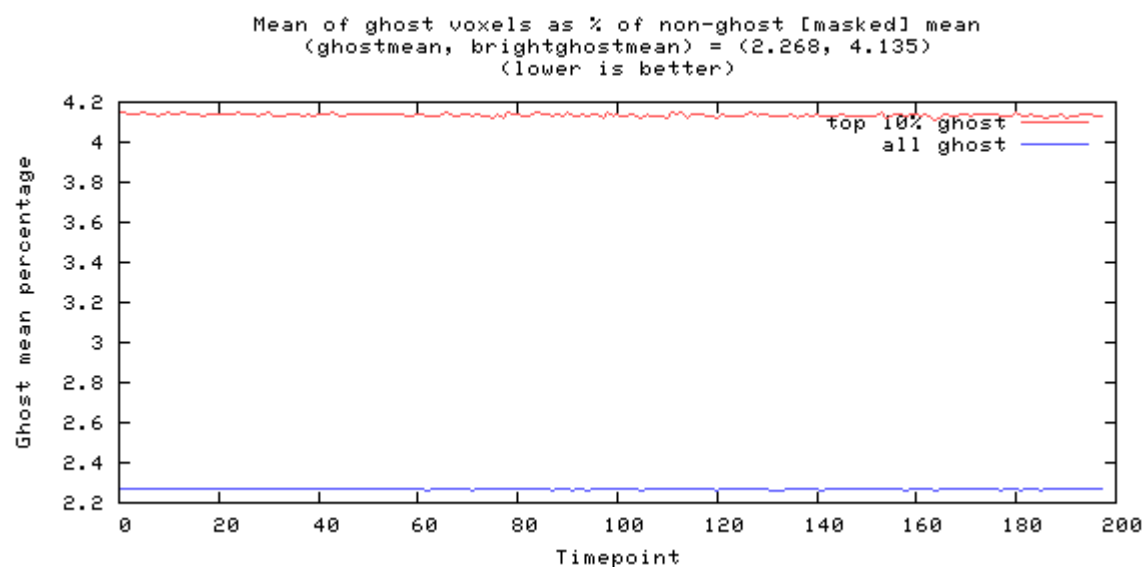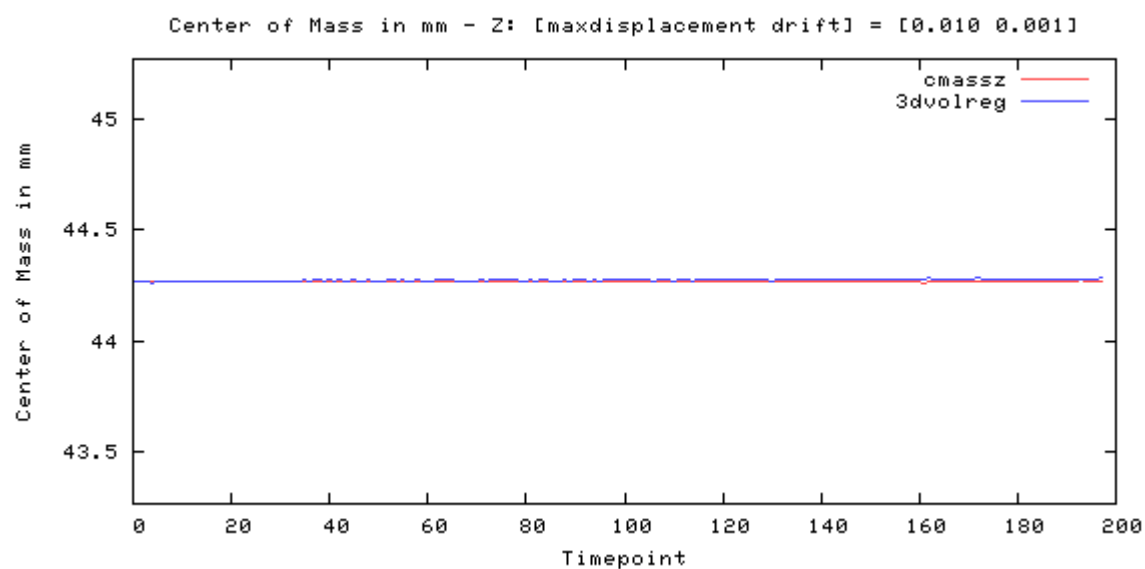


---

## Odd-even difference image



image min: -277, image max: 268

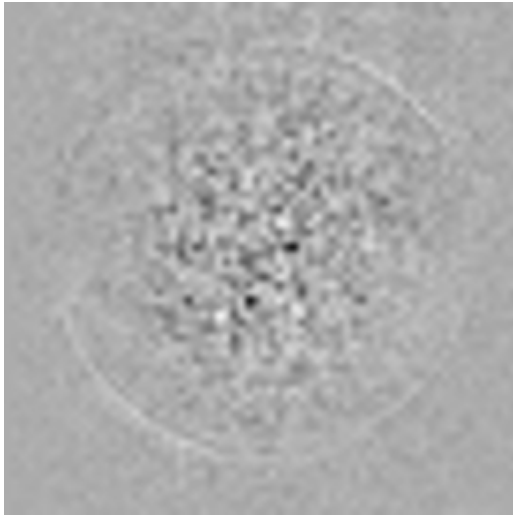

---

### Mean image



image min: 0, image max: 1349.77




---

### Standard Deviation image



image min: 0, image max: 9.70123

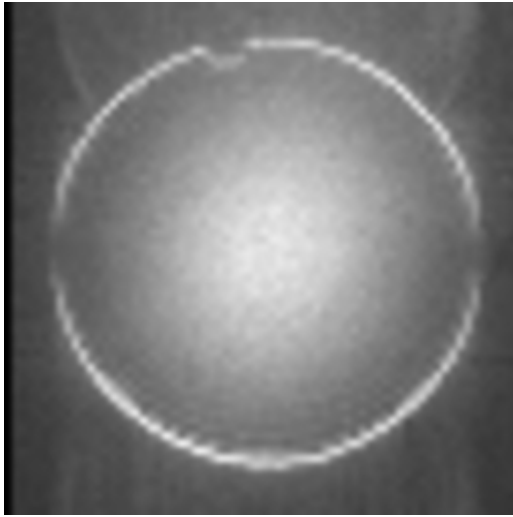


---

### SFNR image



image min: 0, image max: 898.462

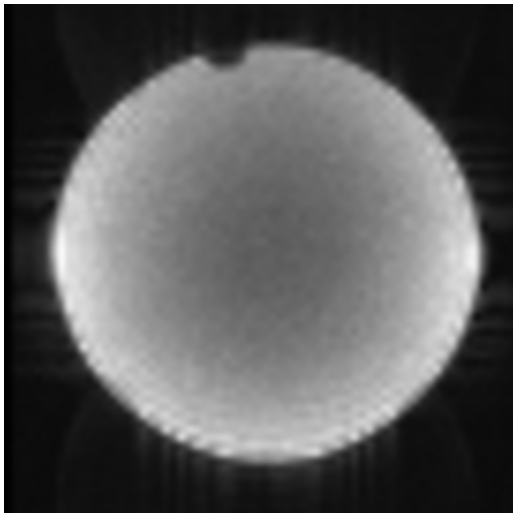

---

### Acquisition parameters

|                  |                       |
|------------------|-----------------------|
| dimensions       | 68x68x30x200          |
| spacing          | 3mmx3mmx3mmx1[units?] |
| gap              | 0mmx0mmx0mmx0[units?] |
| scanner          | AWP213018             |
| psdname          |                       |
| examnumber       | 4061a312-2d94-4d      |
| studyid          | Stability_20250818    |
| seriesnumber     | 8                     |
| runnumber        | 1-200                 |
| scanningsequence | GR/EP                 |

|                             |                                                        |
|-----------------------------|--------------------------------------------------------|
| sequencevariant             | SK                                                     |
| magneticfield               | 3                                                      |
| description                 | ses-64ch_func-bold_acq-siemensbasic_task-tsnr          |
| scandate                    | 2025-08-18                                             |
| scantime                    | 09:50:41                                               |
| tr                          | 2000                                                   |
| te                          | 30                                                     |
| operator                    |                                                        |
| flipangle                   | 79                                                     |
| prescribedslicespacing      | 3                                                      |
| frequencydirection          | 1                                                      |
| acquisitionmatrix           | 68 68                                                  |
| bandwidth                   | 166.668                                                |
| receivecoilname             | HeadNeck_64_CS                                         |
| transmitcoilname            | Body                                                   |
| institution                 | THE UNIVERSITY OF PENNSYLVANIA                         |
| examdescription             | ACTIVE KIRWAN                                          |
| scanner                     | AWP213018                                              |
| scannermanufacturer         | Siemens Healthineers                                   |
| scannermodelname            | MAGNETOM Cima.X                                        |
| ti                          |                                                        |
| protocolname                | ses-64ch_func-bold_acq-siemensbasic_task-tsnr          |
| sliceorder                  | 2,4,6,8,10,12,14,16,18,20,22,24,26,28,30,1,3,5,7,9,11, |
| institutionaddress          | Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,6  |
| institutionaldepartmentname | MindCORE Neuroimaging Suite 140                        |
| acquisitiondatetime         | 20250818095040.030000                                  |
| volumetricproperties        | DISTORTED                                              |
| scanoptions                 | FS/PER                                                 |
| mracquisitiontype           | 2D                                                     |
| sequencename                | *epfid2d1_68                                           |
| angioflag                   | N                                                      |
| numaverages                 | 1                                                      |
| imagingfrequency            | 123.253297                                             |
| imagednucleus               | 1H                                                     |
| spacingbetweenslices        | 3                                                      |
| numphaseencodingsteps       | 68                                                     |
| echotrainlength             | 68                                                     |
| percentsampling             | 100                                                    |
| percentphasefieldofview     | 100                                                    |
| scannerserialnumber         | 213018                                                 |
| softwareversions            | syngo MR XA61                                          |

|                                  |                  |
|----------------------------------|------------------|
| triggertime                      |                  |
| beatrejectionflag                | N                |
| cardiacnumberofimages            | 1                |
| inplane phase encoding direction | COL              |
| variableflipangleflag            | N                |
| sar                              | 0.09313652477200 |
| contentqualification             | PRODUCT          |
| acquisitionnumber                | 1                |
| instancenumber                   | 1                |
| numtemporalpositions             | 200              |
| samplesperpixel                  | 1                |
| photometricinterpretation        | MONOCHROME2      |
| rows                             | 68               |
| columns                          | 68               |
| bitsallocated                    | 16               |
| bitsstored                       | 12               |
| highbit                          | 11               |
| lossyimagecompression            | 00               |
| echonumber                       | 0                |
| pixelbandwidth                   | 2451             |