






Phantom QA results for / fbirn_out/ ses20ch_funcbold_acqsiemensstd_task slice 17




1d acqsiemensstd tasktsnr 24.xml [percent fluct (trend removed), drift, driftfit]


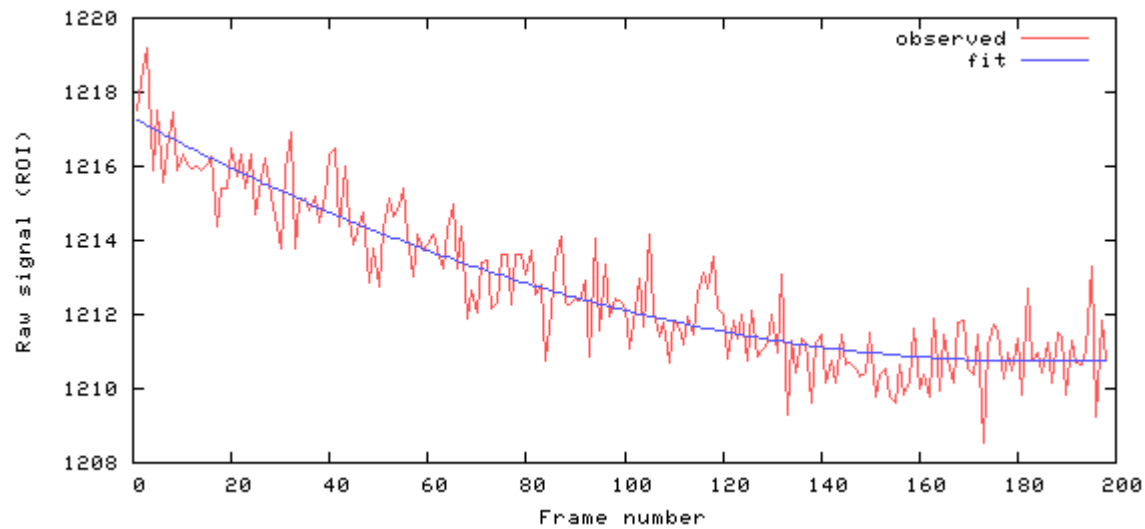


[mean, SNR, SFNR] = [1212.8 138.3 154.9]

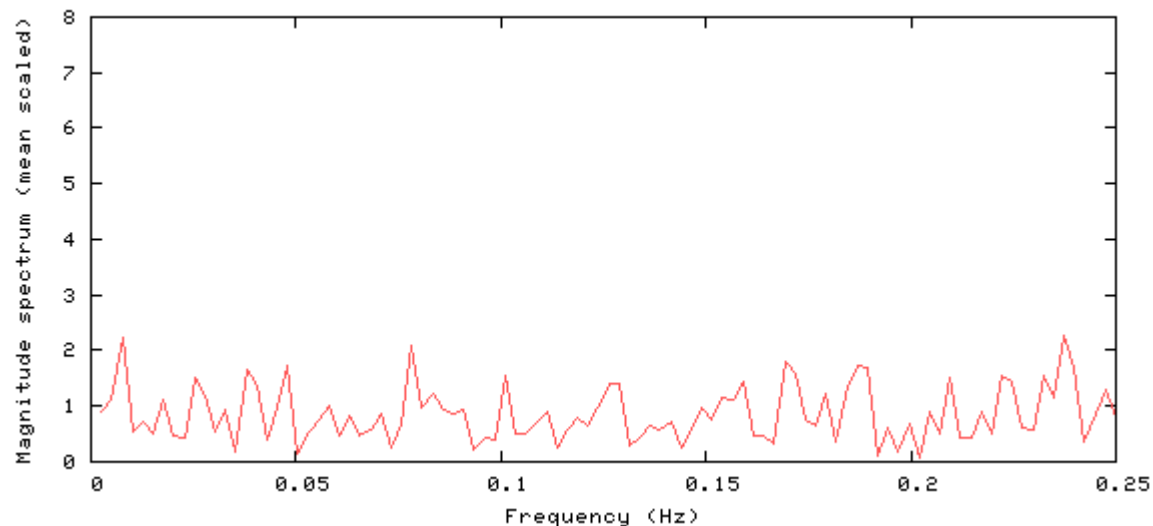

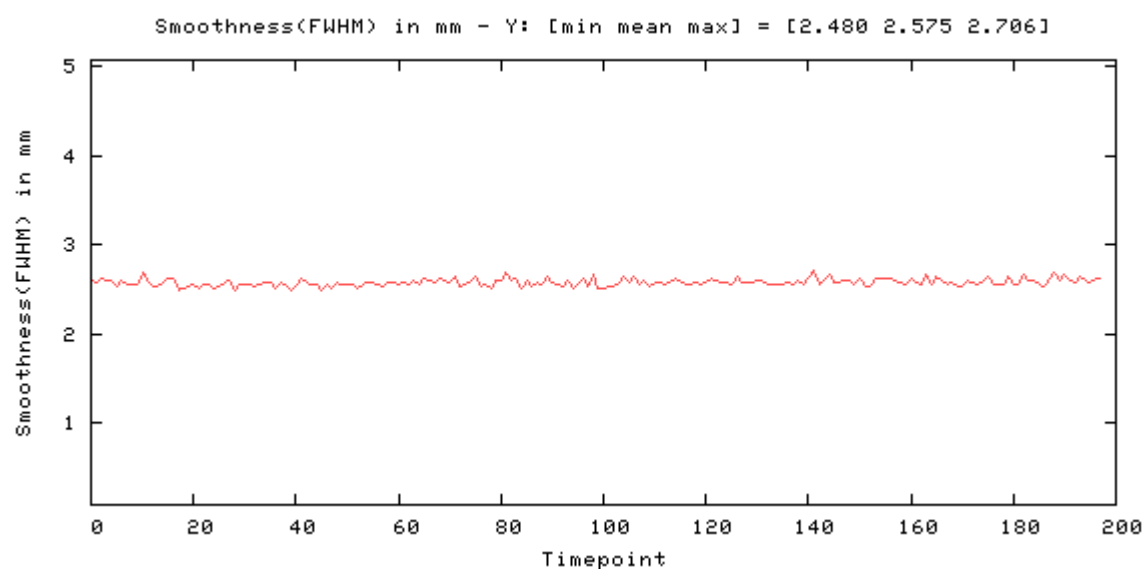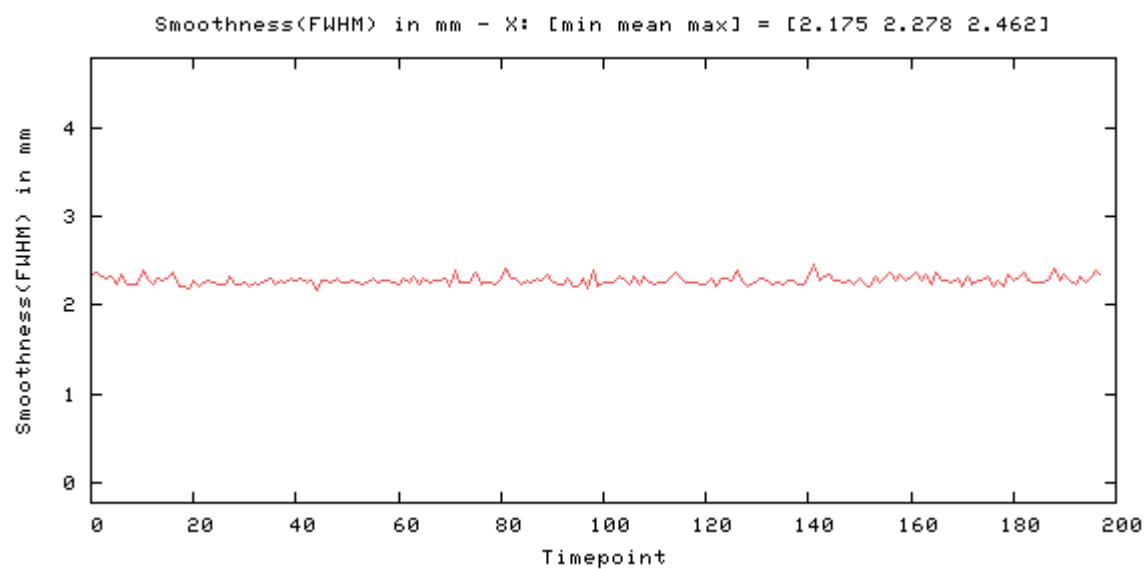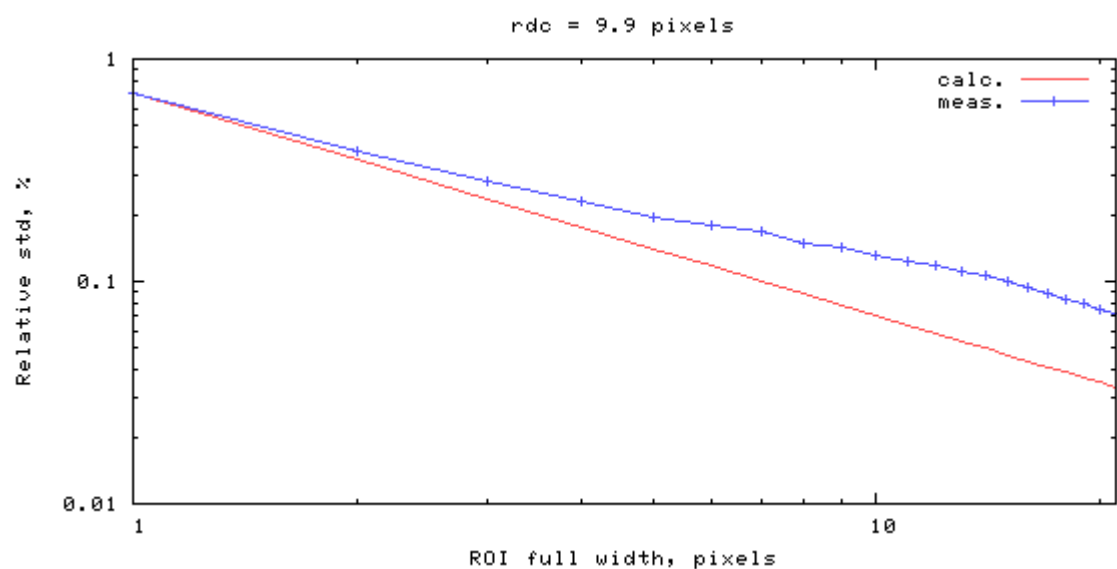

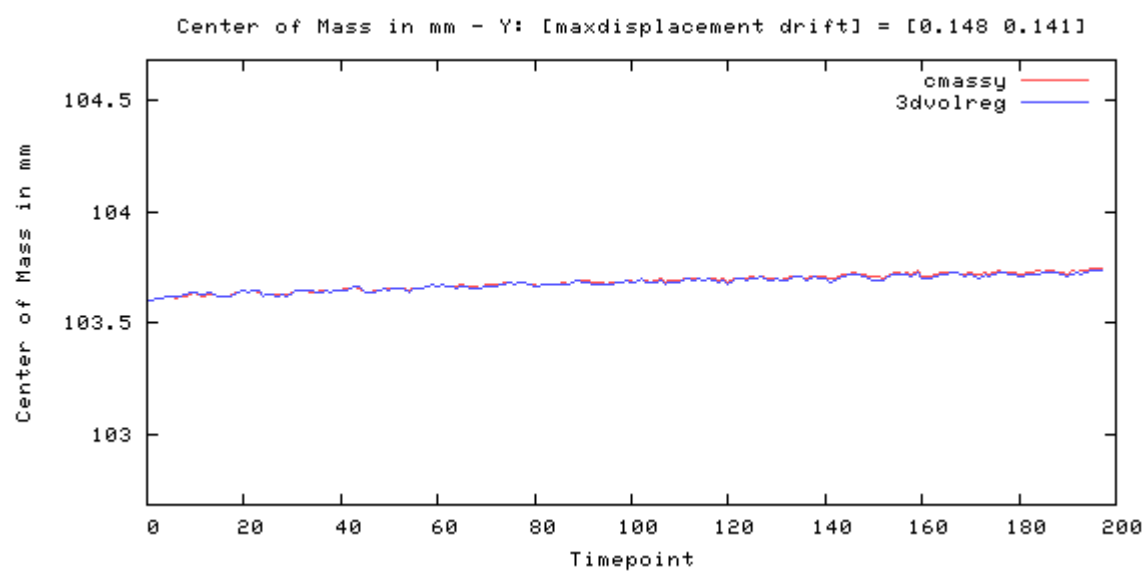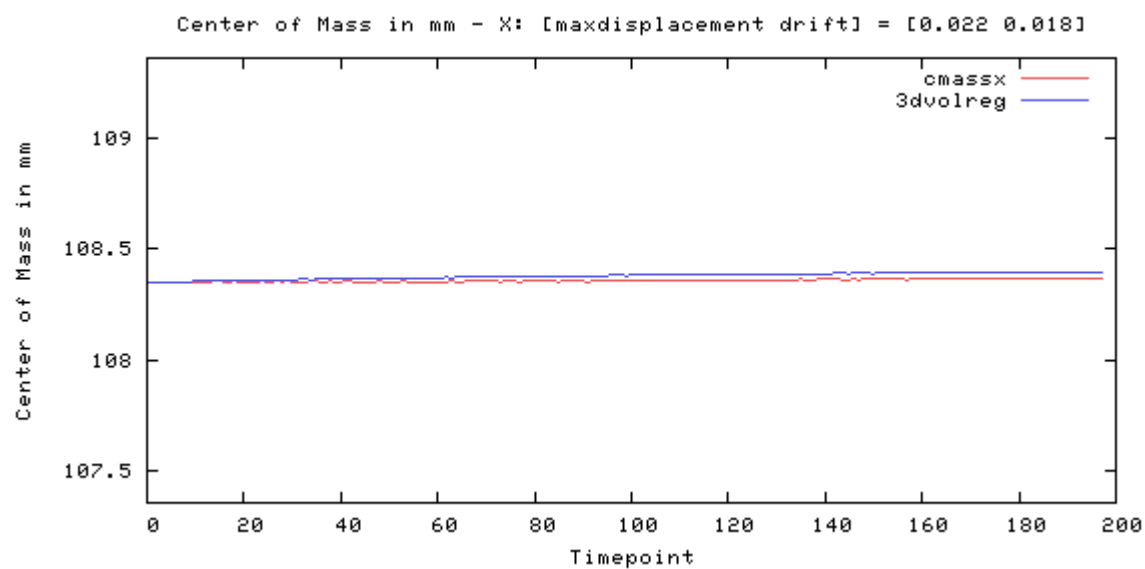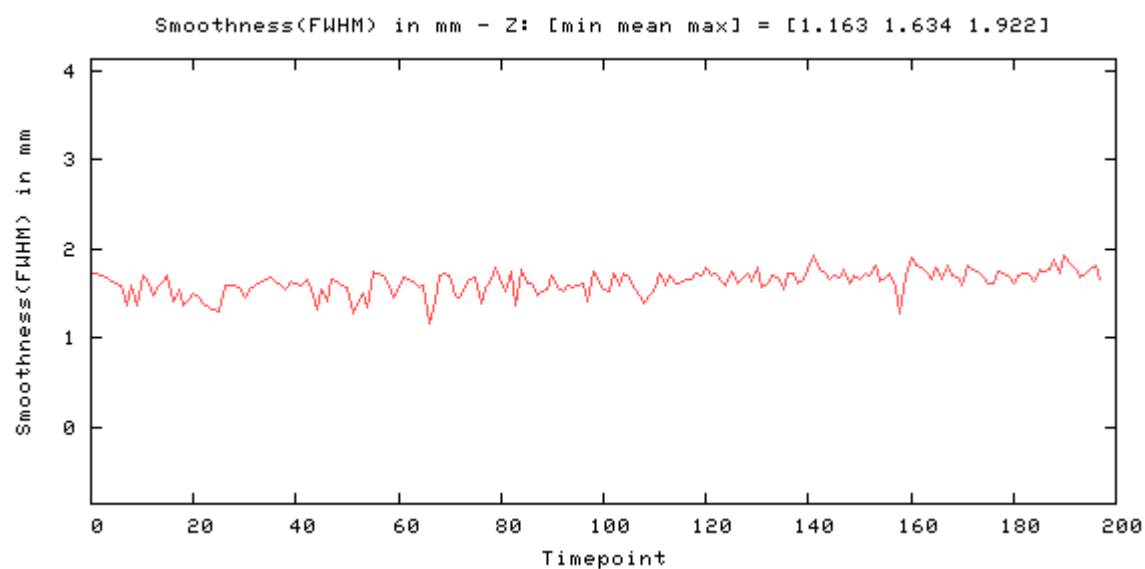


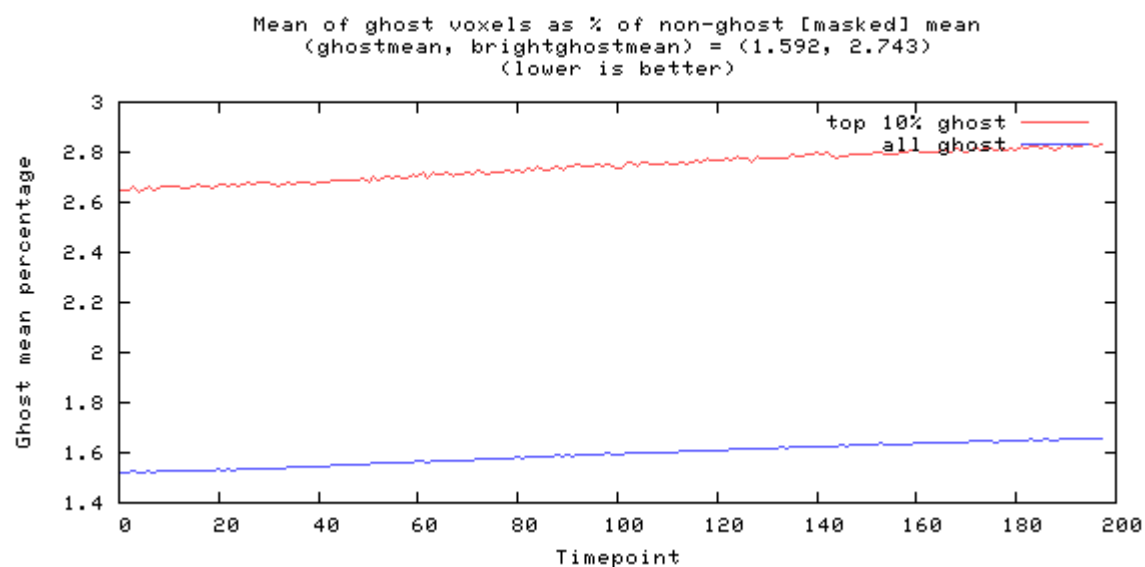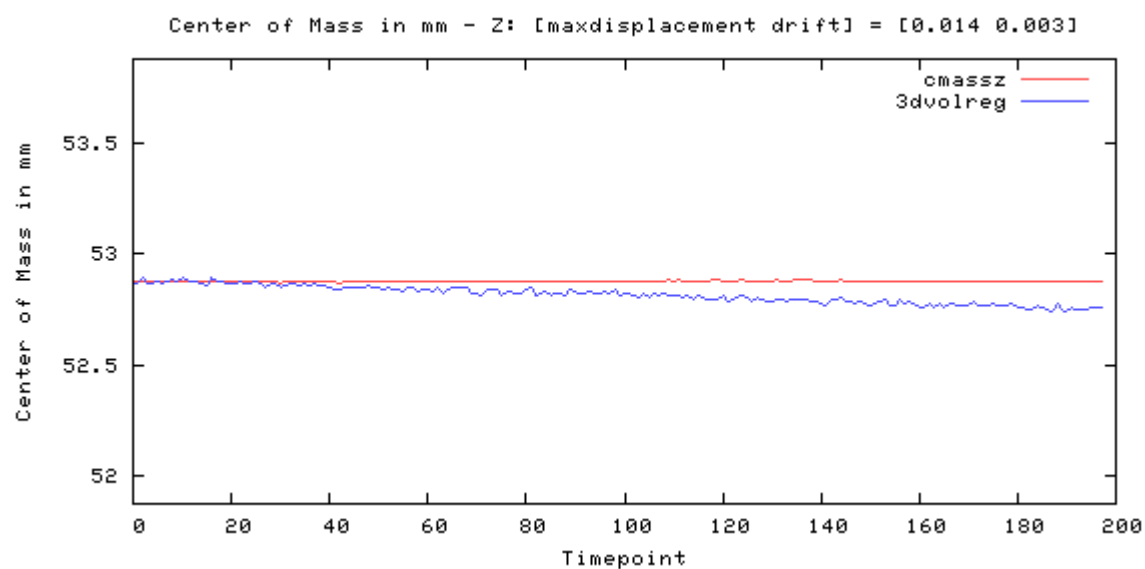

Odd-even difference image



image min: -432, image max: 351

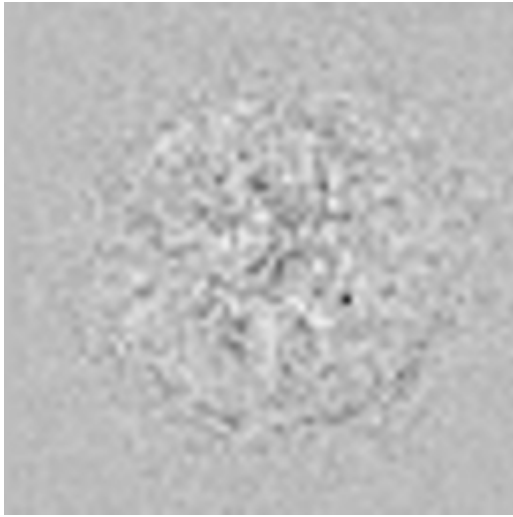

Mean image




image min: 0, image max: 1482.81



Standard Deviation image




image min: 0, image max: 10.3527

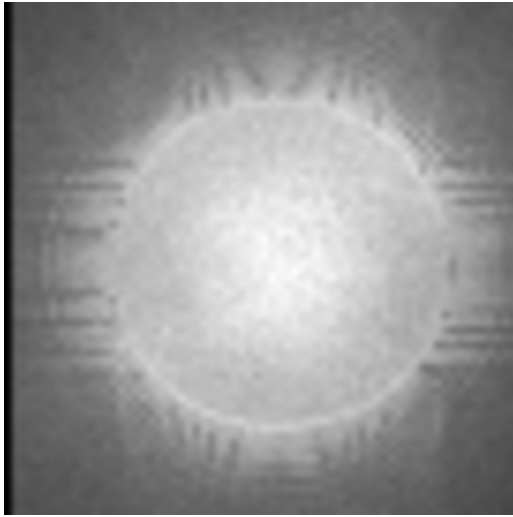

SFNR image



image min: 0, image max: 310.432

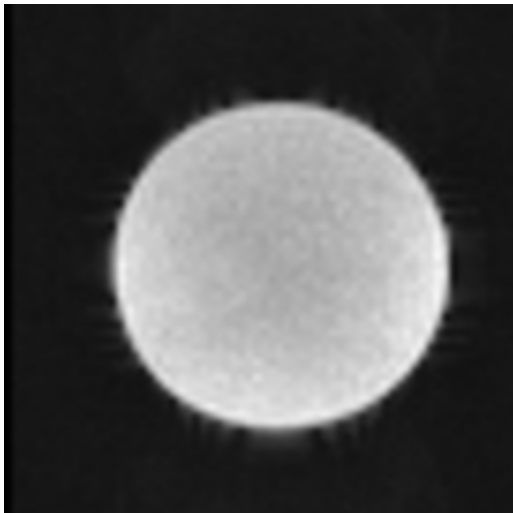

Acquisition parameters

dimensions	68x68x36x200
spacing	3mmx3mmx3mmx1[units?]
gap	0mmx0mmx0mmx0[units?]
scanner	AWP213018
psdname	
examnumber	0dd057a4-d800-49
studyid	Stability_20250902
seriesnumber	24
runnumber	1-200
scanningsequence	GR/EP

sequencevariant	SK
magneticfield	3
description	ses-20ch_func-bold_acq-siemensstd_task-tsnr
scandate	2025-09-02
scantime	10:31:41
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	3
frequencydirection	1
acquisitionmatrix	68 68
bandwidth	166.668
receivecoilname	HeadNeck_20_TCS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
ti	
protocolname	ses-20ch_func-bold_acq-siemensstd_task-tsnr
sliceorder	2,4,6,8,10,12,14,16,18,20,22,24,26,28,30,32,34,36,1,3
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,6
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20250902103140.987500
volumetricproperties	DISTORTED
scanoptions	FS/PER
mracquisitiontype	2D
sequencename	*epfid2d1_68
angioflag	N
numaverages	1
imagingfrequency	123.252975
imagednucleus	1H
spacingbetweenslices	3
numphaseencodingsteps	68
echoTRAINlength	33
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61

triggertime	
beatrejectionflag	N
cardiacnumberofimages	1
inplane phase encoding direction	COL
variableflipangleflag	N
sar	0.10939145174042
contentqualification	PRODUCT
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
samplesperpixel	1
photometricinterpretation	MONOCHROME2
rows	68
columns	68
bitsallocated	16
bitsstored	12
highbit	11
lossyimagecompression	00
echonumber	0
pixelbandwidth	2451