






Phantom QA results for / fbirn_out/ ses64ch_funcbold_acqcmrrstd_taskwar slice 34




old acqcmrrstd taskwarmup 2.xml [percent fluct (trend removed), drift, driftfit]


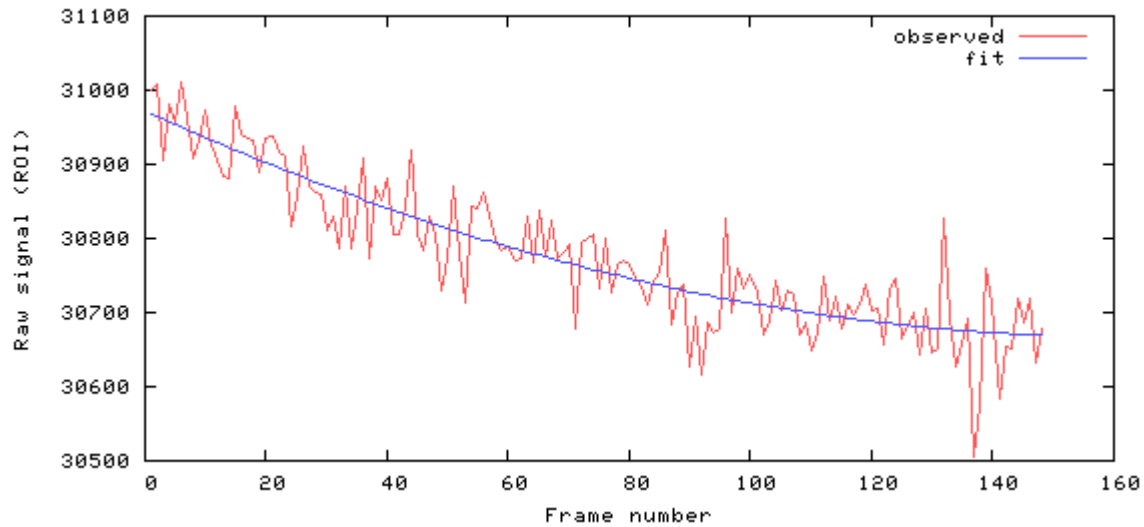


[mean, SNR, SFNR] = [30778.0 113.7 119.9]

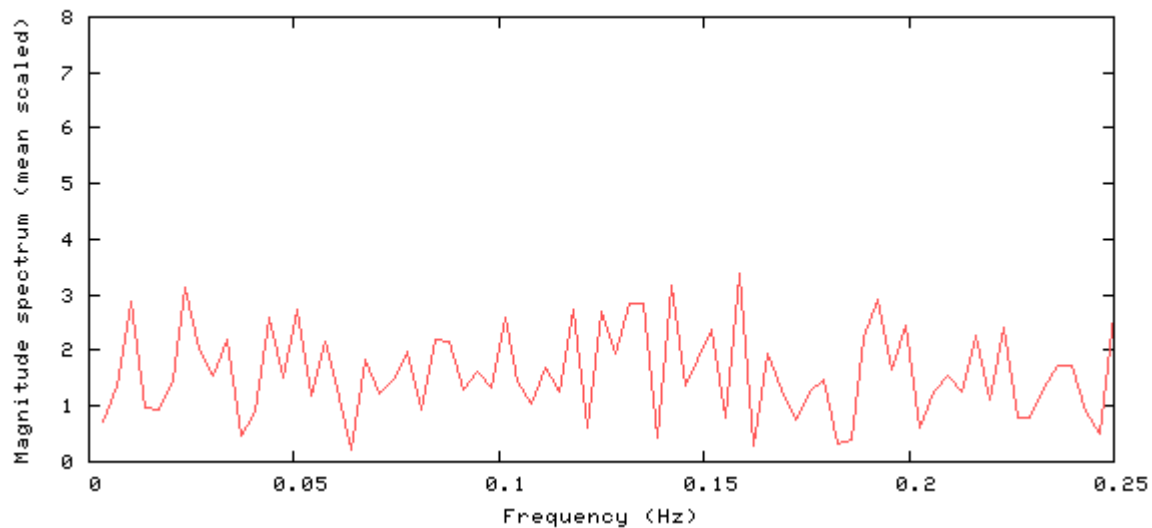

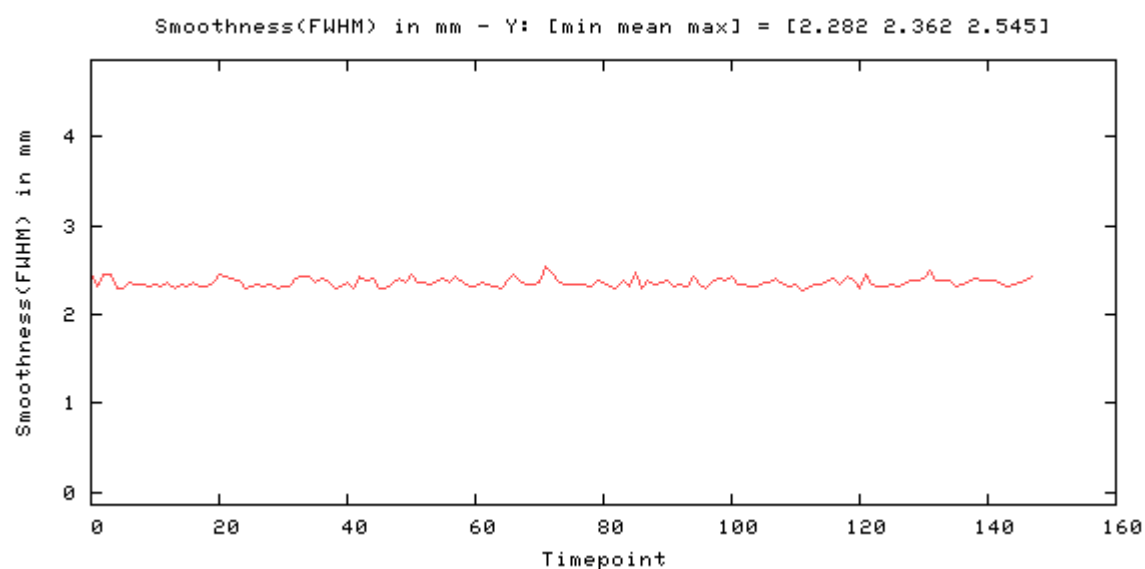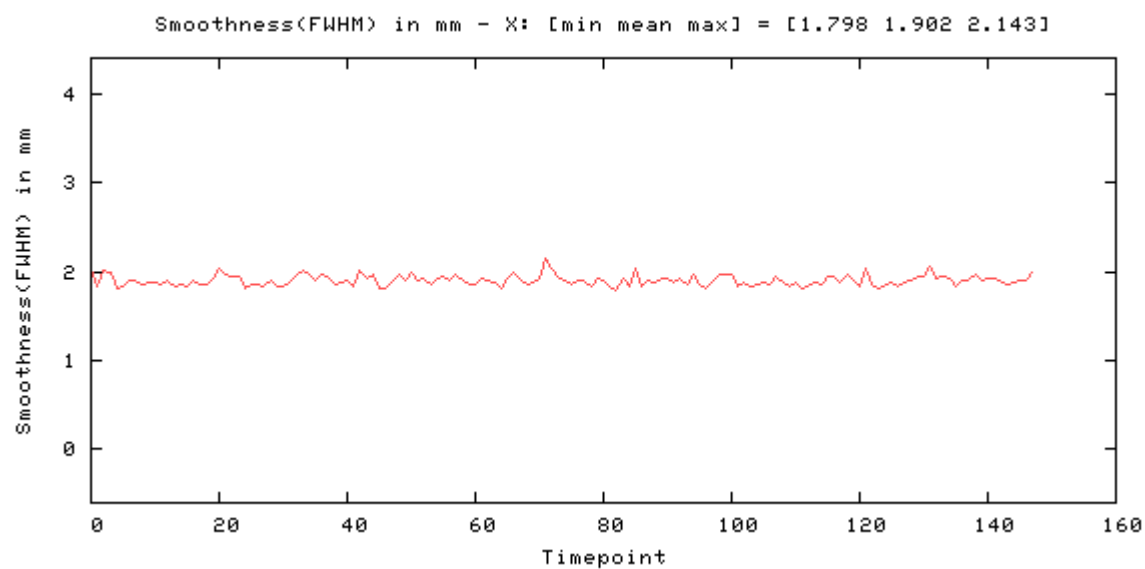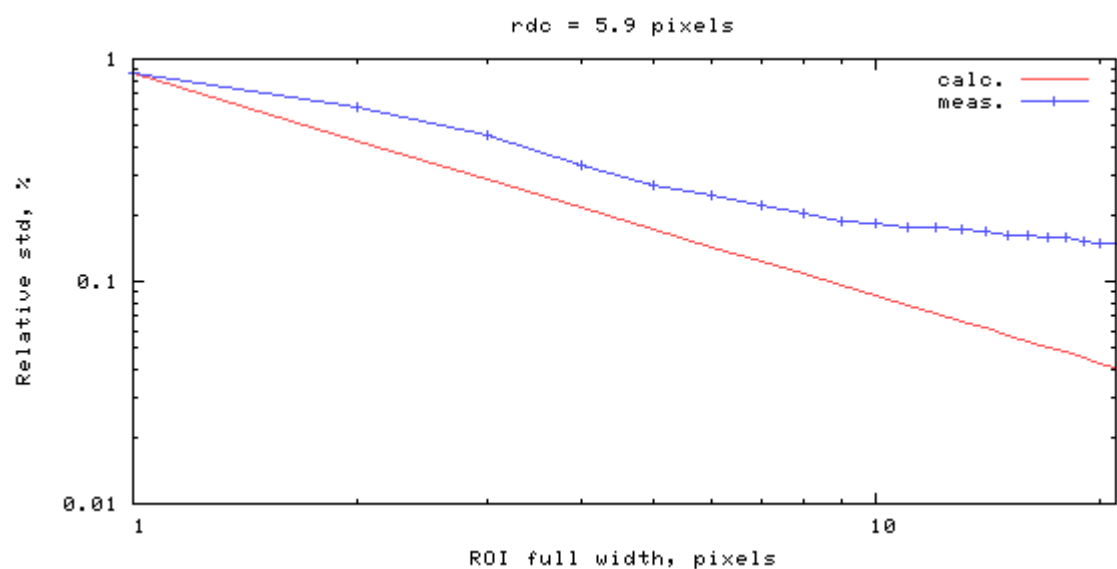

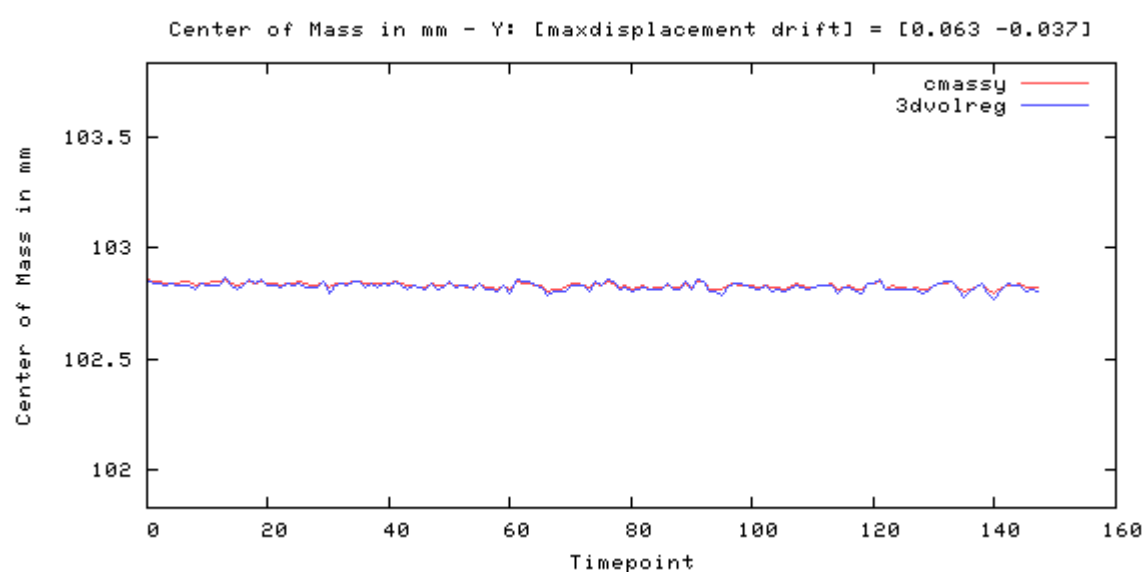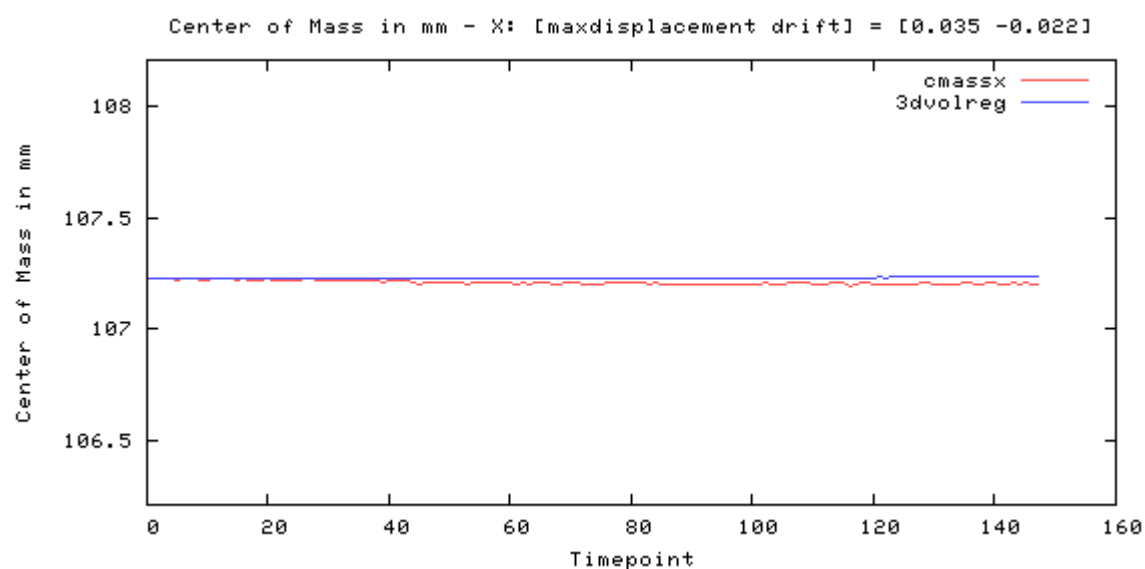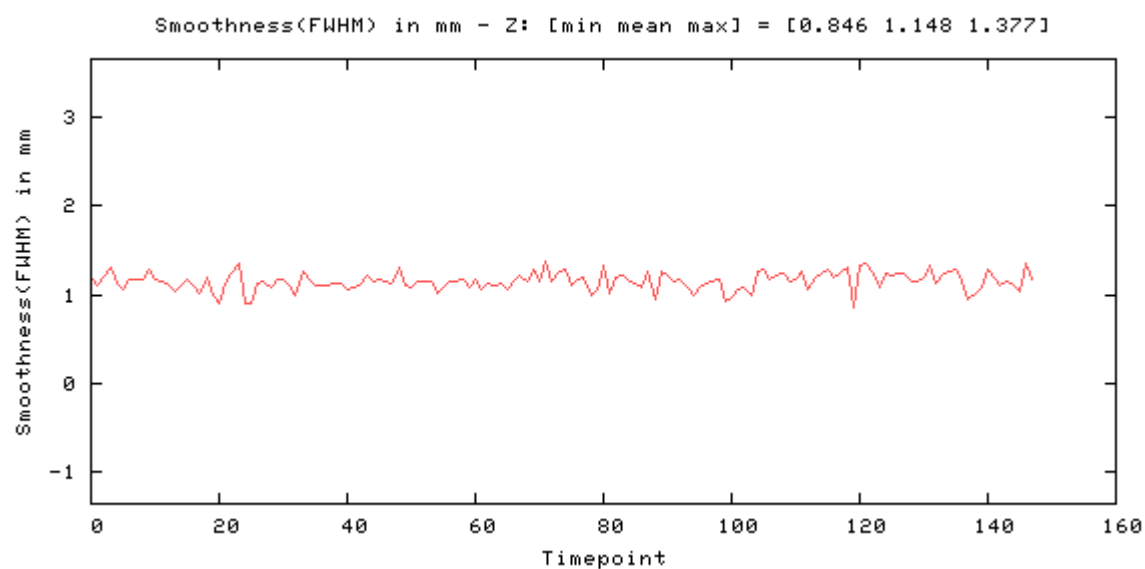


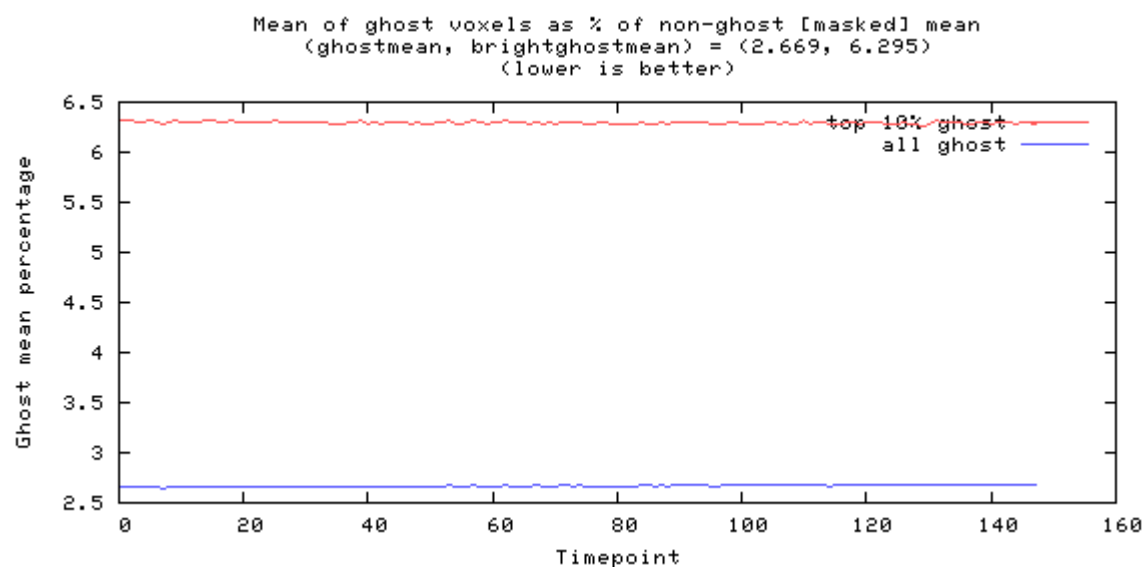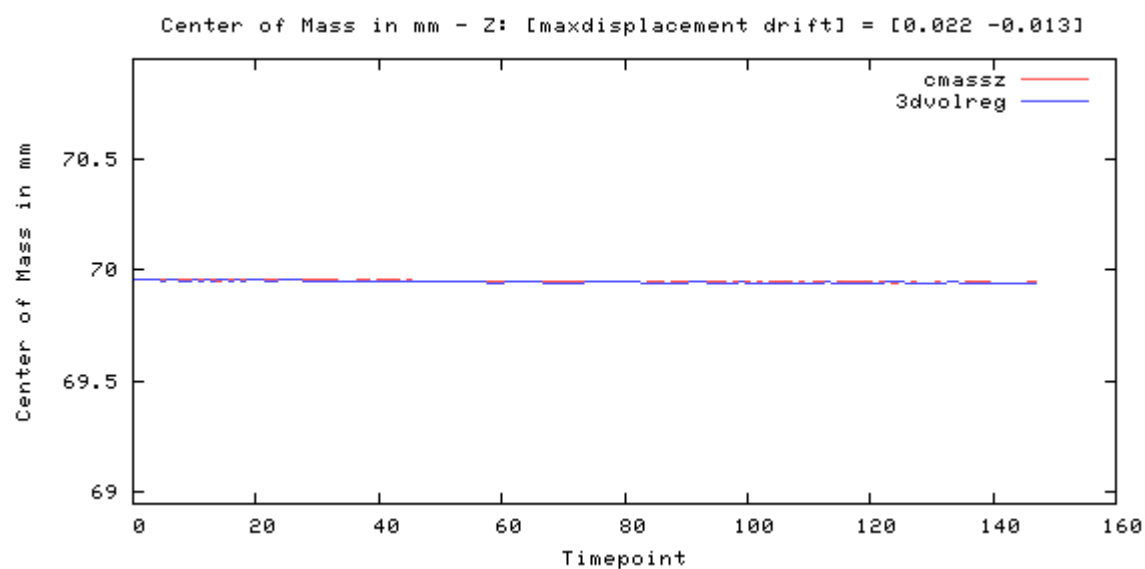

Odd-even difference image




image min: -9891, image max: 8967

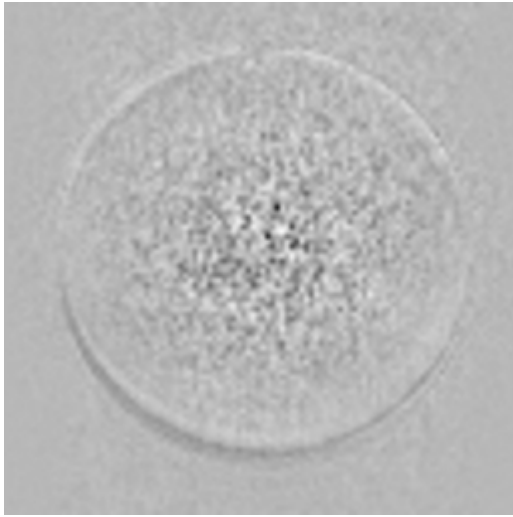

Mean image




image min: 0, image max: 43428.9

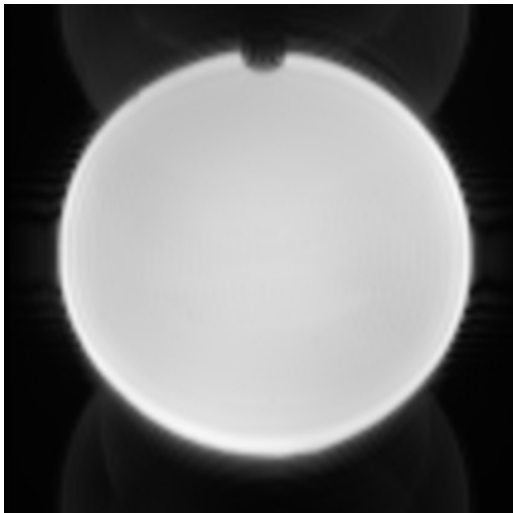

Standard Deviation image




image min: 0, image max: 323.258

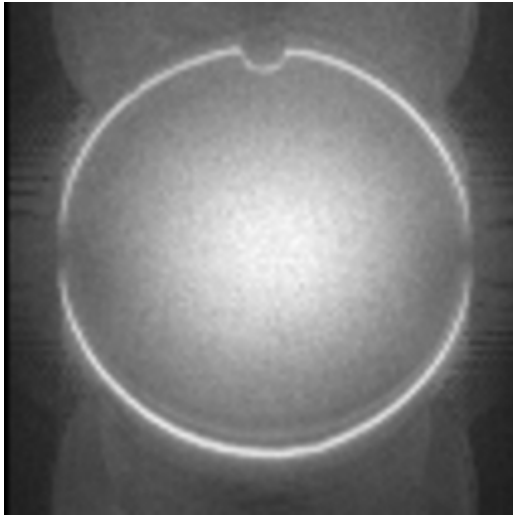

SFNR image



image min: 0, image max: 792.135

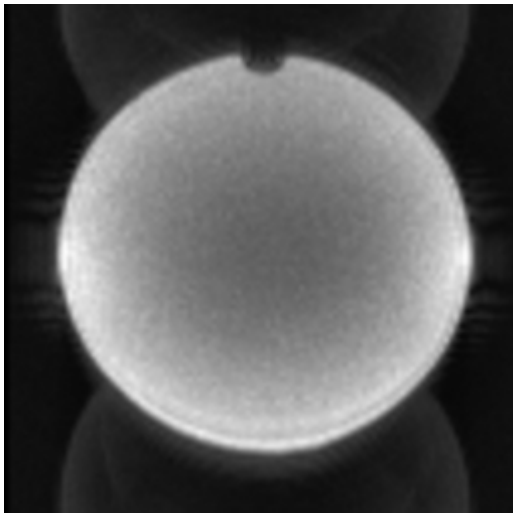

Acquisition parameters

dimensions	106x106x69x150
spacing	2mmx2mmx2mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	bdb969a6-da0b-4d
studyid	2025.12.23-09:07:19- STD-1.3.12.2.1107.5.99.3
seriesnumber	2
runnumber	1

magneticfield	3
description	ses-64ch_func-bold_acq-cmrrstd_task-warmup
scandate	2025-12-23
scantime	09:09:34
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-64ch_func-bold_acq-cmrrstd_task-warmup
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20251223090932.285000
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	2
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	RESEARCH

acquisitionduration	317.89499999999998
resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	150
imagecomments	Not for diagnostic use, Unaliased MB3/ PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	69
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00