






Phantom QA results for /fbirn_out/ses64ch_funcbold_acqsiemensstd_task slice 17




ld acqsiemensstd tasktsnr 4.xml [percent fluct (trend removed), drift, driftfit] :


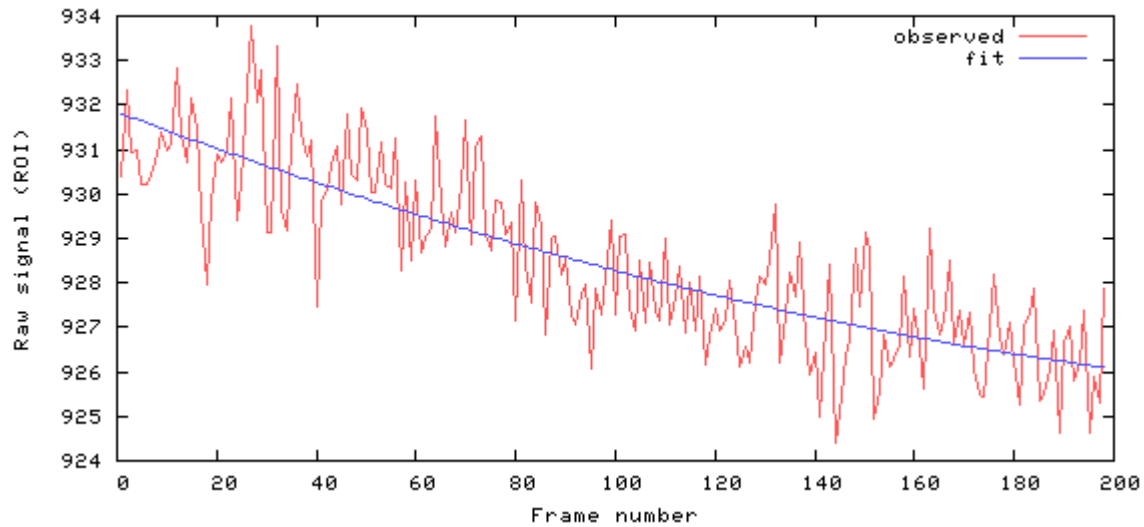


[mean, SNR, SFNR] = [928.5 102.2 104.3]

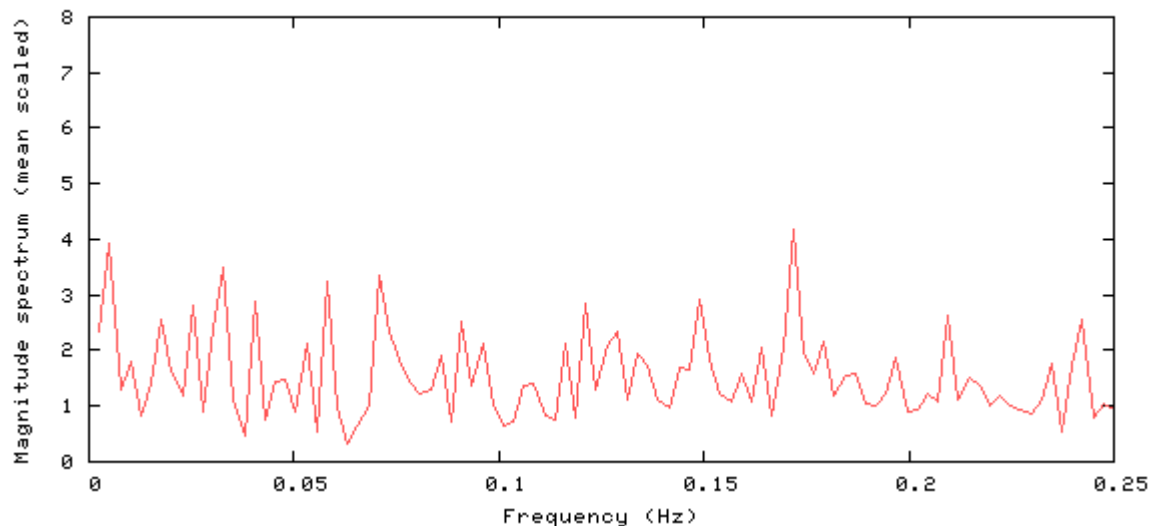

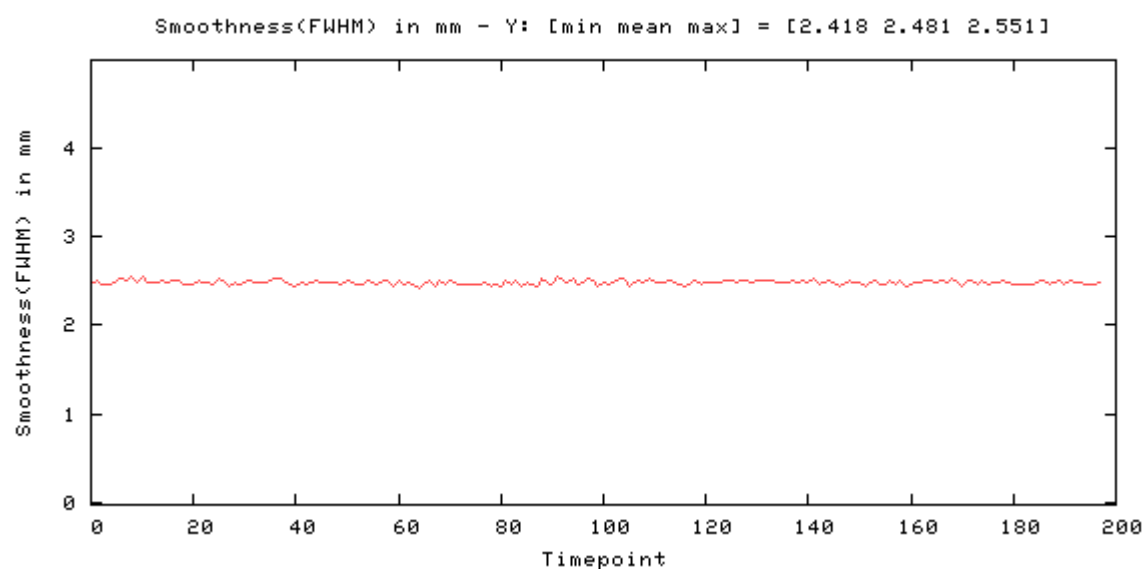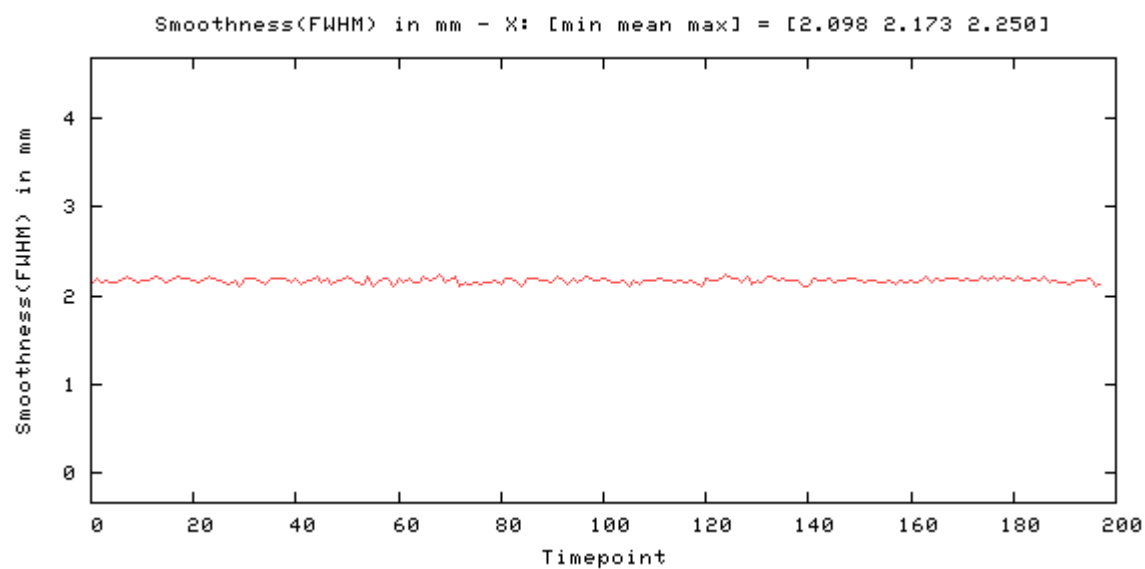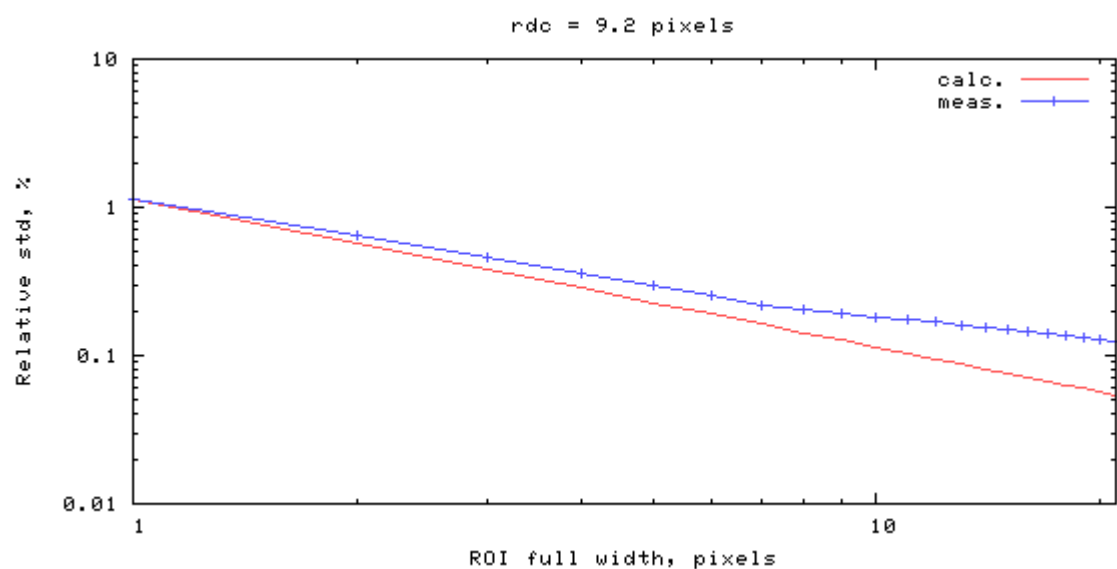

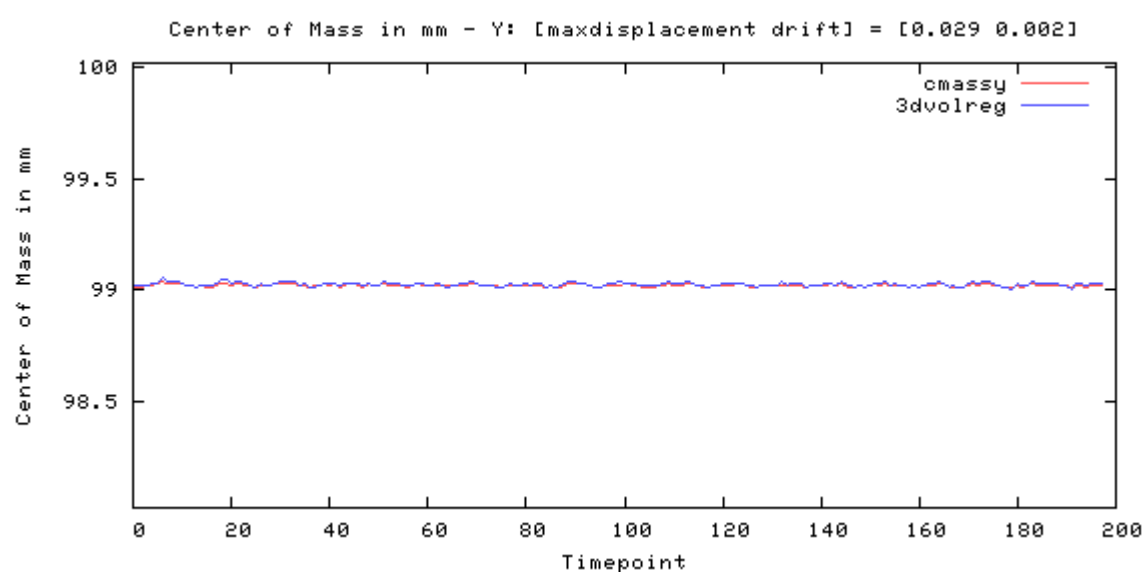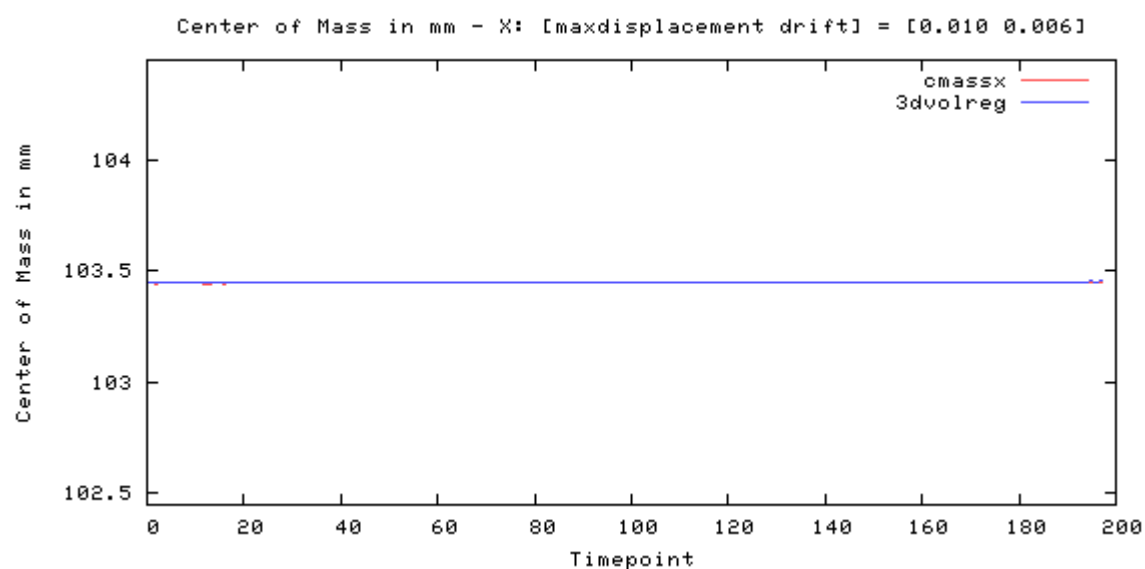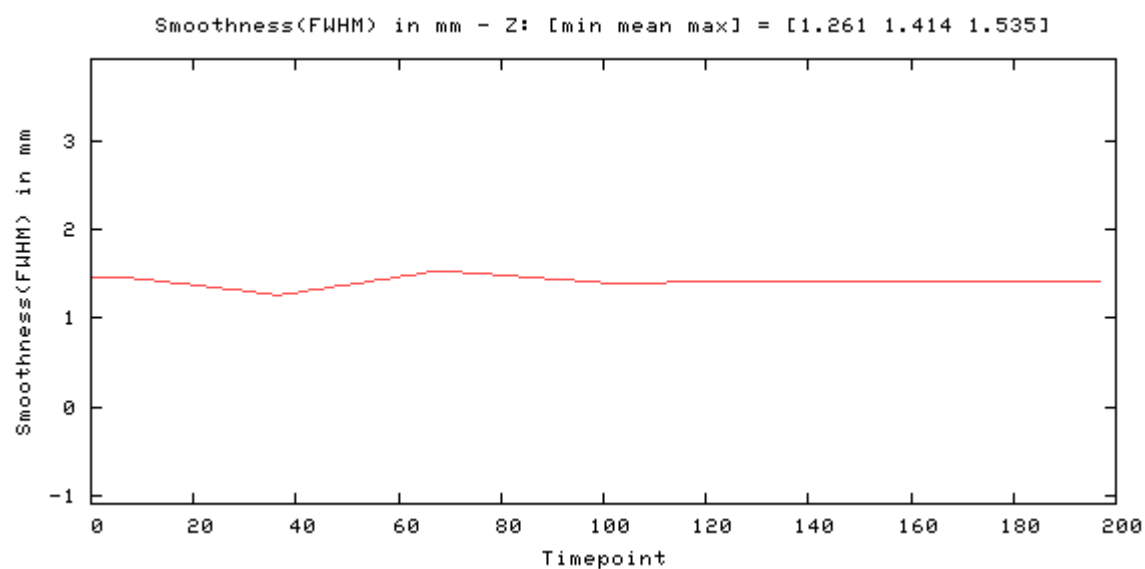


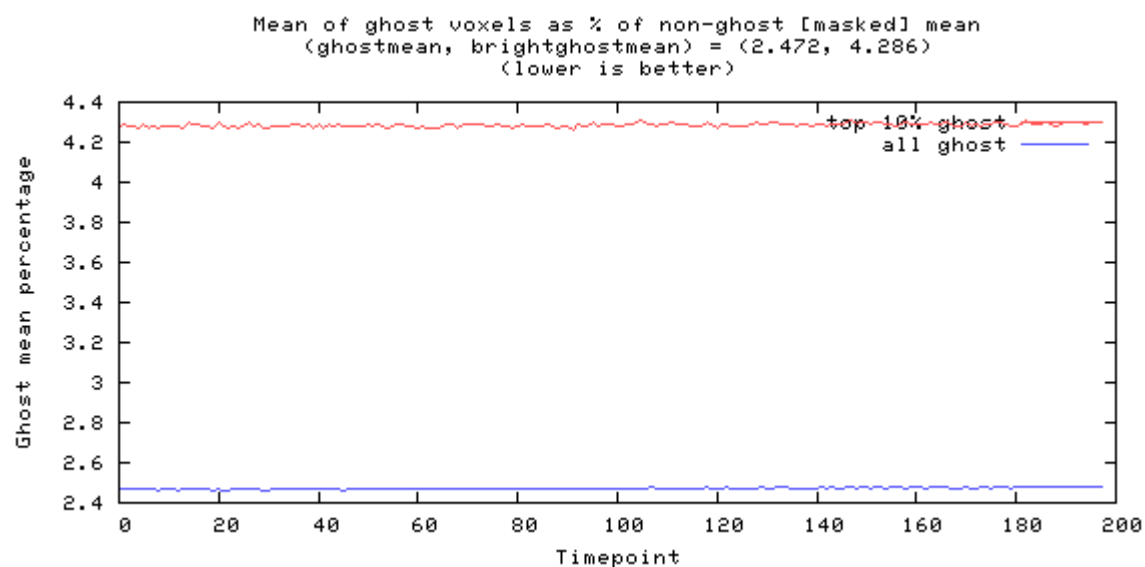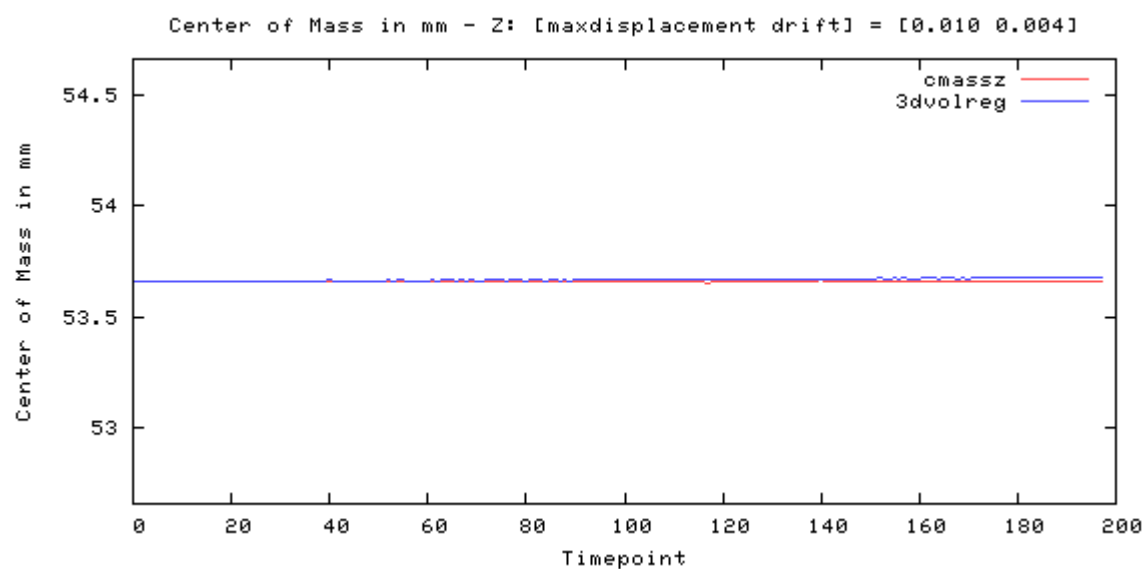

Odd-even difference image



image min: -412, image max: 496

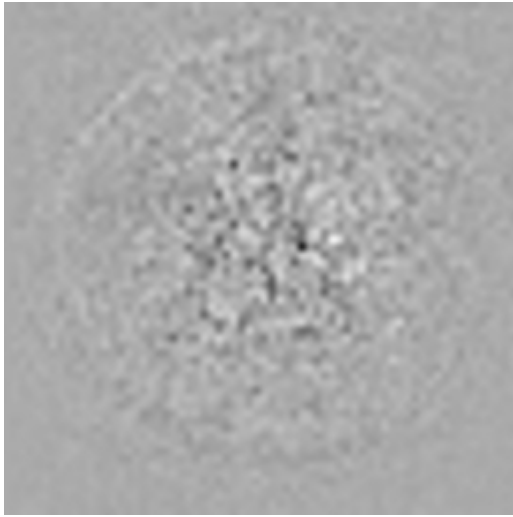

Mean image




image min: 0, image max: 1285.57



Standard Deviation image




image min: 0, image max: 12.5178

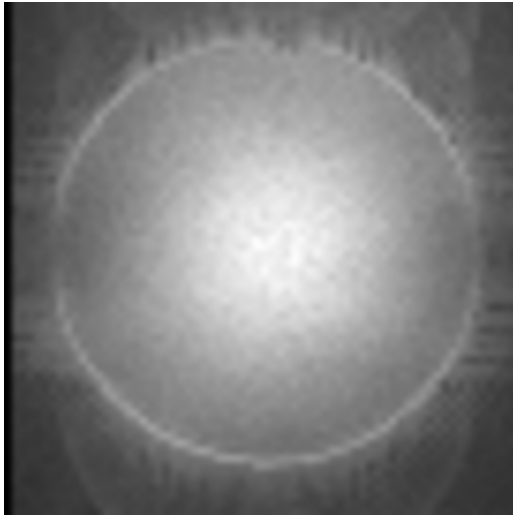

SFNR image



image min: 0, image max: 588.961

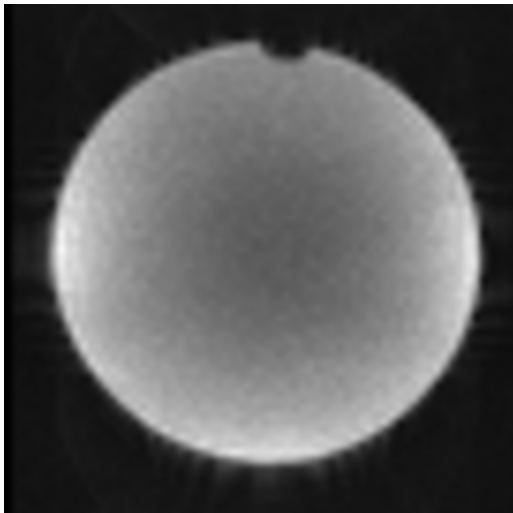

Acquisition parameters

dimensions	68x68x36x200
spacing	3mmx3mmx3mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	a82539b1-a6bd-44
studyid	Stability_20260112
seriesnumber	4
runnumber	1
magneticfield	3

description	ses-64ch_func-bold_acq-siemensstd_task-tsnr
scandate	2026-01-12
scantime	09:16:31
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	3
frequencydirection	1
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-64ch_func-bold_acq-siemensstd_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260112091629.620000
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	3
echotrainlength	33
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	PRODUCT
acquisitionduration	407.92500000000001

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	36
rows	68
columns	68
bitsallocated	16
bitsstored	12
highbit	11
lossyimagecompression	00