






Phantom QA results for / fbirn_out/ ses64ch_funcbold_acqcmrrstd_taskwar slice 34




old acqcmrrstd taskwarmup 2.xml [percent fluct (trend removed), drift, driftfit]


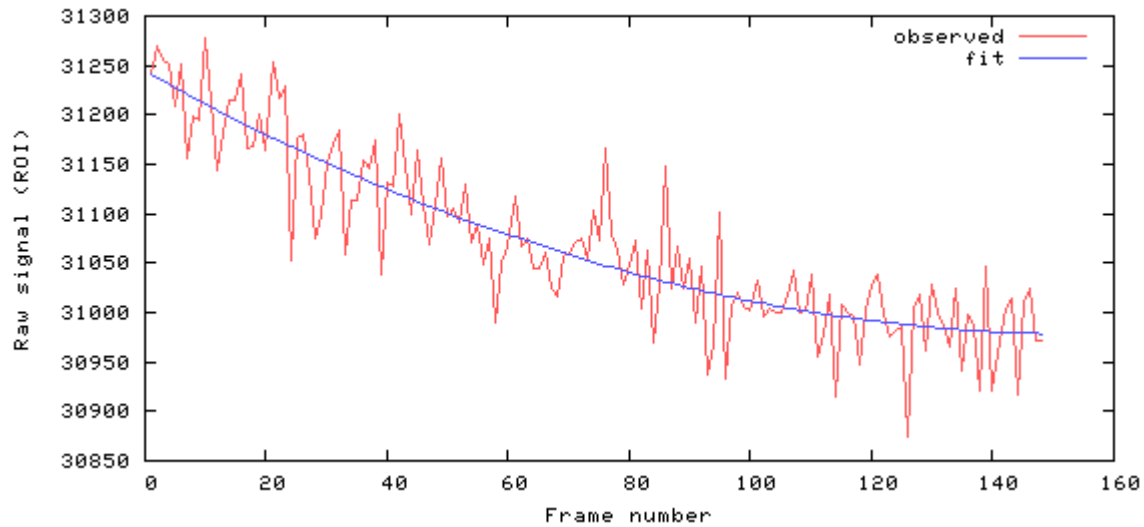


[mean, SNR, SFNR] = [31070.3 118.3 121.2]

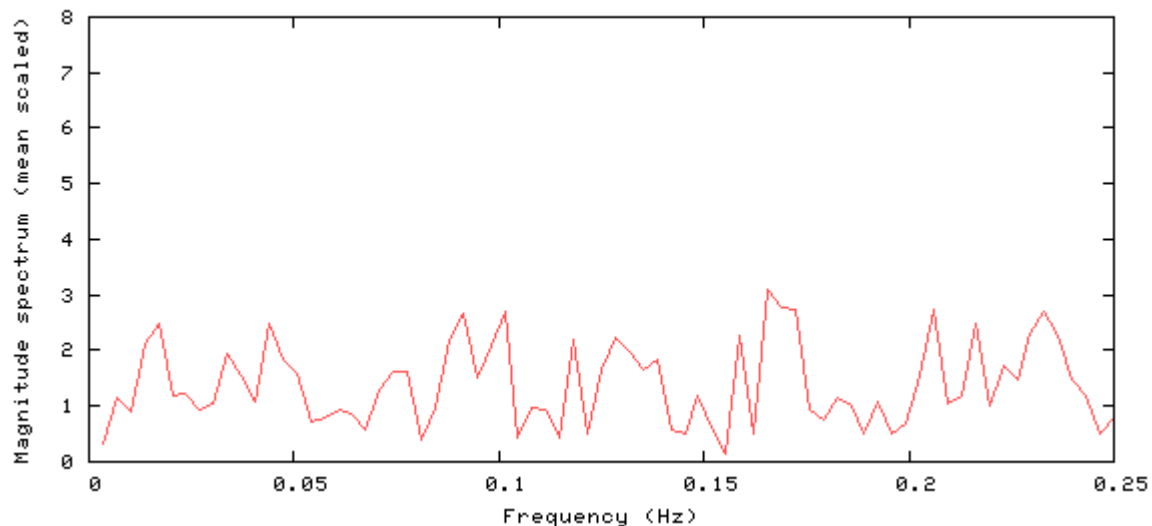

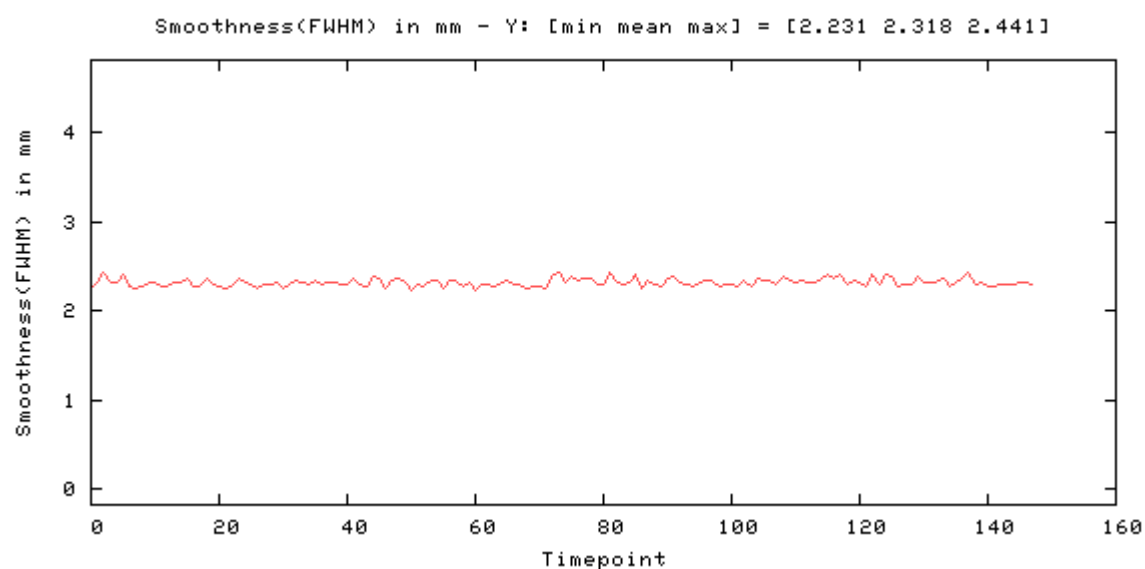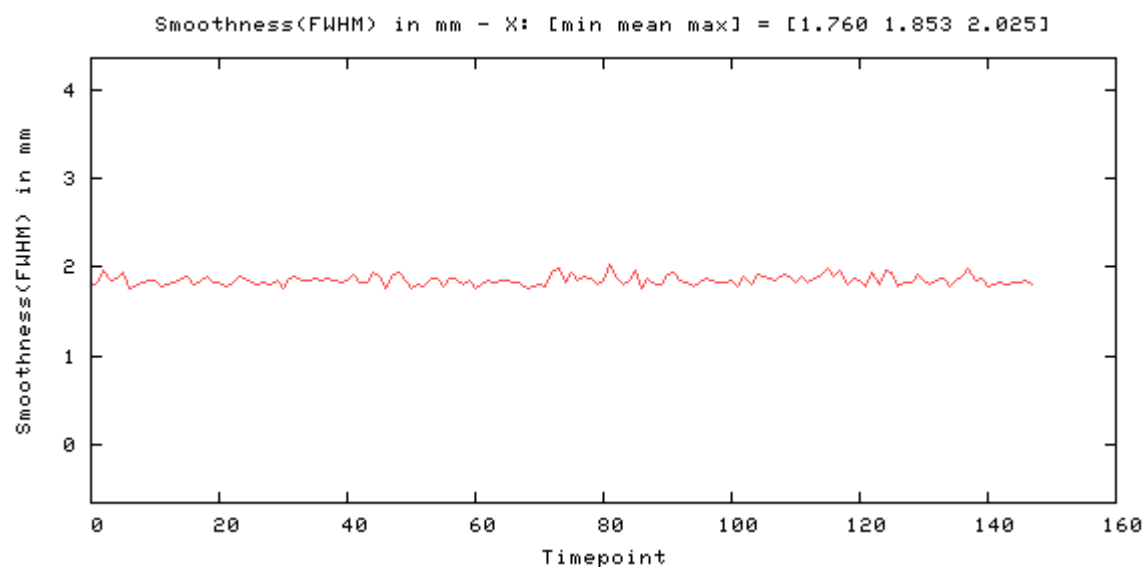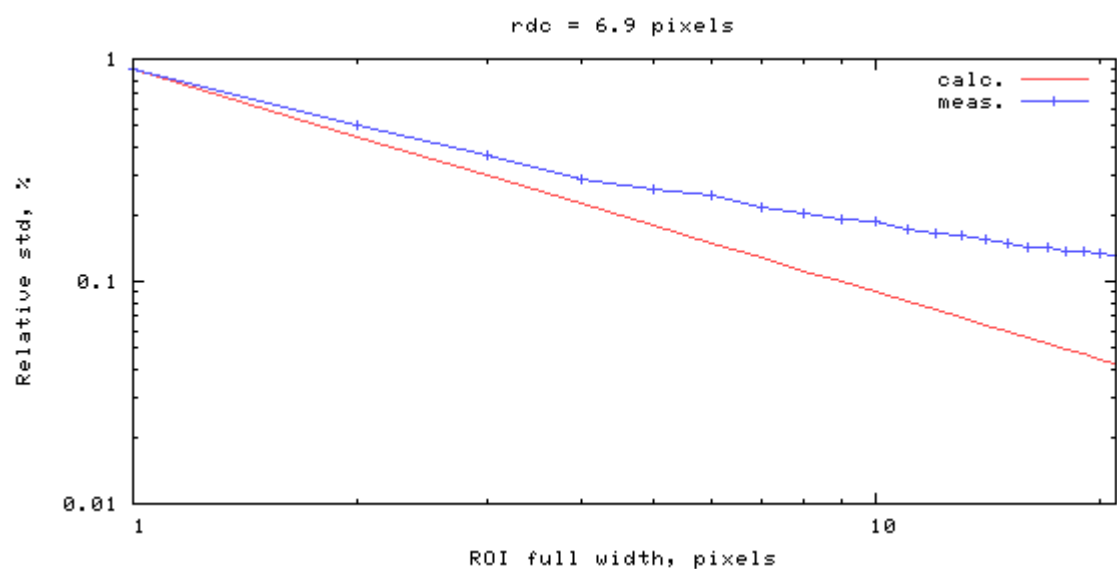

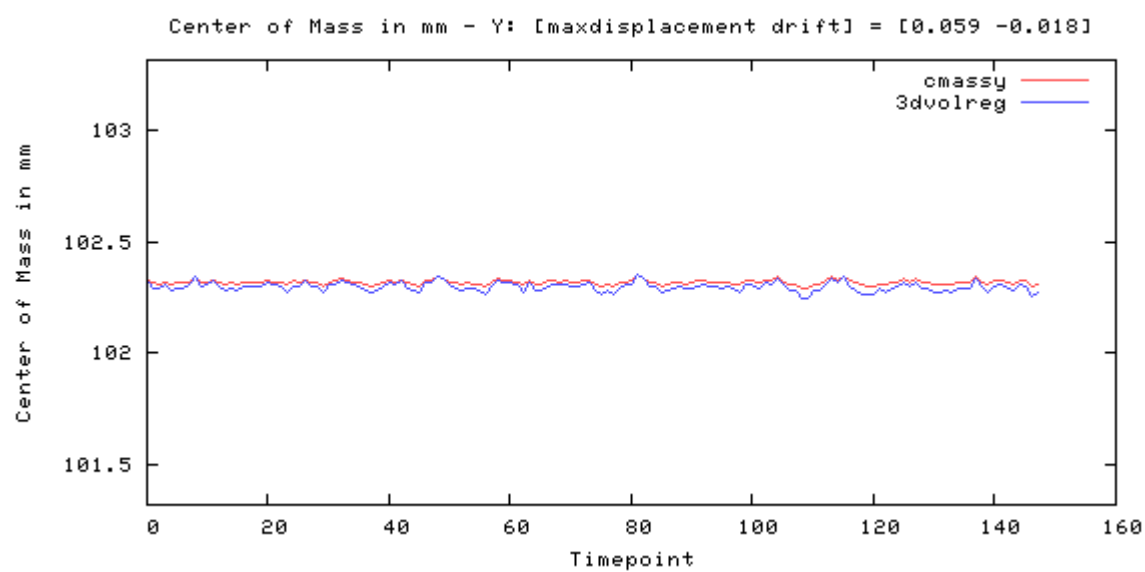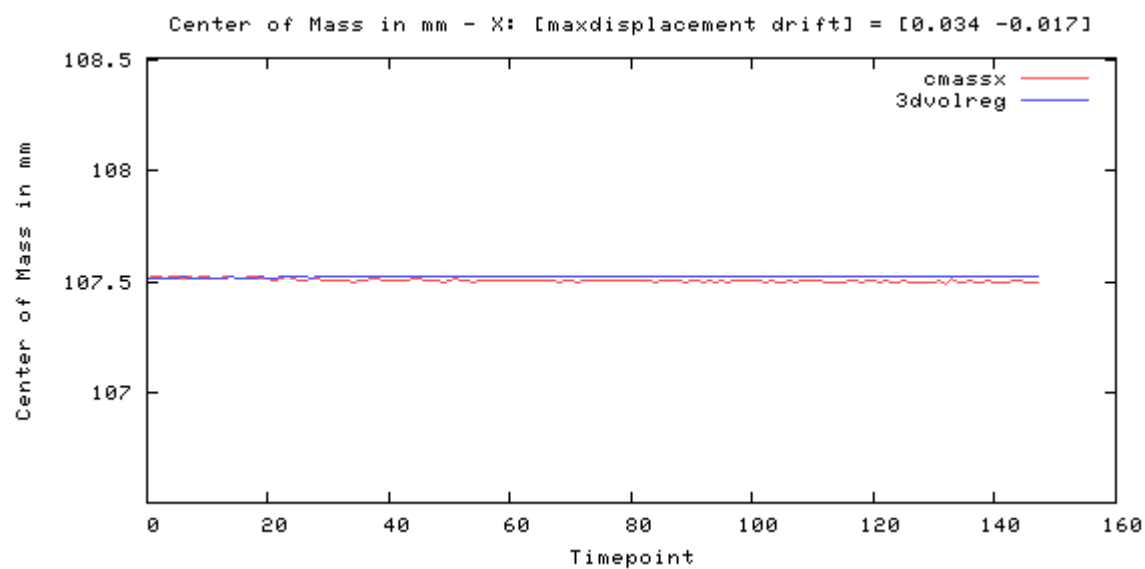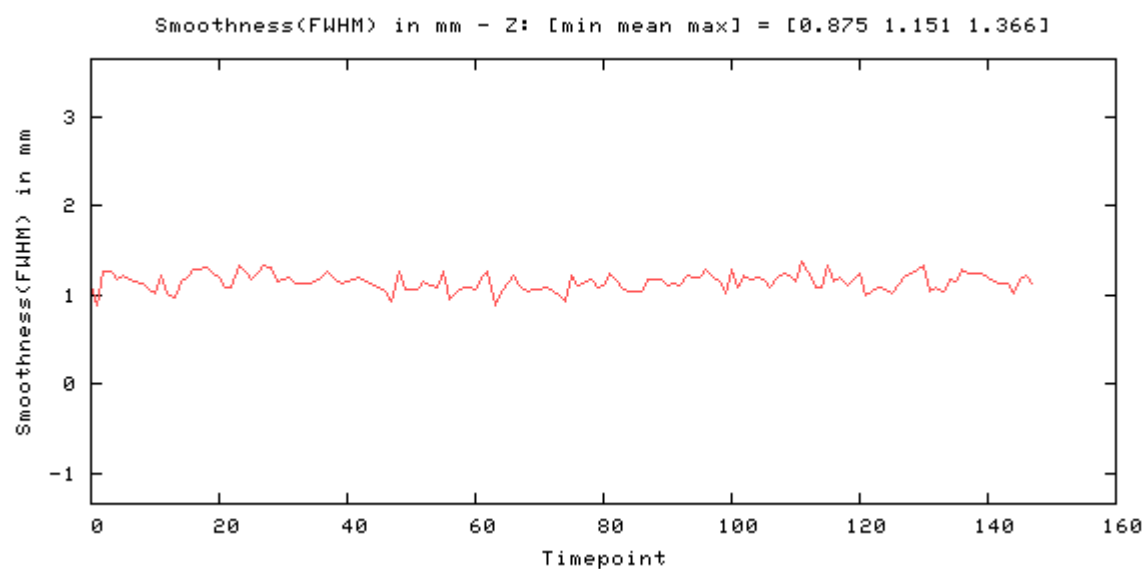

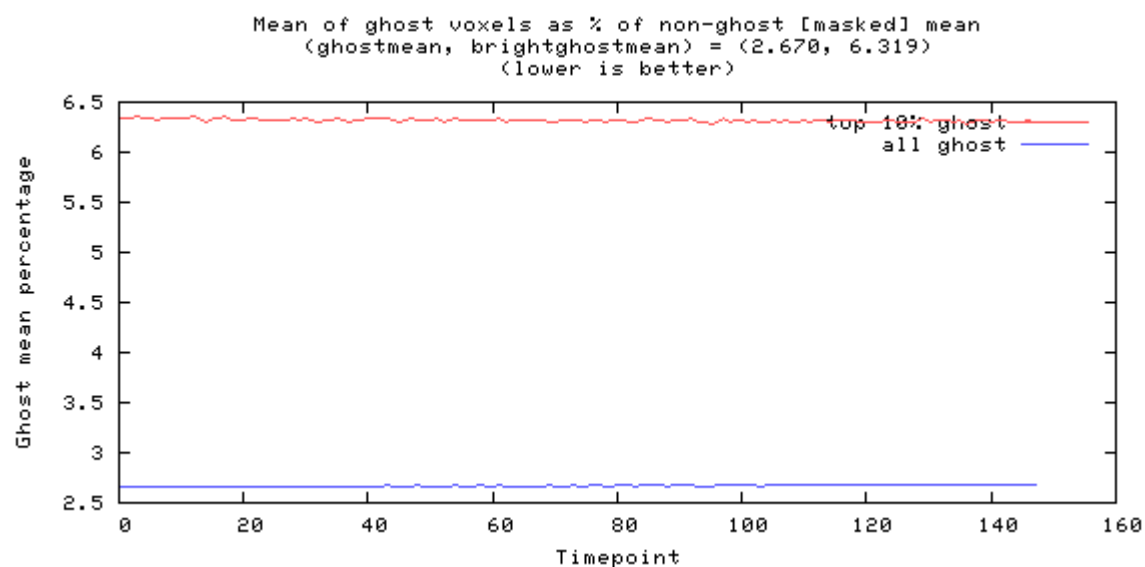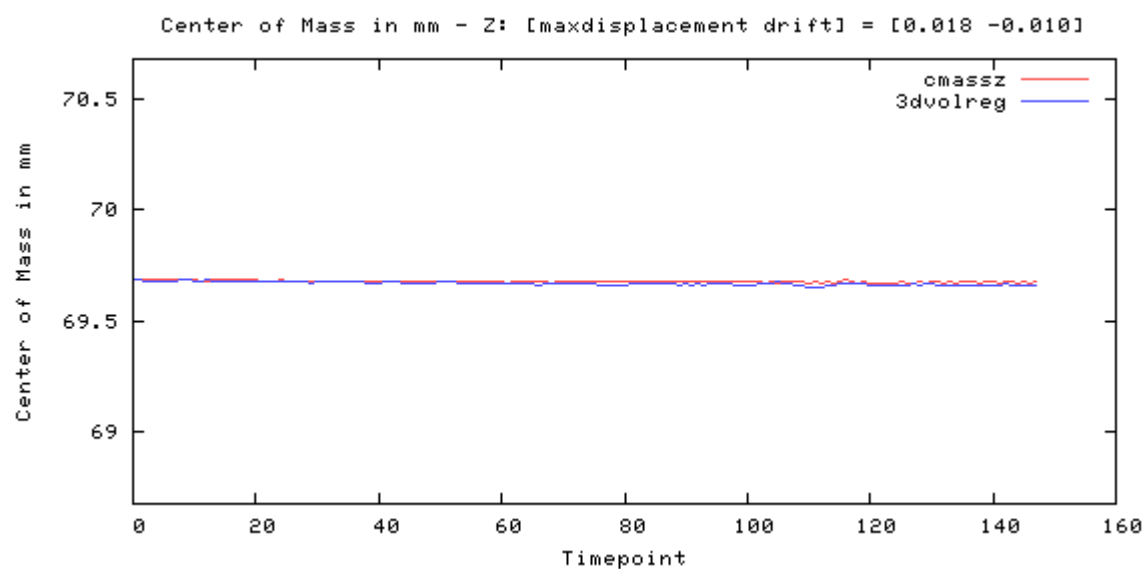


Odd-even difference image

-9187 10721




image min: -12463, image max: 10721

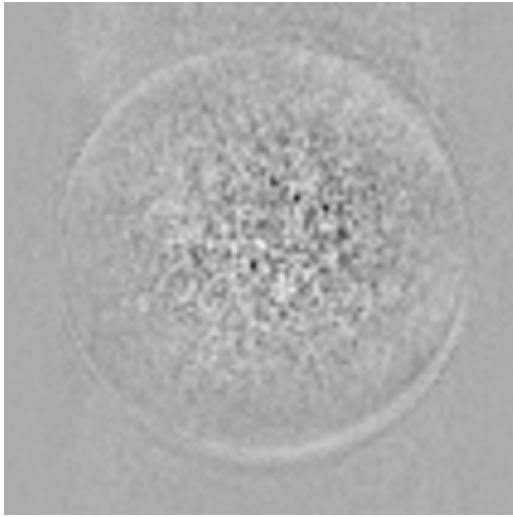

Mean image




image min: 0, image max: 42980.1

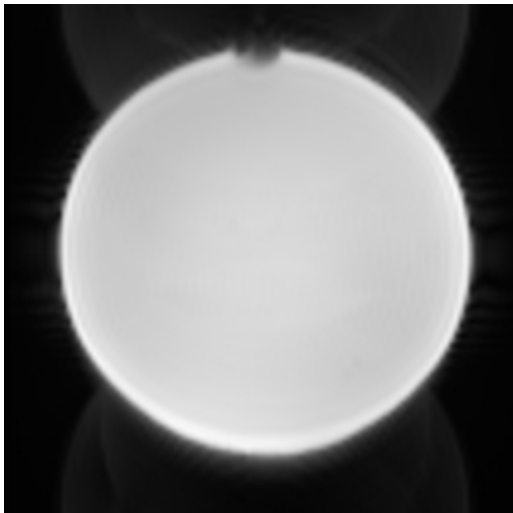

Standard Deviation image




image min: 0, image max: 354.065

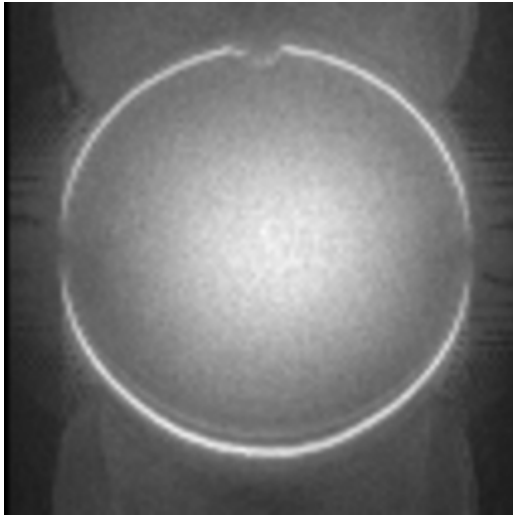

SFNR image



image min: 0, image max: 745.548

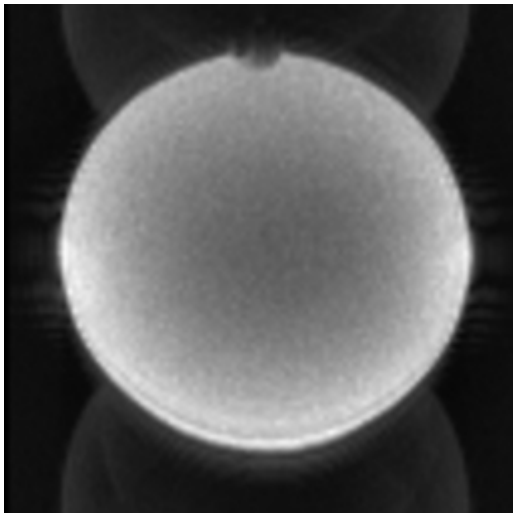

Acquisition parameters

dimensions	106x106x69x150
spacing	2mmx2mmx2mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	c9dbdd47-8558-45
studyid	Stability_20260120
seriesnumber	2
runnumber	1
magneticfield	3

description	ses-64ch_func-bold_acq-cmrrstd_task-warmup
scandate	2026-01-20
scantime	09:28:07
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-64ch_func-bold_acq-cmrrstd_task-warmup
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260120092805.420000
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	2
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	RESEARCH
acquisitionduration	317.89499999999998

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	150
imagecomments	Not for diagnostic use, Unaliased MB3/ PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	69
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00