


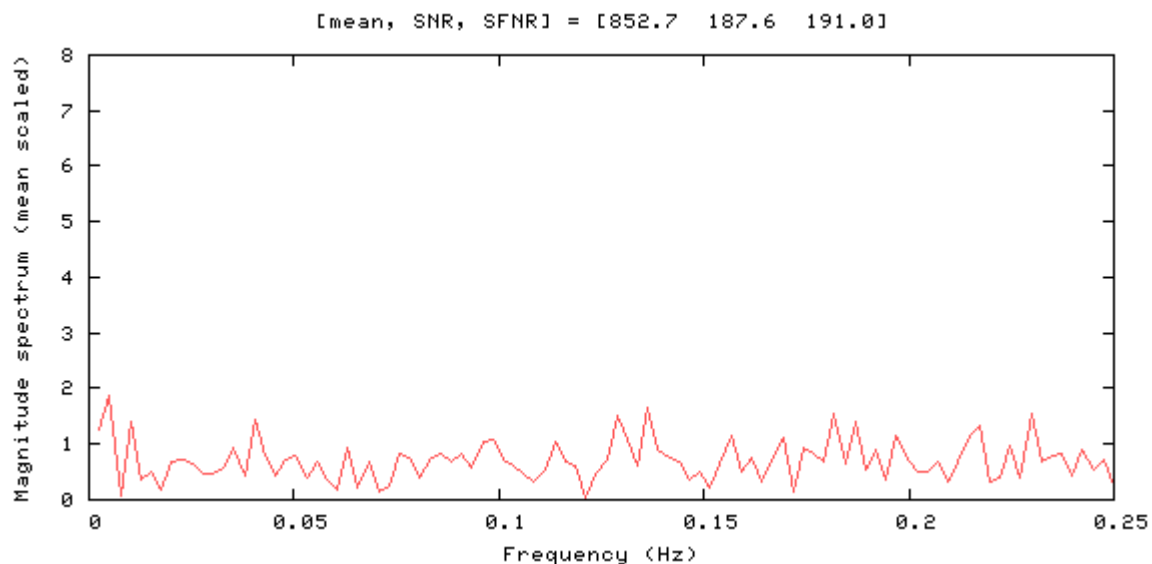
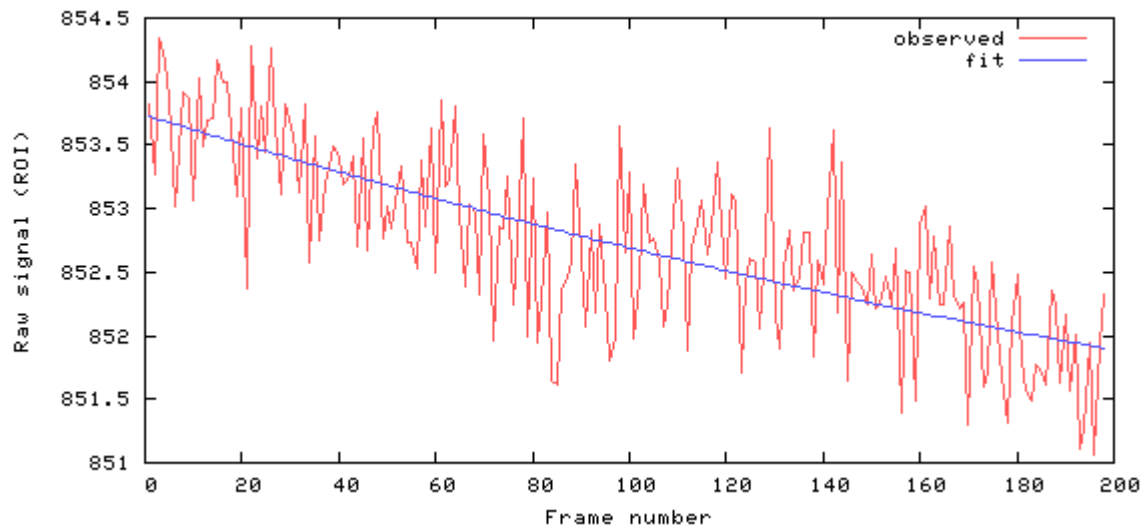
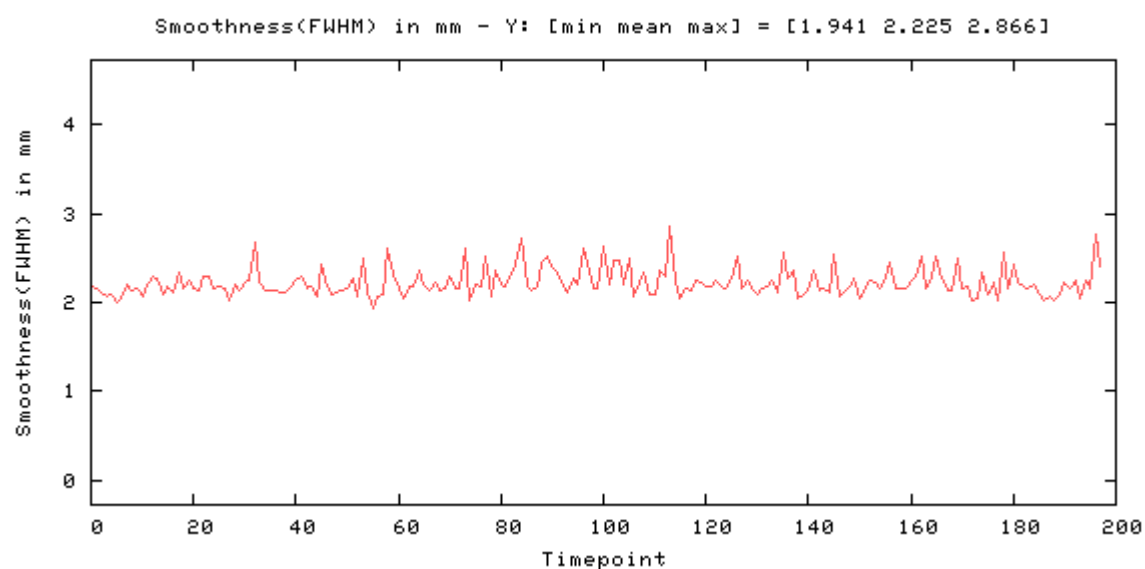
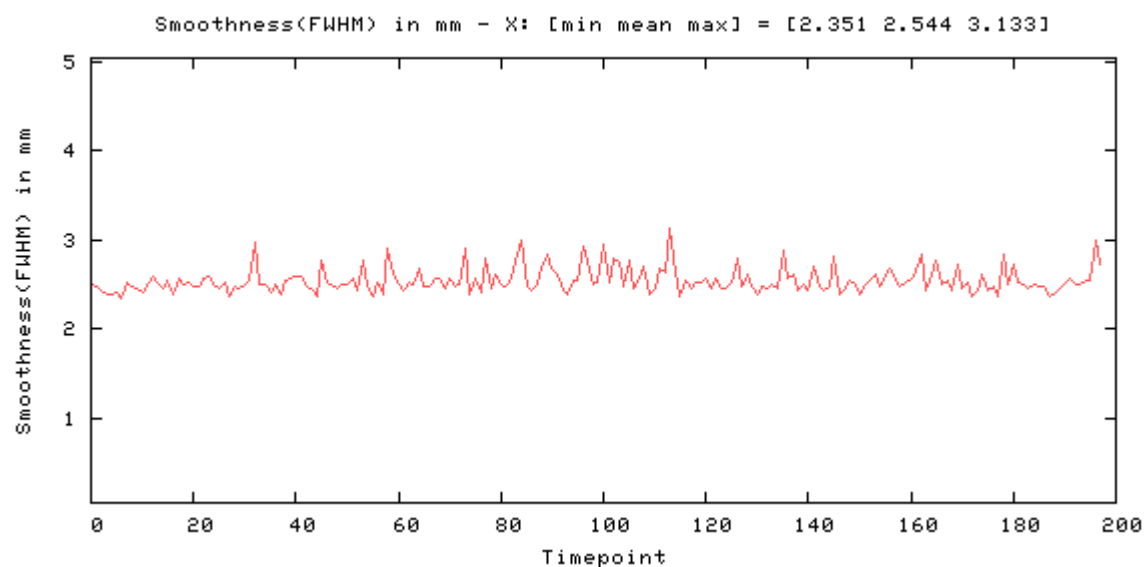
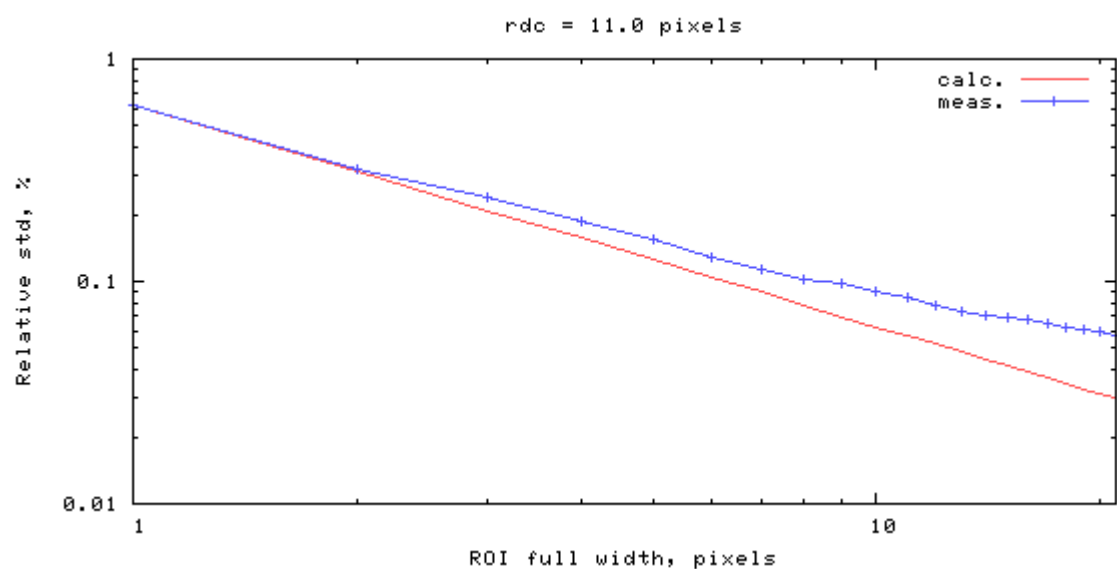




Phantom QA results for /fbirn_out/ses-64ch_func-bold_acq-siemensbasic_task-tsnr.xml, slice 14

old_acq-siemensbasic_task-tsnr.xml [percent fluct (trend removed), drift, driftfi-

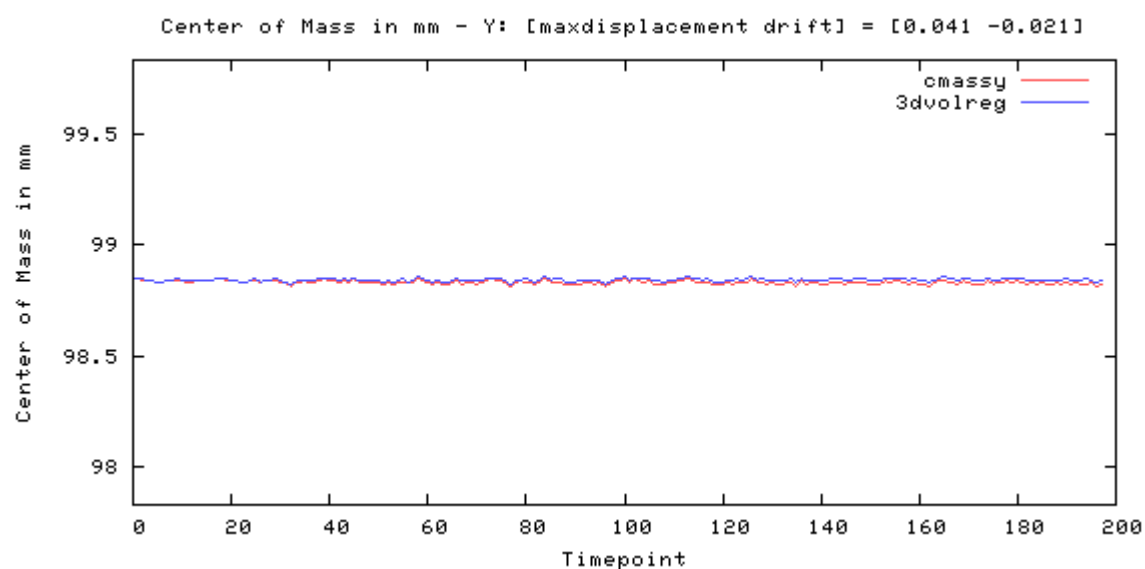

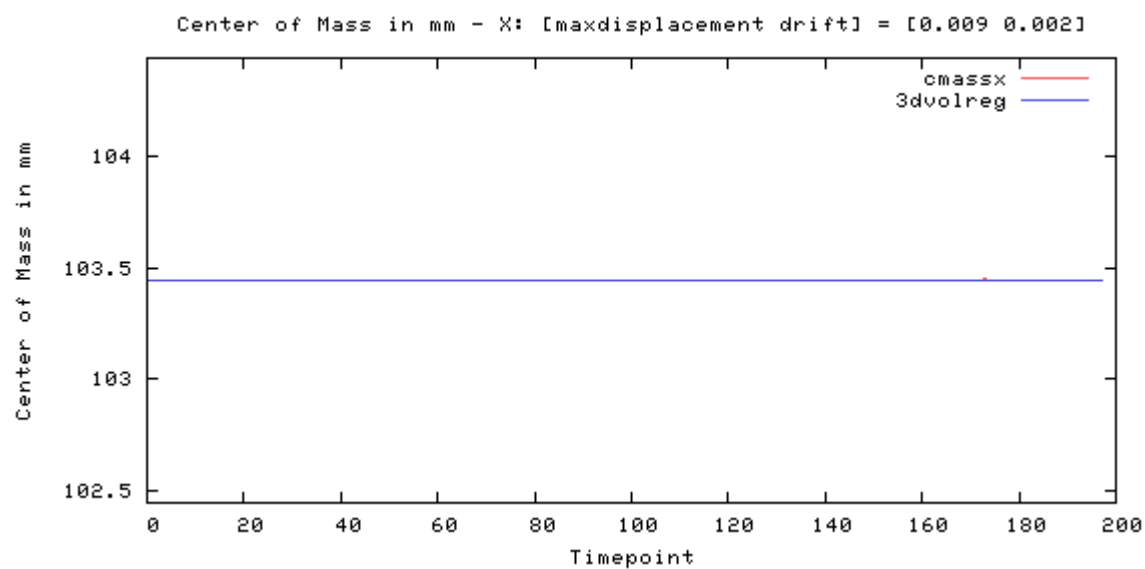

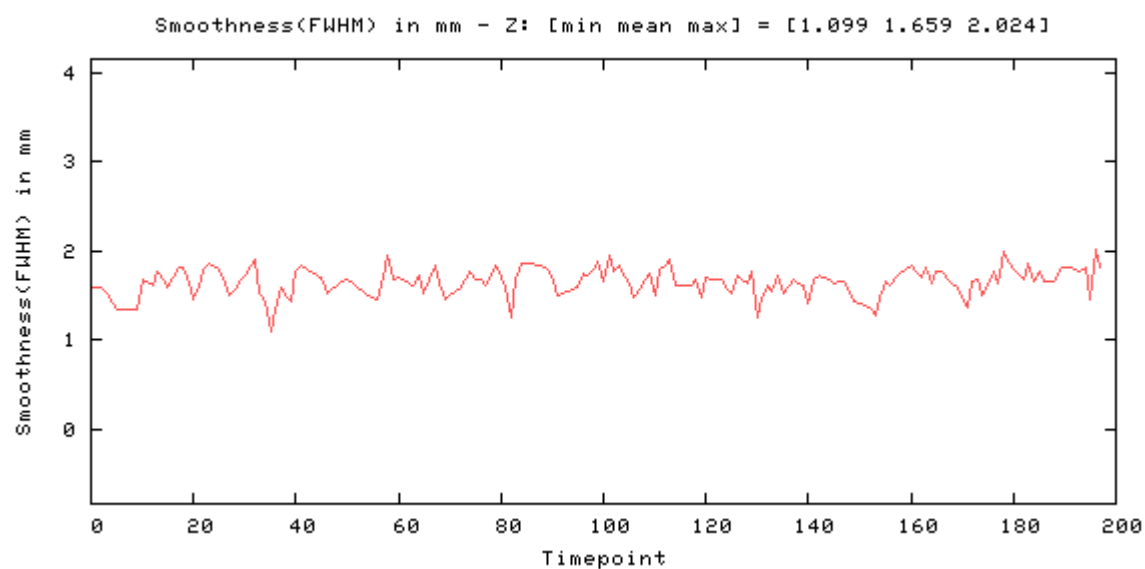


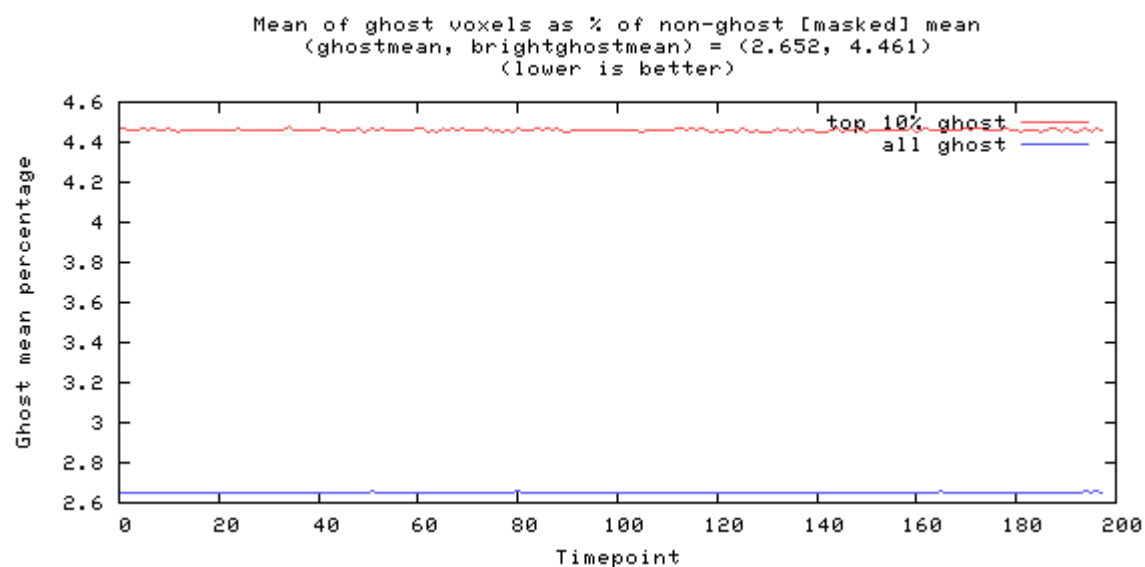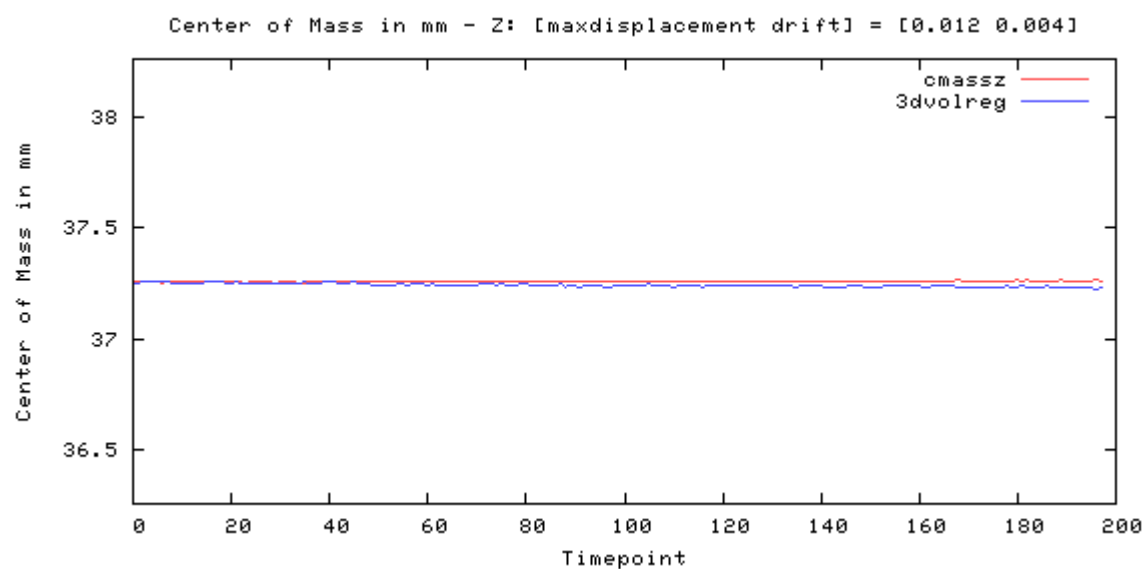

Odd-even difference image



image min: -199, image max: 206

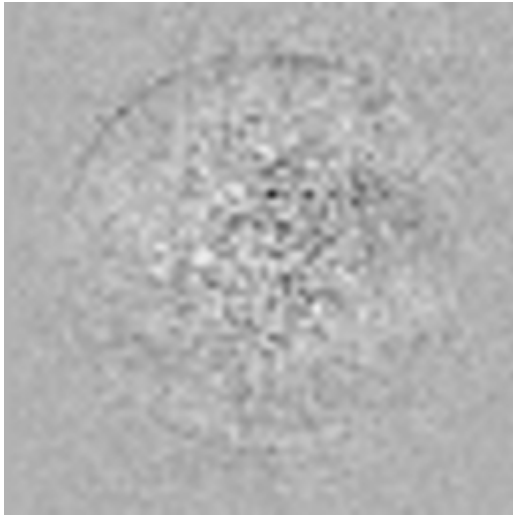

Mean image




image min: 0, image max: 1263.89



Standard Deviation image




image min: 0, image max: 7.35089

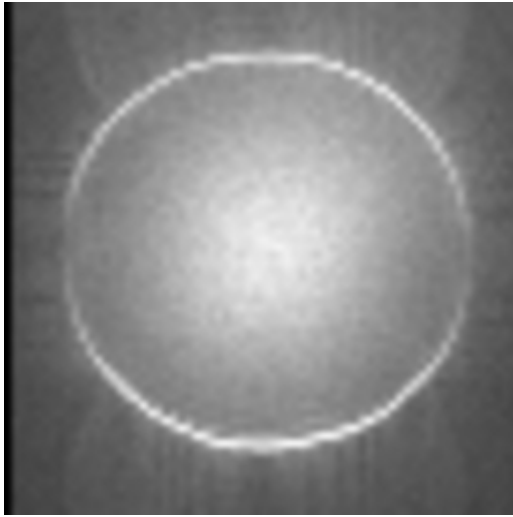

SFNR image



image min: 0, image max: 777.515

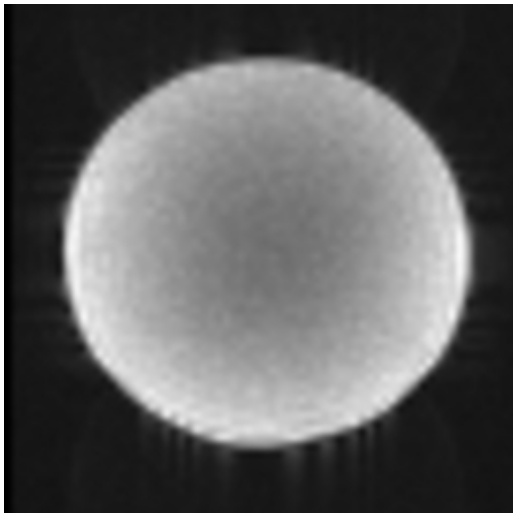

Acquisition parameters

dimensions	68x68x30x200
spacing	3mmx3mmx3mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	310dc38b-b5e3-41
studyid	Stability_20260720
seriesnumber	8
runnumber	1
magneticfield	3

description	ses-64ch_func-bold_acq-siemensbasic_task-tsnr
scandate	2026-07-20
scantime	09:57:13
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	3
frequencydirection	1
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-64ch_func-bold_acq-siemensbasic_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260720095711.097500
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	3
echotrainlength	68
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	PRODUCT
acquisitionduration	405.91500000000002

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	30
rows	68
columns	68
bitsallocated	16
bitsstored	12
highbit	11
lossyimagecompression	00