






Phantom QA results for /fbirn_out/ses64ch_funcbold_acqsiemensbasic_taslice 14




1d acqsiemensbasic tasktsnr 8.xml [percent fluct (trend removed), drift, driftfit:


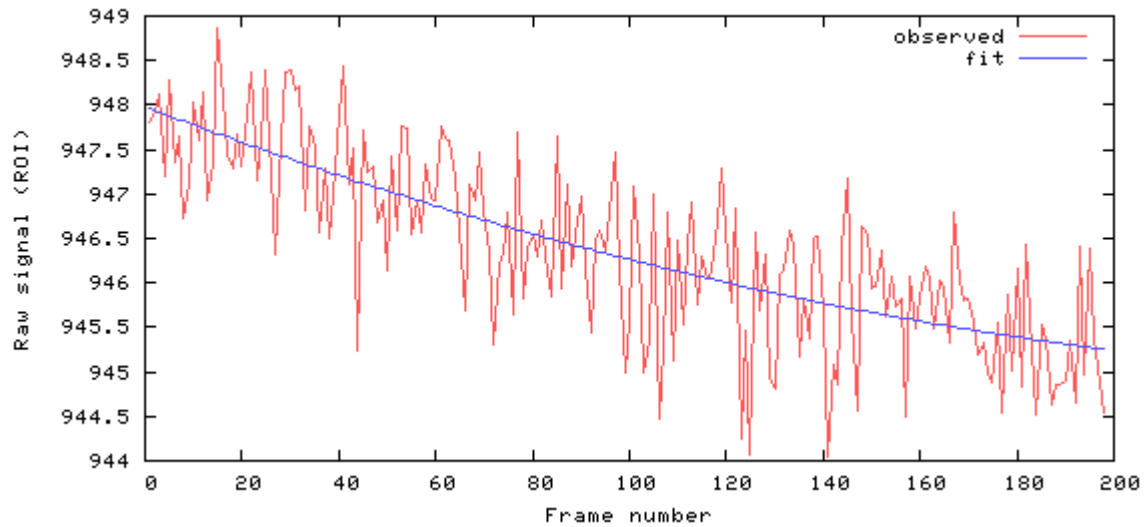


[mean, SNR, SFNR] = [946.4 164.0 168.4]

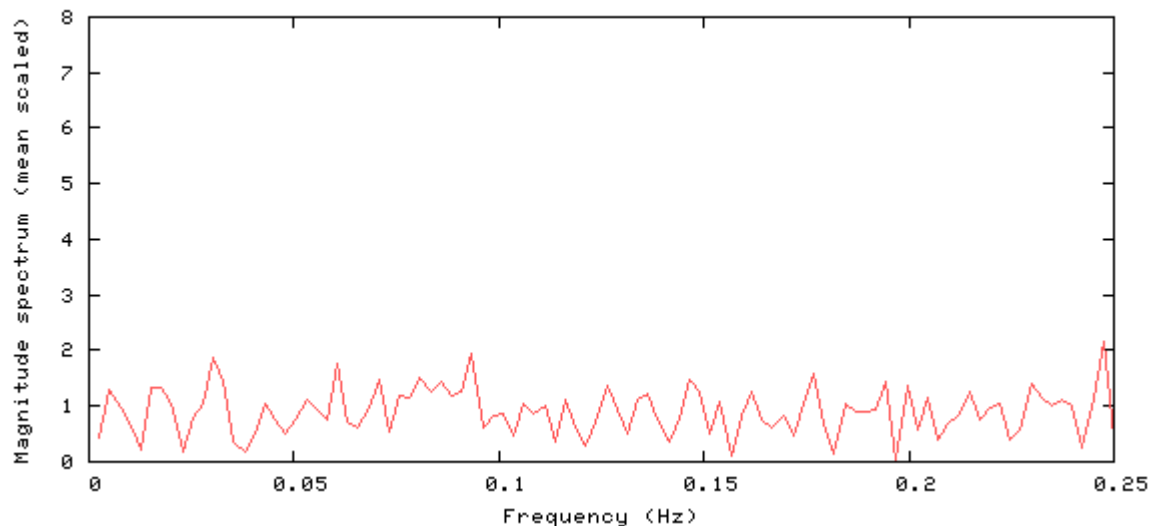

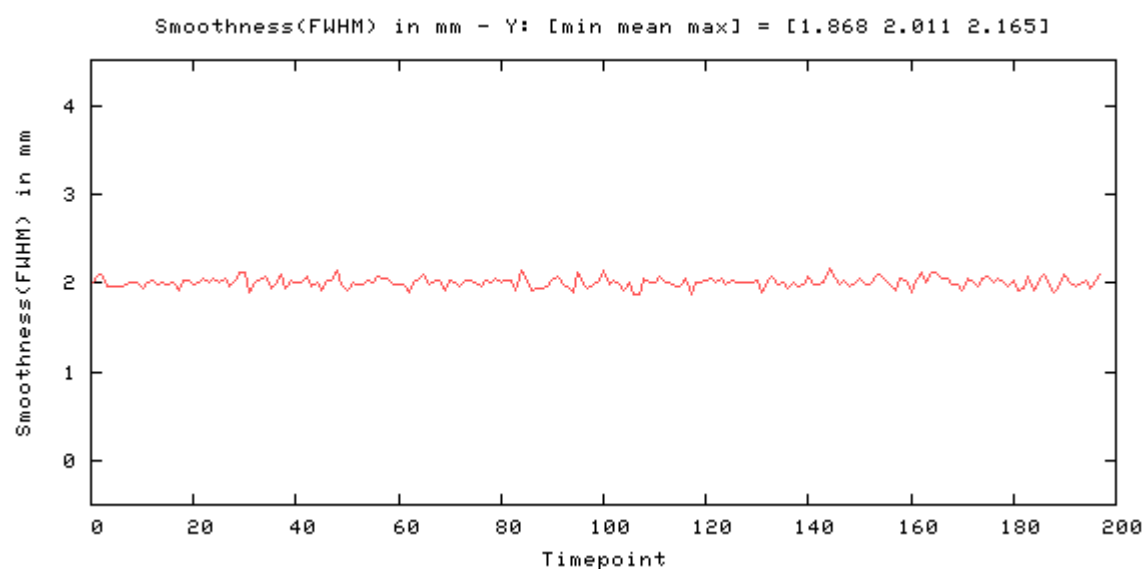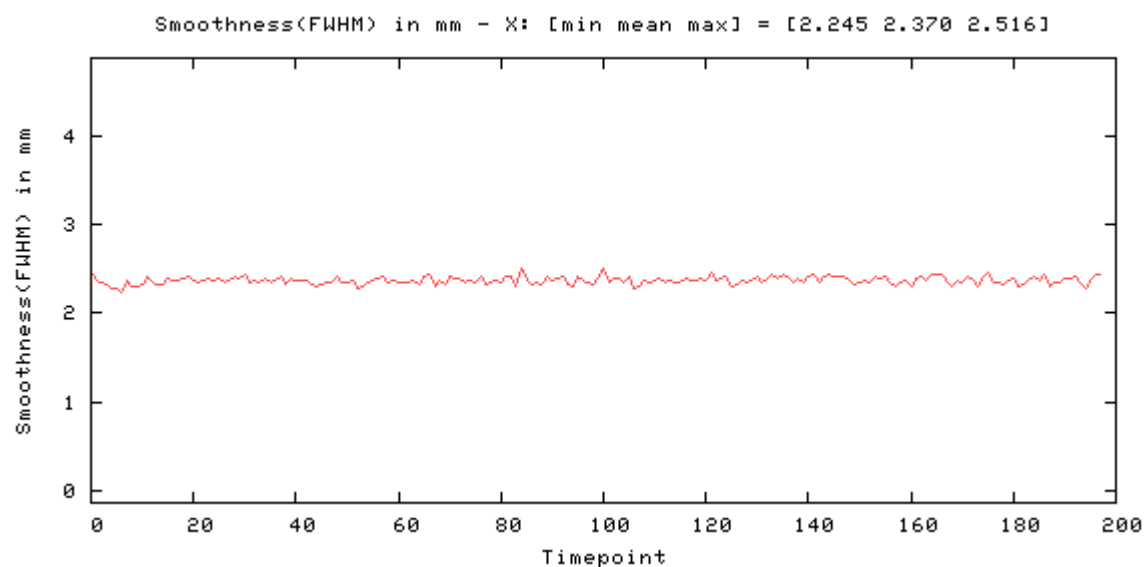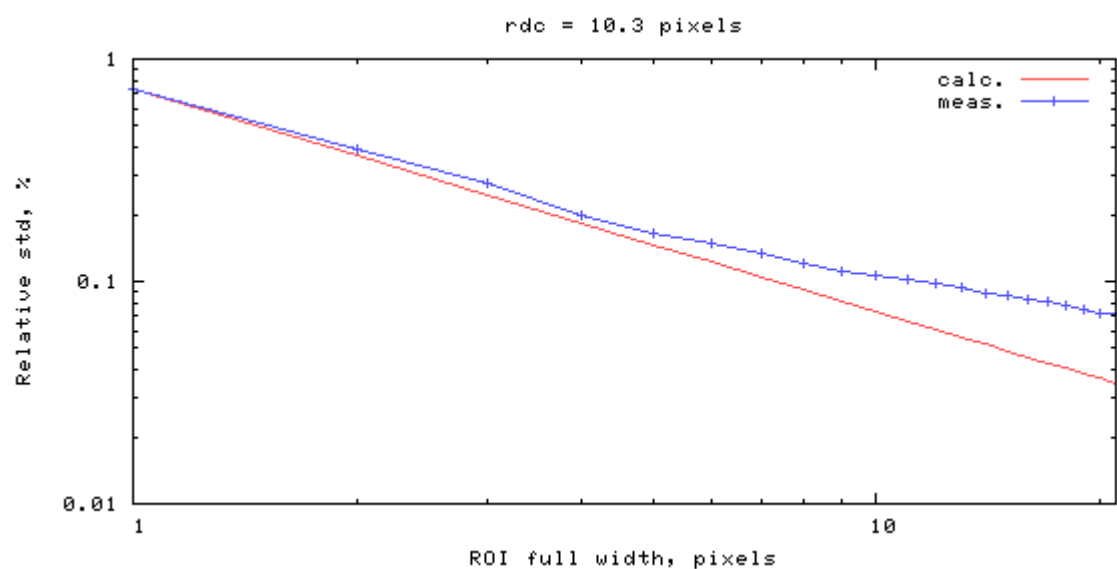

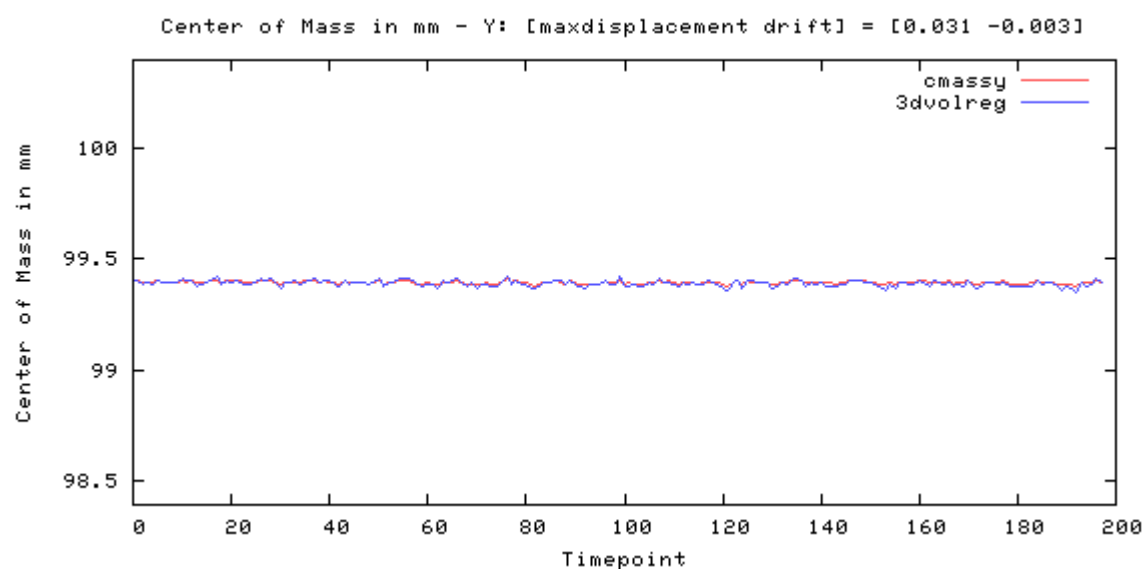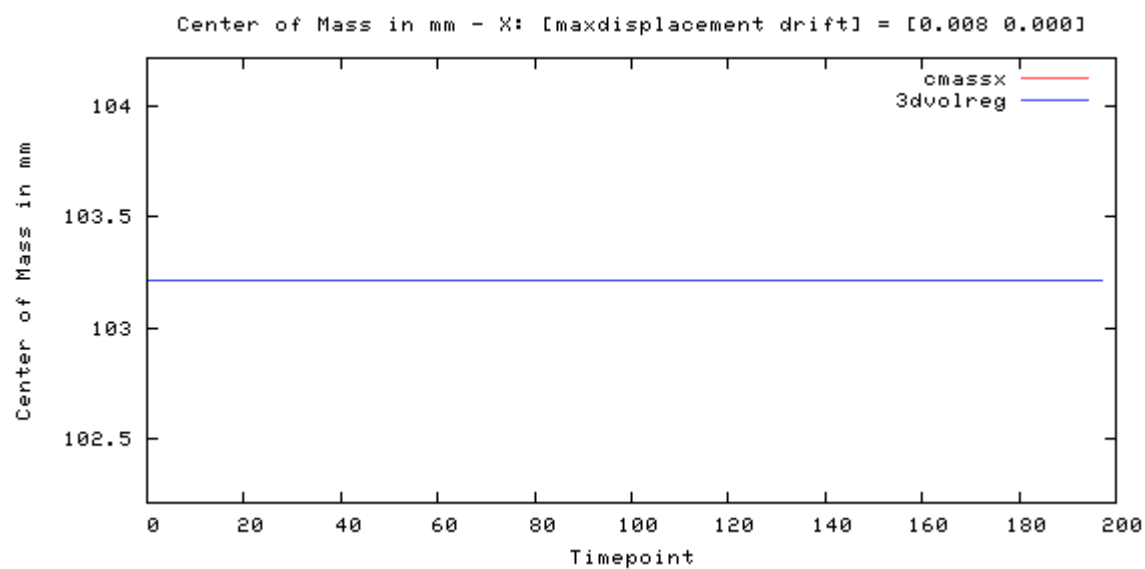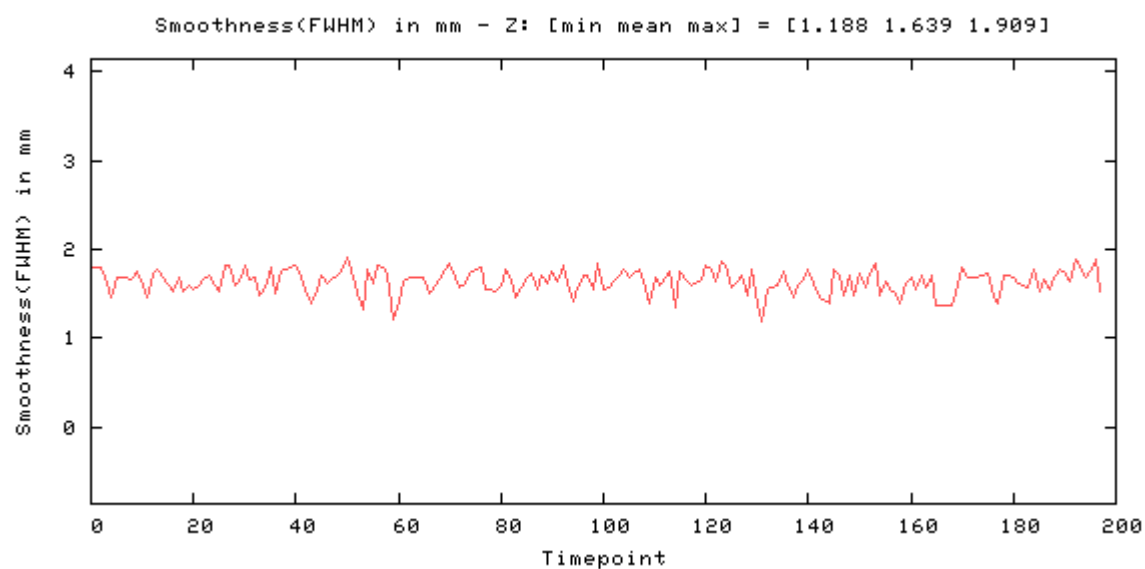


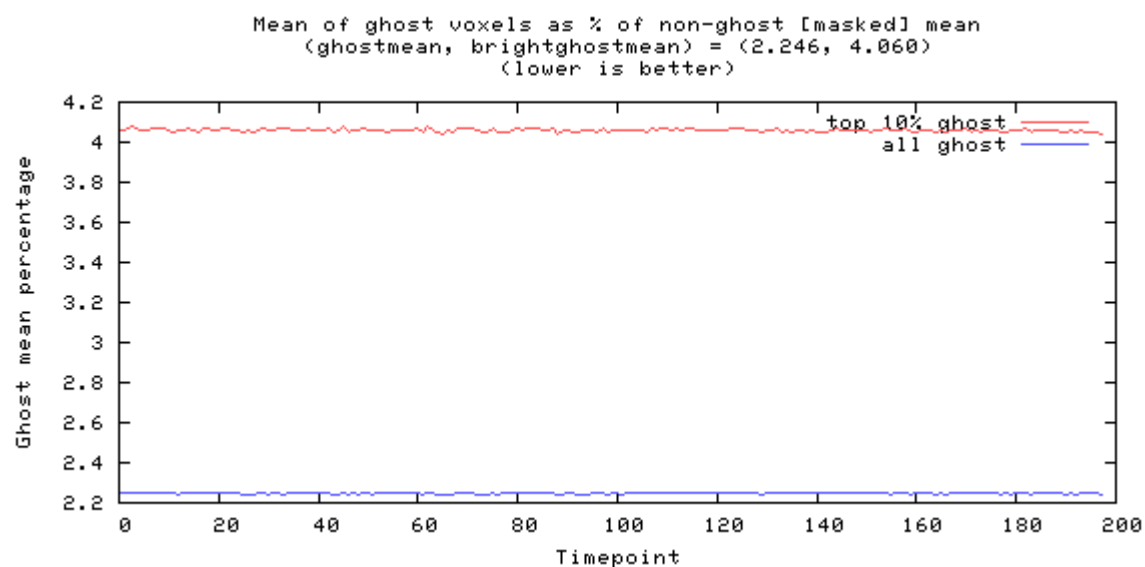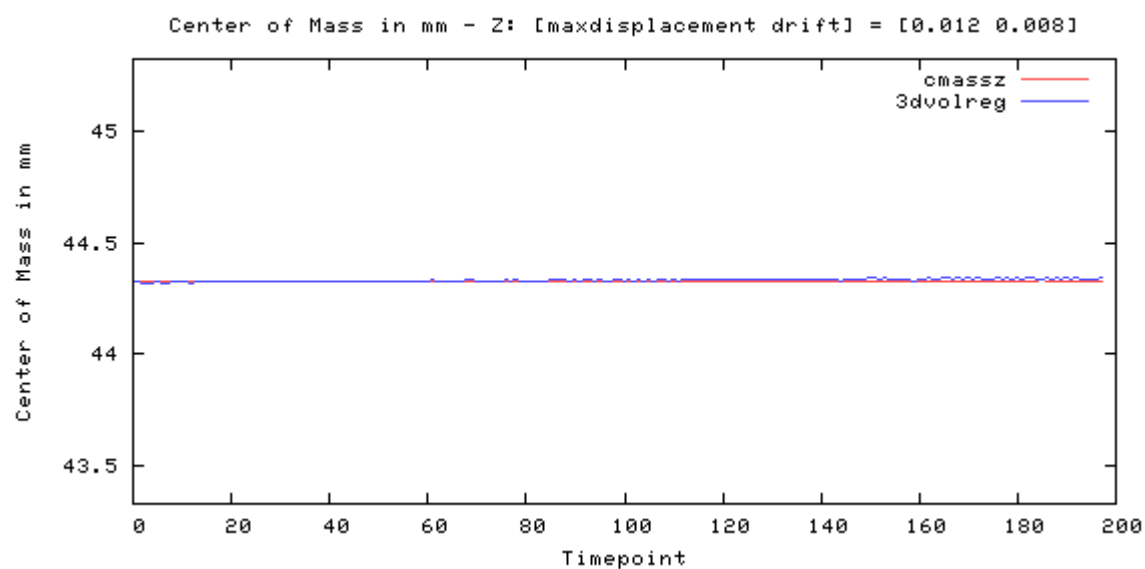

Odd-even difference image



image min: -281, image max: 229

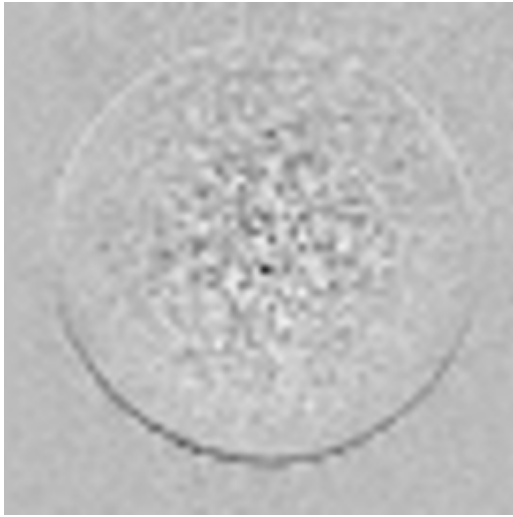

Mean image




image min: 0, image max: 1326.92



Standard Deviation image




image min: 0, image max: 7.86193

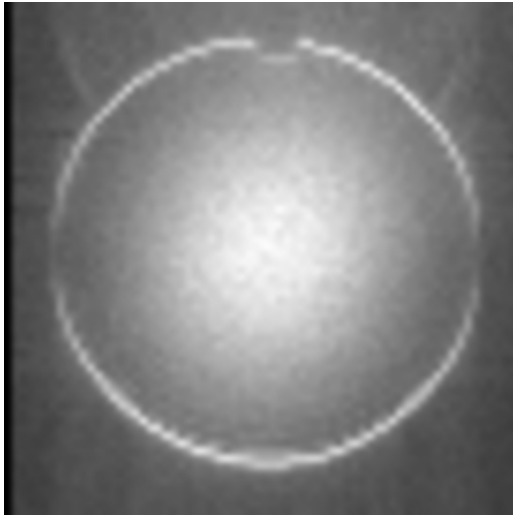

SFNR image



image min: 0, image max: 973.847

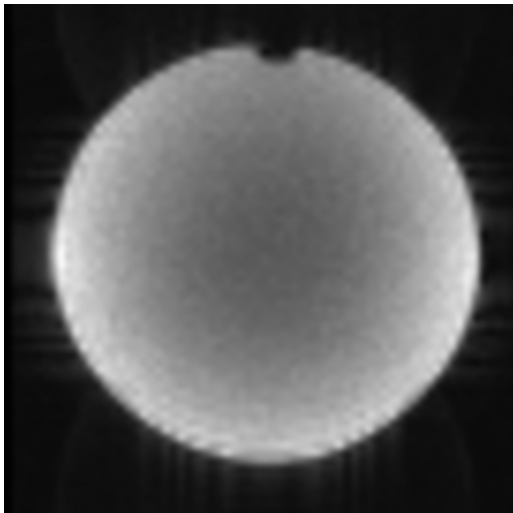

Acquisition parameters

dimensions	68x68x30x200
spacing	3mmx3mmx3mmx1[units?]
gap	0mmx0mmx0mmx0[units?]
scanner	AWP213018
psdname	
examnumber	d9a9a0b2-2cfa-4d
studyid	Stability_20250408
seriesnumber	8
runnumber	1-200
scanningsequence	GR/EP

sequencevariant	SK
magneticfield	3
description	ses-64ch_func-bold_acq-siemensbasic_task-tsnr
scandate	2025-04-08
scantime	09:55:25
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	3
frequencydirection	1
acquisitionmatrix	68 68
bandwidth	166.668
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
ti	
protocolname	ses-64ch_func-bold_acq-siemensbasic_task-tsnr
sliceorder	2,4,6,8,10,12,14,16,18,20,22,24,26,28,30,1,3,5,7,9,11,
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,6
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20250408095524.460000
volumetricproperties	DISTORTED
scanoptions	FS/PER
mracquisitiontype	2D
sequencename	*epfid2d1_68
angioflag	N
numaverages	1
imagingfrequency	123.253285
imagednucleus	1H
spacingbetweenslices	3
numphaseencodingsteps	68
echoTRAINlength	68
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61

triggertime	
beatrejectionflag	N
cardiacnumberofimages	1
inplane phase encoding direction	COL
variableflipangleflag	N
sar	0.08194738917667
contentqualification	PRODUCT
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
samplesperpixel	1
photometricinterpretation	MONOCHROME2
rows	68
columns	68
bitsallocated	16
bitsstored	12
highbit	11
lossyimagecompression	00
echonumber	0
pixelbandwidth	2451