






Phantom QA results for / fbirn_out/ ses64ch_funcbold_acqcmrrstd_taskwar slice 34




old acqcmrrstd taskwarmup 2.xml [percent fluct (trend removed), drift, driftfit]


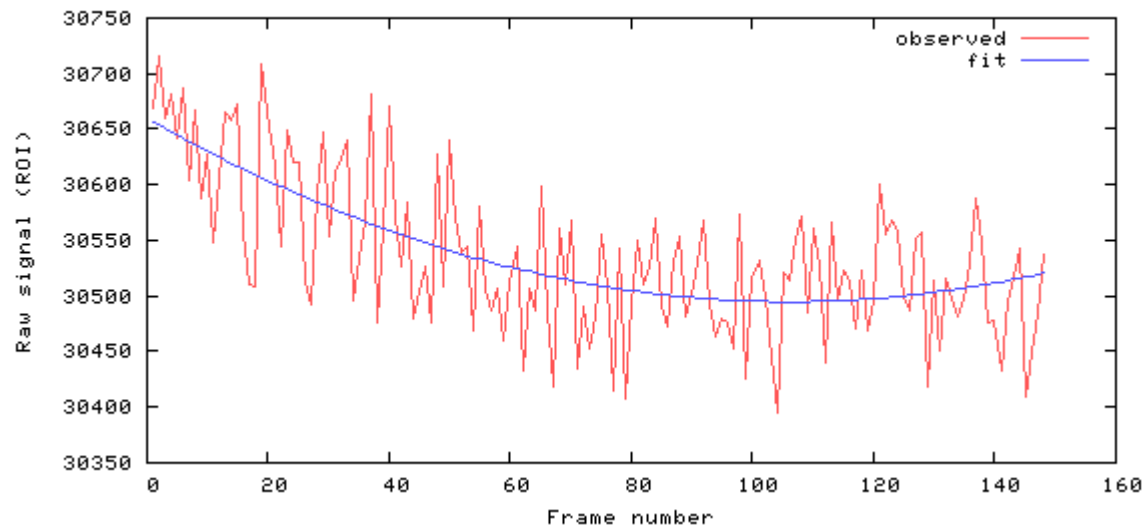


[mean, SNR, SFNR] = [30536.3 85.1 93.7]

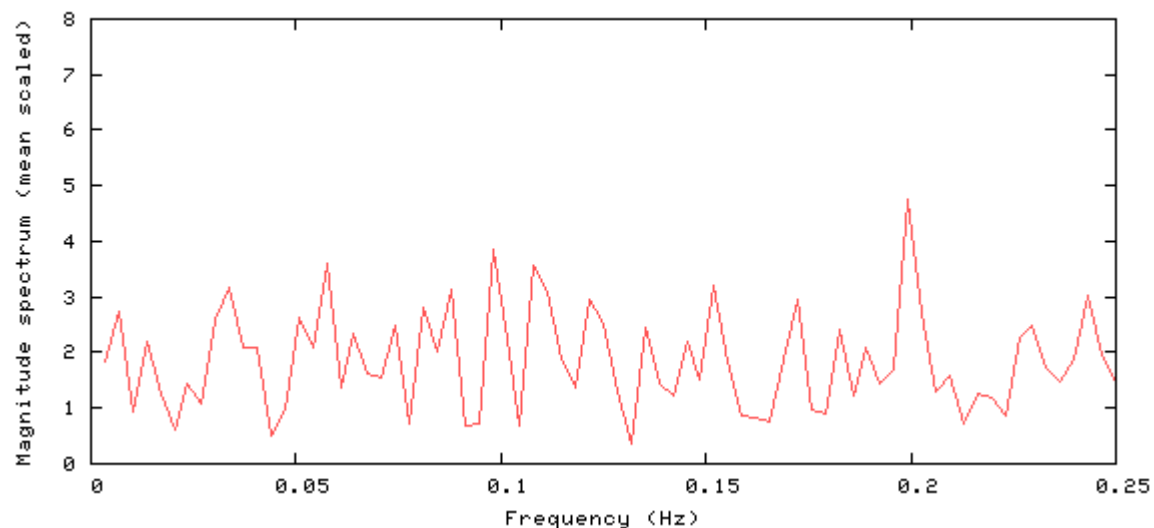

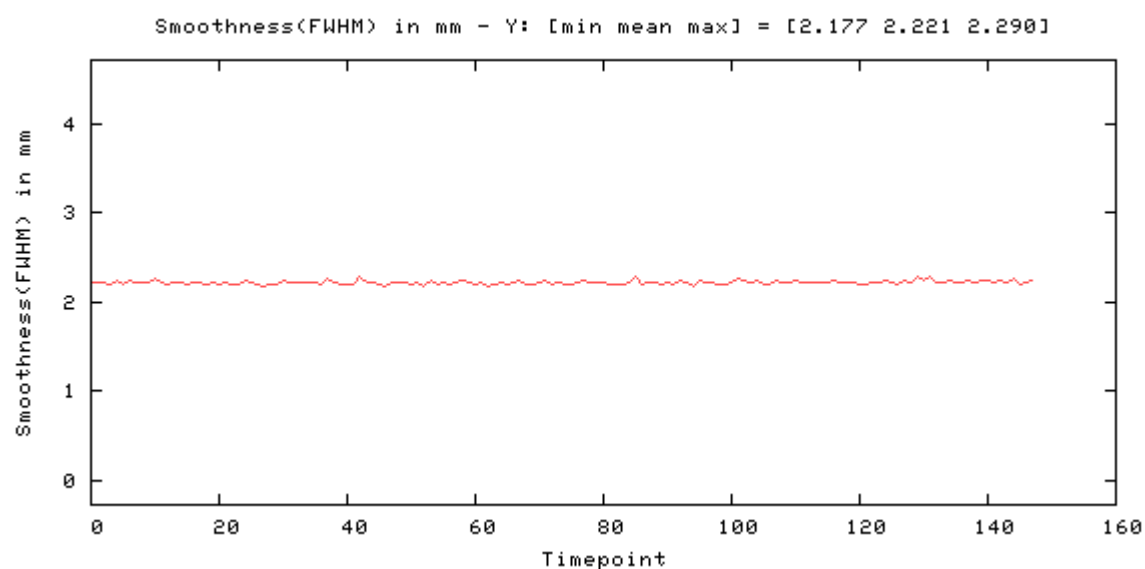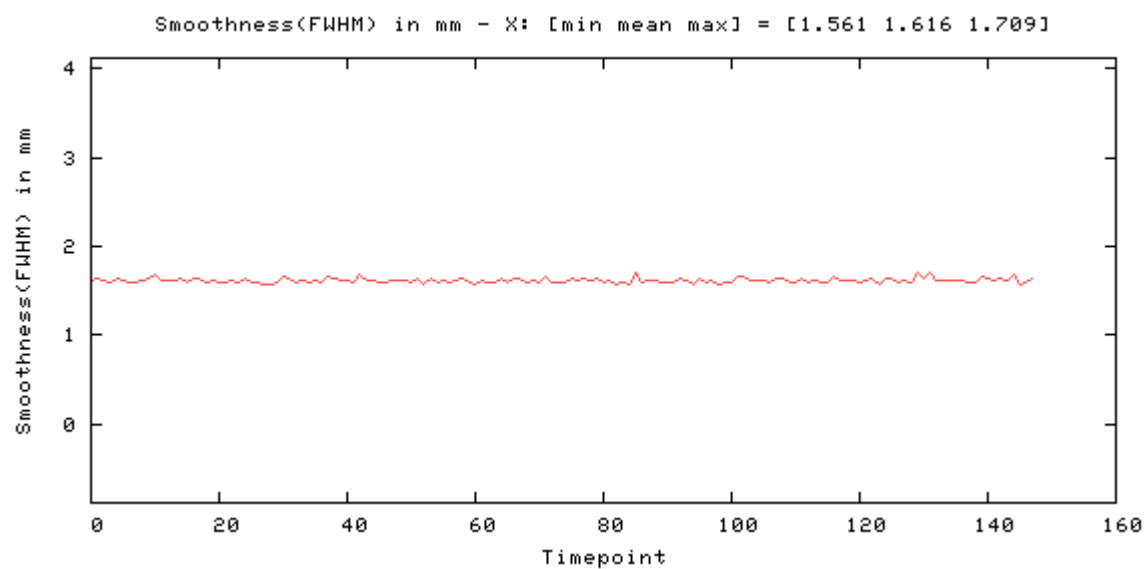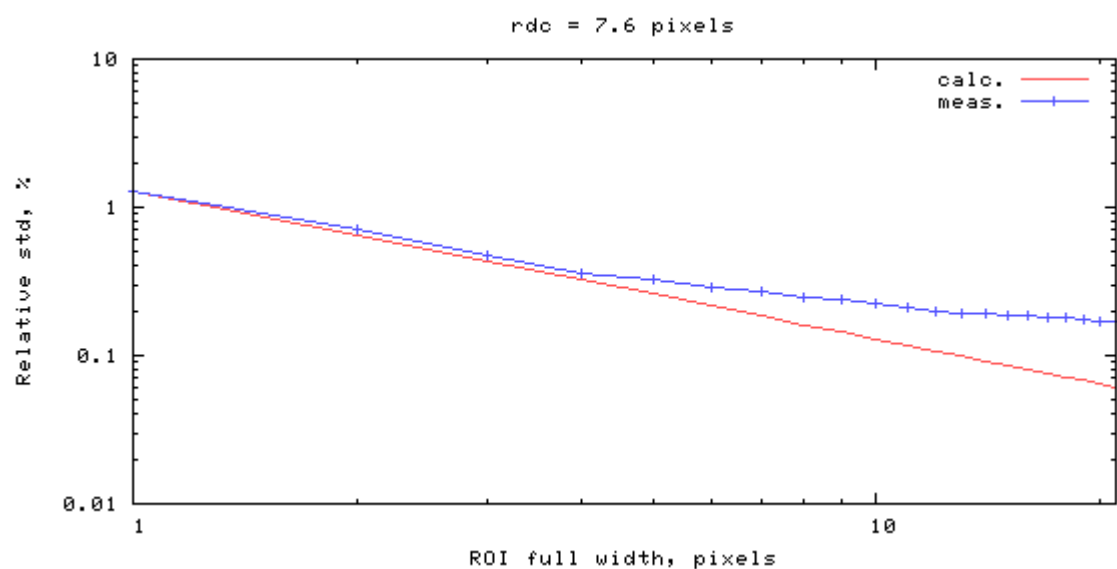

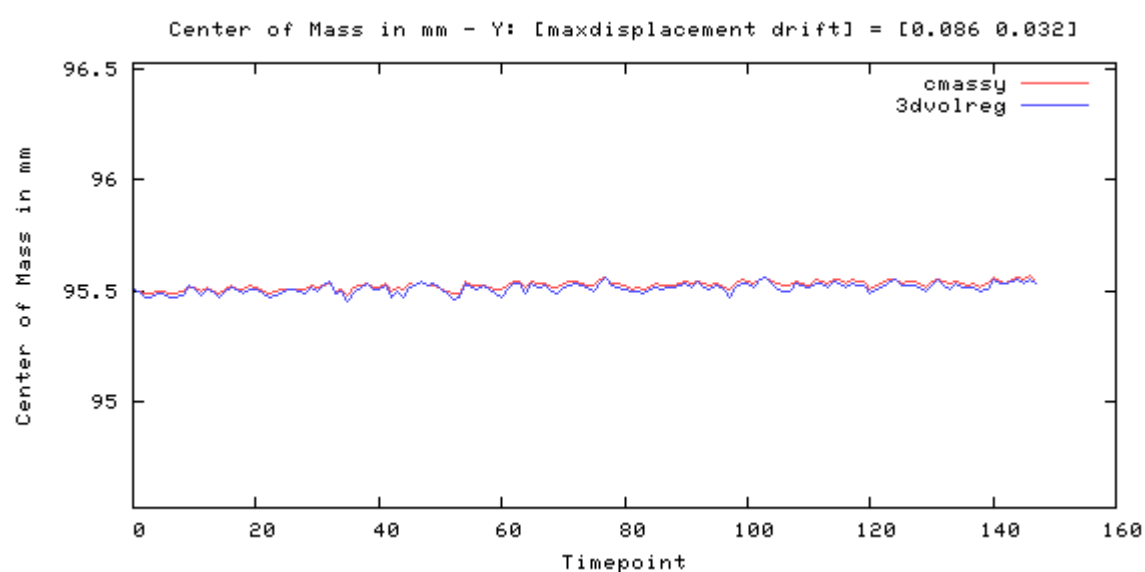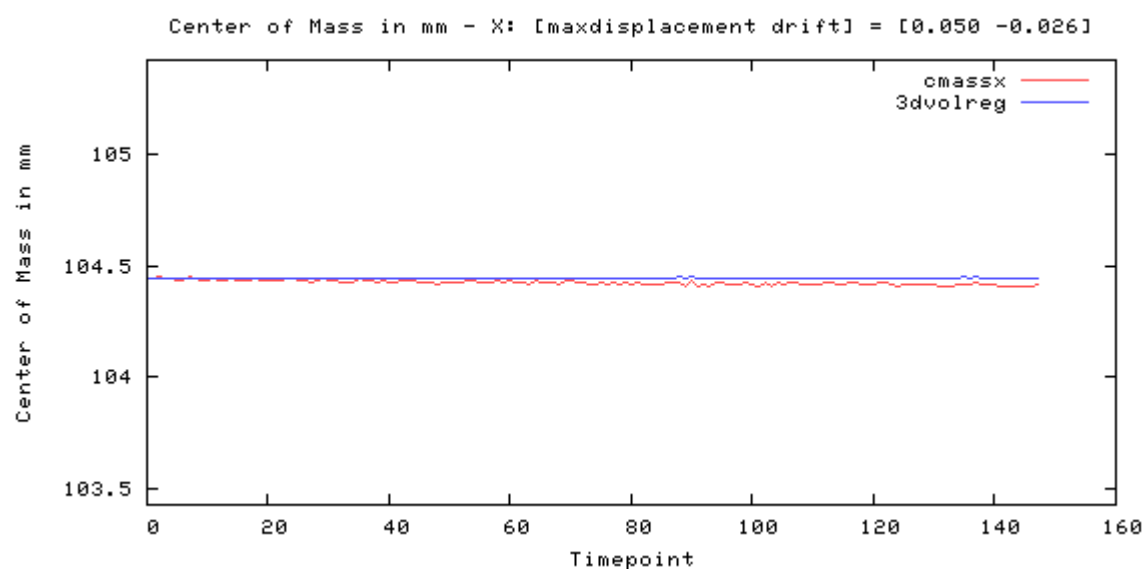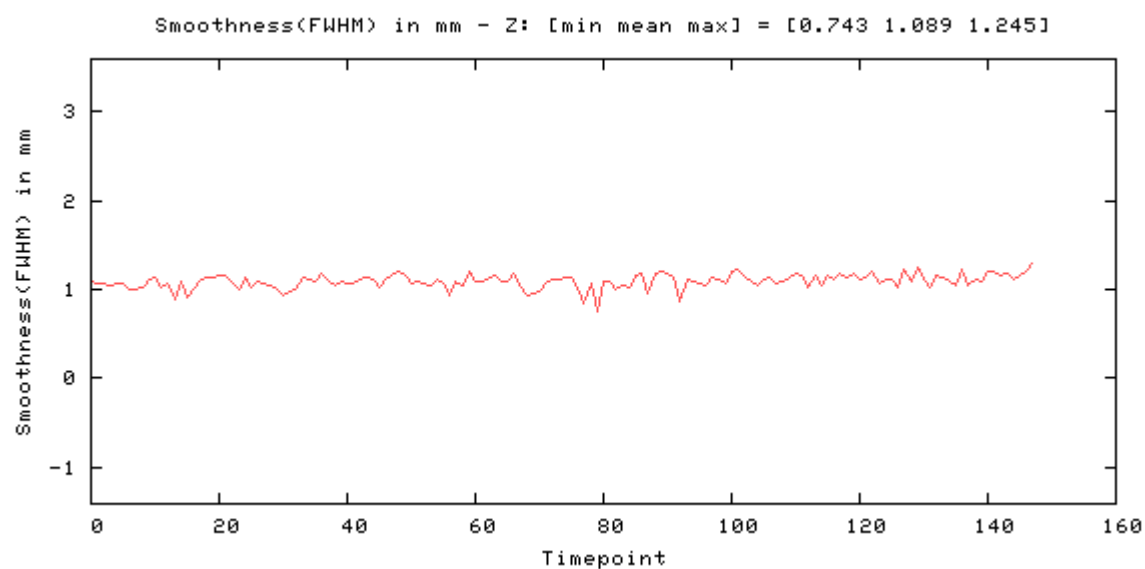

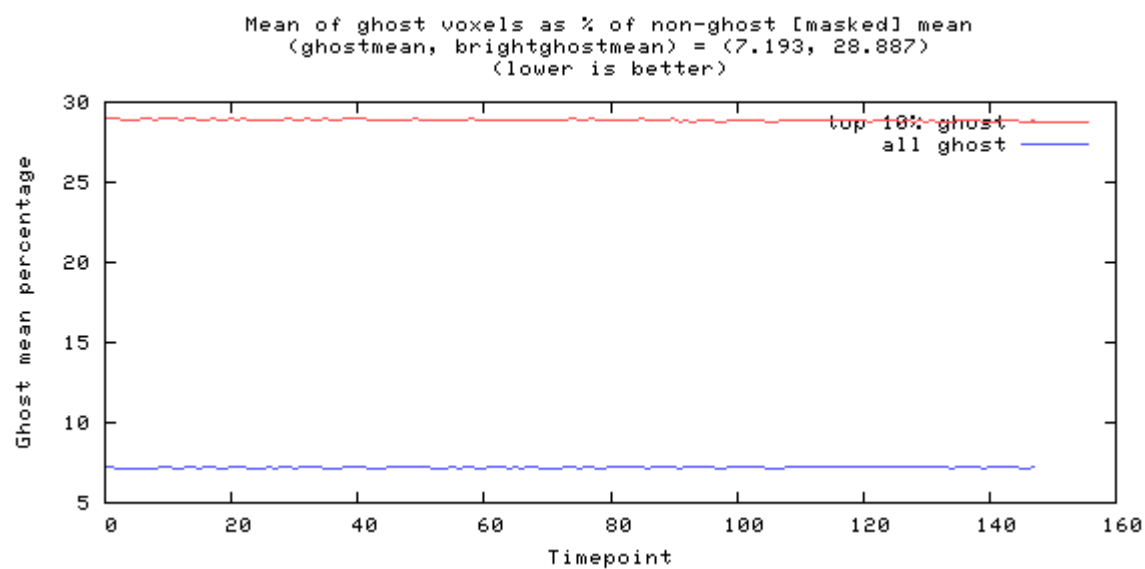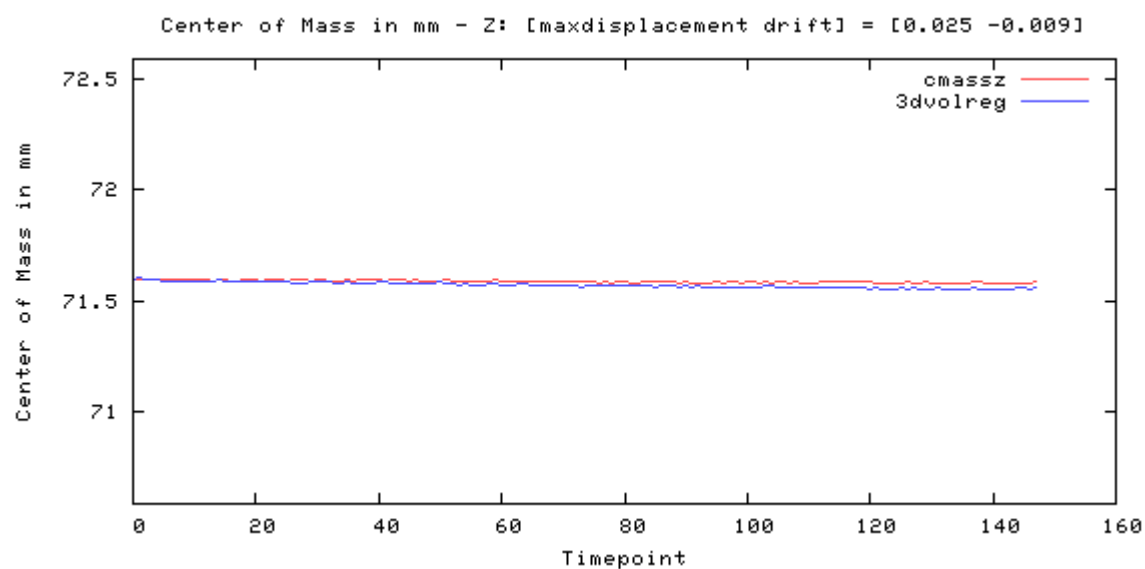


Odd-even difference image

-20408 23715




image min: -20408, image max: 23715

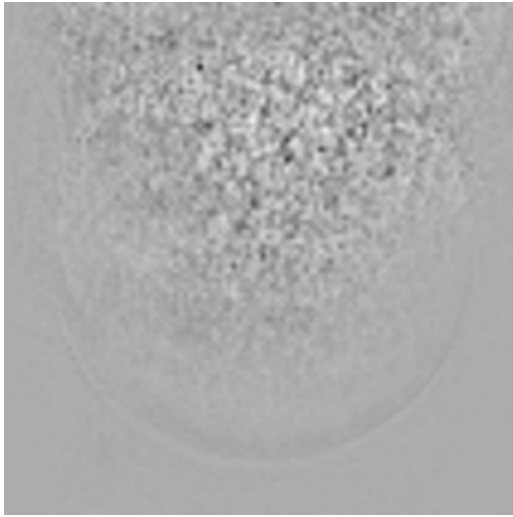

Mean image




image min: 0, image max: 60539.3

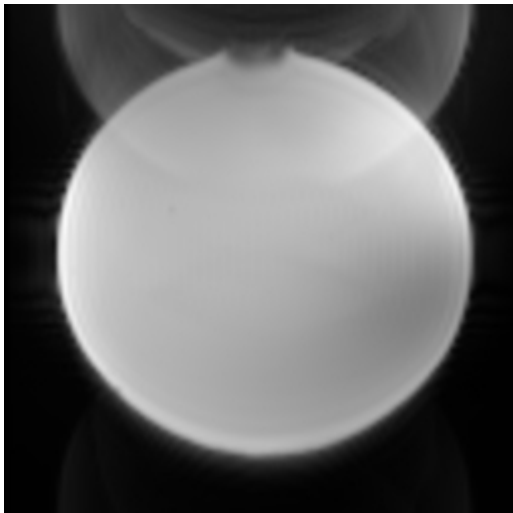

Standard Deviation image




image min: 0, image max: 699.82

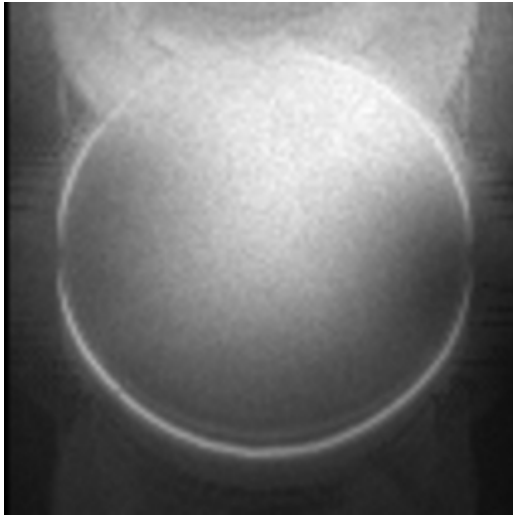

SFNR image



image min: 0, image max: 648.169

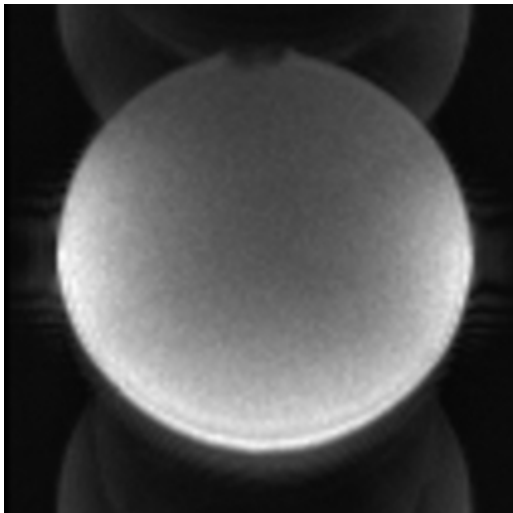

Acquisition parameters

dimensions	106x106x69x150
spacing	2mmx2mmx2mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	18e745c6-8069-43
studyid	Stability_20260622
seriesnumber	2
runnumber	1
magneticfield	3

description	ses-64ch_func-bold_acq-cmrrstd_task-warmup
scandate	2026-06-22
scantime	10:54:11
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	2
frequencydirection	1
receivecoilname	HeadNeck_64_CS
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-64ch_func-bold_acq-cmrrstd_task-warmup
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260622105409.210000
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	2
echotrainlength	80
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	RESEARCH
acquisitionduration	317.89499999999998

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	150
imagecomments	Not for diagnostic use, Unaliased MB3/ PE3/LB
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	69
rows	106
columns	106
bitsallocated	16
bitsstored	16
highbit	15
lossyimagecompression	00