






Phantom QA results for /fbirn_out/ses-32ch_func-bold_acq-siemensstd_task-tsnr.xml, slice 17




old acq-siemensstd task-tsnr.xml [percent fluct (trend removed), drift, driftfit]


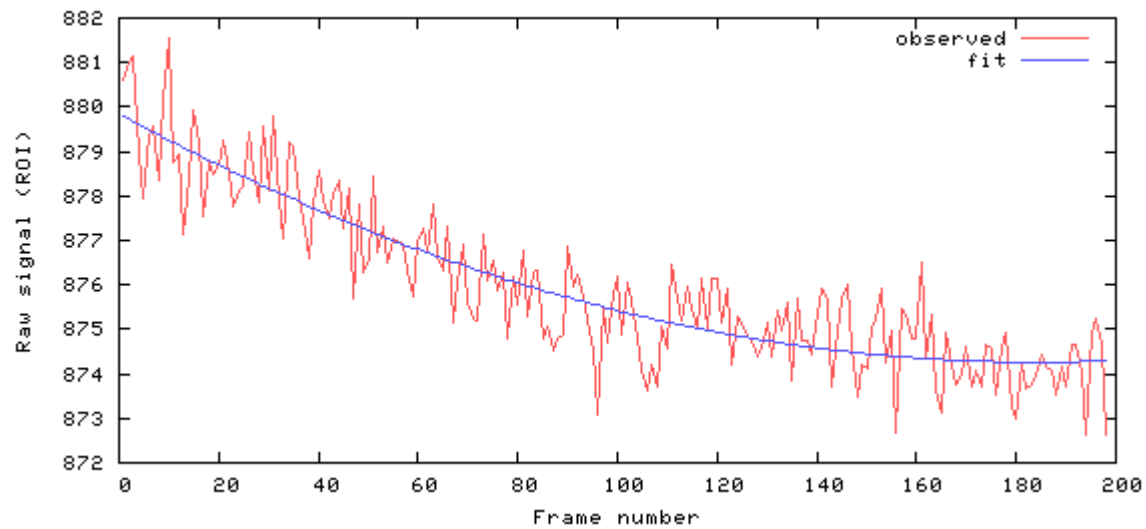


[mean, SNR, SFNR] = [876.0 119.0 124.2]

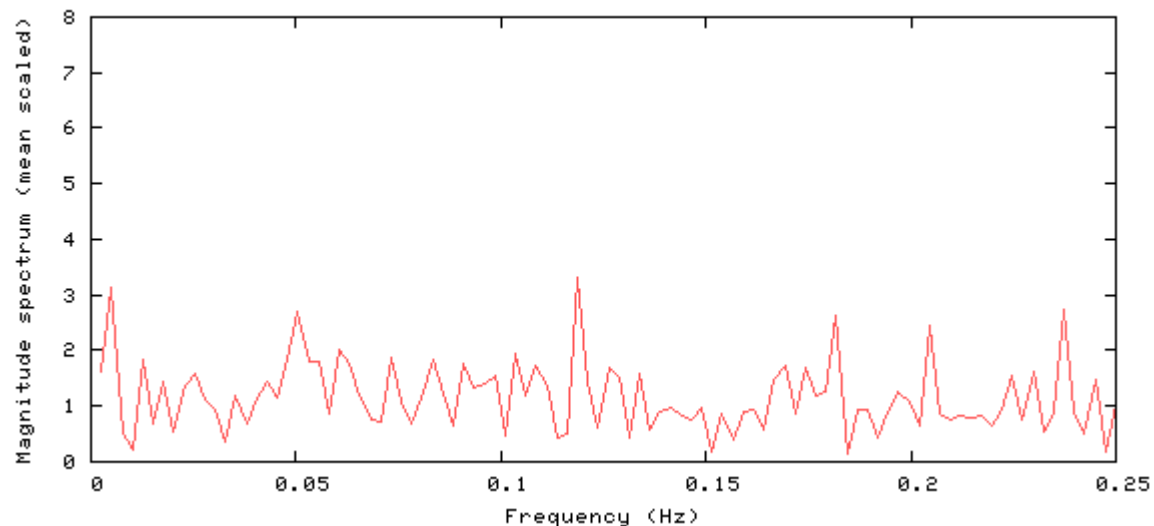

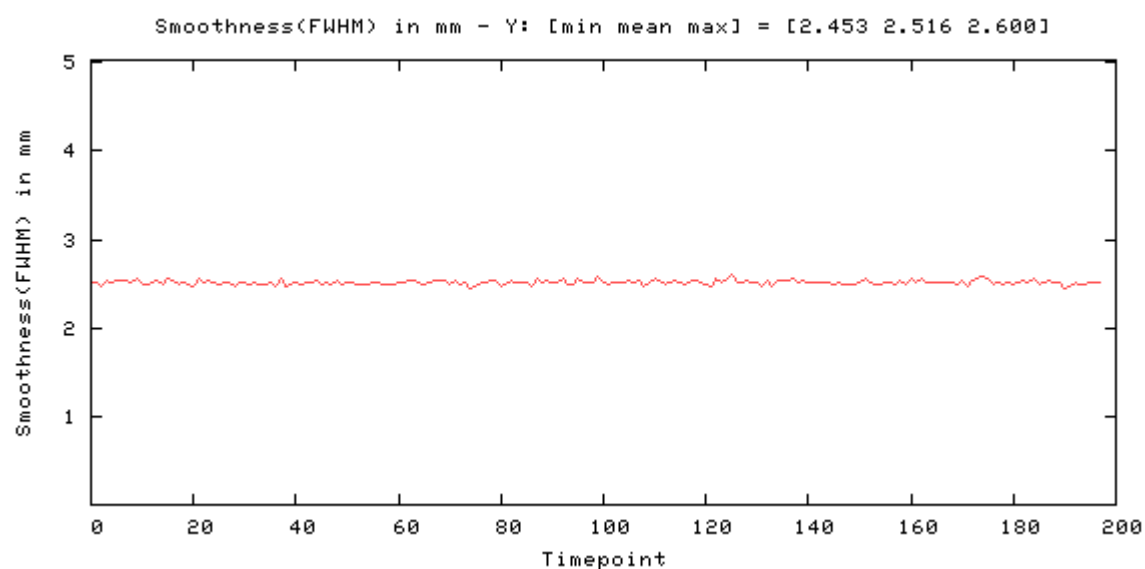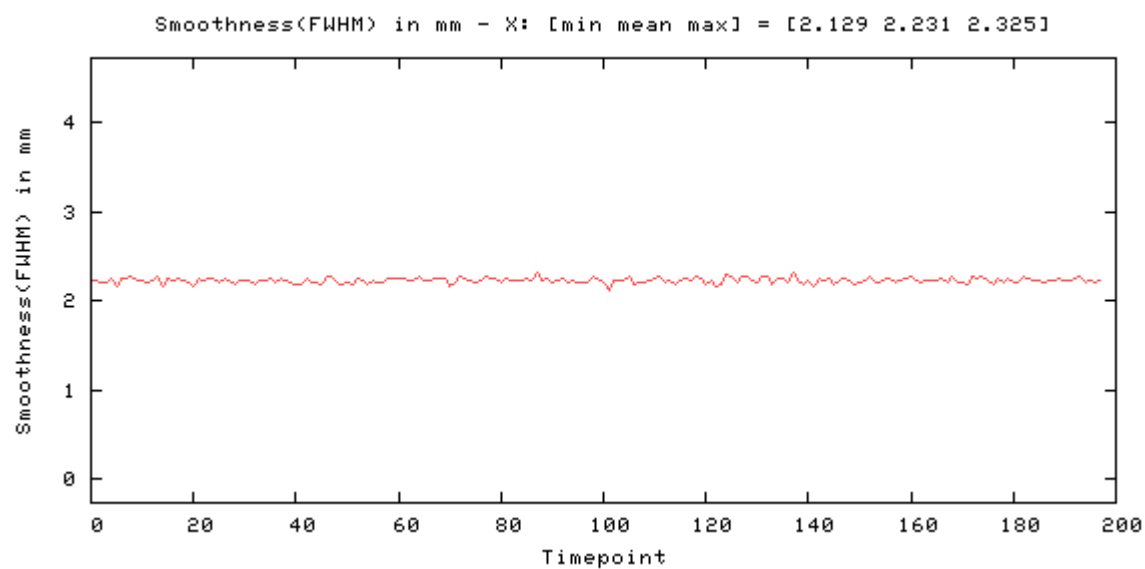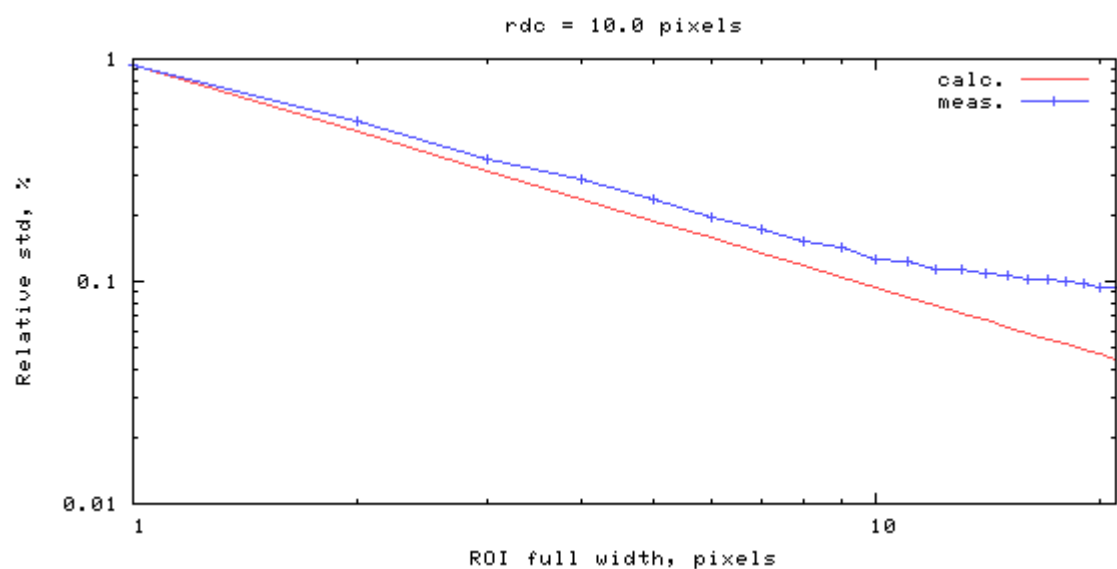

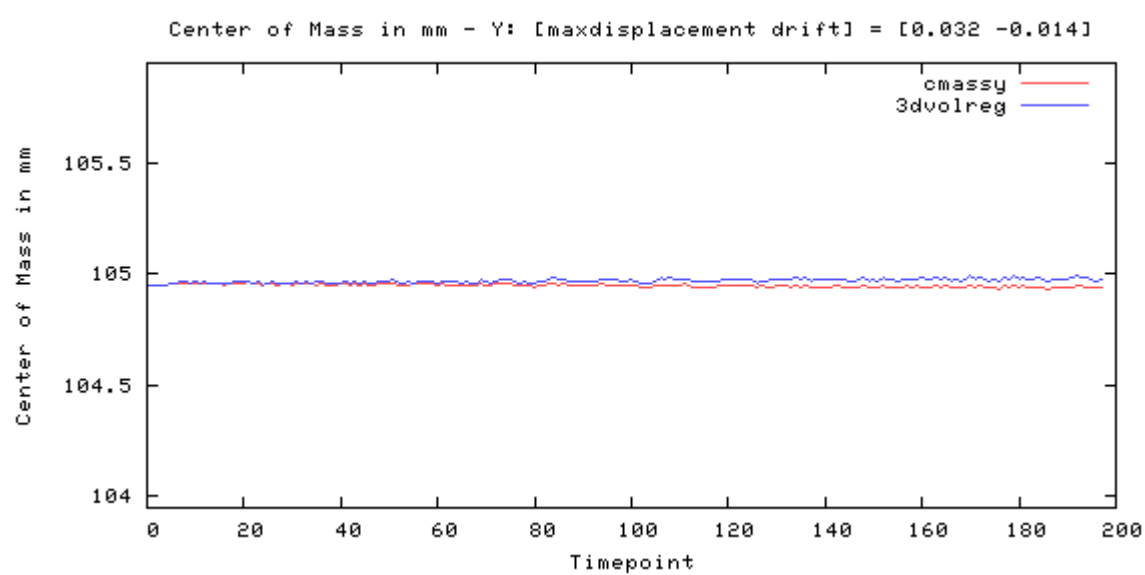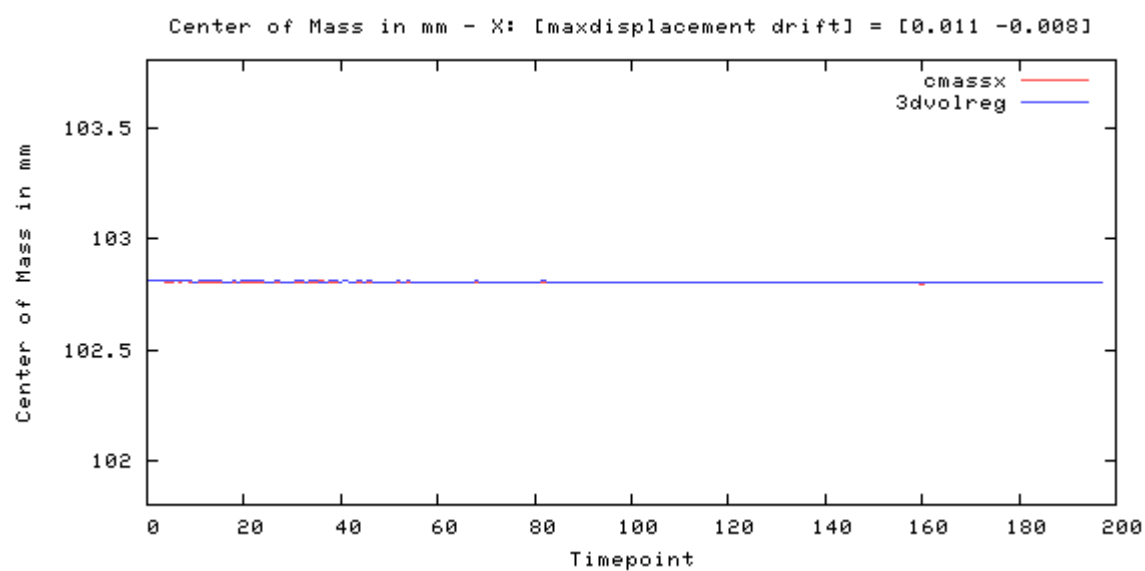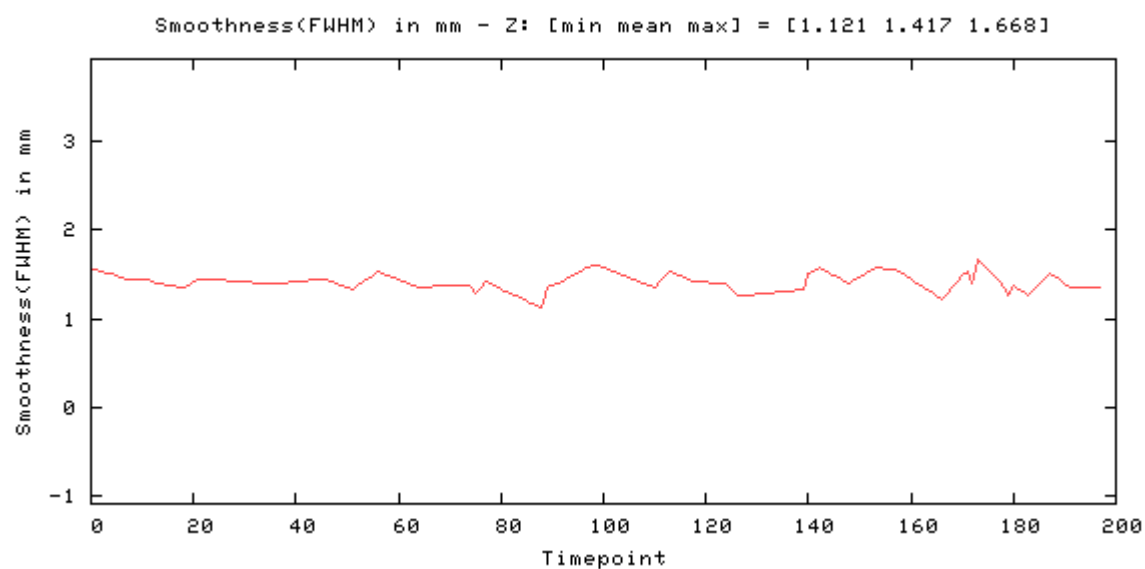


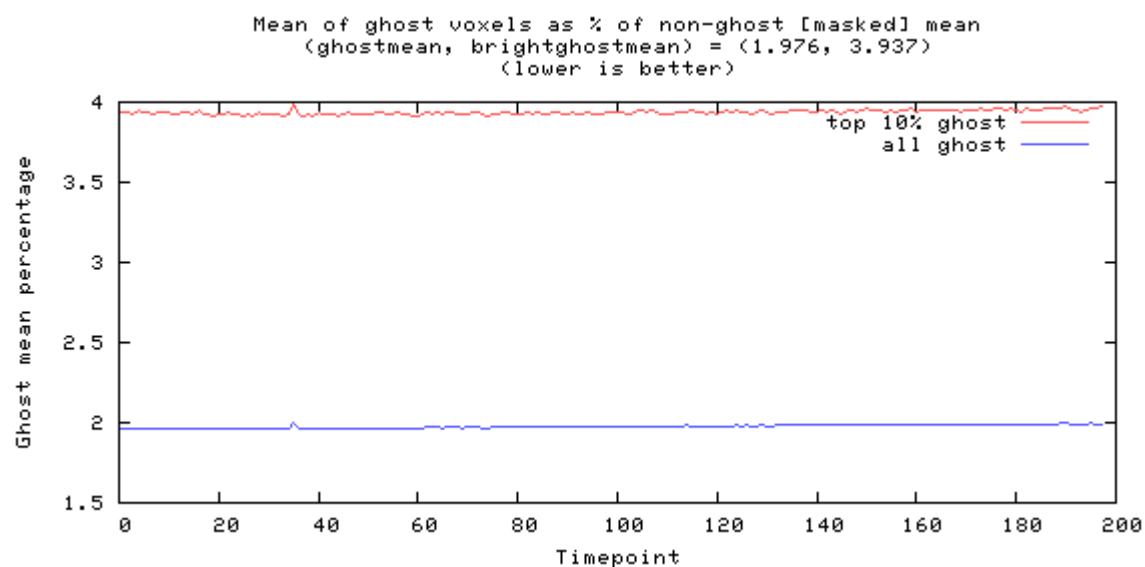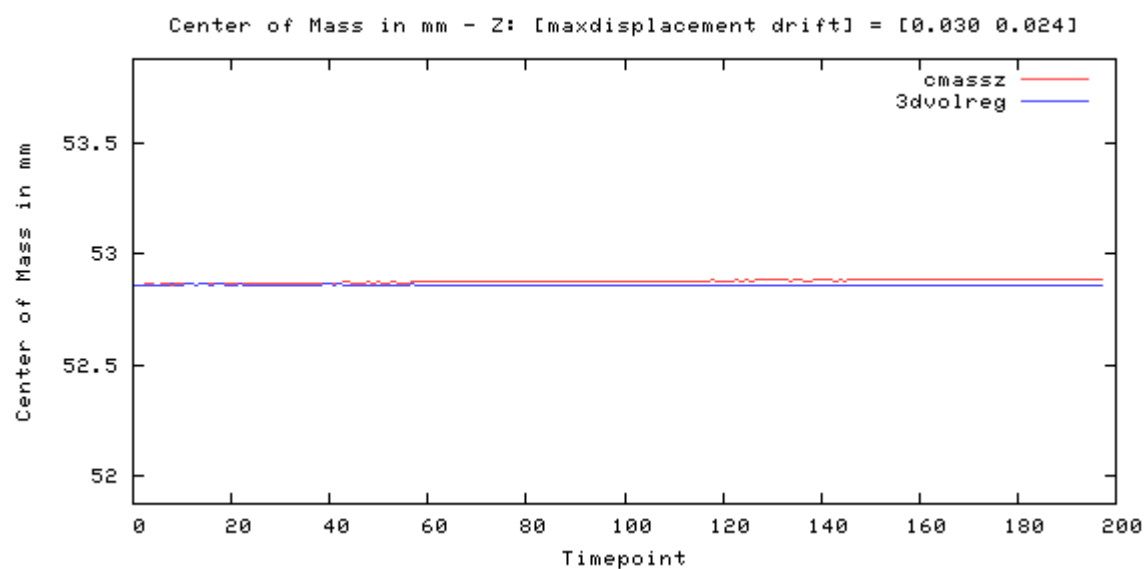

Odd-even difference image



image min: -337, image max: 411

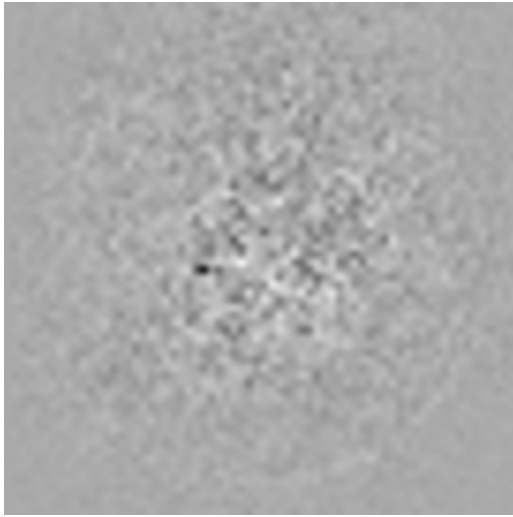

Mean image




image min: 0, image max: 1212.54



Standard Deviation image




image min: 0, image max: 9.36754

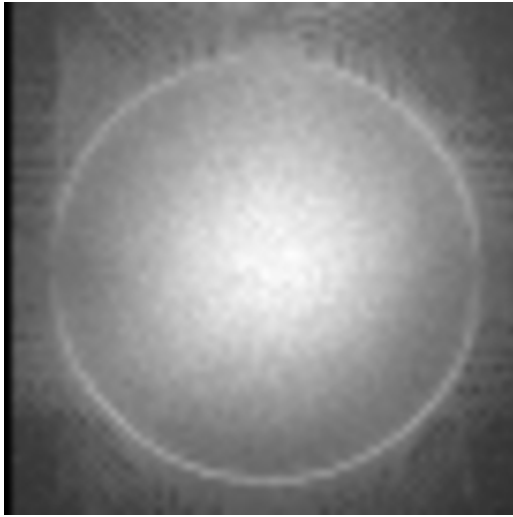

SFNR image



image min: 0, image max: 605.099

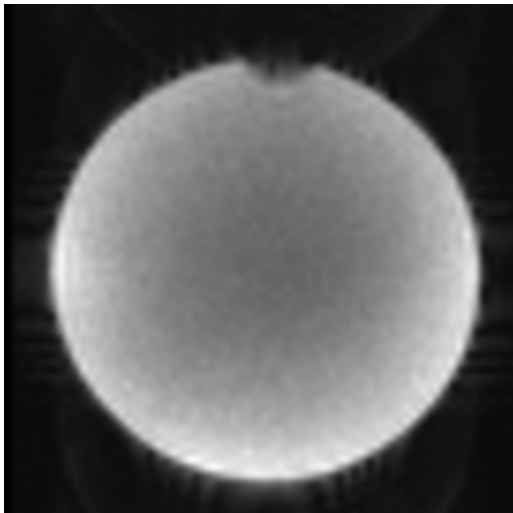

Acquisition parameters

dimensions	68x68x36x200
spacing	3mmx3mmx3mmx2000ms
gap	0mmx0mmx0mmx0ms
scanner	AWP213018
psdname	
examnumber	883e1130-f719-43
studyid	Stability_20260817
seriesnumber	16
runnumber	1
magneticfield	3

description	ses-32ch_func-bold_acq-siemensstd_task-tsnr
scandate	2026-08-17
scantime	09:57:20
tr	2000
te	30
operator	
flipangle	79
prescribedslicespacing	3
frequencydirection	1
receivecoilname	Head_32
transmitcoilname	Body
institution	THE UNIVERSITY OF PENNSYLVANIA
examdescription	ACTIVE KIRWAN
scanner	AWP213018
scannermanufacturer	Siemens Healthineers
scannermodelname	MAGNETOM Cima.X
protocolname	ses-32ch_func-bold_acq-siemensstd_task-tsnr
institutionaddress	Grays Ferry Ave. BLDG 250 STE 140 3401,Philadelphia,69 - KEYSTONE,US,19146
institutionaldepartmentname	MindCORE Neuroimaging Suite 140
acquisitiondatetime	20260817095718.272500
frametype	ORIGINAL
pixelpresentation	MONOCHROME
volumetricproperties	DISTORTED
volumebasedcalculationtechnique	NONE
compleximagecomponent	MAGNITUDE
acquisitioncontrast	UNKNOWN
mracquisitiontype	2D
numaverages	1
spacingbetweenslices	3
echotrainlength	33
percentsampling	100
percentphasefieldofview	100
scannerserialnumber	213018
softwareversions	syngo MR XA61
inplanephaseencodingdirection	COLUMN
variableflipangleflag	N
contentqualification	PRODUCT
acquisitionduration	407.92500000000001

resonantnucleus	1H
applicablesafetystandardagency	IEC
acquisitionnumber	1
instancenumber	1
numtemporalpositions	200
samplesperpixel	1
photometricinterpretation	MONOCHROME2
numberofframes	36
rows	68
columns	68
bitsallocated	16
bitsstored	12
highbit	11
lossyimagecompression	00