






# Phantom QA results for / fbirn\_out/ ses32ch\_funcbold\_acqsiemensbasic\_ta slice 14




.d acqsiemensbasic tasktsnr 17.xml [percent fluct (trend removed), drift, driftfit


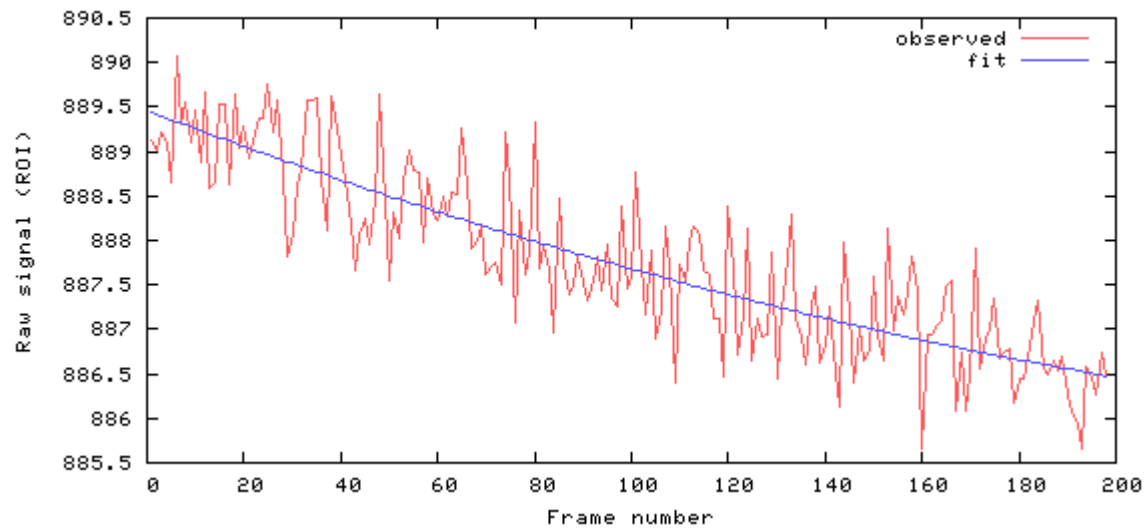


[mean, SNR, SFNR] = [887.8 184.1 183.6]

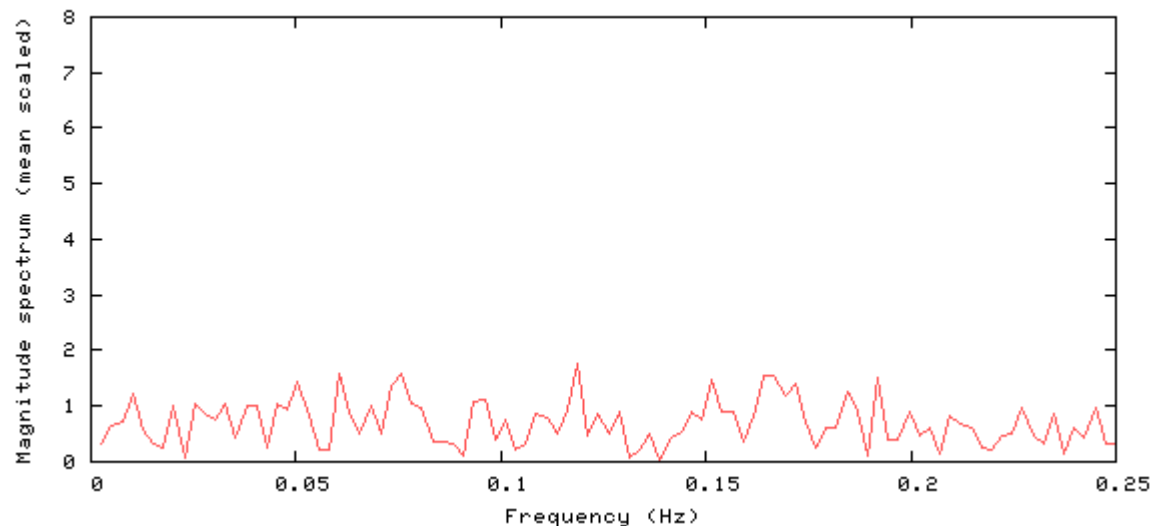

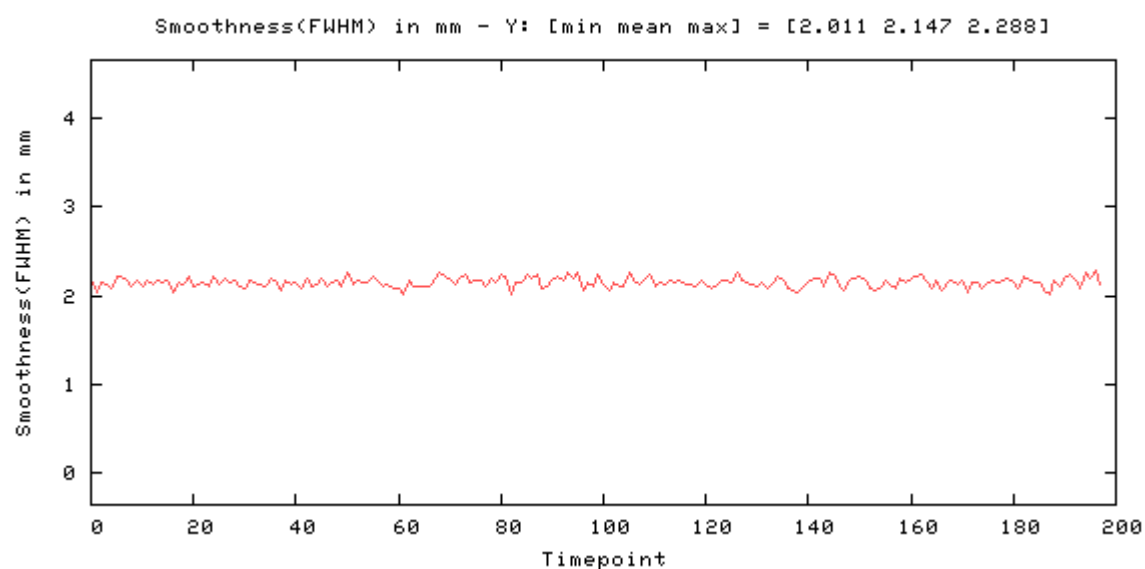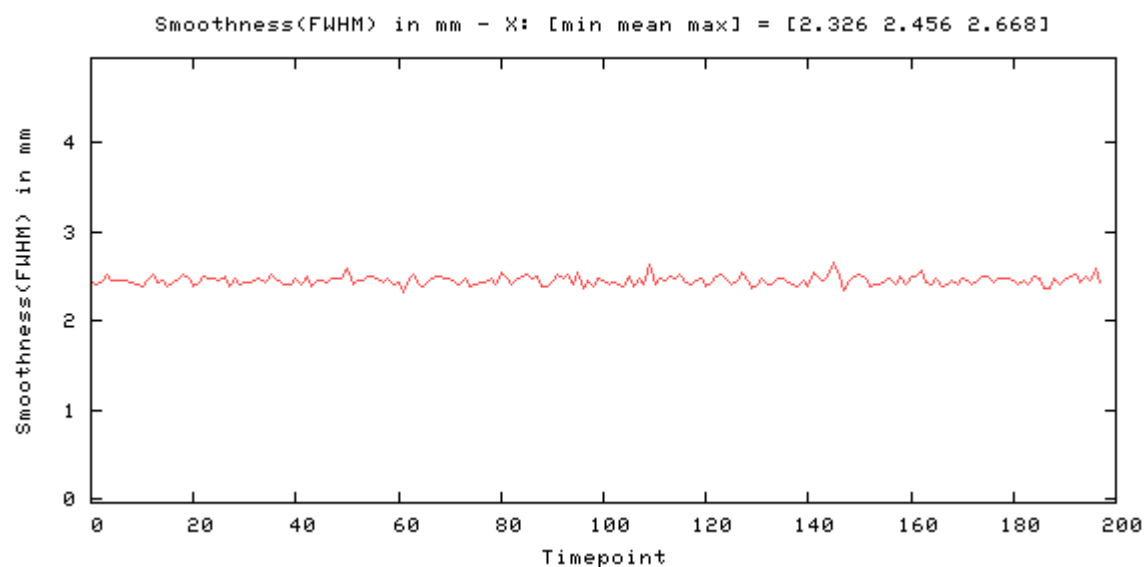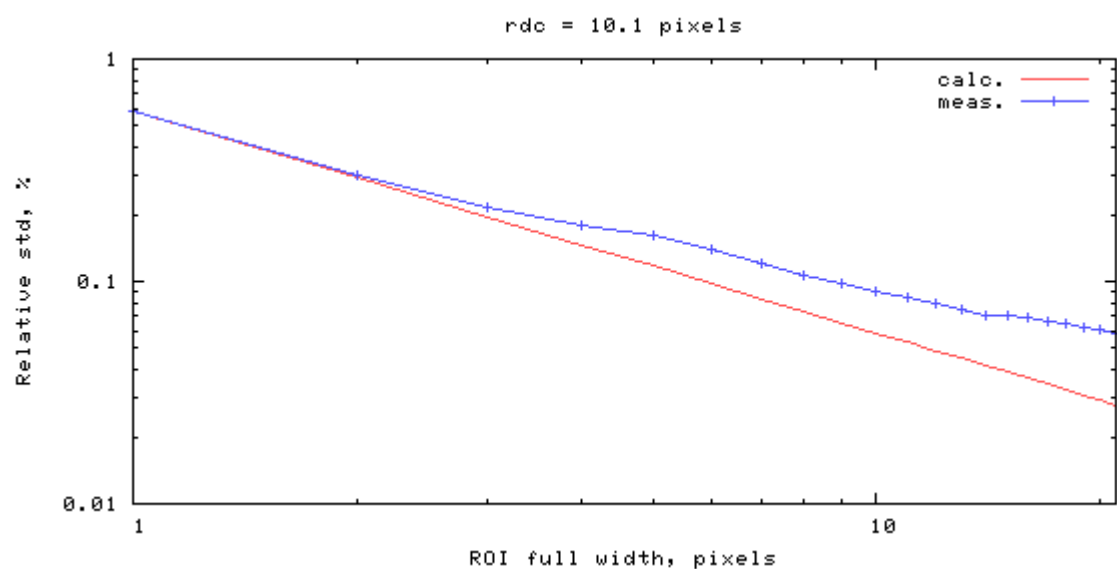

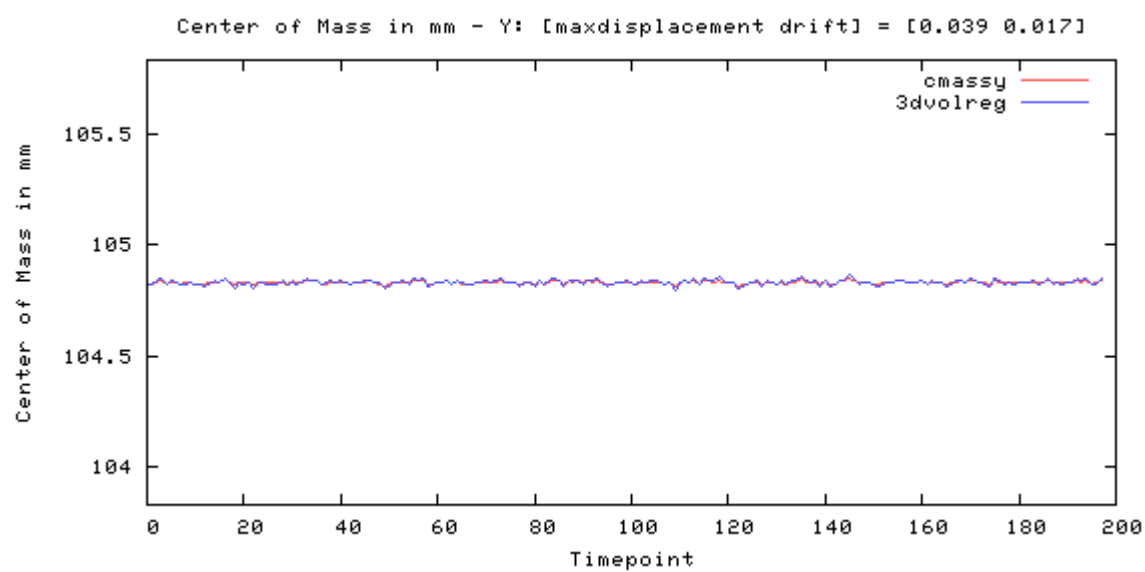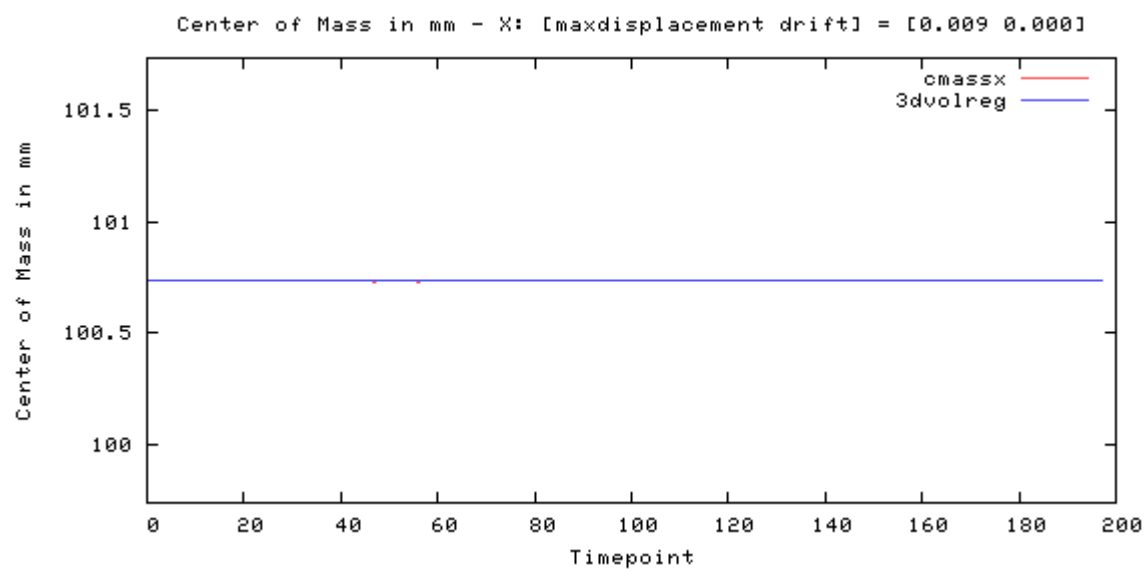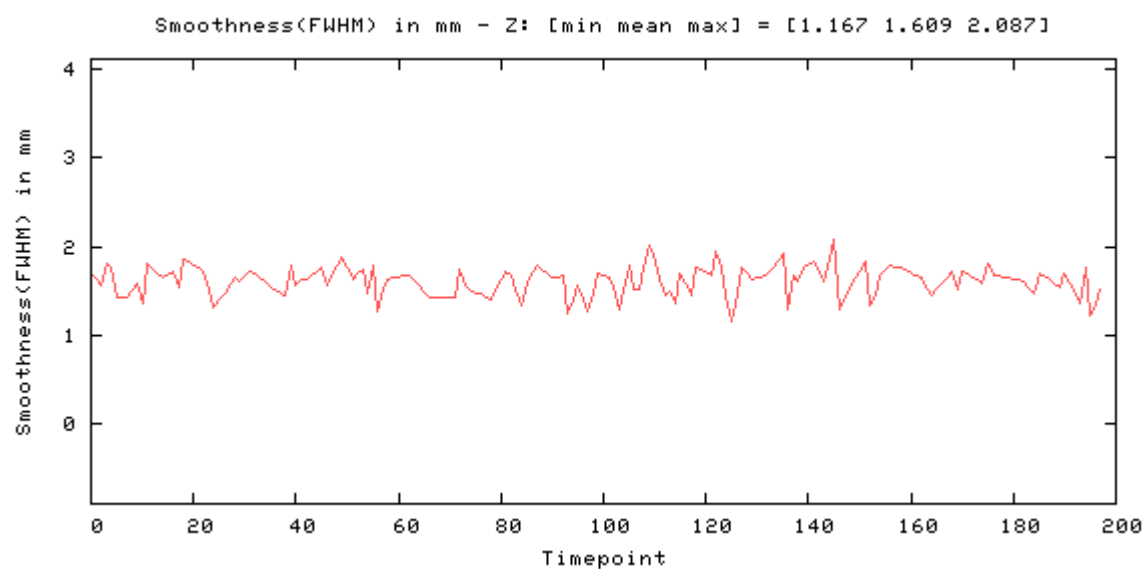

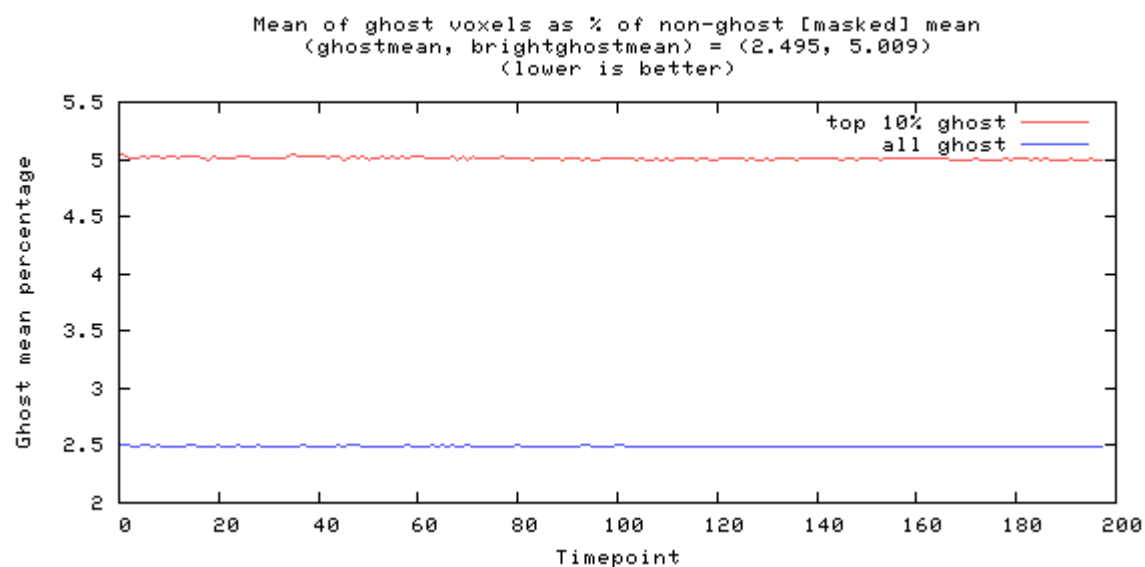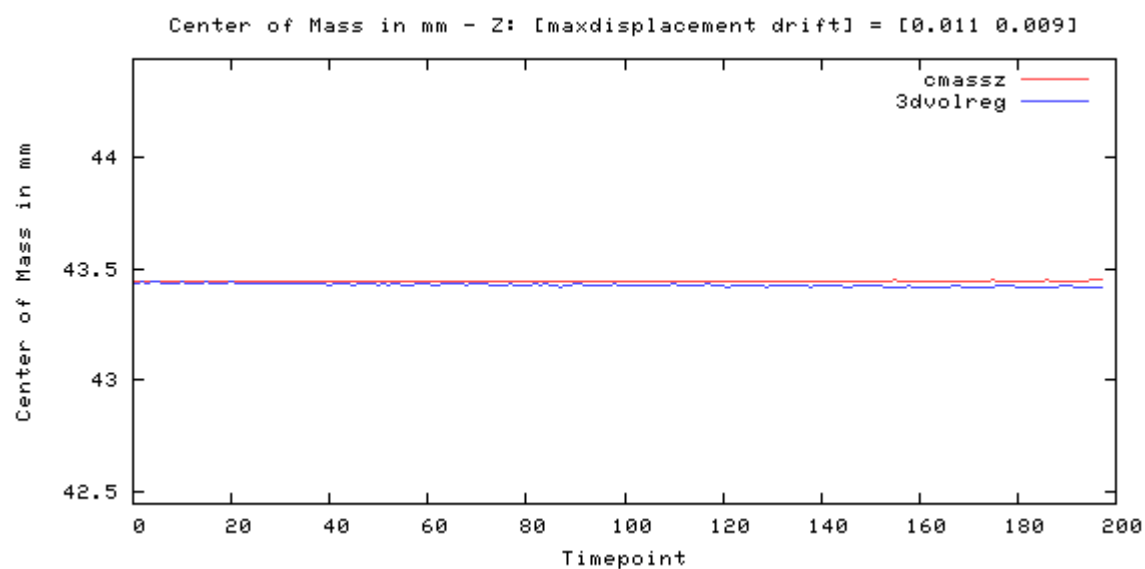


---

## Odd-even difference image



image min: -194, image max: 176

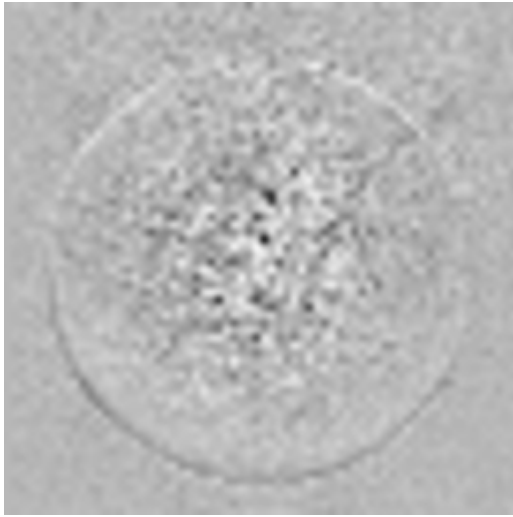

---

### Mean image



image min: 0, image max: 1264.32




---

### Standard Deviation image



image min: 0, image max: 6.45367

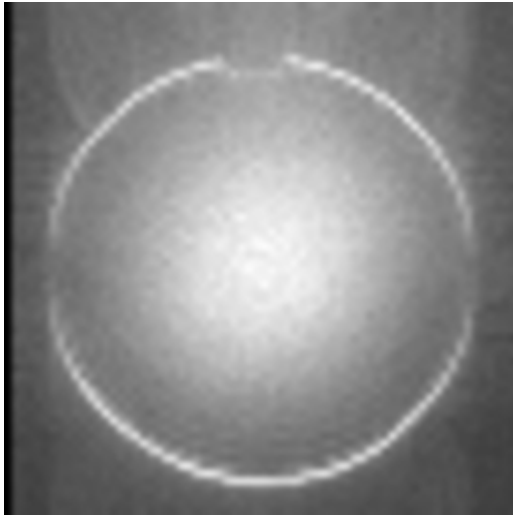


---

### SFNR image



image min: 0, image max: 830.048

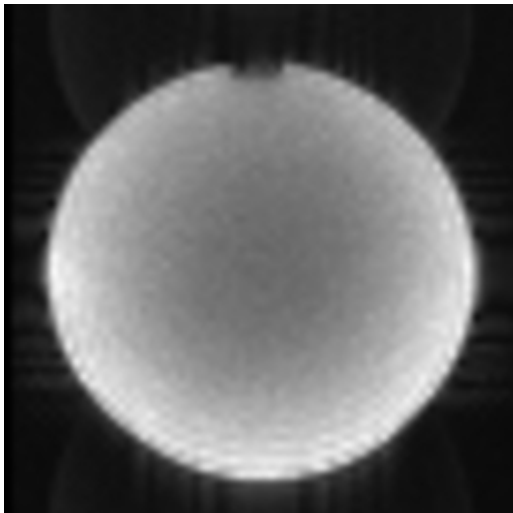

---

### Acquisition parameters

|               |                    |
|---------------|--------------------|
| dimensions    | 68x68x30x200       |
| spacing       | 3mmx3mmx3mmx2000ms |
| gap           | 0mmx0mmx0mmx0ms    |
| scanner       | AWP213018          |
| psdname       |                    |
| examnumber    | a821de8a-18fb-43   |
| studyid       | Stability_20260518 |
| seriesnumber  | 17                 |
| runnumber     | 1                  |
| magneticfield | 3                  |

|                                 |                                                                                  |
|---------------------------------|----------------------------------------------------------------------------------|
| description                     | ses-32ch_func-bold_acq-siemensbasic_task-tsnr                                    |
| scandate                        | 2026-05-18                                                                       |
| scantime                        | 10:19:44                                                                         |
| tr                              | 2000                                                                             |
| te                              | 30                                                                               |
| operator                        |                                                                                  |
| flipangle                       | 79                                                                               |
| prescribedslicespacing          | 3                                                                                |
| frequencydirection              | 1                                                                                |
| receivecoilname                 | Head_32                                                                          |
| transmitcoilname                | Body                                                                             |
| institution                     | THE UNIVERSITY OF PENNSYLVANIA                                                   |
| examdescription                 | ACTIVE KIRWAN                                                                    |
| scanner                         | AWP213018                                                                        |
| scannermanufacturer             | Siemens Healthineers                                                             |
| scannermodelname                | MAGNETOM Cima.X                                                                  |
| protocolname                    | ses-32ch_func-bold_acq-siemensbasic_task-tsnr                                    |
| institutionaddress              | Grays Ferry Ave. BLDG 250 STE 140<br>3401,Philadelphia,69 -<br>KEYSTONE,US,19146 |
| institutionaldepartmentname     | MindCORE Neuroimaging Suite 140                                                  |
| acquisitiondatetime             | 20260518101942.160000                                                            |
| frametype                       | ORIGINAL                                                                         |
| pixelpresentation               | MONOCHROME                                                                       |
| volumetricproperties            | DISTORTED                                                                        |
| volumebasedcalculationtechnique | NONE                                                                             |
| compleximagecomponent           | MAGNITUDE                                                                        |
| acquisitioncontrast             | UNKNOWN                                                                          |
| mracquisitiontype               | 2D                                                                               |
| numaverages                     | 1                                                                                |
| spacingbetweenslices            | 3                                                                                |
| echotrainlength                 | 68                                                                               |
| percentsampling                 | 100                                                                              |
| percentphasefieldofview         | 100                                                                              |
| scannerserialnumber             | 213018                                                                           |
| softwareversions                | syngo MR XA61                                                                    |
| inplanephaseencodingdirection   | COLUMN                                                                           |
| variableflipangleflag           | N                                                                                |
| contentqualification            | PRODUCT                                                                          |
| acquisitionduration             | 405.91200000000003                                                               |

|                                |             |
|--------------------------------|-------------|
| resonantnucleus                | 1H          |
| applicablesafetystandardagency | IEC         |
| acquisitionnumber              | 1           |
| instancenumber                 | 1           |
| numtemporalpositions           | 200         |
| samplesperpixel                | 1           |
| photometricinterpretation      | MONOCHROME2 |
| numberofframes                 | 30          |
| rows                           | 68          |
| columns                        | 68          |
| bitsallocated                  | 16          |
| bitsstored                     | 12          |
| highbit                        | 11          |
| lossyimagecompression          | 00          |